

KNQZE KF9

KNQZE Foot Massager with Heat and Remote (Model KF9) - Instruction Manual

Welcome to your new KNQZE Foot Massager. Please read this manual carefully before use to ensure safe and optimal operation.

IMPORTANT SAFETY INFORMATION

WARNING: May cause burns if used for extended periods at high heat settings. Avoid use if you have a blood clotting disorder, deep vein thrombosis, or are pregnant. Do not use on broken skin or open wounds. Keep away from water. Unplug the massager when not in use. Always follow the instructions in the user manual.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

PACKAGE CONTENTS

- KNQZE Foot Massager Machine
- Remote Control
- UL Listed Power Adapter
- User Manual



Image: The KNQZE Foot Massager Machine, power adapter, remote control, and user manual laid out.

Video: An unboxing video demonstrating the contents of the KNQZE Foot Massager package, including the massager unit, power adapter, remote, and user manual.

PRODUCT FEATURES


- **3D Kneading & Shiatsu Rollers:** Provides deep tissue massage for forefoot, arch, and heel.
- **360° Full-Wrapped Air Compression:** Adjustable intensities (OFF/L/M/H) for instep and ankle.
- **3 Heat Levels:** Select from 113°F, 122°F, or 131°F for soothing warmth.
- **3 Timer Settings:** Choose massage durations of 15, 20, or 30 minutes with auto shut-off.
- **Dual Control Options:** Features an intuitive LCD touchscreen and a convenient wireless remote control.
- **Removable & Washable Foot Covers:** Ensures easy maintenance and hygiene.
- **Compact & Lightweight Design:** Easy to carry and store for home or office use.

3 Level Heating Therapy



Image: A diagram illustrating the 3-section comprehensive massage mechanism, highlighting kneading, rolling, scraping, and pushing actions on the forefoot, arch, and heel.

3-Section Comprehensive Massage for Ultimate Relaxation

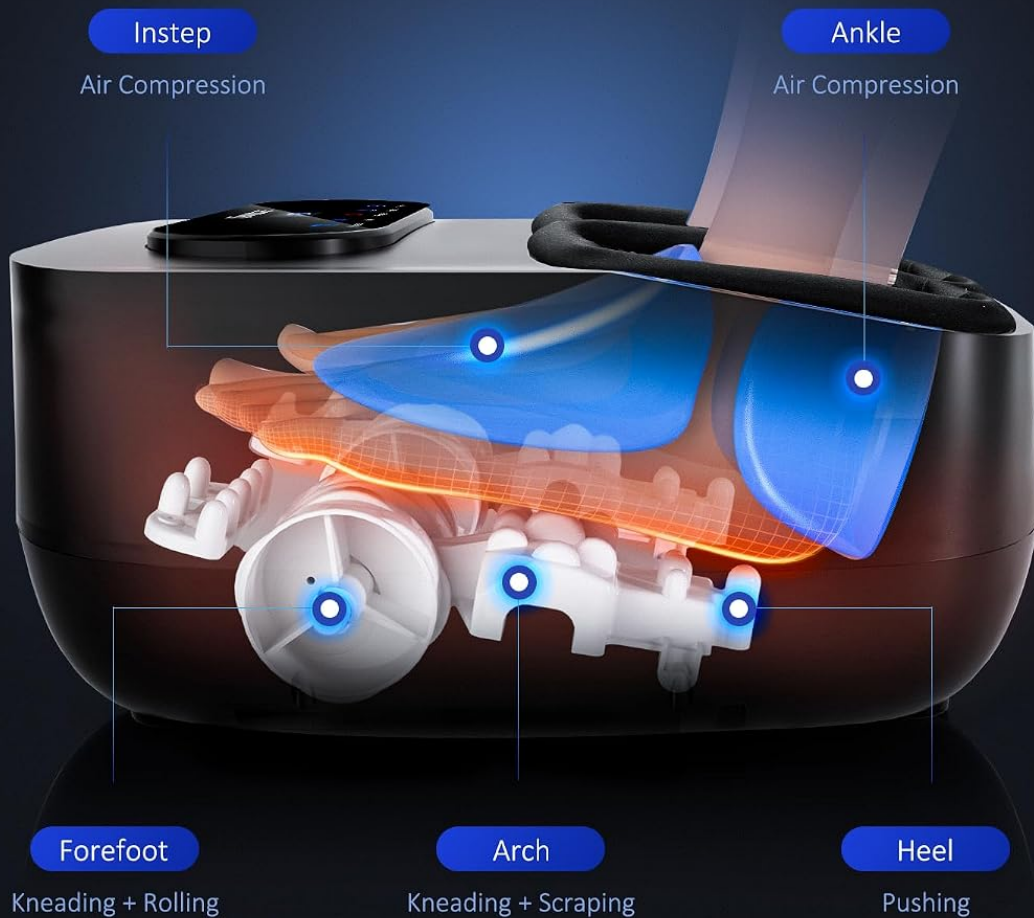

Kneading



Rolling



Acupressure



Scraping



Firm Pressure

Image: A woman using the foot massager, demonstrating the 3-level heating therapy with visual representations of improved circulation (blue to red feet).

2 Control Options

Larger LCD Touch Panel x Remote Control



Image: A detailed view of the foot massager's LCD touch panel and the remote control, showing buttons for Knead, Air, Power, Heat, and Timer.

SETUP

Unboxing

- Carefully remove all components from the product packaging.
- Inspect the foot massager and all accessories for any signs of damage.

Power Connection

- Connect the provided UL Listed Power Adapter to the power input port located on the back of the foot massager.
- Plug the power adapter into a standard electrical outlet.

Remote Control Preparation

- Before first use, remove the clear plastic battery tab from the remote control to activate the battery.

OPERATING INSTRUCTIONS

Placing Your Feet

- Sit comfortably in a chair or on a sofa.
- Place your feet into the foot chambers of the massager. Ensure your feet are centered for optimal massage contact.

Using the Control Panel (LCD Touchscreen)

- **Power Button:** Press to turn the massager ON. The red light will flash, and the green light will stay on. Press again to turn OFF all functions.
- **Knead Button:** Press to cycle through kneading intensities: Green light (low speed), Blue light (medium speed), Red light (high speed). Press again to turn off kneading.
- **Air Button:** Press to cycle through air compression intensities: Green light (low intensity), Blue light (medium intensity), Red light (high intensity). Press again to turn off air compression.
- **Heat Button:** Press to cycle through heat levels: Green light (113°F), Blue light (122°F), Red light (131°F). Press again to turn off heating.
- **Timer Button:** Press to select massage duration: Green light (15 minutes), Blue light (20 minutes), Red light (30 minutes). The massager will automatically shut off after the selected time.

Using the Remote Control

- The remote control provides convenient access to all functions available on the LCD touchscreen.
- Use the corresponding buttons on the remote to power on/off, adjust kneading, air compression, heat, and timer settings without bending over.

Video: A demonstration of the foot massager's LCD touch panel and remote control in action, showing how to adjust various settings.

MAINTENANCE

Removable Foot Covers

- To maintain hygiene, the foot covers are removable and washable.
- Unzip the foot covers from the massager.
- Machine wash the covers on a gentle cycle with mild detergent.
- Allow the covers to air dry completely before reattaching them to the massager.

Unit Cleaning

- Ensure the massager is unplugged from the power outlet before cleaning.
- Wipe the exterior of the massager with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the unit in water. The unit is made of water-resistant ABS material, but it is not waterproof.



Image: A close-up showing the detachable and washable foot cover being unzipped and removed from the foot massager for cleaning.

TROUBLESHOOTING

- **Unit not powering on:** Check if the power adapter is securely connected to both the massager and the electrical outlet. Ensure the remote control's battery tab has been removed.
- **Massage feels too intense or too weak:** Adjust the kneading and air compression intensity levels using either the control panel or the remote control. Start with lower settings and gradually increase to your comfort level.
- **Heat function not felt:** Ensure the heat function is activated. Allow a few minutes for the heating elements to warm up.
- **Unusual noises or malfunctions:** Discontinue use immediately and unplug the device. Contact KNQZE customer support for assistance.

SPECIFICATIONS

Feature	Detail
Product Dimensions	15.7 x 13.38 x 7.48 inches
Item Weight	7 Pounds
Model Name	KF9
Power Source	Corded Electric

Feature	Detail
Voltage Input	AC 100-240V ~50/60Hz
Output	24VDC 1.5A
Nominal Power	36W
Material Type	Plastic, Polyester
Massage Therapy Technique	Compression, Kneading, Rolling, Shiatsu
Heat Levels	3 (113°F, 122°F, 131°F)
Timer Settings	3 (15, 20, 30 minutes)

WARRANTY AND SUPPORT

The KNQZE foot massager comes with a one-year warranty from the date of purchase.

If you encounter any issues with your product or require technical assistance, please contact KNQZE customer support. For more information, you can visit the [KNQZE Store on Amazon](#).