



Manuals.plus /

› **NaturaleBio** /

› NaturaleBio Neck Massager Instruction Manual - Model 8057949208556

NaturaleBio 8057949208556

NaturaleBio Neck Massager Instruction Manual

Model: 8057949208556

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the NaturaleBio Neck Massager. Failure to follow these safety guidelines may result in injury or damage to the device.

- Always follow the instructions provided in this user manual for safe and effective use.
- Do not use the massager on broken skin or areas with inflammation.
- Keep the massager away from water to prevent electrical shock.
- Do not leave the massager unattended while in use.
- Allow the massager to cool down before storing it.
- Consult a healthcare professional before use if you have any medical conditions or are pregnant.
- This device is not intended to diagnose, treat, cure, or prevent any disease or health condition.

2. PRODUCT OVERVIEW

The NaturaleBio Neck Massager is a portable electric device designed to provide deep tissue massage and heat therapy for the neck and surrounding areas. It utilizes electric pulse technology to help relieve muscle stiffness and promote relaxation.

Key Features:

- **Electric Pulse Technology:** Effectively relieves neck soreness and stiffness.
- **Heat Function:** Provides constant temperature hot compress for enhanced comfort.
- **Multi-Mode Regulation:** Features 10 massage modes and 16 intensity levels for personalized massage.
- **Ergonomic Design:** Conforms to various neck curves for a comfortable fit.
- **Portable and Lightweight:** Compact design (130g ultra-light weight) for use at home, office, or while traveling.
- **Versatile Use:** Suitable for the neck, and can also be applied to the back, waist, and legs for general

muscle relaxation.

Package Contents:

- 1 x Neck Massager Device
- 1 x Remote Control
- 1 x USB Charging Cable
- 1 x User Manual
- Hydrogel Pads



The package includes the neck massage device, remote control, USB charging cable, user manual, and hydrogel pads.

3. SETUP

Charging the Device:

Before first use, fully charge the neck massager using the provided USB charging cable. Connect the cable to the massager and a compatible USB power source. The indicator light will show charging status and turn off or change color when fully charged.

Preparation for Use:

1. **Moisten Skin:** Ensure the skin on your neck is clean and slightly moist. You may use a small amount of water, aloe vera gel, or a non-oily cream/lotion on the contact areas for better conductivity and comfort.
2. **Apply Hydrogel Pads:** If using, tear off the blue protective film from the transparent hydrogel pads and apply them to the metal electrodes of the massager.
3. **Position the Massager:** Place the massager around your neck, ensuring the electrodes make firm contact with your skin. Adjust for a comfortable and secure fit.

Apply Gel, Cream Or Lotion To Moisten The Skin Before Wearing



Follow these steps for optimal use: moisten skin, apply hydrogel, power on, and adjust intensity.

4. OPERATING INSTRUCTIONS

Powering On/Off:

- **To Turn On:** Press and hold the power button on the device or remote control for a few seconds until the display illuminates.
- **To Turn Off:** Press and hold the power button again until the display turns off.

Selecting Massage Mode:

The massager offers 10 distinct massage modes. Use the 'Mode' button on the device or remote control to cycle

through the available modes. The current mode will be displayed on the screen.

Adjusting Intensity Levels:

There are 16 intensity levels available. Use the '+' and '-' buttons on the device or remote control to increase or decrease the massage intensity. It is recommended to start from the lowest level and gradually increase to a comfortable setting.

16 Levels of Intensity

Meet your individual massage needs

- STAA Mode
- PLS Mode
- Tui Na Mode
- Acupuncture Mode
- Hammer Mode
- Mild Mode
- Kneading Mode
- Auto Mode

— [1-5 Levels: Falling Pillow Pain | 6-10 Levels: Daily Cervical Care | 11-16 Levels: Chronic Neck Strain] +

*It is recommended to start from the lowest levels.

Choose from 10 massage modes and 16 intensity levels to customize your massage experience.

Heat Function:

The massager features a hot compress function. Activate or deactivate this function using the dedicated heat button on the device or remote control. The heat is maintained at a constant temperature for comfort.




ENHANCE LITERARY JOURNEY

EXPERIENCE COMFORT
AND BOOST FOCUS

The NatureBio neck massager is designed for comfortable use during various activities, such as reading.

ELEVATE WORK COMFORT

WITH ERGONOMIC
NECK MASSAGER

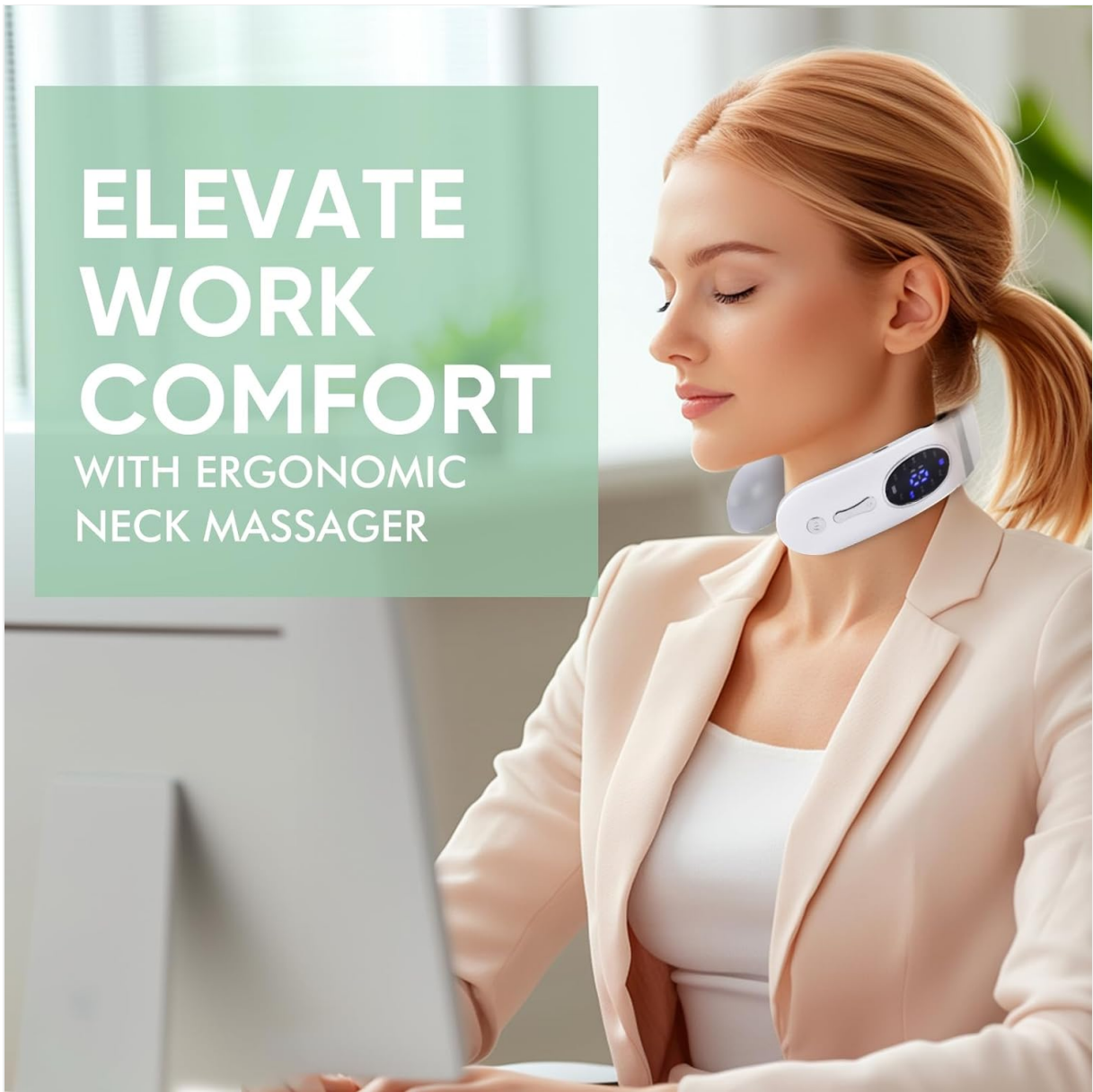

The ergonomic design allows for comfortable wear, enhancing comfort during work or study.

5. MAINTENANCE

Cleaning:

- Ensure the device is turned off and unplugged before cleaning.
- Wipe the surface of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- Clean the electrode pads gently after each use to remove any residue.

Storage:

- Store the massager in a cool, dry place away from direct sunlight and extreme temperatures.
- Keep out of reach of children and pets.
- If hydrogel pads are used, ensure they are stored properly to maintain their adhesive quality.

6. TROUBLESHOOTING

If you encounter any issues with your NaturaleBio Neck Massager, please refer to the following common problems and solutions:

- **Device does not turn on:** Ensure the device is fully charged. Press and hold the power button for a few seconds.
- **No sensation or weak pulse:** Ensure the electrode pads are making firm contact with moist skin. Increase the intensity level gradually. Check if hydrogel pads are properly applied or need replacement.
- **Discomfort during use:** Reduce the intensity level. Ensure skin is adequately moistened. Reposition the massager for a better fit. Discontinue use if discomfort persists.
- **Remote control not working:** Check the battery in the remote control. Ensure there are no obstructions between the remote and the massager.

If the problem persists after attempting these solutions, please contact customer support.

7. SPECIFICATIONS

Product Dimensions	0.39 x 0.39 x 0.39 inches
Item Weight	9.14 ounces (approximately 259g)
Item Model Number	8057949208556
Manufacturer	NaturaleBio
Power Source	Battery Powered (Rechargeable via USB)
Use For	Neck (also suitable for back, waist, legs)
Color	White
Massager Form	Wrap Around Neck Massager

NECK SOOTHING PROFESSIONAL MASSAGER



DIZZINESS & HEADACHE



SORENESS AND PAIN



NECK SORENESS AND STIFFNESS

This massager is designed to alleviate symptoms associated with neck soreness, stiffness, and related discomforts.

8. WARRANTY INFORMATION

NaturaleBio products are manufactured to high-quality standards. This product is covered by a standard manufacturer's warranty against defects in materials and workmanship from the date of purchase. Please retain your proof of purchase for warranty claims. The warranty does not cover damage caused by misuse, accident, unauthorized modification, or normal wear and tear.

9. CUSTOMER SUPPORT

For any questions, concerns, or technical assistance regarding your NaturaleBio Neck Massager, please contact our customer support team. Refer to the contact information provided on the product packaging or our official website.

Thoughtful Love, Warm Gifts

Neck Soothing Magic



The NaturelBio neck massager is a suitable gift for promoting relaxation and well-being.