

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [PRITECH](#) /

› [PRITECH Electric Foot Callus Remover \(BCM-1852\) User Manual](#)

PRITECH BCM-1852

PRITECH Electric Foot Callus Remover (Model BCM-1852)

User Manual

INTRODUCTION

The PRITECH Electric Foot Callus Remover (Model BCM-1852) is designed to effectively remove dead skin and calluses from your feet, providing a smooth and healthy appearance. This rechargeable, portable electronic foot file features multiple roller heads and adjustable speeds for personalized foot care.

SAFETY INFORMATION

- Always use on dry skin.
- Avoid applying excessive pressure to prevent injury.
- Do not use on broken or irritated skin.
- Keep the product away from water when not in use.

PRODUCT OVERVIEW

The PRITECH Electric Foot Callus Remover comes with various components designed for optimal performance and user convenience.



The PRITECH Electric Foot Callus Remover, shown with hands and feet in the background, highlighting its purpose for foot care.



The device includes three interchangeable quartz sand grinding heads: Extra Coarse for tough calluses, Regular Coarse for dead skin, and Gentle Coarse for daily maintenance.

SETUP

Charging the Device

The device is equipped with a 1200mAh lithium battery. Before first use, fully charge the device.

1. Locate the USB charging port at the bottom of the device.
2. Connect the provided USB charging cable to the device and a suitable USB power source (e.g., laptop, power bank, wall adapter).
3. The LED power display will indicate the charging status. A full charge takes approximately 130 minutes and

provides about 120 minutes of continuous use.

FAST USB CHARGING



Laptop



Socket



Power Bank



Car Charger



1200mah Battery
Charger Time: **130** mins
Using Time: **120** mins

The device features fast USB charging, compatible with various USB power sources like laptops, power banks, and wall sockets.

Attaching Roller Heads

To change or attach a roller head:

1. Ensure the device is turned off.
2. Gently press the release button on the side of the roller head and pull the head upwards to remove it.
3. Align the desired roller head with the device and push it down until it clicks securely into place.



Illustration of how to easily attach and detach the different roller heads for various exfoliation needs.

OPERATING INSTRUCTIONS

Follow these steps for effective and safe use of your callus remover:

1. **Preparation:** Start with dry feet for safety and better results.
2. **Power On:** Press the power button to turn on the device. The LED display will show the current speed setting and battery level.
3. **Speed Adjustment:** The device offers three speed modes (1700 r/pm, 1900 r/pm, 2100 r/pm). Press the power button repeatedly to cycle through the speeds to select your preferred intensity.
4. **Application:** Gently apply the foot callus remover to the affected areas, avoiding broken or irritated skin. Use light pressure and do not go too deep into the skin to avoid pain or injury.
5. **Targeting Calluses:** For tough calluses, you may need to repeat the process a few times. Move the device slowly over the area.
6. **Post-Treatment:** After use, clean your feet and moisturize with foot cream.

THREE SPEEDS ADJUSTABLE



HIGH SPEED: 2100 R/PM



NORMEL SPEED: 1900 R/PM



LOW SPEED: 1700 R/PM

ON/OFF
Speed Switch

20
%
POWERLEFT
2

The device features three adjustable speeds (1700, 1900, and 2100 RPM) to cater to different callus removal needs, from gentle to powerful.

1. WASH

Soaked your feet in the warm water for about 10 minutes then dry out with towel.



2. GRIND

Turn on the electric callus remover and slide over the callused feet gently to achieve the desired result.



3. PROTECT

Clean out your feet and moisturize with foot cream.



A visual guide to the three-step foot care process: 1. Wash feet, 2. Grind calluses with the device, 3. Protect and moisturize feet.

MAINTENANCE


Cleaning the Device

The PRITECH Electric Foot Callus Remover has a waterproof design (IPX6 rated), allowing for easy cleaning.

1. Ensure the device is turned off and disconnected from any charging cable.
2. Remove the roller head as described in the Setup section.
3. Rinse the roller head and the main body of the device under running water. You may use mild soap if needed.
4. Use the provided cleaning brush to remove any remaining dead skin particles from the roller head and device.
5. Dry the foot callus remover thoroughly before storing it. Avoid frequent submersion of the entire device to prevent potential issues.

PRITECH CALLUS REMOVERS

Waterproof Level: IPX6

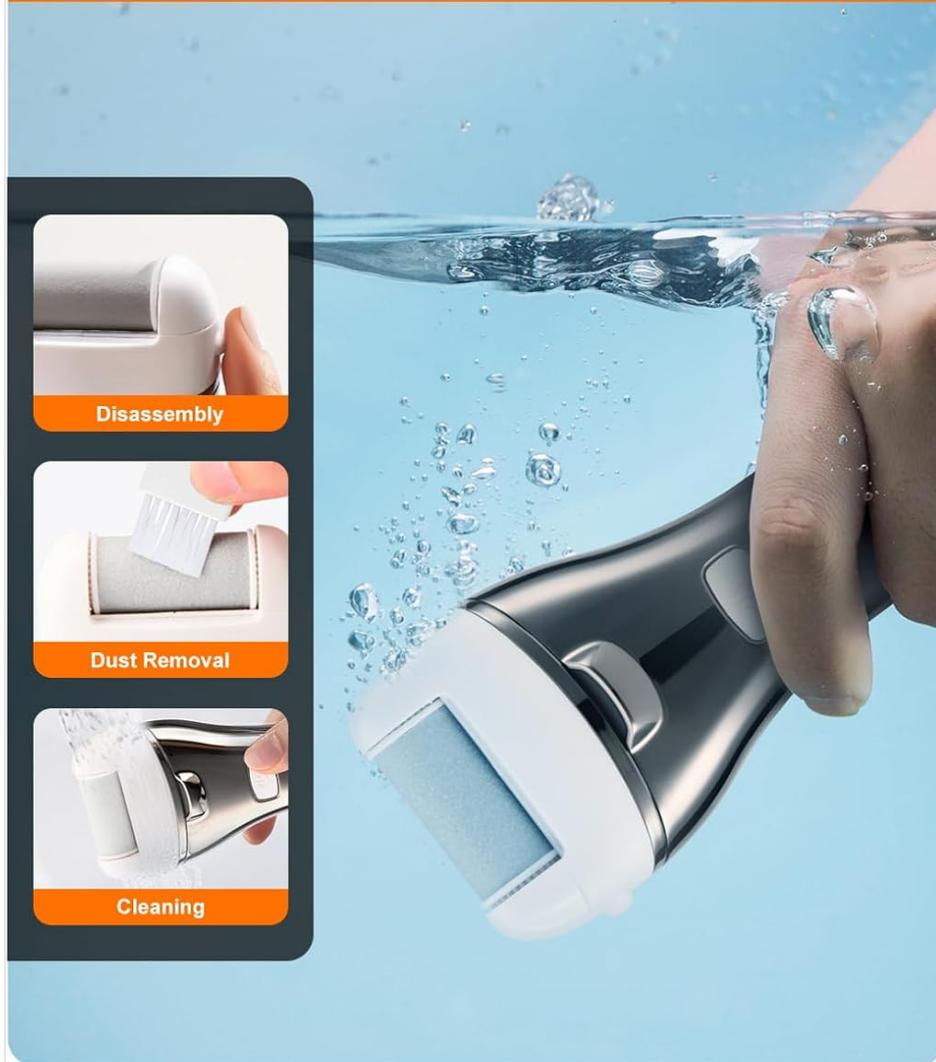

Other
IPX4

The device is IPX6 waterproof, allowing for easy cleaning under running water. The image shows disassembly, dust removal with a brush, and rinsing.

Replacing Roller Heads

Roller heads should be replaced when they show signs of wear or become less effective. Replacement heads can be purchased separately.

1. Follow the steps for 'Attaching Roller Heads' to remove the old head and install a new one.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on	Low battery or no charge	Charge the device fully using the provided USB cable.

Problem	Possible Cause	Solution
Roller head stops spinning during use	Excessive pressure applied	Reduce pressure on the skin. The device has a safety feature to prevent injury from too much force.
Reduced effectiveness of roller head	Roller head is worn out or clogged	Clean the roller head or replace it with a new one.
Device not charging	Faulty cable or power source	Try a different USB cable or power adapter. Ensure the charging port is clean.

SPECIFICATIONS

Feature	Detail
Model Number	BCM-1852
Product Dimensions	6.93 x 5.12 x 1.97 inches
Item Weight	5.29 ounces (0.15 Kilograms)
Battery	1200mAh Lithium (1 AA battery required, included)
Charging Time	Approx. 130 minutes
Operating Time	Approx. 120 minutes
Speed Settings	3 (1700 r/pm, 1900 r/pm, 2100 r/pm)
Waterproof Rating	IPX6
Material	Plastic
Color	White

WARRANTY AND SUPPORT

PRITECH products are manufactured to high-quality standards. For any issues, please contact our customer support. If you receive the product with any problems such as wear or stains, please contact us in the first instance and we will provide you with a solution.

For support, please visit the [PRITECH Store on Amazon](#).

LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

