

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [SALLOOT](#) /

› [SALLOOT TRIPLE-Pulse IPL Laser Hair Removal Device \(Model DT015\) User Manual](#)

SALLOOT DT015

SALLOOT TRIPLE-Pulse IPL Laser Hair Removal Device (Model DT015) User Manual

Model: DT015 | Brand: SALLOOT

INTRODUCTION

This manual provides essential instructions for the safe and effective use of your SALLOOT TRIPLE-Pulse IPL Laser Hair Removal Device, Model DT015. Please read this manual thoroughly before first use and retain it for future reference.

SAFETY INFORMATION

Important Safeguards

- Always wear the provided protective eyewear during operation to shield your eyes from the intense pulsed light.
- Ensure your skin is clean, dry, and free from any lotions, creams, or deodorants before treatment.
- Do not use the device on irritated, sunburned, or broken skin.
- Consult a medical professional if you have any pre-existing skin conditions or concerns.
- Keep the device out of reach of children.

Skin Tone and Hair Color Suitability

The SALLOOT IPL device works by targeting melanin in the hair follicle. Therefore, its effectiveness and safety depend on your natural skin tone and hair color. IPL is most effective on lighter skin tones with darker hair. It is **not suitable** for very dark skin tones or very light hair colors (red, blonde, grey, white) as it may cause skin damage or be ineffective.

35 Adjustable Sensitivity Levels for the Whole Body

Safe and Efficient full-body hair removal can be achieved within 2-8 minutes



Figure 1: IPL Hair Removal Skin Tone and Hair Color Suitability Chart. This chart illustrates which skin tones and hair colors are safe and effective for IPL treatment. Darker skin tones and lighter hair colors are generally not suitable.

 **PCOS Friendly**

Can I Use IPL?

IPL primarily targets melanin, which means darker skin tones with higher melanin levels increase the risk of burns. Light hair color lacks sufficient melanin for IPL to be effective.

Step 1 Identify your skin tone



Step 2 Identify the color of your unwanted hair



✓ Skin Tone + ✓ Hair Color = Applicable

Either Skin Tone or Hair Color ✗ = Not Applicable



Figure 2: Detailed suitability guide for IPL hair removal, including considerations for PCOS. It shows that IPL primarily targets melanin, making it less effective or unsafe for light hair colors and very dark skin tones.

If you are unsure about your suitability, please refer to the skin tone and hair color chart in the manual or consult customer support.

PRODUCT OVERVIEW

The SALLOOT TRIPLE-Pulse IPL Laser Hair Removal Device (Model DT015) is designed for effective and comfortable at-home hair reduction. Key features include:

- **TRIPLE-Pulse Technology:** Delivers 25J energy density (7-8J/cm²) to target hair follicles for up to 97% hair reduction in 2 weeks.
- **12cm² Sapphire Ice-Cooling:** Instantly cools the skin surface to 37°F (3°C) and maintains temperatures below 59°F (15°C) for comfortable, pain-free treatments, even on sensitive areas.
- **Fast Full-Body Sessions:** Auto-glide mode allows for full-body treatment in 2-8 minutes.
- **AI Coaching App:** Analyzes skin tone and hair color to provide personalized treatment plans, intensity recommendations, and session reminders.
- **Ergonomic Design:** 30-degree angled design for comfortable handling and easy reach on various body parts.
- **LCD Screen:** Displays intensity level and 7 modes in real-time.

Unrivaled Triple-Pulses IPL Technology for Unmatched Efficiency

Triple Pulses™ Light & IGBT FlasherPro Technology

200 flashes per minute deliver stable, effective pulses, enhancing hair removal efficiency

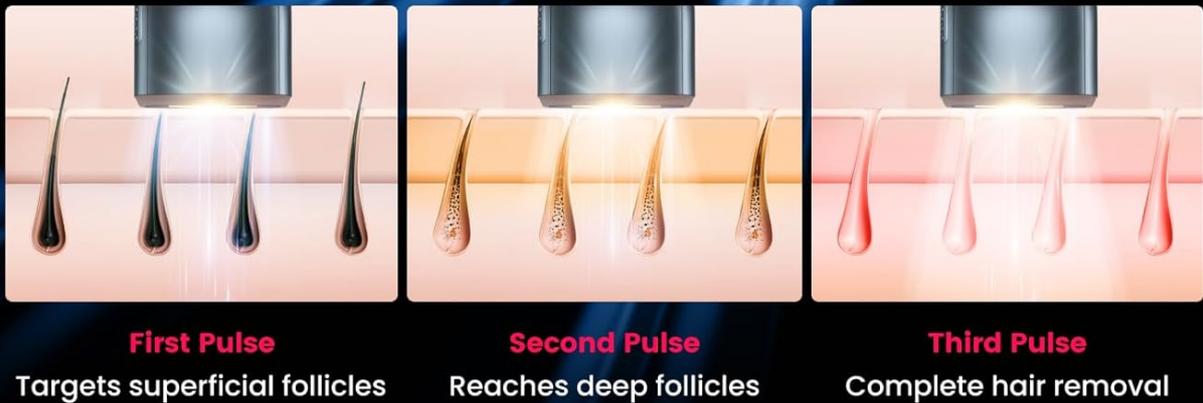


Figure 3: Illustration of Triple-Pulse IPL Technology. This image shows how the device delivers three pulses to target superficial and deep hair follicles for comprehensive hair removal.

Quadruple Sapphire Frost Cube™ Technology



Speedy Cooling

Instantly cools to 59°F upon startup



Extreme Comfy Cooling

Lowers to 37°F



Long-lasting Cooling

Maintains 59°F for 35 minutes of continuous use



12cm² Biggest full-screen sapphire

Unparalleled Comfort



Figure 4: Features of Quadruple Sapphire Frost Cube Technology. This image highlights the speedy cooling, extreme comfy cooling, long-lasting cooling, and the large 12cm² full-screen sapphire window for enhanced comfort during treatment.

BURNING VS PAINLESS ZAPS



Figure 5: Comparison of traditional hair removal methods versus the SALLOOT IPL device with 37°F Ice-cooling. This image emphasizes the pain-free experience provided by the cooling technology.

SETUP

1. **Prepare Your Skin:** Shave the area to be treated to ensure the IPL light targets the hair follicle directly, not the surface hair.
2. **Wear Protective Eyewear:** Always put on the provided safety goggles to protect your eyes from the intense light flashes.
3. **Connect to Power:** Plug the device into a power outlet.
4. **Download the App (Optional but Recommended):** Download the "CTONE" app from your app store. The app provides personalized treatment plans and guidance.

AI Adaptive™ Personalized Treatments for Optimal Results

- **Personalized** Full Body Treatment
- Precision **Reminder**
- **Step-by-Step** Tutorials
- Results **Tracking**



Figure 6: The Smart AI APP provides personalized treatment plans, tracks treated areas, and sends reminders for optimal results.

Your browser does not support the video tag.

Video 1: This video demonstrates the initial setup process, including shaving the treatment area, wearing protective eyewear, and connecting the device to the CTONE app for personalized settings. It also highlights the triple-pulse technology and sapphire cooling.

OPERATING INSTRUCTIONS

1. **Power On:** Long press the power button to turn on the device.
2. **Select Mode and Intensity:** The device features 7 modes for different body areas and 5 intensity levels. Use the LCD screen to select the appropriate mode (e.g., legs, underarms, bikini, face) and intensity level based on your skin tone and hair color. The AI app can assist with this selection.
3. **Apply to Skin:** Place the device's treatment window firmly against your skin at a 90-degree angle. The device will detect skin contact.
4. **Initiate Flash:** Press the flash button. For continuous treatment, use the auto-glide mode by holding down the flash button and gliding the device over the skin.

5. Treatment Schedule:

- Weeks 1-2: 3-4 times per week.
- Weeks 3-4: 2-3 times per week.
- After 4 weeks: 2 times per month for long-lasting results.

6. Post-Treatment Care: After each session, you may apply a soothing lotion or gel to the treated area.

Full Body Treatment in Just 2-8 Minutes



- 

Step 1 1
Prep your skin by shaving
- 

Step 2 2
Wear sunglasses to protect your eyes
- 

Step 3 3
Choose from 5 intensity levels
- 

Step 4 4
Press the flash button

Figure 7: Step-by-step guide for full body IPL treatment, including skin preparation, eye protection, intensity selection, and activating the flash button.

Less Hair Growth After Just 1 Use

- ✓ Permanent Results in 12 Weeks
- ✓ No More Razor Burn Ingrown Hairs
- ✓ Dolphin-smooth, touchably soft skin

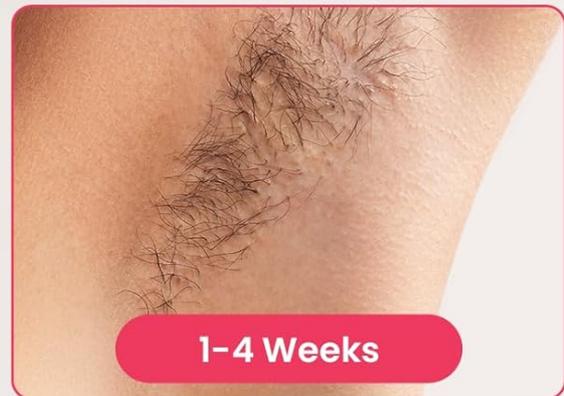


Figure 8: Visual representation of hair growth reduction over 12 weeks, showing less hair growth after just one use and permanent results over time.

Instant AI Analysis Hair Removal Reports

Like a home beautician, smart personalized hair removal scheduling and tracking tailored to you


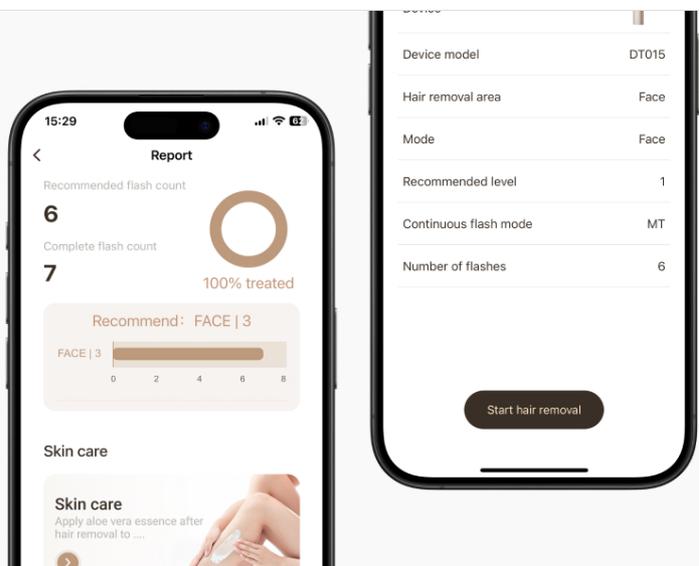


Figure 9: Screenshot of the AI app providing instant hair removal reports and personalized scheduling.

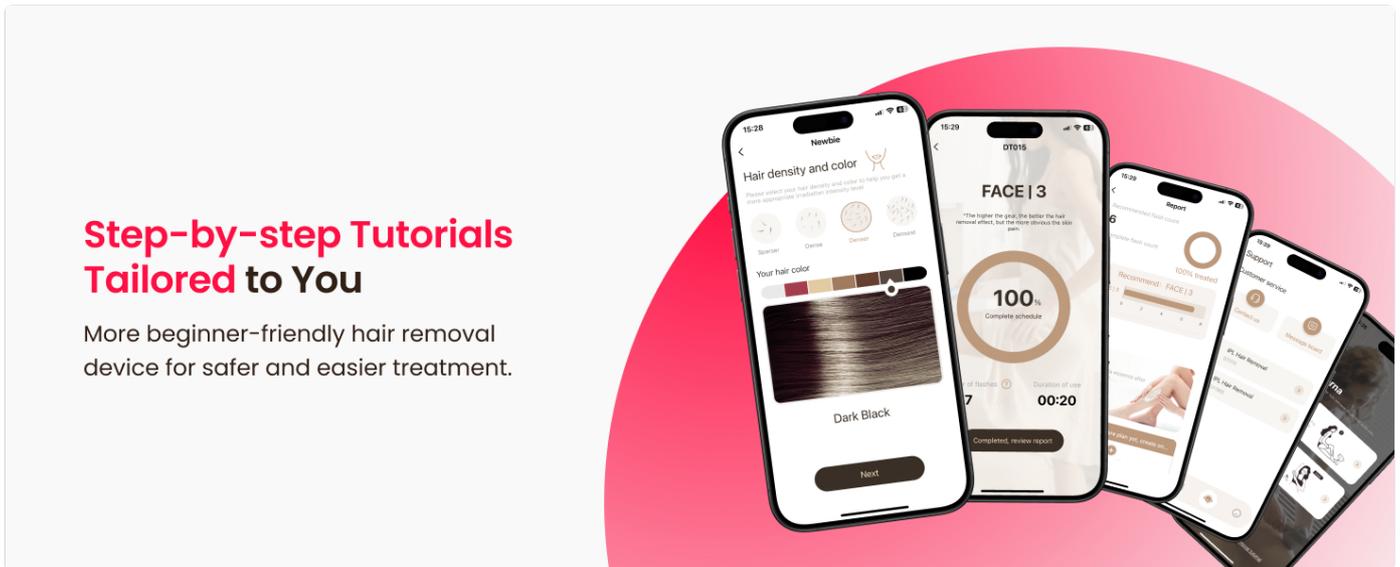

Step-by-step Tutorials Tailored to You

More beginner-friendly hair removal device for safer and easier treatment.

Figure 10: The app offers step-by-step tutorials tailored for easier and safer treatment.

MAINTENANCE

- **Cleaning:** After each use, ensure the treatment window is clean. Use a soft, dry cloth to gently wipe away any residue. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.
- **Inspection:** Regularly inspect the device and its power cord for any signs of damage. If any damage is found, discontinue use and contact customer support.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not properly connected to power, power button not pressed long enough.	Ensure the power cord is securely plugged in. Press and hold the power button for a few seconds.
No flash emitted.	Device not in full contact with skin, incorrect skin tone/hair color, device malfunction.	Ensure the treatment window is fully pressed against the skin. Verify suitability with the skin tone chart. If issues persist, contact support.
Skin irritation or discomfort.	Intensity level too high, unsuitable skin type, improper use.	Reduce the intensity level. Review skin tone suitability. Ensure skin is properly prepared. Apply a soothing gel after treatment.
Hair regrowth is not reducing.	Inconsistent use, incorrect intensity, unsuitable hair color.	Follow the recommended treatment schedule. Increase intensity if comfortable and suitable. Verify hair color suitability.

SPECIFICATIONS

Model Number: DT015

Brand: SALLOOT

Technology: TRIPLE-Pulse IPL with Sapphire Ice-Cooling

Energy Density: Up to 7-8J/cm² (25J total)

Cooling Temperature: Down to 37°F (3°C)

Treatment Window Size: 12 cm²

Modes: 7 body modes

Intensity Levels: 5 adjustable levels

App Compatibility: Smart AI APP (CTONE)

UPC: 642380741281

ASIN: B0DPM2XGRC

Package Dimensions: 9.65 x 6.14 x 4.41 inches

Item Weight: 2.8 Pounds

WARRANTY & SUPPORT

SALLOOT offers an extended 3-year support for this device, including a free trial and no-reason return policy. For any questions, technical assistance, or warranty claims, please contact SALLOOT customer support through their official website or the contact information provided with your purchase.

You can also visit the [SALLOOT Store on Amazon](#) for more information and support resources.