



Manuals.plus /

› Norelie /

› Norelie RedRevive Cold Laser Therapy Device User Manual

Norelie RedRevive

Norelie RedRevive Cold Laser Therapy Device User Manual

Model: RedRevive

1. INTRODUCTION

The Norelie RedRevive Cold Laser Therapy Device is a handheld unit designed to provide dual-wavelength (650nm and 808nm) light therapy for deeper tissue support. It aims to help ease discomfort, reduce tension, and promote natural recovery through consistent daily use. This device is intended for daily wellness support and is not designed for instant relief. It features customizable settings for intensity and session duration, and its portable design allows for convenient use at home or on the go.

2. SAFETY INFORMATION

- Discontinue use if irritation, rash, or adverse reaction occurs.
- Not recommended for use if pregnant, or if you have pacemakers, electronic implants, or metastatic conditions.
- Avoid direct eye exposure to the laser. Always use the provided protective glasses during operation.
- Keep the device away from flammable substances during operation.

3. PRODUCT OVERVIEW

The RedRevive device utilizes dual wavelengths (650nm and 808nm) to target both superficial and deeper tissues. The 650nm red light is ideal for surface-level treatment, supporting skin, shallow muscles, and tendons. The 808nm infrared light penetrates deeper to reach joints, ligaments, and tissues, helping to alleviate stiffness and discomfort.

The device offers 4 intensity levels and an adjustable timer (5/10/15/20 minutes) for personalized therapy sessions. Its lightweight, cordless design with a 2600mAh rechargeable battery ensures portability and extended use.



A Gift of Wellness for the Ones You Love

Figure 1: Visual representation of the three-step operation process for the Norelie RedRevive device.

4. PACKAGE CONTENTS

Your Norelie RedRevive Cold Laser Therapy Device package includes:

- Norelie RedRevive Therapy Wand
- Storage Box
- Protective Glasses
- User Guide
- Charging Adapter
- Charging Cable
- Adjustable Stand (for hands-free use)

Consistent Use Daily, Emerge Real Benefits



Week 1-2 – Building the Habit

Relief may be subtle or unnoticeable as your body adapts.



Week 3-4: Gradual Improvements

You'll start moving comfortably



Week 5+: Confident & Active

Daily consistent use delivers lasting results

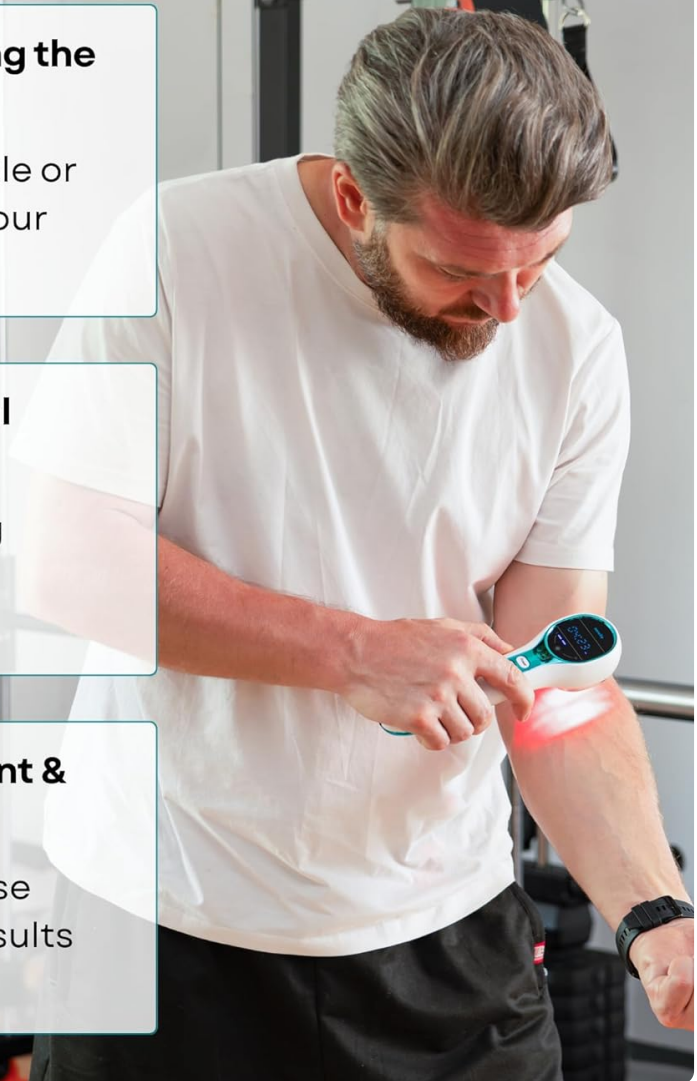


Figure 2: All components included with the Norelie RedRevive Cold Laser Therapy Device.

5. SETUP

1. **Charge the Device:** Before first use, fully charge the RedRevive wand using the provided charging adapter and cable. The charging port is typically located at the base of the device.
2. **Assemble the Stand (Optional):** If using the adjustable stand for hands-free operation, assemble it according to the instructions provided in the user guide. Securely attach the device to the stand's clip.

6. OPERATING INSTRUCTIONS

1. **Wear Protective Glasses:** Always put on the provided protective glasses before turning on the device to shield your eyes from the laser light.
2. **Power On:** Press and hold the power button (usually the bottom button) for 2 seconds until the display illuminates.
3. **Adjust Timer:** Press the 'Timer' button to cycle through the available session durations (5, 10, 15, or 20 minutes).
4. **Select Intensity:** Press the 'Power' button to choose from 4 intensity levels (L1-L4).

5. **Begin Therapy:** Once the desired timer and intensity are set, press the power button again to start the therapy. The device will emit red and infrared light.
6. **Apply to Target Area:** Hold the device 0.5-1.2 inches (1-3 cm) away from the skin on the target area. You can gently glide the wand or hold it stationary.
7. **End Session:** The device will automatically turn off once the timer expires. You can also manually turn it off by pressing and holding the power button.

Your browser does not support the video tag.

Video 1: Official Norelie RedRevive red light laser therapy demonstration, showing device operation and application.

7. USAGE GUIDELINES


- **Session Duration:** Use for 15 minutes per session.
- **Consistency:** For optimal results, use consistently on a daily basis.
- **Expected Results:** This device is designed for gradual improvement. Some users may notice changes within the first week, while others may require 2–4 weeks of daily use depending on individual needs and condition.
- **Application:** The device can be used on various body parts for discomfort, tingling, numbness, and to support circulation and relaxation. It is also suitable for pets.

PAIN RELIEF AND EASY OPERATION



**Natural healing, non-invasive treatment,
no cross-infection**

Relieve Pain. Restore Motion.

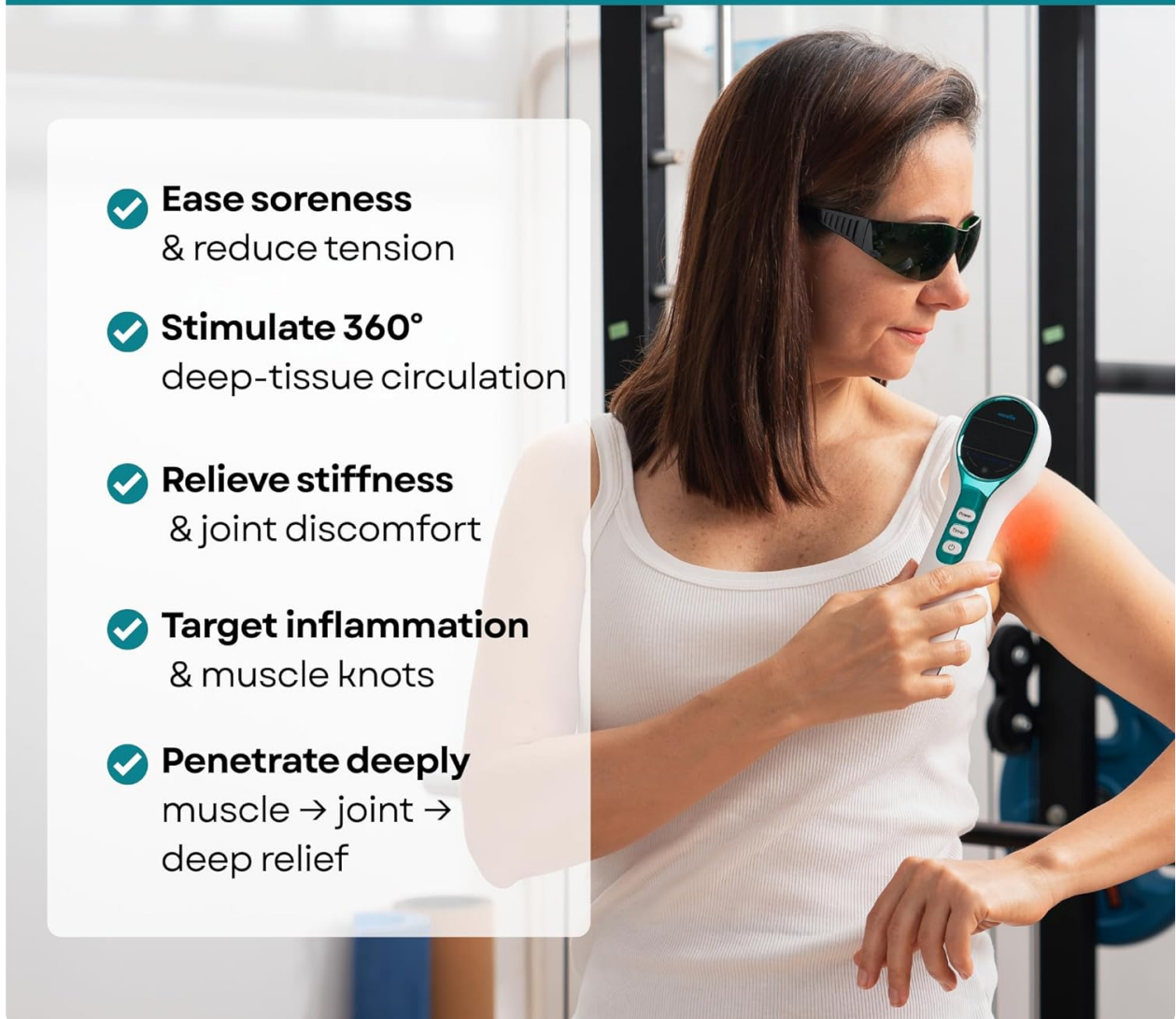

- ✓ **Ease soreness**
& reduce tension
- ✓ **Stimulate 360°**
deep-tissue circulation
- ✓ **Relieve stiffness**
& joint discomfort
- ✓ **Target inflammation**
& muscle knots
- ✓ **Penetrate deeply**
muscle → joint →
deep relief

Figure 4: The Norelie RedRevive device applied to a human shoulder for targeted therapy.

8. MAINTENANCE AND CARE

- **Cleaning:** Gently wipe the device with a soft, dry cloth after each use. Do not immerse the device in water or use near a water source.
- **Storage:** Store the device in its provided storage box in a cool, dry place when not in use.
- **Charging:** Recharge the device as needed. The 2600mAh battery provides extended use.

9. TROUBLESHOOTING

If you encounter any issues with your Norelie RedRevive device, please refer to the user guide for basic troubleshooting steps. If the problem persists, contact customer support for assistance.

10. SPECIFICATIONS

| Feature | Detail |
|--------------------|------------------------------------|
| Product Dimensions | 3 x 2.5 x 9.5 inches |
| Item Weight | 15.52 ounces |
| Model Number | RedRevive |
| Batteries | 1 9V batteries required (included) |
| Power Source | Battery Powered |
| Material | Lightweight, Durable Plastic |

11. WARRANTY AND SUPPORT

Your satisfaction is our priority. Register your RedRevive Wand to activate a 90-day extended warranty. If your device arrives defective or doesn't meet expectations, we offer replacement or support for 90 days. For assistance, contact us via the "Contact Seller" button on the product page or visit the Norelie Store.

12. LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.