

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [COMFIER](#) /

› [COMFIER CF-6419M Cordless Shiatsu Neck and Back Massager with Heat Instruction Manual](#)

## COMFIER CF-6419

# COMFIER CF-6419M Cordless Shiatsu Neck and Back Massager with Heat Instruction Manual

Model: CF-6419

## INTRODUCTION

Thank you for choosing the COMFIER CF-6419M Cordless Shiatsu Neck and Back Massager with Heat. This manual provides essential information for the safe and effective use of your new massager. Please read it thoroughly before operation and retain it for future reference.



Image: The COMFIER CF-6419M massager pillow, showcasing its ergonomic design and compact size.

## IMPORTANT SAFETY INFORMATION

To reduce the risk of electric shock, burns, fire, or injury, always follow these basic precautions:

- Always unplug the massager when not in use.
- Do not use on injured or sensitive areas of the body.
- Keep away from water or damp environments.
- Do not leave unattended while in use.
- Allow the massager to cool down completely before storing.

- Do not use if the charging cord or any part of the device is damaged.
- Consult a doctor before use if you have any medical conditions, are pregnant, or have a pacemaker.
- Keep out of reach of children and pets.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## PACKAGE CONTENTS

---

- COMFIER CF-6419M Shiatsu Neck and Back Massager with Heat
- USB Type-C Charging Cable
- User Manual

## PRODUCT OVERVIEW

---

The COMFIER CF-6419M is a versatile, cordless massager designed for deep tissue relief. It features 3 adjustable shiatsu modes, 2 heat levels, and flexible massage nodes. The device is compact and lightweight, making it suitable for various environments.



Image: A person relaxing with the COMFIER massager positioned on their neck, illustrating its ergonomic fit.

## Controls and Features:

- **Power/Mode Button:** Turns the device on/off and cycles through massage modes (M1, M2, M3).
- **Direction Button:** Changes the rotation direction of the massage nodes.
- **Heat Button:** Activates and adjusts the heat levels (Low, High, Off).
- **Type-C Charging Port:** For convenient recharging.
- **Adjustable Velcro Straps:** For securing the massager to chairs or headrests.
- **Flexible & Gentle Massage Nodes:** Designed for comfortable deep tissue relief.
- **Premium Anti-Fray Fabric:** Durable and skin-friendly exterior.

## Easy to Use

**M1:** 1 min auto reverse - 2 beeps

**M2:** Change direction - 1 beep

**M3:** Change direction - 1 beep



**Heat**  
High/Low/Off



**Direction**  
M1/M2/M3



**Power/Mode**  
On / Off  
Low / High / Only heating

Image: Close-up of the control panel on the massager, showing the power, direction, and heat buttons, along with the Type-C charging port.

## SETUP

1. **Initial Charge:** Before first use, fully charge the massager using the provided USB Type-C charging cable. Connect the cable to the Type-C port on the massager and to a compatible USB power source. The battery indicator will show charging status. A full charge provides up to 1.5 hours of use.
2. **Positioning:** The massager can be used on various body parts. For hands-free use, utilize the adjustable Velcro straps to secure it to an office chair, car seat headrest, or other suitable seating surfaces.

# Velceo Straps

Hands-Free



Image: A person using the massager secured to an office chair with its adjustable Velcro straps, demonstrating hands-free operation.

## OPERATING INSTRUCTIONS

1. **Power On:** Press and hold the Power/Mode button to turn on the massager. The massage nodes will begin rotating.
2. **Select Massage Mode:**
  - **M1 (Auto Reverse):** The massage nodes will automatically reverse direction every minute. (Default mode upon power-on)
  - **M2 (Forward Rotation):** Press the Power/Mode button once to activate continuous forward rotation.

- **M3 (Reverse Rotation):** Press the Power/Mode button a second time to activate continuous reverse rotation.

3. **Adjust Heat Function:** Press the Heat button to cycle through heat settings:

- **Low Heat:** Approximately 107.6°F (42°C).
- **High Heat:** Approximately 113°F (45°C).
- Press again to turn heat off. Heat reaches desired temperature within 30-60 seconds.



Image: Visual representation of the two heat levels available on the massager, Low Heat (107.6°F) and High Heat (113°F).

4. **Auto Shut-Off:** The massager is equipped with an auto shut-off feature after 10 minutes to prevent overuse and ensure safety. You can restart the device after a brief rest period.
5. **Power Off:** Press and hold the Power/Mode button to turn off the massager.

### Recommended Usage Areas:

The versatile design allows for use on various body parts:

- Neck and Shoulders
- Upper and Lower Back
- Waist
- Arms and Legs
- Feet and Hands
- Abdomen (for gentle warmth and relaxation)

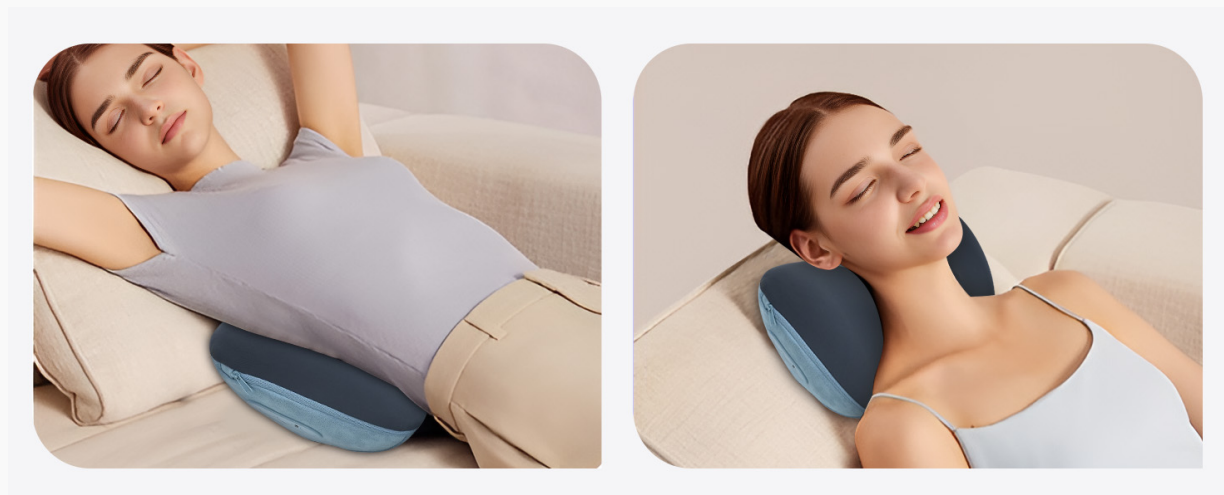


Image: Demonstrations of the massager being used on both the neck and lower back for targeted relief.



## Official Product Video:

Your browser does not support the video tag.

Video: An official product demonstration of the COMFIER Shiatsu Neck and Back Massager with Heat, showcasing its features and various uses.

## MAINTENANCE

---

- **Cleaning:** Unplug the massager before cleaning. Wipe the surface with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the device in water. The exterior fabric is designed to be durable and resistant to wear.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Ensure the device is fully cooled down before storage.

# Premium Anti-Fray Fabric

Resists Wear & Tear



Image: Close-up of the massager's premium anti-fray fabric, highlighting its high elasticity, supportability, and breathable qualities.

## TROUBLESHOOTING GUIDE

---

# TROUBLESHOOTING GUIDE

Device switches off suddenly	The device will switch off when the timer runs out. Press the "Power" button to continue using it. If the device has run for 10 mins, restart after resting the unit for 15 mins.
	Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 mins before restarting
Device fail to start	Check if device has been used for more than 10 mins. If so, let device rest for 15 mins before continuing. Device may require longer time to cool in warmer environments
	The battery is low, please charge in time
Was too small	The Size of this product is 14.6 x 4.3 x 7.5 inch, this neck massage pillow is relatively small in size, please choose to purchase according to your actual situation.
Massage feels too hard	Add the towel to the neck massager pillow before using it

Image: A detailed troubleshooting guide table for common issues with the massager.

Issue	Solution
Device switches off suddenly	The device will switch off when the timer runs out (10 minutes). Press the "Power" button to continue using it. If the device has run for 10 minutes, restart after resting the unit for 15 minutes. Over-heat protection function activated: Device is switched off for your safety. Please rest device for 15-30 minutes before restarting.
Device fails to start	Check if device has been used for more than 10 minutes. If so, let device rest for 15 minutes before continuing. Device may require longer time to cool in warmer environments. The battery is low, please charge in time.
Massage feels too hard	Add a towel to the neck massager pillow before using it to soften the massage intensity.

## SPECIFICATIONS

- **Product Dimensions:** 14.69 x 4.72 x 7.87 inches
- **Item Weight:** 3.99 Pounds
- **Item Model Number:** CF-6419
- **Batteries:** 1 AA batteries required (internal rechargeable)
- **Power Source:** Battery Powered (Cordless)
- **Material:** Exterior material 100% Polyester
- **Brand:** COMFIER
- **Color:** Royal Blue

## WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided on the product packaging or visit the official COMFIER website. Keep your purchase receipt for warranty claims.



## Related Documents - CF-6419

 <p>User Manual</p>  <p>CF-6701 DEEP-KNEADING SHIATSU NECK &amp; SHOULDER MASSAGER</p>	<p><a href="#">Comfier CF-6701 Deep-Kneading Shiatsu Neck &amp; Shoulder Massager User Manual</a></p> <p>User manual for the Comfier CF-6701 Deep-Kneading Shiatsu Neck &amp; Shoulder Massager, covering product features, technical data, setup, operation, safety instructions, care, troubleshooting, and warranty information.</p>
 <p>Shiatsu Neck and Back Massager</p>  <p>User Manual CF-2309A</p>	<p><a href="#">Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</a></p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
 <p>Shiatsu neck-and-back-massage DE Shiatsu neck and back massager GB Massager de cou et de dos Shiatsu FR Massaggiatore shiatsu per collo e schiena IT Masajeador shiatsu de para cuello y espalda ES Shiatsu neck- en rugmassagetoerust. NL</p>  <p>CE User Manual CF-2309A</p>	<p><a href="#">Comfier CF-2309A Shiatsu Neck and Back Massager User Manual</a></p> <p>Comprehensive user manual for the Comfier CF-2309A Shiatsu Neck and Back Massager, covering safety instructions, features, operation, maintenance, technical specifications, and warranty information.</p>
 <p>Shiatsu massage pillow with heat</p>  <p>CE User Manual CF-1228</p>	<p><a href="#">Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</a></p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>
 <p>shiatsu neck and shoulder massager with heat GB neck &amp; shoulder shiatsu massager with heat GB masajeador shiatsu para el cou et les épaules avec chaleur FR massaggiatore shiatsu per collo e spalle con calore IT masajeador de cuello y hombros shiatsu con calor ES</p>  <p>CE User Manual CF-6302</p>	<p><a href="#">COMFIER CF-6302 Shiatsu Nacken- und Schultermassagegerät mit Wärmefunktion – Bedienungsanleitung</a></p> <p>Entdecken Sie die COMFIER CF-6302 Bedienungsanleitung für das Shiatsu Nacken- und Schultermassagegerät mit Wärmefunktion. Erfahren Sie mehr über die Funktionen, Anwendung und Pflege dieses entspannenden Massagegeräts für Nacken und Schultern.</p>
 <p>Cordless Neck Massager</p>  <p>User Manual CF-6026</p>	<p><a href="#">Comfier Cordless Neck Massager CF-6026 User Manual</a></p> <p>Comprehensive user manual for the Comfier Cordless Neck Massager (Model CF-6026), detailing safety instructions, features, technical specifications, usage guidelines, troubleshooting tips, care and maintenance procedures, and warranty information.</p>

