

Taylor 755841932CH

Taylor Digital Body Weight Scale Model 755841932CH Instruction Manual

INTRODUCTION

Welcome to your new Taylor Digital Body Weight Scale, Model 755841932CH. This scale is designed to provide accurate and consistent body weight measurements with ease of use. Featuring a sleek design, instant-on technology, and a large display, it is a reliable tool for monitoring your weight.

IMPORTANT SAFETY INFORMATION

- Always place the scale on a hard, flat, and stable surface. Uneven or soft surfaces like carpets can affect accuracy.
- The scale surface is glass and can be slippery when wet. Ensure the surface is dry before stepping on.
- Do not jump on the scale or use it for purposes other than body weight measurement.
- Keep the scale away from extreme temperatures and humidity.
- This scale is intended for household use only.

PACKAGE CONTENTS

- Taylor Digital Body Weight Scale (Model 755841932CH)
- 2 x AAA Batteries

SETUP

1. Battery Installation

Before first use, open the battery compartment located on the underside of the scale. Insert the two included AAA batteries, ensuring correct polarity (+/-). Close the compartment securely.



This image displays the underside of the scale, highlighting the battery compartment and the four anti-skid foot pads.

2. Initial Placement

Place the scale on a hard, flat, and stable surface. Optimal surfaces include tile, hardwood, or concrete. Avoid placing the scale on carpets, rugs, or uneven flooring, as this can lead to inaccurate readings.



The scale is shown positioned on a hard, flat wooden floor, demonstrating proper placement for accurate readings.

3. Unit Selection

On the underside of the scale, locate the unit switch. Slide the switch to select your preferred unit of measurement: pounds (LB), kilograms (KG), or stones (ST).

OPERATING INSTRUCTIONS

1. Taking a Measurement


Gently step onto the center of the scale platform. The scale features instant-on technology, so it will automatically activate and begin measuring. Remain still until your weight is displayed and stabilizes. The reading will flash briefly, then lock.



This image shows the top surface of the scale with a digital display indicating a weight of 134.6 pounds.

2. Understanding the Display

The scale is equipped with a large 3.2" x 1.5" LCD display for clear readability. After displaying your weight, the scale will automatically turn off after a few seconds of inactivity to conserve battery life.

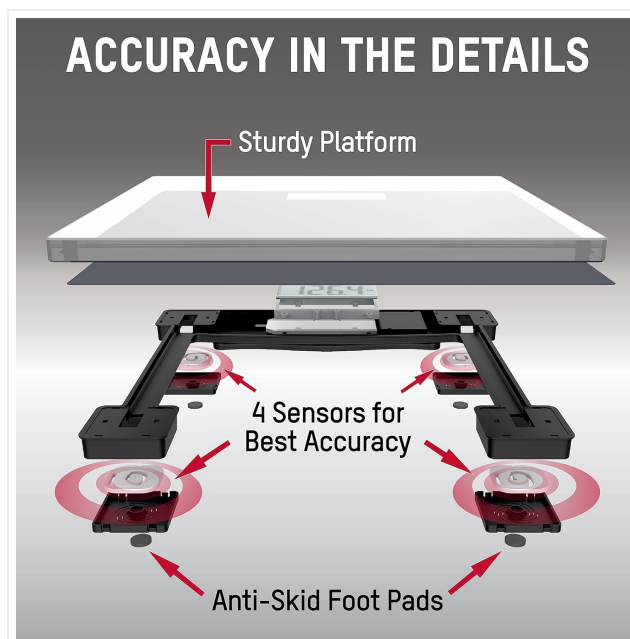

A graphic highlights the dimensions of the scale's large LCD display, measuring 3.2 inches by 1.5 inches, designed for easy readability.

MAINTENANCE

1. Cleaning

To clean your scale, wipe the surface with a damp cloth and a mild glass cleaner or detergent. Avoid using abrasive cleaners, harsh chemicals, or submerging the scale in water, as this can damage the electronic components.

2. Battery Replacement

If the display shows a "Lo" indicator or becomes dim, it is time to replace the batteries. Open the battery compartment on the underside of the scale and replace the old batteries with two new AAA batteries,

observing correct polarity.

3. Storage

Store the scale in a dry, cool place away from direct sunlight and extreme temperatures. Do not store anything on top of the scale when not in use, as this can affect its calibration over time.

TROUBLESHOOTING

1. Inaccurate or Inconsistent Readings

- **Surface:** Ensure the scale is placed on a hard, flat, and stable surface. Soft surfaces like carpet or uneven tiles can cause significant inaccuracies.
- **Recalibration:** If you move the scale, or if readings seem inconsistent, perform a quick recalibration. Step onto the scale briefly, step off, wait for the display to show "0.0", then step on again for your measurement.
- **Batteries:** Low batteries can affect accuracy. Replace batteries if the "Lo" indicator appears or if readings are erratic.
- **Movement:** Ensure you stand still on the scale during measurement. Any movement can cause fluctuations.

2. Display Not Working or Showing "Lo"

- **Batteries:** Check if the batteries are correctly installed with the proper polarity. If the display shows "Lo", replace the 2 AAA batteries with new ones.
- **Contact:** Ensure there is no obstruction preventing proper contact between the battery and the terminals.

SPECIFICATIONS

Feature	Detail
Model Number	755841932CH
Capacity	400 lb (180 kg / 28 st)
Increments	0.2 lb (0.1 kg)
Display	3.2" x 1.5" LCD
Power	2 x AAA Batteries (included)
Dimensions	11.8 x 11.8 x 1.04 inches
Special Features	Instant-On, Auto-Off, Anti-Skid Foot Pads



This diagram illustrates the internal components of the scale, specifically highlighting the four sensors for accurate measurement and the anti-skid foot pads for stability.

WARRANTY AND SUPPORT

Warranty Information

This Taylor Digital Body Weight Scale, Model 755841932CH, is covered by a **5-Year Limited Warranty**. Please retain your proof of purchase for any warranty claims. The warranty covers defects in materials and workmanship under normal use.

Customer Support

For further assistance, troubleshooting tips, or to inquire about warranty service, please visit the official Taylor website or contact their customer service department. You can also find more information and products at the official Taylor store on Amazon:

[Visit the Taylor Store on Amazon](#)