



Manuals.plus /

› TYHGSF /

› IPL Laser Hair Removal Device User Manual

## TYHGSF 04B-1

# IPL Laser Hair Removal Device User Manual

Brand: TYHGSF | Model: 04B-1

## 1. INTRODUCTION

Thank you for choosing the TYHGSF IPL Laser Hair Removal Device. This device utilizes advanced Intense Pulsed Light (IPL) technology to provide a safe, effective, and convenient solution for long-term hair reduction in the comfort of your home. Please read this manual thoroughly before first use to ensure proper operation and safety.



The TYHGSF IPL Laser Hair Removal Device, shown with included protective glasses and a razor.

## 2. SAFETY INFORMATION

**IMPORTANT: This device is NOT suitable for all skin and hair types. Please refer to the skin and hair color applicability chart before use.**

- Always wear the provided protective glasses during operation to shield your eyes from the intense light flashes.
- Do not use on tattooed areas, moles, birthmarks, or areas with skin diseases.
- Avoid direct exposure to the eyes.
- Do not use on sensitive areas such as nipples, genitals, or around the eyes.

- Keep out of reach of children.
- Consult a doctor before use if you have any pre-existing medical conditions or are taking medications that may increase light sensitivity.
- Do not use if you are pregnant or breastfeeding.



Key steps for safe and effective use of the IPL device.

### 3. PRODUCT FEATURES

- **Permanent Hair Reduction:** Utilizes Intense Pulsed Light (IPL) technology to target hair follicles and inhibit regrowth.
- **High Flash Capacity:** Offers 999,999 flashes for long-term use, eliminating the need for replacement cartridges.
- **5 Adjustable Power Levels:** Allows customization of intensity to suit different skin sensitivities and treatment areas. Start at level 1 and gradually increase.
- **Dual Operating Modes:**
  - *Manual Mode:* Ideal for small, precise areas like the bikini line, armpits, and upper lip.
  - *Automatic Mode:* Suitable for larger body areas such as legs, arms, and back, providing continuous flashes.
- **Ergonomic Design:** Comfortable to hold and easy to maneuver for effective treatment.
- **Economical Home Use:** Provides salon-quality hair removal results from the convenience of your home, saving time and money.

# NEWEST IPL HAIR REMOVAL TECHNOLOGY



Manual &  
Auto Modes



5 Energy  
Levels



Permanent Hair  
Removal



Eyes  
Protection



Safe &  
Painless



Overview of the IPL device's core features.

# FREE SWITCHING BETWEEN TWO FLASH MODES

## Auto Mode



legs



Arms



Back

Suitable for large areas like arms legs and back

## Manual Mode



Armpit



Lips



Bikini line

Suitable for small areas like facial armpits and bikini area



Switching between Automatic and Manual flash modes for different body areas.

## 4. PACKAGE CONTENTS

Upon opening your TYHGSF IPL Laser Hair Removal Device package, you should find the following items:

- 1 x TYHGSF IPL Laser Hair Removal Device

- 1 x Power Adapter
- 1 x Protective Glasses
- 1 x Razor
- 1 x User Manual (this document)

## 5. GETTING STARTED (SETUP)

1. **Prepare the Skin:** Ensure the treatment area is clean, dry, and free of any lotions, oils, or deodorants. Shave the hair in the area you intend to treat. Do not wax or pluck, as IPL targets the hair follicle.
2. **Connect Power:** Plug the power adapter into the device and then into a power outlet. The device will power on, and the display will light up.
3. **Wear Protective Glasses:** Always put on the provided protective glasses before operating the device to protect your eyes from the intense light flashes.
4. **Select Energy Level:** Press the power button (usually also the On/Off button) to adjust the energy level. Start with level 1 for your first treatment and gradually increase to a comfortable yet effective level as your skin adapts.

## 6. OPERATING INSTRUCTIONS

### How to Use:

1. After preparing the skin and selecting the energy level, place the flash window of the device firmly against your skin, ensuring full contact.
2. For Manual Mode (small areas): Press the flash button once to emit a single flash. Lift the device and move to the next adjacent area.
3. For Automatic Mode (large areas): Press and hold the flash button until the device starts flashing continuously. Glide the device slowly over the skin, ensuring each area is treated once.
4. Treat each area 1-2 times per session. Avoid overlapping flashes excessively.


### Recommended Treatment Schedule:

Consistency is key for optimal results. The hair growth cycle typically requires several sessions.

- **Weeks 1-4:** Use 2-3 times per week. Hair growth will begin to slow down.
- **Weeks 5-12:** Use 2 times per week. Hair will become thinner and sparser.
- **After 12 Weeks:** Use once per week, then transition to once a month for maintenance as needed. Visible results are typically seen after 4 weeks, with permanent hair reduction achieved in 8 weeks.

# LASER HAIR REMOVAL

Remove unwanted hair safely and effectively

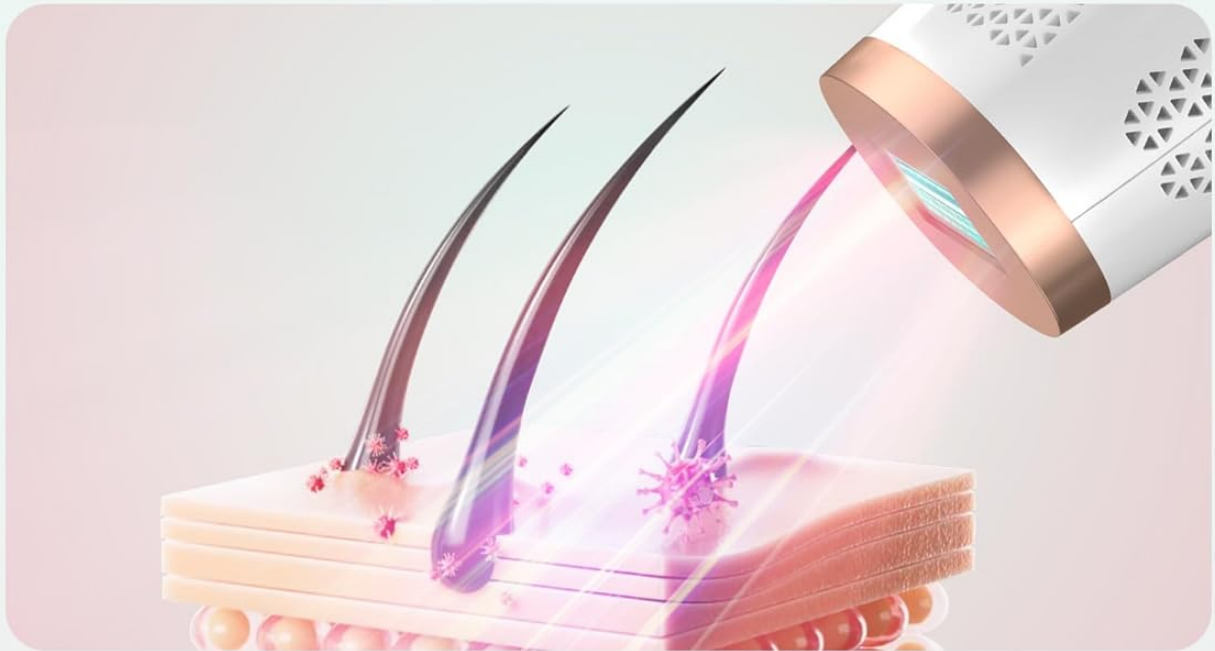

**1-4 Weeks**  
2-3 times per week

Partially dormant



**5-12 Weeks**  
2 times per week

Hair growth slowly



**12 Weeks later**  
Once per week

Skin becomes smoother,  
Continued monthly  
treatment as needed.

The hair removal process and expected results timeline.

## Applicable Body Areas:

The device is suitable for various body parts:

- Face (e.g., upper lip, chin)
- Armpits
- Arms
- Back

- Legs
- Bikini Line

# WHOLE BODY HAIR REMOVAL



Recommended treatment areas and approximate times.

## 7. SKIN AND HAIR SUITABILITY

IPL technology works by targeting melanin in the hair follicle. Therefore, it is most effective on darker hair colors and lighter skin tones. It is **not effective** on red, light blonde, grey, or white hair, as these hair colors contain insufficient melanin to absorb the light energy.

The device is also **not suitable** for very dark skin tones (Type VI on the Fitzpatrick scale) as it may cause skin damage due to high melanin content in the skin.

# SKIN AND HAIR COLOR APPLICABLE

Comprehensive care, adapt to different skin defects



## Skin Color



White


Rice White

Cream Color

Light Brown

Brown

## Hair Color

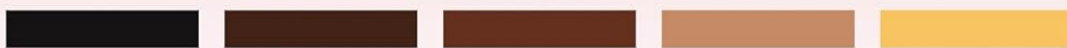

Black

Dark Brown

Brown

Light Brown

Yellow

Skin and hair color applicability chart. Consult this chart before use.

Hair Color	Skin Color					
	White	Lvory	Beige	Light Brown	Brown	Black
White/Grey	X	X	X	X	X	X
Red	X	X	X	X	X	X
Light Blonde	X	X	X	X	X	X
Dark Blonde	4/5	4/5	4/5	1/2/3	X	X
Brown	4/5	4/5	4/5	1/2/3	X	X
Dark Brown	4/5	4/5	4/5	1/2/3	X	X
Black	4/5	4/5	4/5	1/2/3	X	X

## IPL HAIR REMOVAL

Leaves you with silky smooth skin

(Note: X is not applicable)

Note: If you are not sure your skin types please refer to our skin tone table before using

Note: Your skin may react differently in different exposure times and occasions

Detailed compatibility chart for skin and hair colors. 'X' means not applicable.

## 8. MAINTENANCE

- **Cleaning:** After each use, ensure the device is unplugged and cool. Gently wipe the flash window with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.

## 9. TROUBLESHOOTING

If you encounter any issues with your device, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power outlet issue.	Ensure the power adapter is securely connected to the device and the wall outlet. Try a different outlet.
No flash emitted.	Flash window not in full contact with skin; device in standby mode.	Ensure the flash window is pressed firmly and flat against the skin. Press the power button to wake the device from standby.
Skin irritation after use.	Energy level too high for your skin type; improper use.	Reduce the energy level for subsequent treatments. Ensure you are not treating the same area multiple times in one session. Apply a soothing lotion if needed.
No visible hair reduction after several sessions.	Inconsistent use; hair/skin type not suitable; energy level too low.	Ensure consistent use as per the recommended schedule. Verify your skin and hair type against the suitability chart. Gradually increase the energy level if comfortable.

If the problem persists, please contact customer support.

## 10. SPECIFICATIONS

- **Manufacturer:** TYHGSF
- **Model Reference:** 04B-1
- **ASIN:** B0DQHFSMTG
- **EAN:** 0786463416663
- **Flash Capacity:** 999,999 Flashes
- **Power Levels:** 5
- **Power Source:** Corded Electric (Does not operate on battery)
- **Package Dimensions:** 35.9 x 27.7 x 10.3 cm
- **Item Weight:** 770 g

## 11. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact the seller directly through the platform where the product was purchased. Keep your proof of purchase for any warranty claims.