

[Manuals.plus](#) /

› [KIWIBIRD](#) /

› KIWIBIRD Water Dental Flosser H1-HB Instruction Manual

KIWIBIRD H1-HB

KIWIBIRD Water Dental Flosser H1-HB Instruction Manual

Model: H1-HB | Brand: KIWIBIRD

1. IMPORTANT SAFETY INFORMATION

WARNING: To reduce the risk of electric shock, do not immerse the unit in water or other liquid. Do not use while bathing. Do not reach for a unit that has fallen into water. Unplug immediately.

This product is not a toy; keep out of reach of children. Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Always unplug the product before cleaning or if it will not be used for an extended period.

2. PRODUCT OVERVIEW

The KIWIBIRD Water Dental Flosser H1-HB is a cordless oral irrigator designed for enhanced oral hygiene. It features advanced UVC sanitization, multiple cleaning modes, and a large water tank for convenient use.



Image: KIWIBIRD Water Dental Flosser H1-HB, showcasing its sleek design and main components.

Key Features:

- **Advanced 360° UV-C Sanitization:** Enhances cleanliness of components with a 90-second timer.
- **Deep Cleaning:** Delivers 1800 ultra-fine water pulses per minute (40-120 PSI) for plaque removal and gum

health.

- **Personalized Flossing:** 3 pressure modes (gentle, normal, strong) and 4 interchangeable 360° rotating tips.
- **Upgraded BPA-Free 300ML Water Tank:** Supports uninterrupted flossing with leakproof design and dual water inlets.
- **Portable Design & Long Battery Life:** 1800mAh rechargeable battery provides up to 70 days of use on a single charge, with USB compatibility.
- **IPX7 Waterproof:** Safe for use in the shower, with an ergonomic, textured handle and built-in filter.

3. SETUP

1. **Charge the Device:** Before first use, fully charge the flosser using the provided USB cable. The charging indicator will show progress.
2. **Attach the Nozzle:** Select the desired irrigator tip from the 4 interchangeable tips. Insert the tip into the hole on the top of the device until it clicks into place.
3. **Fill the Water Tank:** Twist the water tank counter-clockwise to detach it from the main body. Fill the 300ML water tank with fresh water. Ensure the tank is securely attached to the device by twisting clockwise until it locks.



Image: Visual guide for adding water, inserting the tip, selecting mode, and turning on the device.

4. OPERATING INSTRUCTIONS

1. **Select Mode:** Press the 'Mode' button to cycle through the 3 pressure modes: Gentle, Normal, and Strong. Choose the mode suitable for your comfort and cleaning needs.
2. **Position the Tip:** Lean over a sink. Place the irrigator tip into your mouth, directing it towards your teeth and gumline at a 90-degree angle.
3. **Start Flossing:** Press the 'Power' button to turn on the device. The water stream will begin.
4. **Clean Your Teeth:** Guide the water stream along your gumline, pausing briefly between teeth. Move the tip in a circular motion around each tooth, ensuring to cover all areas.
5. **Turn Off:** Press the 'Power' button again to turn off the device when finished.

Your browser does not support the video tag.

Video: A demonstration of the KIWIBIRD Water Dental Flosser H1 in use, showing how to operate the device and its

5. MODES AND TIPS

Pressure Modes:

- **Gentle Mode:** Ideal for beginners or those with sensitive gums.
- **Normal Mode:** For regular daily cleaning.
- **Strong Mode:** Provides deep cleaning with high water pressure for effective plaque removal.




Image: Illustration of the four professional modes available on the KIWIBIRD Water Dental Flosser, including Sensitive, Standard, Strong, and UV Sanitization Mode.

Interchangeable Tips:

- **Standard Tip:** For daily dental hygiene needs.
- **Sensitive Tip:** Designed for those with sensitive teeth.
- **Orthodontic Tip:** Effectively cleans around braces, bridges, and other dental work.

- **Tongue Scraper:** Helps remove tongue coating for fresher breath.

4 Different Flossing Tips

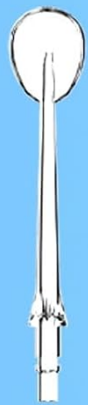

Tongue Scraper
Removes tongue coating



Sensitive Tip
For those have sensitive teeth



Standard Tip
For daily dental hygiene needs



Orthodontic Tip
Cleans braces, bridges, and remnant



Image: Depiction of the four different flossing tips: Tongue Scraper, Sensitive Tip, Standard Tip, and Orthodontic Tip.

6. UVC SANITIZATION

The KIWIBIRD Water Dental Flosser features an integrated UVC light for enhanced hygiene. This technology sanitizes the water tank, internal pipes, and nozzles.

1. **Activate Sanitization:** After use, ensure the water tank is empty or contains clean water. Press the 'Sterilization' button.
2. **Automatic Cycle:** The UVC light will activate for a 90-second cycle, effectively eradicating 99.9% of bacteria.
3. **Completion:** The device will automatically turn off once the sanitization cycle is complete.

Sanitization UV-C Light

Water tank + Jet tips, with 360° all-round sanitization!



Image: Close-up of the water flosser's tank with a blue UVC light active, illustrating the sanitization process.

Your browser does not support the video tag.

Video: A detailed look at the UVC light feature of the KIWIBIRD Water Dental Flosser, explaining its sanitization benefits.

7. CHARGING AND BATTERY LIFE

The KIWIBIRD Water Dental Flosser is equipped with an 1800mAh rechargeable battery, offering extended use on a single charge.

- **Long-lasting Power:** Enjoy up to 70 days of use on a single charge (based on sensitive mode, once a day, consuming one tank of water each time).
- **Convenient Charging:** The device is USB compatible, allowing you to charge it from various sources such as power banks, computers, or power adapters.
- **Charging Indicator:** The charging indicator lights will show the battery status and charging progress.

70 Days Runtime

Use in a single charge



Fast Charging



Standby Time



POWERBANK



COMPUTER



POWER ADAPTER



Image: Visual representation of the KIWI BIRD Water Dental Flosser's 70-day runtime and various USB charging options.

8. MAINTENANCE AND CLEANING

Regular cleaning and maintenance ensure the longevity and hygienic operation of your water flosser.

- **After Each Use:** After use, rinse the irrigator tip and the water tank with fresh water.
- **Tank Cleaning:** The extra-large water tank features a leakproof drainage design and dual water inlets to prevent water buildup, keeping it fresh and clean. It is made with BPA-free materials for safety.
- **Nozzle Storage:** Store the device and nozzles in a dry place.
- **Filter Maintenance:** The built-in filter prevents clogging and ensures long-lasting performance. Periodically check and clean the filter if necessary.

300ML Larger

Designed for easy filling & cleaning

Bpa Accepted safe and secure



Image: Close-up of the KIWIBIRD Water Dental Flosser's 300ML water tank, highlighting its design for easy filling and cleaning.

9. TROUBLESHOOTING

If you encounter any issues with your KIWIBIRD Water Dental Flosser, please refer to the following common solutions:

- **No Water Flow:** Ensure the water tank is properly filled and securely attached. Check if the nozzle is correctly inserted.
- **Weak Water Pressure:** Verify that the selected mode is appropriate. Ensure the battery is sufficiently charged.
- **Device Not Turning On:** Check the battery level and charge the device if needed. Ensure the power button is pressed firmly.
- **UVC Light Not Activating:** Confirm the 'Sterilization' button is pressed. Ensure the device has sufficient power.

If problems persist, please contact KIWIBIRD customer support for assistance.

10. SPECIFICATIONS

Feature	Detail
Model Number	H1-HB
Product Dimensions	1.57 x 3.15 x 8.66 inches; 9.35 ounces
Item Weight	265 Grams
Battery	1 P76 battery required (included)
Power Source	Battery Powered
Tank Volume	300 Milliliters
Water Pressure Range	40-120 PSI
Water Pulses	1800 per minute
Waterproof Rating	IPX7
Material Type	Plastic
Color	Black & White
Certifications	CE, FCC, FDA, RoHS

11. WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the product packaging or contact KIWIBIRD customer service directly. Keep your purchase receipt for any warranty claims.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question