



Manuals.plus /

› INIA /

› INIA X3B002 3-in-1 Skin Tightening Facial Device Instruction Manual

## INIA X3B002

# INIA X3B002 3-in-1 Skin Tightening Facial Device Instruction Manual

Model: X3B002

## 1. INTRODUCTION

The INIA X3B002 3-in-1 Skin Tightening Device is designed to enhance your skincare routine with its multi-mode functionality. This device combines advanced technologies to help improve skin appearance, promote natural facial contours, and enhance the absorption of skincare products. It is suitable for use on the face and neck.

Key features include:

- **3-in-1 Multi-Mode Design:** Features LIFT, TIGHTEN, and INFUSE modes with adjustable intensity.
- **LIFT Mode:** Utilizes microcurrents and red light therapy to help reshape and tone facial muscles, promoting natural contours.
- **TIGHTEN Mode:** Employs red light technology and thermal-frequency technology to help reduce the appearance of fine lines and improve skin firmness and elasticity.
- **INFUSE Mode:** Incorporates negative ion technology to enhance the absorption and effectiveness of your skincare products.
- **Ergonomic Design:** Features a flexible 3D head for consistent application across various facial contours.

## 2. SAFETY PRECAUTIONS

Please read all instructions carefully before using the INIA X3B002 device. Keep this manual for future

reference.

- Do not use the device on broken, irritated, or infected skin.
- Avoid direct contact with eyes.
- Consult a physician before use if you have any medical conditions, are pregnant, or have electronic implants.
- Keep out of reach of children.
- Do not immerse the device in water. Clean with a damp cloth as instructed.
- Use only with the recommended conductive gel or essence to ensure proper function and avoid skin irritation.
- Discontinue use if you experience any discomfort or adverse reactions.

### 3. PRODUCT COMPONENTS

Your INIA X3B002 package includes:

- INIA X3B002 3-in-1 Skin Tightening Device
- Hydrating Conductive Gel (300g)
- Charging Cable (USB-C)
- User Manual

# Full Skin Rejuvenation From Within

**95%** Reported a Noticeable Lift



**97%** Experienced Fuller, More Defined Cheeks



**98%** Achieved Radiant, Luminous Skin



Image: The INIA X3B002 device and its accompanying hydrating conductive gel.

## 4. DEVICE OVERVIEW

The INIA X3B002 features an intuitive design for ease of use. It includes a multi-functional head with 9 advanced zinc-alloy electrodes and a flexible 3D head for optimal skin contact.

### Modes and Indicators:

- **LIFT Mode:** Indicated by purple light. Focuses on contouring and muscle toning.
- **TIGHTEN Mode:** Indicated by red light. Targets skin firmness and elasticity.
- **INFUSE Mode:** Indicated by blue light. Enhances product absorption.
- **Intensity Levels:** Three adjustable intensity levels for each mode, indicated by LED lights.

### TIGHTEN Mode

Enhance Firmness:  
Unlock firmer, more elastic skin  
with our Thermal-Frequency Technology  
designed to tighten and  
refine your complexion.

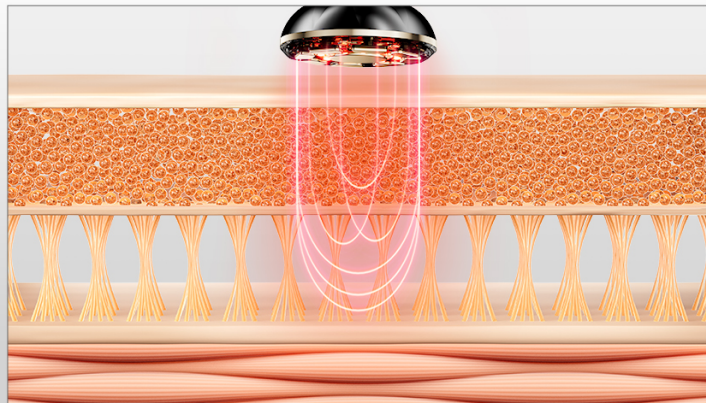


Image: The device head with 9 advanced zinc-alloy electrodes and flexible 3D design for efficient penetration and perfect fit.

## 5. GETTING STARTED (SETUP)

Before your first use, ensure the device is fully charged. A full charge typically takes 2-3 hours.

### Preparation Steps:

1. **Cleanse Skin:** Thoroughly clean and dry your face and hands. Ensure your skin is free from any makeup, dirt, or other products.
2. **Apply Conductive Gel:** Squeeze a sufficient amount of INIA Hydrating Conductive Gel onto your palm. Apply the gel evenly to the area of your face or neck you intend to treat. The gel is essential for optimal conductivity and comfort during treatment.



Image: Visual guide for preparing your skin and applying the conductive gel before using the device.

## 6. OPERATING INSTRUCTIONS

The INIA X3B002 device offers three distinct modes to target different skincare needs. Press and hold the Power/Mode button for 1 second to switch on the device and select your desired treatment mode.

### 6.1. LIFT Mode (Purple Light)

This mode helps to elevate tired-looking skin and promote the appearance of natural facial contours through microcurrents and red light therapy.

- **Duration:** 10-minute smart partition (4 minutes for left cheek, 4 minutes for right cheek, 2 minutes for forehead).
- **Application:** Gently glide the device upwards along your jawline, cheekbones, and forehead. Follow the natural contours of your face.



Image: Detailed illustration of LIFT Mode's action on skin layers.

### 6.2. TIGHTEN Mode (Red Light)

This mode employs red light technology and thermal-frequency technology to help reduce the appearance of fine lines and wrinkles, promoting a smoother and firmer complexion.

- **Duration:** 10-minute smart partition (4 minutes for left cheek, 4 minutes for right cheek, 2 minutes for forehead).
- **Application:** Move the device in small, circular motions or short, upward strokes over areas with fine lines and areas needing firmness.



Image: Detailed illustration of TIGHTEN Mode's action on skin layers.

### 6.3. INFUSE Mode (Blue Light)

Features negative ion technology to help enhance the absorption of your favorite serums and moisturizers, maximizing their effectiveness.

- **Duration:** 4-minute symmetrical treatment (2 minutes on each cheek).
- **Application:** After applying your essence or serum, gently glide the device over the treated areas to aid absorption.

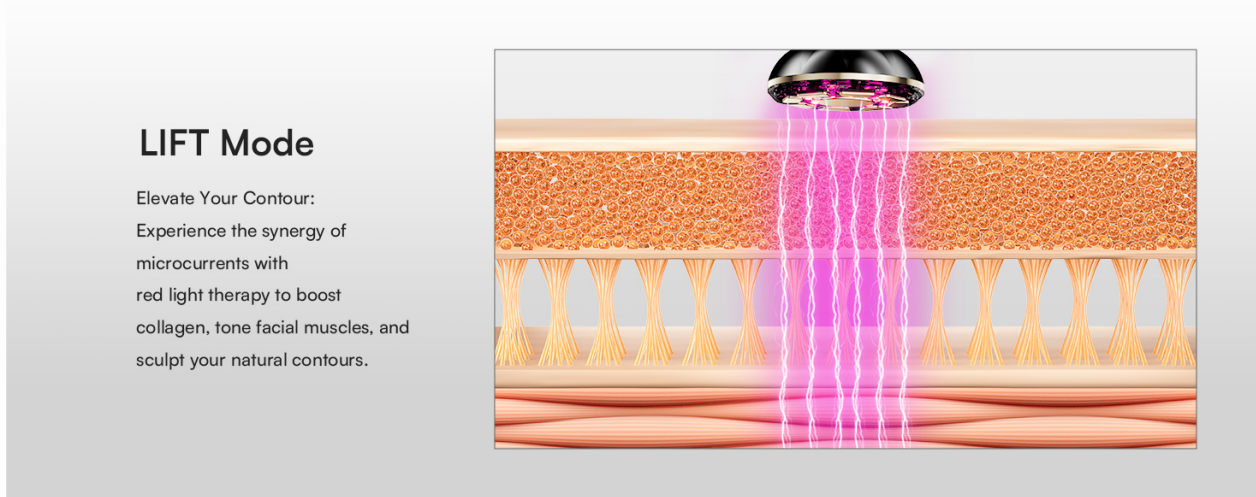


Image: Detailed illustration of INFUSE Mode's action on skin layers.

### 6.4. Skincare Routine Suggestion

It is recommended to dynamically adjust the frequency of use according to your skin condition (2-5 times/week). The total duration of a single treatment should not exceed 15 minutes (including preparation time).

A suggested weekly routine:

Day	Recommended Mode(s)
Monday	LIFT + INFUSE
Tuesday	TIGHTEN + INFUSE
Wednesday	LIFT + INFUSE
Thursday	TIGHTEN + INFUSE
Friday	LIFT + INFUSE
Saturday	TIGHTEN + INFUSE
Sunday	Rest Day

Note: It is advised not to use LIFT and TIGHTEN modes on the same day to prevent dry skin.

## Your Everyday At-Home Beauty Salon

Three Times a Week  
Boosting Collagen and Elastin, Restoring Skin Elasticity

**Visible Results:**

▲ **+102.87%**

Cheek Muscle  
Volume Increased

▲ **+133.13%**

Mandibular Line  
Clarity Enhanced

▲ **+225.62%**

Cheeks Tightened

▲ **-92.14%**

Wrinkles Reduced

\*Data Source: Scientific trials by INIA Lab, 127 participants, 4 weeks.



Image: A visual representation of the suggested weekly skincare routine.

## Official Product Video:

Your browser does not support the video tag.

Video: Official INIA FLARE demonstration, showcasing the device and its use with the 300g gel.

## 7. MAINTENANCE AND CARE

Proper cleaning and storage will ensure the longevity and effectiveness of your INIA X3B002 device.

- **After Each Use:** Thoroughly clean the electrode head with a soft, damp cloth or a wet paper towel after turning off the device.
- **Drying:** Ensure the device is completely dry before storing.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture.
- **Avoid Harsh Chemicals:** Do not use alcohol, abrasive cleaners, or harsh chemicals to clean the device, as this may damage its surface or internal components.

## 8. TROUBLESHOOTING

If you encounter any issues with your INIA X3B002 device, please refer to the following common solutions:

- **Device not turning on:** Ensure the device is fully charged. Press and hold the Power/Mode button for at least 1 second.
- **No sensation during use:** Ensure sufficient conductive gel is applied to the skin. Increase the intensity level if needed.
- **Skin irritation:** Discontinue use immediately. Ensure you are using the recommended conductive gel and that your skin is clean before treatment. Consult a dermatologist if irritation persists.
- **Lights not working:** Check if the device is powered on and fully charged. If the issue persists, contact customer support.

## 9. SPECIFICATIONS

Feature	Detail
Model Number	X3B002
Brand	INIA
Package Dimensions	10.28 x 7.4 x 3.43 inches
Item Weight	2.23 Pounds
Color	Black
Material	Plastic
Style	Modern, Sleek, Luxury
Product Benefits	Skin Care (Lifting, Tightening, Infusion)

## 10. WARRANTY AND SUPPORT

The INIA X3B002 device comes with a manufacturer's warranty. For detailed warranty information, please refer to the product packaging or contact INIA customer support.

For any questions, technical support, or service inquiries, please visit the official INIA website or contact our customer service team. Keep your purchase receipt as proof of purchase for warranty claims.

