

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [COSLUS](#) /

› [COSLUS Advanced Water Flosser Model E2 Instruction Manual](#)

COSLUS E2

COSLUS Advanced Water Flosser Model E2 Instruction Manual

Model: E2

INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your COSLUS Advanced Water Flosser Model E2. Designed for advanced oral hygiene, this device offers customizable cleaning experiences for various dental needs, including orthodontic care, braces, and sensitive gums. Please read this manual thoroughly before operating the device.

SAFETY INFORMATION

Keep out of reach of children.

PACKAGE CONTENTS

- COSLUS Advanced Water Flosser Main Unit
- 5 Jet Tips (Standard, Orthodontic, Tongue Cleaner, Periodontal)
- Patented Nozzle Storage Box
- Power Adapter
- User Manual

PRODUCT OVERVIEW

The COSLUS Advanced Water Flosser Model E2 is engineered for superior oral cleaning. It features advanced high-frequency pulses (10-160 PSI, 1400 pulses/min) to effectively remove dental plaque. With 2 distinct modes (Floss and Massage) and 12 adjustable pressure settings, it caters to a wide range of users from ages 6 to 60, including those with braces or sensitive gums. The device includes a large 600ml detachable water tank for extended use and a patented nozzle storage box for hygienic organization of the included jet tips.



Figure 1: COSLUS Advanced Water Flosser Model E2 with its accessories.

SETUP INSTRUCTIONS

1. **Plug in the Device:** Connect the power adapter to the device and plug it into a suitable electrical outlet. The device supports 110-220V.
2. **Fill the Water Tank:** Detach the 600ml water tank from the base. Fill it with warm water (below 40°C / 104°F) for comfort. Reattach the tank securely to the base.
3. **Install the Nozzle:** Select the desired jet tip. Insert the nozzle into the hole at the top of the handle until you hear a "click," ensuring it is securely in place without any gaps.
4. **Select Mode and Pressure:** Press the "ON/OFF" button on the main unit to power on. Choose between 'Floss' or 'Massage' mode using the 'M' button. Adjust the water pressure to your desired level using the dial (1-12). For first-time users or those with sensitive teeth, it is recommended to start with a lower setting (gear 1-2).

Home Comforts, Travel Essentials



Figure 2: Step-by-step setup and initial adjustment guide.

Video 1: Demonstrates the setup process, including plugging in, filling the tank, and installing a nozzle.

OPERATING INSTRUCTIONS

1. **Position the Nozzle:** Lean over the sink. Place the nozzle at a 90-degree angle to your gum line, ensuring your lips are slightly closed around the nozzle to prevent splashing while allowing water to drain into the sink.
2. **Start Water Flow:** Slide the ON/OFF switch on the handle to the 'ON' position to begin water flow.
3. **Guide Along Gumline:** Guide the water jet along your gumline and between teeth, pausing briefly between each tooth.
4. **Turn Off:** After cleaning, slide the ON/OFF switch on the handle to the 'OFF' position. Press the main unit's power button to turn off the device completely.

Video 2: Overview of the water flosser's features and basic operation.

MAINTENANCE

- **Empty and Clean Tank:** After each use, empty any remaining water from the tank. The 600ml detachable tank is designed for easy cleaning.
- **Clean Nozzles:** Rinse the jet tips under running water. The patented nozzle storage box features quick-drain ventilation to ensure tips remain dry and clean.
- **Store Properly:** Store the main unit and nozzles in a dry, clean place. The compact design with a sturdy upright stand helps keep your bathroom organized.

Easy to Use, Stays Dry & Clean

Removable and washable vertical storage design



Figure 3: Patented nozzle storage box for dry and clean accessories.

TROUBLESHOOTING

Phenomenon: You can hear the machine running, but no water comes out.

Possible Cause: The machine is on, water pressure level is above 6th level, and the handle is turned off.

Reason: This is the automatic protection state of the device under certain conditions. This is completely normal and not a fault. Please use it without worry!

SPECIFICATIONS






- **Product Dimensions:** 8.07 x 5.7 x 4.92 inches
- **Item Model Number:** E2
- **Weight:** 1.92 Pounds
- **Brand:** COSLUS


- **Power Source:** Corded Electric
- **Special Features:** 24 Personalized Cleaning Levels, Auto Shut-Off, Multiple Operation Modes, Multiple Pressure Settings, Multiple Tips, Lightweight, Large Water Tank, Non-Slip Grip, Rotatable Tip
- **Product Benefits:** Deep Cleansing, Gentle Cleansing, Gum Health, Prevents Bad Breath, Prevents Cavities, Removes Plaque
- **Unit Count:** 1 Count

WARRANTY & SUPPORT

For any questions, support, or warranty inquiries, please refer to the contact information provided in your product packaging or visit the official COSLUS website.

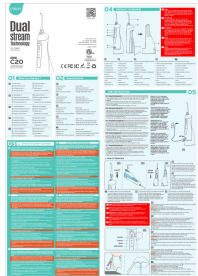Related Documents - E2

	<p>COSLUS Water Dental Flosser F5020E User Manual</p> <p>User guide for the COSLUS Water Dental Flosser F5020E, a portable cordless oral irrigator. Includes safety precautions, product overview, and operating instructions for charging and use.</p>
	<p>COSLUS AOW10 Oral Irrigator User Manual - Features, Setup, and Usage Guide</p> <p>Comprehensive user manual for the COSLUS AOW10 Oral Irrigator. Learn about product features, setup instructions, operating procedures, tip usage, maintenance, and troubleshooting for optimal oral hygiene.</p>
	<p>COSLUS C20 Dual Stream Water Dental Flosser User Manual & Guide</p> <p>Comprehensive user manual for the COSLUS C20 Dual Stream portable cordless water dental flosser. Learn about features, specifications, usage, maintenance, and troubleshooting for effective oral hygiene.</p>
	<p>COSLUS F5020E Oral Irrigator User Manual and Guide</p> <p>Comprehensive user manual for the COSLUS F5020E Oral Irrigator, covering safety precautions, product features, operating instructions, cleaning, maintenance, troubleshooting, and disposal. Learn how to effectively use and care for your COSLUS oral irrigator.</p>
	<p>COSLUS F5020E Oral Irrigator User Manual</p> <p>This user manual provides instructions for the COSLUS F5020E Oral Irrigator, covering product overview, charging, operation, cleaning, maintenance, troubleshooting, and disposal.</p>

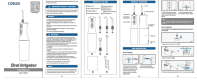

[COSLUS Oral Irrigator F5029A User Manual - Guide to Oral Hygiene](#)

Comprehensive user manual for the COSLUS Oral Irrigator (Model F5029A). Learn how to use, charge, maintain, and troubleshoot your device for effective oral hygiene. Includes safety precautions and disposal information.