

Burrki Cordless Heating Pad with Massager

Burrki Cordless Heating Pad with Massager User Manual

Model: Cordless Heating Pad with Massager for Back Pain Relief

Brand: Burrki

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your Burrki Cordless Heating Pad with Massager. Designed for targeted heat therapy and soothing vibration massage, this device offers relief for discomfort in the lower back, hips, waist, and abdomen, including menstrual cramps, abdominal pain, stomach upset, leg arthritis, and muscle tension.

Please read this manual thoroughly before first use and retain it for future reference.

SAFETY INFORMATION

WARNING: This is an electric heating pad. To prevent burns and fire hazards, do not leave the pad unattended while in use. Always follow the instructions for proper use and storage. Keep out of reach of children and pets. Do not use if damaged. Consult a doctor before use if you have any medical conditions or are pregnant. In case of accidental contact with skin, immediately cool the affected area with cool water and seek medical attention if necessary.

Indications for Use:

- Menstrual cramps
- Abdominal pain
- Stomach upset
- Leg arthritis
- Muscle tension
- Lower back pain
- Hip pain
- Waist pain

PRODUCT COMPONENTS

Your Burrki Cordless Heating Pad with Massager package includes the following items:

- Rechargeable Cordless Heating Pad with Massager
- Extension Strap (for waists up to 60 inches)
- 1 ft Type-C Charge Cable
- 6.6 ft Type-C Charge Cable
- User Manual



Image: The main product, a blue cordless heating pad with a digital display and adjustable straps.

LOCK SCREEN FUNCTION



15 20 30 Minutes Timer

Short Press the power button again after turn on

Image: Contents of the product package, including the heating pad, extension strap, two USB-C cables, and user manual.

SETUP

1. Charging the Device:

The heating pad is designed with a built-in battery for cordless use. Before first use, fully charge the device using one of the provided USB-C cables. The device can also be used while plugged in.

CORDLESS RECHARGEABLE LARGE CAPACITY BATTERY

Support External Power Bank (not included)



**Battery
3000mAh**

Image: A person wearing the heating pad, highlighting its cordless rechargeable battery with 3000mAh capacity.

CAN BE USED WHILE CHARGING

(NO ADAPTER INCLUDED)



Image: A person using the heating pad while it is plugged into a power outlet, demonstrating its ability to be used during charging.

2. Adjusting the Fit:

The heating pad features two-way adjustable straps and an included extension band to fit waists up to 47-60 inches. Adjust the straps to customize the fit to your comfort. The front Velcro ensures a secure attachment.

ADJUSTABLE SUPPORT STRAPS WITH EXTENSION BAND-UP TO 60" WAIST



Fits Most People

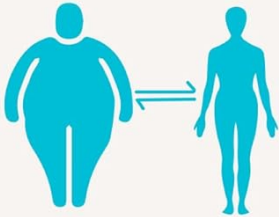


Image: Diagram showing the adjustable support straps and extension band, indicating a fit for waists up to 60 inches.

UP TO 60-INCH WAIST
WITH EXTENSION BELT

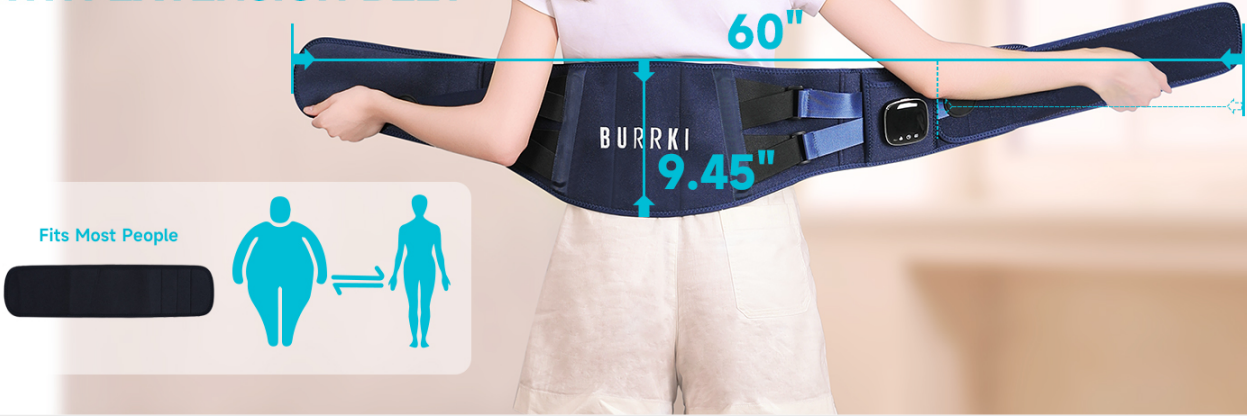


Image: A sequence of five images demonstrating how to adjust and secure the elastic straps of the heating pad for a comfortable fit.

OPERATING INSTRUCTIONS

1. Power On/Off:

Press and hold the power button on the control panel to turn the device on or off.

2. Adjusting Heat Levels:

The heating pad offers three adjustable heat levels:

- Low: 113°F (45°C)
- Medium: 131°F (55°C)
- High: 149°F (65°C)

Press the heat button on the control panel to cycle through the heat settings. The current setting will be displayed on the LCD screen.



Image: A person wearing the heating pad, with an overlay showing the three heating temperature levels and five vibration massage

intensities.

3. Adjusting Vibration Massage:

The device features five adjustable vibration intensities. Press the massage button on the control panel to cycle through the vibration modes. You can use heating and massage functions simultaneously or independently.

4. Setting the Timer:

Choose from three timing options: 15, 20, or 30 minutes. Press the timer button to select your desired duration. The device will automatically shut off after the selected time.

5. Lock Screen Function:

Short press the power button after turning the device on to lock the screen. This prevents accidental changes to your settings during use.

LOWER BACK AND LUMBAR SUPPORT HEATING PAD WITH MASSAGER

- *Targeted Heat Therapy*
- *Soothing Vibration Massage*
- *Effective Pain Relief*
- *Lumbar Supports Pad*



Image: Close-up of the control panel showing the digital timer (30 minutes) and an illustration of how to short press the power button to activate the lock screen function.

Usage Guidelines:

- Use for a maximum of 30 minutes at a time.
- Allow at least a 30-minute break between uses.
- Avoid using while sleeping or under blankets.
- Keep away from direct contact with skin to prevent burns.

MAINTENANCE

Cleaning:

Clean the pad with a damp cloth and mild soap. Air dry completely before storing or next use. Do not immerse the device in water.

Storage:

Store the heating pad in a cool, dry place away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Ensure the device is fully charged. Connect to power using the provided USB-C cable.
Heat is not felt or is too low.	Incorrect heat setting or environmental factors.	Increase the heat level. Note that actual warmth may vary slightly with environment or clothing. Ensure direct contact with clothing, not bare skin.
Vibration is not working.	Vibration function is off or intensity is too low.	Press the massage button to activate vibration and cycle through intensities.
Device turns off unexpectedly.	Timer has expired or low battery.	Check the timer setting. Recharge the device if the battery is low.
Velcro does not hold securely.	Improper adjustment or worn Velcro.	Ensure the belt is adjusted snugly. The front Velcro uses a super-dense adhesive for secure grip. If worn, contact customer support.

SPECIFICATIONS

- **Product Dimensions:** 7.87 x 1.57 x 7.87 inches
- **Weight:** 1.81 Pounds

- **Batteries:** 1 A battery required (included)
- **Heat Levels:** 3 (113°F-149°F / 45°C-65°C)
- **Vibration Intensities:** 5
- **Timer Settings:** 15, 20, 30 minutes
- **Waist Fit:** Up to 47-60 inches (with extension band)
- **Material:** Faux Leather
- **Special Features:** Automatic Shut Off, Lightweight, Overheat Protection, Portable, Timer
- **Manufacturer:** Burrki

WARRANTY AND SUPPORT

Burrki offers an 18-month warranty for this product, providing peace of mind for your purchase. For any product inquiries, technical support, or warranty claims, please contact Burrki customer service directly.

You can visit the official Burrki store on Amazon for more information and support: [Burrki Amazon Store](#)

Please refer to the contact information provided with your product packaging for the most up-to-date customer service details.