

KNQZE KF9

KNQZE KF9 Shiatsu Foot Massager Machine with Heat

Instruction Manual
Model: KF9 | Brand: KNQZE

INTRODUCTION

Thank you for choosing the KNQZE KF9 Shiatsu Foot Massager Machine with Heat. This manual provides essential information for the safe and effective use of your new foot massager. Please read it thoroughly before operation and retain it for future reference.

IMPORTANT SAFETY INSTRUCTIONS

- Always unplug the massager before use and when not in use.
- Do not use if you have any foot injuries or medical conditions that may be adversely affected by heat or massage. Consult a physician if unsure.
- Keep the massager dry and do not submerge it in water or any other liquid.
- Do not use while sleeping or with bare feet if you have sensitive skin.
- Keep the massager out of reach of children and pets.
- Always follow the manufacturer's instructions for use and maintenance.
- Do not rely solely on the information presented here; always read product labels, warnings, and directions before using or consuming a product.

PACKAGE CONTENTS

The KNQZE KF9 Foot Massager package includes:

- KNQZE KF9 Foot Massager Machine
- Remote Control (requires 1 AAA battery, not included)
- Power Adapter

- User Manual (this document)

PRODUCT FEATURES

The KNQZE KF9 Foot Massager is designed to provide a comprehensive and customizable foot massage experience. Key features include:

- **360° Air Squeezing Technology:** Gently compresses the entire foot for enhanced relaxation.
- **3-Section Rolling Kneading System:** Ergonomic massage heads target the forefoot, arch, and heel with Shiatsu-inspired techniques.
- **Adjustable Intensity Levels:** Choose from 3 levels for kneading and air compression.
- **Soothing Heat Function:** Gentle heat can be activated to further relax muscles.
- **Timer Settings:** Select massage durations of 15, 20, or 30 minutes with auto shut-off.
- **Dual Control Options:** Operate via the integrated LCD touch panel or the convenient remote control.
- **Removable and Washable Foot Covers:** Ensures hygiene and easy maintenance.
- **Compact and Portable Design:** Easy to store and move between locations.


2 Control Options

Larger LCD Touch Panel x Remote Control



Image: The KNOZE KF9 Foot Massager showing its LCD touch panel and accompanying remote control for easy operation.

3-Section Comprehensive Massage for Ultimate Relaxation

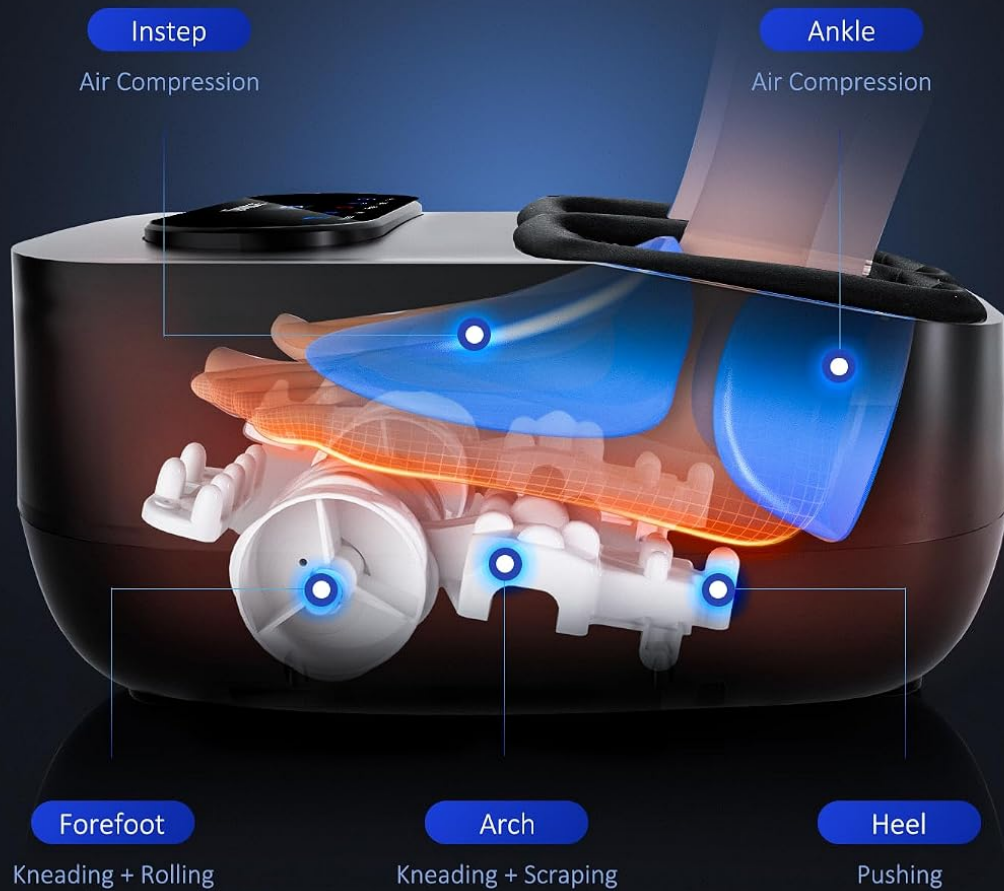

Kneading



Rolling



Acupressure



Scraping



Firm Pressure

Image: An internal diagram of the KNQZE KF9 Foot Massager illustrating the 3-section comprehensive massage system targeting the instep, ankle, forefoot, arch, and heel.

SETUP

1. Unpack the foot massager and all accessories from the box.
2. Place the massager on a flat, stable surface near a power outlet.
3. Connect the power adapter to the massager's power port.
4. Plug the power adapter into a suitable electrical outlet.
5. Ensure the foot covers are properly installed inside the foot chambers.

Your browser does not support the video tag.

Video: A demonstration of unboxing and setting up a foot massager, including plugging it into a power source. This video is provided by Alyno, a Seller.

OPERATING INSTRUCTIONS

1. **Power On:** Sit comfortably and place your feet into the foot chambers. Press the **Power** button on the LCD touch panel or remote control to turn on the massager.
2. **Select Massage Mode:**
 - **Kneading:** Press the **Knead** button to activate the deep kneading and rolling massage. Adjust intensity (L/M/H) by pressing the button repeatedly.
 - **Air Compression:** Press the **Air** button to activate the air compression massage. Adjust intensity (L/M/H) by pressing the button repeatedly.
 - **Heat:** Press the **Heat** button to turn on the soothing heat function.
3. **Set Timer:** Press the **Timer** button to select your desired massage duration (15, 20, or 30 minutes). The massager will automatically shut off after the selected time.
4. **Adjust Settings:** Use the respective buttons on the touch panel or remote to customize kneading intensity, air compression intensity, and heat level to your preference.
5. **Power Off:** Press the **Power** button again to turn off the massager. Unplug the device from the power outlet after each use.



Image: A user operating the KNQZE KF9 Foot Massager using the remote control while relaxing on a couch.



Image: A woman using the KNQZE KF9 Foot Massager, highlighting the heat therapy feature with visual indicators of temperature levels.

Adjustable Air Compression

Thoroughly relax swollen insteps & sore heels
with different squeeze modes



Image: A man using the KNQZE KF9 Foot Massager, illustrating the adjustable air compression feature with different intensity levels.

MAINTENANCE

- **Cleaning the Foot Covers:** The foot covers are removable and washable. Unzip the covers from the foot chambers, remove them, and wash with mild soap and warm water. Allow them to dry completely before reattaching.
- **Cleaning the Unit:** Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or submerge the unit in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight. Its compact design allows for easy storage in a living room, home office, or bedroom.

Max - US Size 12

Fits Whole Family

Larger massage space suitable for all types of feet

40%

↑

Massaging Area

30%

↑

Comfort Experience

Detachable Sleeves

Waterproof Body

Simple and Elegant Appearance

L 15.7 inch

M 12.8 inch

H 7.1 inch

Non-slip Rubber Pads

Image: The KNQZE KF9 Foot Massager showcasing its detachable and washable foot sleeves, waterproof body, and compact dimensions for easy storage.

Your browser does not support the video tag.

Video: A demonstration of the foot massager's features, including the removable and washable foot sleeves for easy cleaning. This video is provided by Alyno, a Seller.

TROUBLESHOOTING

If you encounter any issues with your KNQZE KF9 Foot Massager, please refer to the following common solutions:

Problem	Possible Cause	Solution
Massager does not turn on.	Not plugged in correctly; power outlet issue.	Ensure the power adapter is securely plugged into both the massager and a working electrical outlet.
Massage intensity feels too weak or too strong.	Incorrect intensity setting selected.	Adjust the kneading or air compression intensity using the respective buttons (L/M/H).
Heat function is not noticeable.	Heat function not activated or insufficient warm-up time.	Ensure the heat function is turned on. Allow 5-10 minutes for the heat to become noticeable.
Unusual noise during operation.	Normal operating sound; potential obstruction.	Some operational noise is normal. Check for any obstructions in the foot chambers. If noise is excessive or unusual, discontinue use and contact support.

SPECIFICATIONS

- **Model Name:** KF9
- **Product Dimensions:** 34.6L x 19.5W x 42.2H Centimetres
- **Item Weight:** 3.34 Kilograms
- **Power Source:** Corded Electric
- **Material:** Plastic
- **Special Features:** Auto Shut-Off, Heat, Remote Control, Timer
- **Massage Therapy Technique:** Kneading, Rolling, Scraping, Shiatsu
- **Fits Foot Size:** Up to Men's US Size 12
- **Batteries for Remote:** 1 AAA batteries required (not included)

WARRANTY AND SUPPORT

The KNQZE KF9 Foot Massager comes with a one-year warranty from the date of purchase. If you have any questions, concerns, or require assistance with your product, please do not hesitate to reach out to our friendly customer support team. We are committed to providing excellent service and ensuring your satisfaction.

For support, please visit the official KNQZE store on Amazon or contact us through the platform where you purchased the product.