

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Brookstone](#) /

› [Brookstone Sonic Electric Power Toothbrush BK1034-WHT-STK-12 User Manual](#)

Brookstone BK1034-WHT-STK-12

Brookstone Sonic Electric Power Toothbrush User Manual

MODEL: BK1034-WHT-STK-12

Introduction

Thank you for choosing the Brookstone Sonic Electric Power Toothbrush. This manual provides essential information for the safe and effective use of your new toothbrush, designed to deliver a deep, thorough clean for healthier teeth and gums.


Product Features

- **Advanced Sonic Technology:** Utilizes powerful sonic vibrations for effective plaque removal.
- **2-Minute Smart Timer:** Ensures dentist-recommended brushing time with 30-second quadrant reminders.
- **5 Brushing Modes:** Offers Clean, White, Polish, Gum Care, and Sensitive modes for personalized oral care.
- **Gentle & Effective Cleaning:** Soft, high-frequency bristles are gentle on gums while removing surface stains.
- **Ergonomic Design:** Lightweight handle provides a comfortable and secure grip.

Package Contents

Ensure all items are present in your package:

- 1 x Brookstone Sonic Electric Power Toothbrush Handle
- 4 x Brush Heads (1 pre-installed, 3 replacements)
- 1 x USB-C Charging Cable
- 1 x User Manual

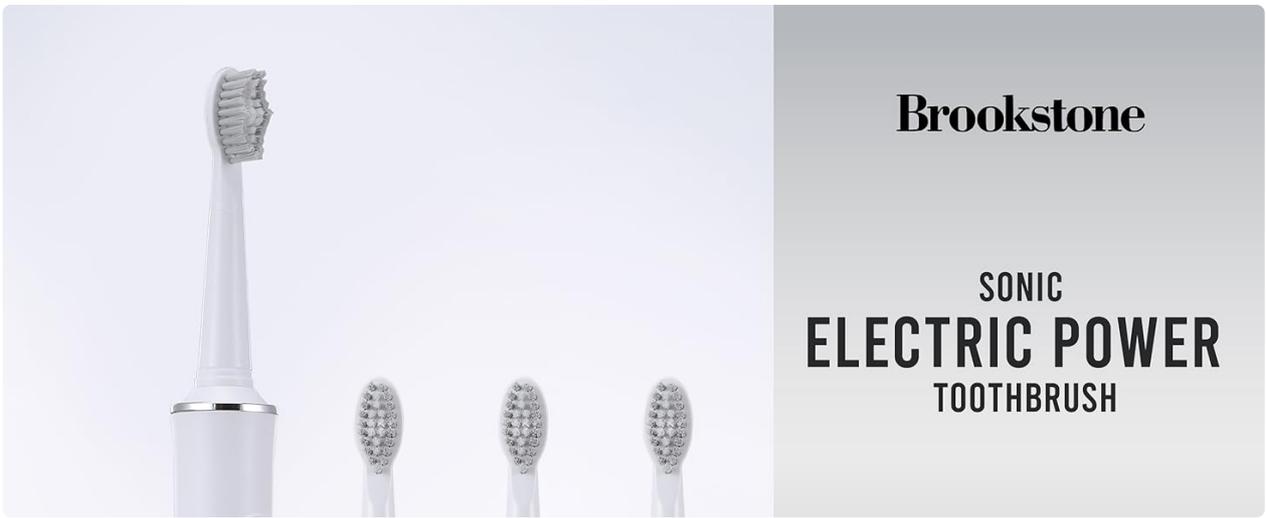

Brookstone

SONIC
ELECTRIC POWER
TOOTHBRUSH

Image: The Brookstone Sonic Electric Power Toothbrush handle shown alongside three replacement brush heads, illustrating the complete package contents.

Product Overview

Familiarize yourself with the components of your toothbrush:

Brookstone



2

MINUTES
BUILT-IN TIMER



TYPE-C
CHARGING



REMOVES
PLAQUE & STAINS

5

OPTIONAL
MODES

Image: A white Brookstone Sonic Electric Power Toothbrush, front view, showing the brush head, power button, and mode indicator lights.

- **Brush Head:** Replaceable head with soft bristles for cleaning.
- **Power Button:** Turns the toothbrush on/off and cycles through modes.
- **Mode Indicator Lights:** Illuminates to show the currently selected brushing mode.
- **Charging Port:** USB-C port located at the bottom of the handle for charging.

Setup

Initial Charging

Before first use, fully charge your toothbrush.

1. Connect the provided USB-C charging cable to the port at the bottom of the toothbrush handle.
2. Connect the other end of the cable to a compatible USB power adapter (not included) or a USB port.
3. The indicator light will show the charging status. Allow the toothbrush to charge completely before its first use.

Brookstone



IPX7
WATERPROOF

Image: A detailed view of the bottom of the toothbrush handle, showing the USB-C charging port for power connection.

Attaching the Brush Head

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks securely into place.

Operating Instructions

Power On/Off

- Press the power button once to turn the toothbrush on.
- Press and hold the power button for approximately 2 seconds to turn the toothbrush off.

Selecting a Brushing Mode

With the toothbrush powered on, press the power button repeatedly to cycle through the available modes:

- **Clean:** Standard daily cleaning.
- **White:** Designed to remove surface stains.
- **Polish:** For brightening and polishing teeth.

- **Gum Care:** Gentle massage for gums.
- **Sensitive:** Extra gentle cleaning for sensitive teeth and gums.

The indicator light next to the selected mode will illuminate.



Image: The toothbrush handle displaying the five different brushing modes, allowing users to select their preferred setting.

Brushing Your Teeth

1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gumline.
3. Press the power button to start the toothbrush.
4. Gently guide the brush head slowly across all surfaces of your teeth. Do not scrub. Allow the sonic vibrations to do the work.
5. The 2-minute smart timer will pause briefly every 30 seconds to indicate when to move to a different quadrant of your mouth (upper right, upper left, lower right, lower left).
6. The toothbrush will automatically turn off after 2 minutes, signaling the completion of the recommended brushing time.



Brookstone

HEALTHIER TEETH & STRONGER GUMS

Image: A woman smiling while brushing her teeth with the Brookstone Sonic Electric Power Toothbrush, demonstrating proper usage.

Maintenance and Care

Cleaning the Brush Head

- After each use, rinse the brush head thoroughly under running water while the toothbrush is still on to remove any remaining toothpaste or debris.
- Turn off the toothbrush, remove the brush head from the handle, and rinse both parts separately.

Cleaning the Handle

- Wipe the toothbrush handle with a damp cloth.
- Do not immerse the handle in water for extended periods, despite its water-resistant design.

Storage

- Store the toothbrush in a dry, upright position.
- Replace brush heads every 3 months, or sooner if the bristles appear frayed or worn.

Brookstone

CHECK OUT THE REST OF OUR WELLNESS MUST-HAVES



Image: The white Brookstone Sonic Electric Power Toothbrush lying on a dark marble surface, highlighting its IPX7 waterproof rating for easy cleaning.

Safety Information

Please read and follow these important safeguards:

- **Choking Hazard:** Keep out of reach of children under 3 years of age.
- Do not use this product for purposes other than oral hygiene.
- This is a battery-powered device; do not ingest chemicals if swallowed. Seek immediate medical attention if this occurs.
- Consult a dentist or physician before use if you have severe periodontal disease, have undergone oral surgery, or have any other oral health concerns.
- Do not attempt to modify, disassemble, or repair the product yourself. Contact customer support for assistance.
- Avoid exposing the charging cable or toothbrush to extreme temperatures or direct sunlight.

Troubleshooting

If you encounter any issues, refer to the following common problems and solutions:

- **Toothbrush not turning on:** Ensure the battery is sufficiently charged. Press the power button firmly.
- **Toothbrush not charging:** Verify that the USB-C cable is securely connected to both the toothbrush and a working power source. Ensure the charging port is clean and free of debris.
- **Reduced brushing effectiveness:** Replace the brush head if the bristles are worn or splayed. Ensure the brush head is securely attached to the handle.
- **Unusual noise or vibration:** Check that the brush head is properly attached. If the issue persists, discontinue use and contact customer support.

Specifications

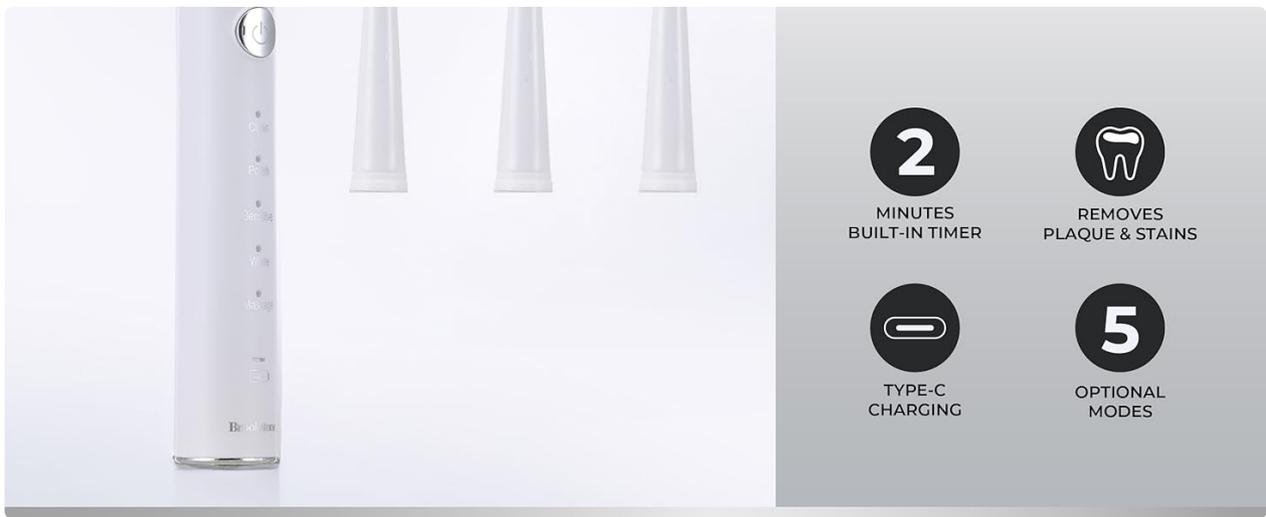


Image: The Brookstone Sonic Electric Power Toothbrush highlighting its operating frequency of 31,000 strokes per minute.

Feature	Detail
Model Number	BK1034-WHT-STK-12
Power Source	Battery Powered
Number of Brushing Modes	5 (Clean, White, Polish, Gum Care, Sensitive)
Smart Timer	2-minute auto-off with 30-second quadrant reminders
Operating Frequency	31,000 strokes/min
Bristle Material	Nylon (Soft)
Handle Material	Plastic
Dimensions (L x W x H)	4.3 x 1.4 x 2 inches
Water Resistance	IPX7
Color	White

Warranty and Support

Warranty Information

For detailed warranty information regarding your Brookstone Sonic Electric Power Toothbrush, please refer to the warranty card included with your product packaging or visit the official Brookstone website.

Customer Support

If you have any questions, require technical assistance, or need to report an issue, please contact Brookstone customer service. You can find contact details and frequently asked questions (FAQs) on the official Brookstone website: www.brookstone.com/support