

Lunix LX10

Lunix LX10 Full Leg Massager User Manual

Model: LX10 (FE-7204B)

INTRODUCTION

The Lunix LX10 Full Leg Massager provides a comprehensive compression massage experience for your feet, calves, and thighs. It features optional heat and cold therapy, an adjustable fit, and a user-friendly handheld controller. This manual will guide you through the setup, operation, and maintenance of your device to ensure optimal performance and a relaxing experience.

LX10 Cordless Full-Leg Massager

With Heat & Cold Therapy



Heat & Cold Therapy



Cordless & Rechargeable



4 Programs + 6 Intensities



Complete Kit with Extras

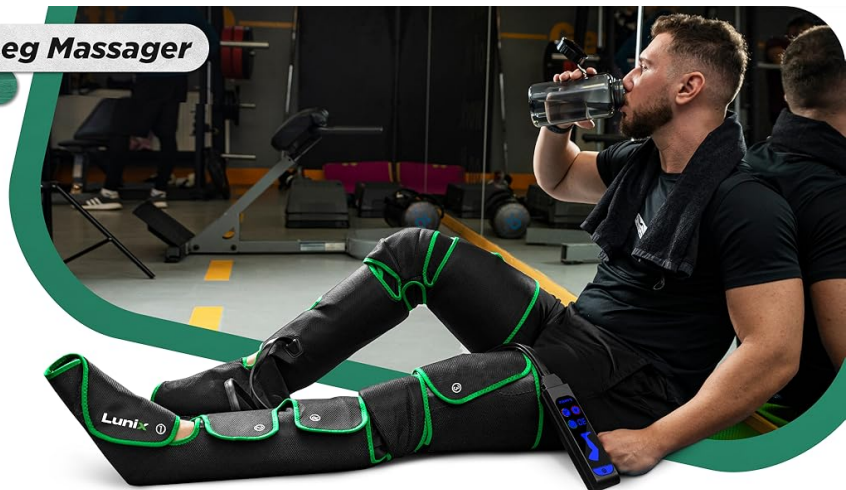


Image: The Lunix LX10 Full Leg Massager illustrating the acupressure technology and compression points across the leg.

SETUP INSTRUCTIONS

1. Unpacking

Carefully remove all components from the packaging. Verify that all items listed below are present:

- Lunix LX10 Leg Massager wraps

- Remote Control
- Plug Adapter
- 4x Extra Straps (Size Extenders)
- 4x Gel Packs
- Travel Bag
- Foot Roller (Bonus Item)

Why choose Lunix

	Lunix	Others
Designed in the US	✓	✗
Customer Support Response < 24hrs	(307) 461-9622	✗
Full Leg Massager	✓	✗
Heat and Cold Therapy	✓	✗
Massage Programs	4 Modes + 6 Intensity Levels	1-2 Modes
Cordless & Rechargeable	✓	✗
Adjustable Fit with Extra Straps	✓ up to 27"	✗
LCD Touchscreen Remote	✓	✗
Accessories Included	Travel Bag + Foot Roller + Gel Packs	✗



Image: All components included in the Lunix LX10 package, including the leg massager, remote, adapter, extra straps, gel packs, travel bag, and foot roller.

2. Assembly and Fit

The LX10 consists of leg wraps designed for full lower body coverage. Unfold the leg wraps and wrap each section around your foot, calf, and thigh, securing them with the Velcro fasteners. Ensure a snug but

comfortable fit. Use the included 4-piece size extenders if additional length or circumference is needed for a comfortable fit. The massager is designed to accommodate thigh areas up to 27 inches and calf areas up to 20 inches.

4 Ways to Enjoy the Leg Massager



Image: The Lunix LX10 Leg Massager demonstrating its flexible design with adjustable straps and size extenders for various leg sizes.

3. Charging the Device

Connect the plug adapter to the massager's charging port and then to a standard power outlet. The battery indicator on the handheld controller will show the charging status. A full charge provides approximately 1 hour of continuous run time.



Image: A user relaxing with the cordless Lunix LX10 Leg Massager, emphasizing its portability and rechargeable battery for use anywhere.

OPERATING INSTRUCTIONS

The Lunix LX10 is controlled via a user-friendly LCD touchscreen remote. Familiarize yourself with the buttons and indicators:

LCD Touchscreen

Remote Control



Image: A detailed view of the Linux LX10's LCD touchscreen remote control, showing buttons for power, program, intensity, timer, and heat.

- **Power Button:** Turns the device on or off.
- **Program Button:** Cycles through massage modes.
- **Intensity Button:** Adjusts compression strength (6 levels).
- **Timer Button:** Sets session duration (15/20/30 minutes).
- **Heat Button:** Activates or deactivates the heat function.
- **Battery Indicator:** Shows current battery level.
- **Time Indicator:** Displays remaining session time.
- **Intensity Level Indicator:** Shows current compression intensity.

1. Power On/Off

Press the **Power Button** on the handheld LCD touchscreen controller to turn the device on or off.

2. Selecting a Program

The LX10 offers 4 targeted programs. Use the **Program Button** to cycle through these options:

- **Full Leg Massage:** Provides compression across the entire leg.
- **Calf Massage:** Focuses compression on the calf area.
- **Thigh Massage:** Focuses compression on the thigh area.
- **Heat Massage:** Activates the heat function in specific regions.



Image: Visual representation of the four available massage modes: Full Leg, Calf, Thigh, and Heat, showing targeted areas of compression or warmth.

3. Adjusting Intensity

The device features 6 intensity levels. Use the **Intensity Button** to increase or decrease the compression strength to your preference. The **Intensity Level Indicator** will display the current setting.

4. Setting Session Time

Choose from 15, 20, or 30-minute session durations using the **Timer Button**. The **Time Indicator** will show the remaining time.

5. Heat Therapy

Press the **Heat Button** to activate the heat function. Heat is typically applied to the knee and foot regions. The heat feature takes approximately 10 minutes to reach optimal temperature (70°F when plugged in, 60°F cordless).



Image: The Lunix LX10 Leg Massager showing the two heat regions located at the knee and foot for targeted warmth.

6. Cold Therapy Option

For cold therapy, pre-cool the 4 gel packs in a freezer. Insert the chilled gel packs into the designated pockets inside the leg wraps. This can be used alone or combined with compression for enhanced relief.

Acupressure Technology

Made by **Lunix**

- Heat Therapy
- Cold Therapy
- 4 Targeted Programs
- 15/20/30' Sessions
- Shut Off Protection

A.Tech
made by
Lunix

Image: The Lunix LX10 Leg Massager with gel packs, demonstrating the option for cold therapy by inserting pre-cooled packs into the wraps.

MAINTENANCE

Cleaning

Disconnect the device from power before cleaning. Wipe the exterior of the leg wraps and controller with a

soft, damp cloth. Do not immerse the device in water or use harsh chemical cleaners.

Storage

When not in use, store the massager and its accessories in the provided travel bag in a cool, dry place, away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Problem	Solution
Device Not Turning On	<ul style="list-style-type: none">• Ensure the battery is charged. Connect the plug adapter and allow it to charge for at least 30 minutes.• Verify all connections between the wraps and the controller are secure.
Weak Compression	<ul style="list-style-type: none">• Check the selected intensity level and increase if necessary.• Ensure the leg wraps are securely fastened and not too loose.
Heat Function Not Working	<ul style="list-style-type: none">• Confirm the Heat Button has been pressed.• Allow sufficient time (approx. 10 minutes) for the heat to reach its operating temperature.• Ensure the device is adequately charged or plugged in, as cordless operation may result in lower heat.
Remote Control Unresponsive	<ul style="list-style-type: none">• Ensure the controller is properly connected to the leg wraps.• If the issue persists, try turning the device off and on again.

SPECIFICATIONS

Product Dimensions	13 x 9.84 x 28.15 inches
Item Weight	6.6 Pounds
Item Model Number	FE-7204B
Power Source	Battery Powered
Batteries	1 x 12V battery (included)
Special Features	Cordless, Heat, LCD Display, Portable, Rechargeable
Target Use Body Part	Calf, Foot, Leg, Thigh
Massage Technique	Acupressure, Compression, Shiatsu
Material Type	Polyester
UPC	850049964821

WARRANTY INFORMATION

The Lunix LX10 Full Leg Massager is backed by a **lifetime warranty**, ensuring long-term peace of mind and consistent performance. For any warranty claims or product issues, please contact our customer support team.

CUSTOMER SUPPORT

For any questions, assistance, or support regarding your Lunix LX10 Full Leg Massager, please contact our US-based customer support team. We are committed to providing reliable support with a guaranteed response time of under 24 hours.

Customer Support Phone: (307) 461-9622



Image: Lunix customer support team members, emphasizing their readiness to assist with product inquiries.