

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Breio](#) /
- › [Breio Scalp 3L Laser Hair Growth Comb Instruction Manual](#)

## Breio Scalp 3L

# Breio Scalp 3L Laser Hair Growth Comb Instruction Manual

Model: Scalp 3L | Brand: Breio

## INTRODUCTION

The Breio Scalp 3L is a specialized device designed for scalp care and hair stimulation. It integrates 650nm infrared laser technology with a 360° scalp massage function. This manual provides essential information for the safe and effective use, maintenance, and troubleshooting of your device.

# LASER - CLASS 3R TECHNOLOGY

The laser light penetrates 3-5mm into the hair root, "awakening" the scalp from the inside out and promoting healthy, shiny hair.

7 x 650NM RED LIGHT

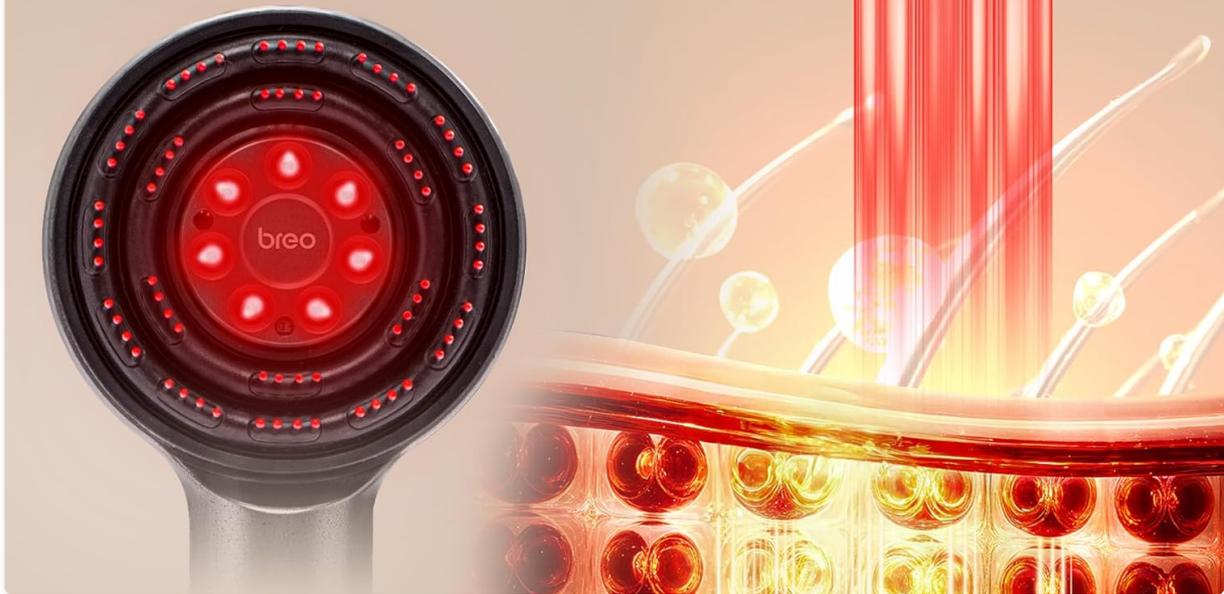


Figure 1: Overview of the Breo Scalp 3L Laser Hair Growth Comb.

## KEY FEATURES

- **650nm Infrared Laser Technology:** Equipped with 7 infrared laser lights that deliver targeted light to hair follicles, stimulating hair growth and revitalizing the scalp.
- **360° Scalp Massage:** Features 72 massage comb teeth covering over 4500mm<sup>2</sup>, providing a gentle, kneading massage to relieve stress and fatigue.
- **Laser Activation Safety Feature:** A skin-detection system activates the laser only when the comb is in safe contact with the scalp, automatically shutting off if the distance is exceeded to protect the eyes.
- **Waterproof Design (IPX7):** Fully waterproof for versatile use in both wet and dry conditions, making it convenient for use in the shower and easy to clean.
- **Compact & Lightweight:** Designed for portability and ease of use in daily scalp care routines.

## SAFETY INFORMATION

- Always keep the device away from your eyes and face.
- Avoid using the device on broken or irritated skin.

- Keep the device away from water while in use, unless specifically designed for wet use (this device is IPX7 waterproof, but caution is always advised).
- Do not use the device while sleeping or lying down.

## SETUP

---

1. **Charging:** Before first use, fully charge the device using the provided Type-C charging cable. A full charge typically takes 2 hours and provides up to 10 days of use (based on 15 minutes/day).
2. **Familiarization:** Locate the power button and observe the massage bristles and laser lights.
3. **Optional: Hair Oil Application:** If using hair growth oils, apply them to your scalp before or during the massage mode for enhanced absorption.

## OPERATING INSTRUCTIONS

---

The Breo Scalp 3L offers two distinct modes: Massage Mode and Laser Mode. These modes operate independently to provide targeted benefits.



Figure 2: The Breo Scalp 3L features separate Massage and Laser modes.

### Usage Recommendations:

- Use 3 to 4 times a week.
- Each session should last approximately 15 minutes.
- For best results, continuous use for 16 weeks or more is recommended.

### Steps for Use:

1. **Power On:** Press and hold the power button to turn on the device.
2. **Select Mode:** Press the power button briefly to cycle between modes:
  - **Massage Mode:** The device will provide a gentle, hand-like massage through its 72 bristles. This mode helps increase blood flow to the scalp and can aid in the absorption of hair growth oils.

- **Laser Mode:** The 650nm infrared laser lights will activate. Ensure the comb is in contact with your scalp for the laser to function due to the skin-detection safety feature.
3. **Application:** Apply gentle pressure while moving the comb through your hair, ensuring the bristles or laser lights make contact with your scalp.
  4. **Power Off:** Press and hold the power button to turn off the device after your session.

Your browser does not support the video tag.

Video 1: Official demonstration of the Breo Scalp 3L Laser Hair Growth Comb in use, highlighting its features and application.

Your browser does not support the video tag.

Video 2: A short clip showcasing the scalp massage and laser functions of the Breo Scalp 3L.

Your browser does not support the video tag.

Video 3: An official video demonstrating the red light therapy for hair growth provided by the Breo Scalp 3L.

## MAINTENANCE

---

### Cleaning:

- Clean the device after each use with a damp cloth and mild soap.
- Ensure the device is completely dry before storing or recharging.

### Storage:

- Store the device in a cool, dry place away from direct sunlight.
- Keep out of reach of children.

## TROUBLESHOOTING

---

### Laser Lights Not Activating:

- The device is equipped with a skin-detection safety feature. The laser lights will only activate when the comb is in direct and safe contact with your scalp.
- Ensure the device is fully charged.
- If the issue persists, contact customer support.

## SPECIFICATIONS

---

Feature	Detail
Product Dimensions	8.35 x 3.35 x 1.97 inches
Item Weight	1.17 Pounds (0.53 Kilograms)
Model Number	Scalp 3L
Power Source	Battery Powered (1 Lithium Metal battery required)
Material	Plastics, Metal

Laser Wavelength	650nm Infrared
Massage Bristles	72
Waterproof Rating	IPX7

## WARRANTY AND SUPPORT

---

Breo products are manufactured to high-quality standards. For specific warranty details and customer support, please refer to the warranty card included with your product or visit the official Breo website. Keep your proof of purchase for any warranty claims.

---

© 2025 Breo. All rights reserved.