

UREVO Urevo-Im01

UREVO Smart Calf Massager (Model Urevo-Im01) User Manual

Your Guide to Operation and Care

1. INTRODUCTION

Thank you for choosing the UREVO Smart Calf Massager. This device is designed to provide dynamic air compression massage, promoting circulation, muscle relaxation, and recovery. Featuring AI massage adaptation, three modes, three intensity levels, and adjustable heat, it offers a personalized experience. This manual provides essential information for safe and effective use of your massager.



Image: The UREVO Smart Calf Massager in use, demonstrating its application on the calves with the accompanying mobile application.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the UREVO Smart Calf Massager. Failure to follow these instructions may result in injury or damage to the device.

- Always follow the instructions provided in this user manual.
- Do not use the massager if you have any medical conditions, are pregnant, or have any injuries or skin conditions in the area to be massaged. Consult a healthcare professional if you have concerns.
- Keep the massager away from water and do not use it while sleeping or taking a bath.
- Do not leave the massager unattended while in use.
- If you experience any discomfort or pain, stop using the massager immediately and consult a healthcare professional.
- Keep the massager clean and dry to prevent any electrical hazards.

3. PACKAGE CONTENTS

Verify that all items are present in your UREVO Smart Calf Massager package:

- 2 x UREVO Relax Massagers (Left and Right)
- 1 x Storage Bag (Tote)
- 2 x Type-C Charging Cables
- 1 x Quick Start Guide (QSG)



Image: Visual representation of the UREVO Smart Calf Massager package contents.

4. PRODUCT COMPONENTS

Familiarize yourself with the main components and controls of your UREVO Smart Calf Massager:

- **Mode Selection / Exit Pause Button (M):** Cycles through massage modes.
- **Intensity Control / Pairing Button:** Adjusts air pressure intensity. Press both left and right buttons simultaneously for 3 seconds to enter/exit pairing mode.
- **Heat Button:** Cycles through 3 heating temperature levels.
- **Power / Vibration / Travel Lock Button:** Press and hold for 3 seconds to power on/off. Press and hold for 5 seconds to enter/exit travel lock mode. In power-on status, press once to adjust vibration intensity (3 levels).
- **Battery Status Light:** Indicates battery level.
- **Charging Port:** Type-C port for charging.

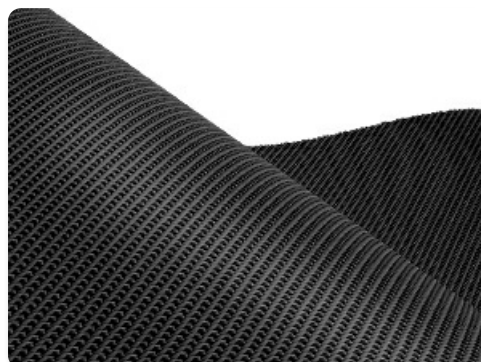


Image: Detailed diagram of the massager's control panel and charging port.

5. SETUP

5.1 Charging the Device

Before first use, fully charge both massagers using the provided Type-C charging cables. The massagers feature a 2500mAh battery, providing up to 150 minutes of use on a full charge.



Image: The massager's long battery life allows for extended relaxation.

5.2 Fitting the Massager

The UREVO Smart Calf Massager is designed with high-elastic materials and an adjustable magic tape for a secure and comfortable fit. It accommodates calves up to 18.9 inches (upper circumference), 14.2 inches (lower circumference), and 13 inches in length. Ensure the massager is snug but not overly tight to allow for proper air compression.



Image: Fitting guide for the UREVO Smart Calf Massager.

5.3 App Download and Connection

For advanced features and personalized control, download the UREVO SmartCoach app from your device's app store. Follow the in-app instructions to connect your massagers via Bluetooth. The app allows for smart comfort detection and custom massage settings.

Smart Leg Messenger with UREVO SmartCoach APP

Provides **10+** smart massage scene,
5 level timing mode (10-30min),
UREVO APP **one-click** control.



Image: The UREVO SmartCoach app provides enhanced control and customization.

6. OPERATING INSTRUCTIONS


6.1 Power On/Off and Travel Lock

- To power on or off, press and hold the Power button for 3 seconds.
- To enter or exit travel lock mode, press and hold the Power button for 5 seconds.

6.2 Selecting Massage Modes

The massager offers three primary modes: Relax, Recovery, and Shaping. You can switch between these modes using the 'M' button on the device or through the UREVO app.

Easy to Use by Dual Wireless Panel/App Control

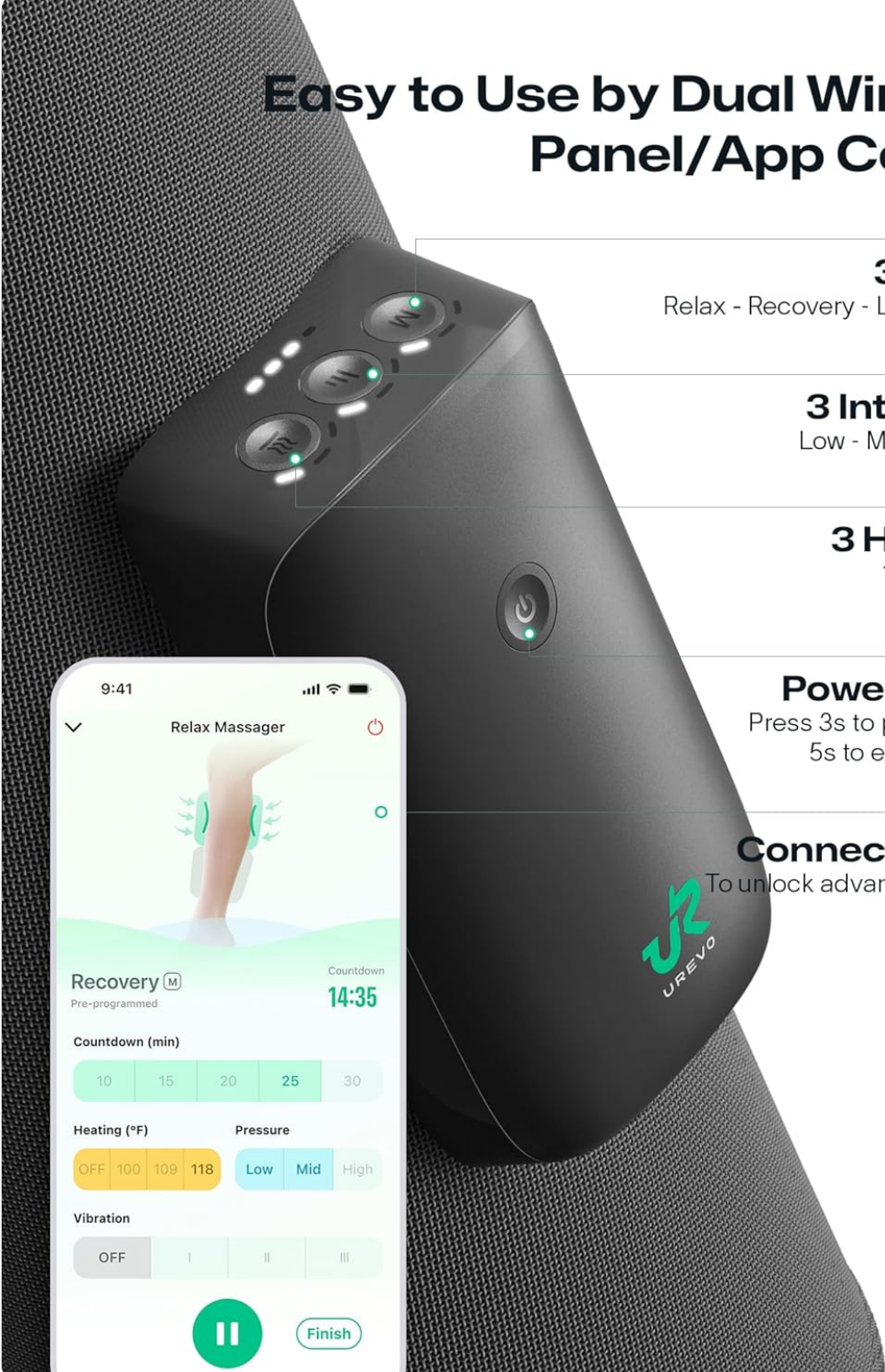

- 3 Modes**
Relax - Recovery - Leg Sculpting
- 3 Intensities**
Low - Medium - High
- 3 Heatings**
100°F - 118°F
- Power Button**
Press 3s to power on/off,
5s to exit travel lock
- Connect to App**
To unlock advanced features

Image: The massager can be controlled via its physical buttons or the UREVO app.

6.3 Adjusting Intensity and Heat

- **Intensity:** There are 3 adjustable air pressure intensity levels. Use the Intensity Control button on the device or the app to adjust.
- **Heat:** Choose from 3 adjustable heat settings (Low, Mid, High) using the Heat button on the device or the app.

3 Heating Levels

Promote Blood Circulation

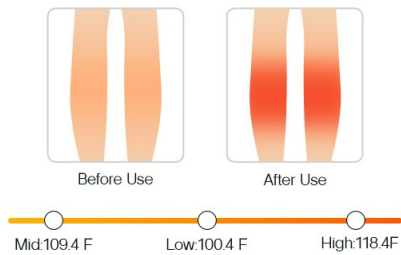


Image: The massager offers three heating levels to enhance comfort and circulation.

6.4 Vibration Modes

The massager includes 3 vibration modes to further relieve muscle tension. In power-on status, press the Power button once to cycle through the vibration intensities. Note: Using the vibration function at intensity level 3 may produce some noise.



Tips: Using the vibration function at intensity level 3 can be a bit noisy.

3 Vibration Modes

Relieve Muscle Tension

Mode 1: Vibrate 1s, pause 3s, cycle



Mode 2: Vibrate 1s, pause 2s, cycle



Mode 3: Vibrate 0.5s, pause 0.5s, cycle



Image: Three vibration modes are available for targeted muscle relief.

6.5 AI Massage Adaptation and Customization

The UREVO SmartCoach app utilizes an advanced algorithm to detect muscle discomfort and assess your leg condition. It recommends the most effective massage mode and intensity for personalized relief. You can also customize settings for each leg independently via the app.

Smart Pain Detection, Personalized Relief

UREVO App built-in precision sensors will analyze your pain tolerance and recommend your personal massage intensity settings.



Image: The UREVO app's Smart Pain Detection feature offers personalized massage settings.

6.6 Pairing and Unpairing

To pair or unpair the left and right massagers for synchronized or independent control, press both Intensity Control buttons simultaneously for 3 seconds. After pairing, the first 'M' level indicator turns blue, enabling One-Click Control. After unpairing, the 'M' indicator turns white, allowing independent control.

6.7 Usage Duration

After 40 minutes of continuous use, allow a 5-minute break before reusing the massager to ensure optimal performance and device longevity.

6.8 Official Product Videos

For visual guidance on using your UREVO Smart Calf Massager, please refer to the official product videos below:

Care Relax Massager Overview

Your browser does not support the video tag.

Video: An overview of the UREVO Care Relax Massager, demonstrating its features and benefits.

Care Relax Massager App Tutorial

Your browser does not support the video tag.

Video: A tutorial on how to use the UREVO Care Relax Massager's mobile application for customized settings.

Unboxing the Care Relax Massager

Your browser does not support the video tag.

Video: An unboxing video showcasing the contents and initial setup of the UREVO Care Relax Massager.

7. MAINTENANCE

To ensure the longevity and hygiene of your UREVO Smart Calf Massager:

- Clean the massager regularly with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the device in water.
- Store the massager in a cool, dry place, preferably in its original tote bag, when not in use.



Image: The included tote bag facilitates easy storage and transport.

8. TROUBLESHOOTING

If you encounter any issues with your UREVO Smart Calf Massager, please refer to the following common solutions:

- **Device not turning on:** Ensure the device is fully charged. Connect it to a power source using the Type-C cable and allow it to charge for at least 30 minutes before attempting to power on again.
- **No compression or weak compression:** Check that the massager is properly fitted and secured around your calf. Ensure the air compression settings are adjusted to your desired intensity level.
- **App connectivity issues:** Make sure Bluetooth is enabled on your smartphone and the massagers are powered on and within range. Try restarting both the app and the massagers.
- **Discomfort during use:** Immediately stop using the massager. Re-adjust the fit, reduce the intensity, or lower the heat setting. If discomfort persists, discontinue use and consult a healthcare professional.

9. SPECIFICATIONS

Feature	Detail
Model Number	Urevo-lm01
Product Dimensions	23.5 x 13 x 18 inches
Item Weight	450 Grams (15.87 ounces)
Power Source	Battery Powered (1 Lithium Ion battery required)
Battery Capacity	2500mAh
Max. Battery Runtime	150 Minutes
Material	Fabric
Calf Fit (Upper)	Up to 18.9 inches
Calf Fit (Lower)	Up to 14.2 inches
Length	13 inches
Manufacturer	Urevo
Date First Available	March 27, 2025

10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the Quick Start Guide included in your

package or visit the official UREVO website. Keep your purchase receipt for warranty claims.

© 2025 UREVO. All rights reserved.