

## VOYOR-HEALTH ET110

# VOYOR-HEALTH Sonic Electric Toothbrush ET110 User Manual

Model: ET110

## 1. INTRODUCTION

---

Thank you for choosing the VOYOR-HEALTH Sonic Electric Toothbrush ET110. This manual provides essential information for the safe and effective use of your new sonic toothbrush. Please read it thoroughly before first use and keep it for future reference.

## 2. IMPORTANT SAFETY INFORMATION

---

Please observe the following safety precautions:

- **WARNING:** Not suitable for children.
- Do not use for purposes other than oral hygiene.
- Keep out of reach of children.
- Swallowing the battery can cause serious injury or death.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## 3. PACKAGE CONTENTS

---

Verify that all items are present in your package:

- 1 x VOYOR-HEALTH Sonic Electric Toothbrush Handle
- 6 x Replacement Brush Heads
- 1 x USB-C Charging Cable
- 1 x User Manual

## 4. PRODUCT OVERVIEW

---



Figure 4.1: VOYOR-HEALTH ET110 Sonic Electric Toothbrush components.

The image above illustrates the main parts of the toothbrush, including the brush head, the power/mode button, indicator lights for modes and charging, and the USB charging port at the base.

## 5. SETUP

### 5.1 Charging the Toothbrush

Before first use, fully charge the toothbrush. A full charge takes approximately 3 hours and provides up to 70 days of normal usage (twice per day).

1. Locate the USB-C charging port at the bottom of the toothbrush handle.
2. Connect the provided USB-C charging cable to the port and plug the other end into a compatible USB power adapter (not included).
3. The charging indicator light will illuminate during charging. It will turn off or change color when fully charged (refer to the indicator light section for specific behavior).

# LONG BATTERY LIFE

Fully charged in 3 hours for 70 days of use



Figure 5.1: Charging the toothbrush.

This image shows the USB-C charging port at the base of the toothbrush handle with the charging cable connected, highlighting the long battery life.

## 5.2 Attaching a Brush Head

Your toothbrush comes with 6 replacement 'W' shape brush heads.

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks into place. Ensure it is securely attached.



Figure 5.2: W-shaped brush head design.

The image displays the W-shaped brush heads, emphasizing their design for deep cleaning between teeth, soft rubber coating for gum protection, and rounded tips for safety.

## 6. OPERATING INSTRUCTIONS

### 6.1 Powering On/Off and Changing Modes

The VOYOR-HEALTH ET110 features three cleaning modes: Clean, White, and Soft.

1. **To Power On:** Press the On/Off and Mode Switch Button once. The toothbrush will start in the last used mode (Memory Function).
2. **To Change Mode:** While the toothbrush is on, press the On/Off and Mode Switch Button repeatedly to cycle through the modes: Clean > White > Soft. The corresponding indicator light will illuminate.
3. **To Power Off:** Press and hold the On/Off and Mode Switch Button for 2 seconds, or cycle through all modes

until the toothbrush turns off.

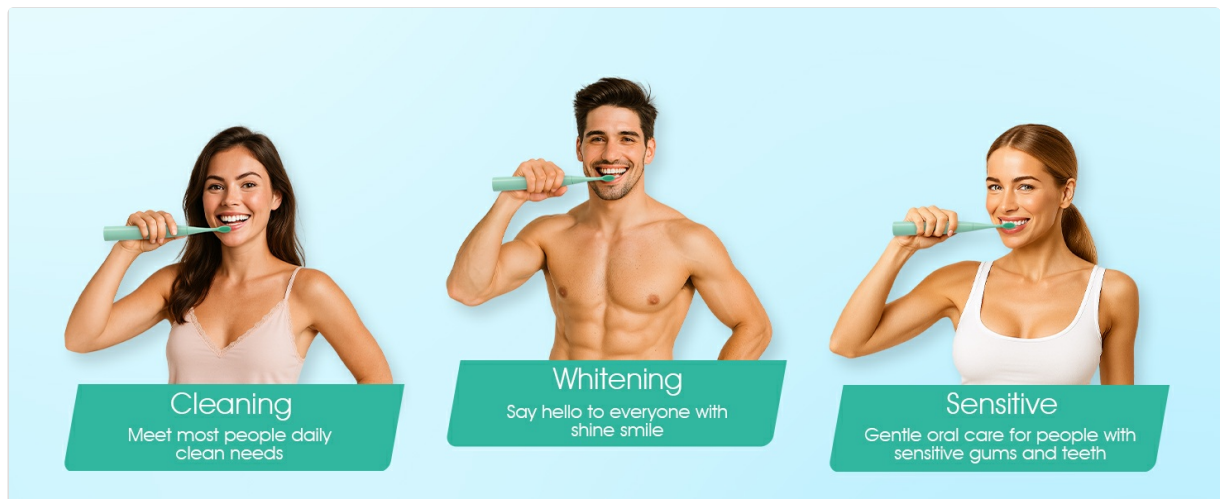


Figure 6.1: Three operating modes.

This image illustrates the three distinct operating modes: Clean for daily oral care, White for deep teeth whitening, and Soft for sensitive teeth and gums.

## 6.2 Smart Timer Function

The toothbrush includes a built-in 2-minute smart timer and a 30-second zone reminder to ensure thorough brushing.

- The toothbrush will pause briefly every 30 seconds to signal you to move to the next quadrant of your mouth.
- After 2 minutes, the toothbrush will automatically turn off, indicating the completion of the dentist-recommended brushing time.



Figure 6.2: Intelligent Timer.

This diagram visually explains the intelligent timer, showing how it prompts users to change brushing areas every 30 seconds for a comprehensive 2-minute clean.

## 6.3 Brushing Technique

For optimal results, follow these brushing guidelines:

1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gum line.
3. Gently move the brush head slowly across each tooth, allowing the sonic vibrations to do the work. Do not

scrub aggressively.

4. Follow the 30-second interval reminders to brush all four quadrants of your mouth (upper right, upper left, lower right, lower left).
5. Rinse your mouth and the brush head thoroughly after brushing.

## 7. MAINTENANCE

### 7.1 Cleaning the Toothbrush

Regular cleaning ensures hygiene and extends the life of your toothbrush.

- After each use, rinse the brush head and bristles thoroughly under running water.
- Wipe the handle with a damp cloth. Do not immerse the entire handle in water for extended periods, although it is IPX7 waterproof.
- Ensure the charging port cover is securely closed when not charging to maintain water resistance.

### 7.2 Replacing Brush Heads

Dentists recommend replacing brush heads every 3 months, or sooner if the bristles become frayed.

1. Grip the brush head firmly and pull it straight off the handle.
2. Align a new brush head with the metal shaft and push it down until it clicks into place.



Figure 7.1: Replacement brush heads.

The image displays the six included replacement brush heads, which provide up to 18 months of use based on a 3-month replacement schedule.

### 7.3 Storage

Store the toothbrush in a clean, dry place when not in use. Ensure it is fully charged before storing for extended periods to maintain battery health.

## 8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Low battery.	Charge the toothbrush fully.
Reduced vibration intensity.	Low battery.	Charge the toothbrush fully.

Problem	Possible Cause	Solution
Brush head feels loose.	Brush head not properly attached.	Remove and re-attach the brush head firmly until it clicks.
Toothbrush stops before 2 minutes.	Low battery.	Charge the toothbrush fully.

## 9. SPECIFICATIONS

Feature	Detail
Model Number	ET110
Brand	VOYOR-HEALTH
Power Source	Battery Powered (Rechargeable)
Battery Capacity	200 mAh
Charging Time	3 hours
Battery Life	Up to 70 days (based on 2 uses/day)
Vibration Frequency	Up to 41,000 VPM
Cleaning Modes	3 (Clean, White, Soft)
Waterproof Rating	IPX7
Smart Timer	2-minute auto-off with 30-second interval reminders
Brush Head Type	'W' shape, soft bristles
Product Dimensions	8.15 x 1.38 x 4.21 inches
Item Weight	8.11 ounces (0.23 Kilograms)
Material	ABS Plastic (Handle), Nylon (Bristles)

## 10. PRODUCT DEMONSTRATION VIDEO

Watch this official video for a visual guide on the VOYOR-HEALTH ET110 Sonic Electric Toothbrush's features, including installation, operating modes, and waterproof capabilities.

Your browser does not support the video tag.

Video 10.1: VOYOR-HEALTH ET110 Sonic Electric Toothbrush Features. This video demonstrates the soft and flexible bristles, easy brush head installation, the three cleaning modes (Clean, White, Soft), the memory function, the long battery life, and the IPX7 waterproof rating of the toothbrush.

## 11. WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official VOYOR-HEALTH website. Keep your proof of purchase for any warranty claims.