

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Nood](#) /
- › [Nood Flasher Pro IPL Laser Hair Removal Device User Manual](#)

Nood T16A

Nood Flasher Pro IPL Laser Hair Removal Device User Manual

Model: T16A | Brand: Nood

[Instructions](#) [Maintenance](#) [Overview](#) [Safety Information](#) [Troubleshooting](#) [Setup](#) [Specifications](#) [Operating](#) [Warranty & Support](#)

1. PRODUCT OVERVIEW

The Nood Flasher Pro is an advanced IPL (Intense Pulsed Light) laser hair removal device designed for at-home use. It offers a near-permanent solution for hair removal and includes an innovative Ageless Lens for anti-aging treatments. This device is engineered for pain-free and effective treatments, suitable for most skin tones and hair types.



Figure 1: Nood Flasher Pro IPL Laser Hair Removal Device and included accessories.

Key Features:

- **Powerful At-Home IPL:** Delivers clinical-grade results for smooth, hair-free skin with visible results after first use.
- **CryoSoothe Cooling Technology:** Ensures pain-free and skin-safe treatments by keeping the skin cool during use. Features 5 comfort levels.
- **Bonus Ageless Lens:** The first IPL device to also offer anti-aging treatments, improving the appearance of fine lines, dark spots, and pores.
- **Long-Term Value:** A one-time purchase with no refills or hidden costs, featuring MaxPulse technology and a lifetime flash replacement guarantee.
- **Versatile Use:** Safe for women, men, and all genders, effective on most skin tones for various body areas.

MEET THE FIRST IPL DEVICE THAT DELIVERS

long lasting hairless & ageless skin



Figure 2: The Nood Flasher Pro with its interchangeable Hairless Lens (for hair removal) and Ageless Lens (for skin rejuvenation).

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the Nood Flasher Pro device. Failure to follow these instructions may result in injury or device malfunction.

Important Precautions:

- **External Use Only:** This device is intended for external use on the skin only.
- **Skin and Hair Type Compatibility:** IPL technology requires a contrast between the melanin in your hair and the melanin in your skin. It is not safe or effective for all hair and skin types. Refer to the skin tone and hair color chart below to determine compatibility.
- **Skin Safety Test:** It is mandatory to complete a Skin Safety Test before your first full treatment. Detailed instructions for this test are included in the packaging.
- **Eye Protection:** Always wear the provided protective eyewear during operation to shield your eyes from the intense pulsed light.
- **Preparation:** Exfoliate and remove any hair from the treatment area before using the device.

- **Avoid Over-Treatment:** Do not treat the same area multiple times in a single session, especially on sensitive areas like the face.

AMERICA'S MOST TRUSTED IPL BRAND

works on most hair and skin types

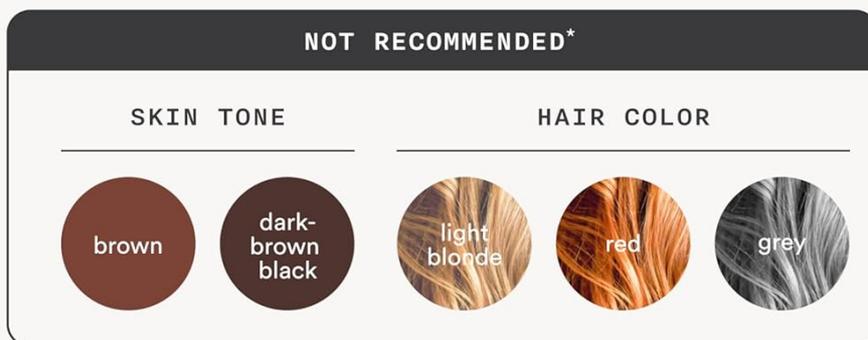
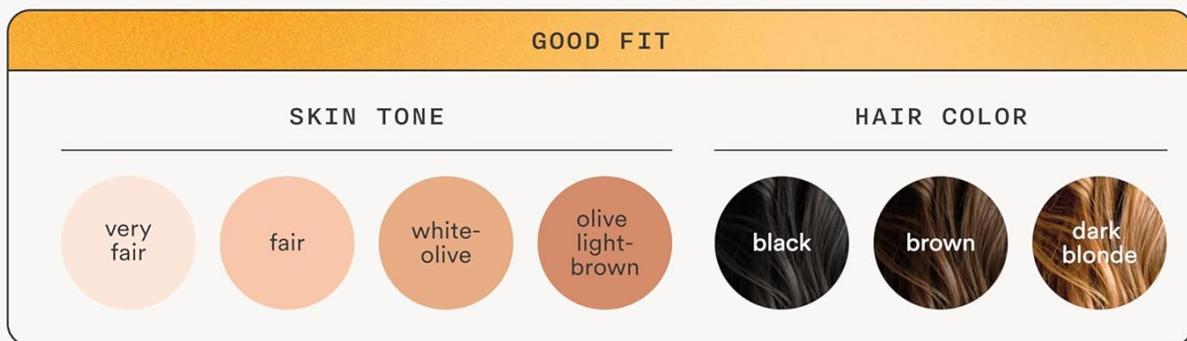


Figure 3: Skin tone and hair color compatibility chart for the Nood Flasher Pro. Not recommended for very dark skin tones or light blonde, red, or grey hair.

3. SETUP

Unboxing and Initial Preparation:

1. Carefully remove all components from the packaging: the Nood Flasher Pro device, power adapter, Hairless Lens, Ageless Lens, and protective eyewear.
2. Connect the power adapter to the device and plug it into a suitable electrical outlet. The device will power on.
3. Select the desired lens (Hairless or Ageless) and gently insert it into the device's treatment window until it clicks securely into place. Ensure the lens is fully seated.
4. Before your first use, perform a Skin Safety Test on a small, inconspicuous area of your skin according to the detailed instructions provided in the separate manual. Wait 24 hours to observe any adverse reactions.

4. OPERATING INSTRUCTIONS

Using Your Nood Flasher Pro:

1. **Prepare the Area:** Ensure the treatment area is clean, dry, and free of hair (shave, do not wax or pluck).
2. **Wear Eye Protection:** Always put on the provided protective eyewear before turning on the device.
3. **Power On:** Long-press the main power button on the device to turn it on. The display will illuminate.
4. **Select Intensity Level:** Use the power level adjustment button (often indicated by a wave or light icon) to select your desired intensity level from 1 to 5. Start with the lowest comfortable setting and gradually increase as your skin tolerates.
5. **Activate CryoSoothe Cooling:** The CryoSoothe cooling technology activates automatically when the device is on, providing a cool sensation to the skin during treatment.
6. **Apply to Skin:** Place the treatment window firmly and flat against the skin in the area you wish to treat. The device will only flash when proper skin contact is detected.
7. **Initiate Flash:** Press the flash button to deliver a pulse of light. For continuous treatment, some models may offer an auto-flash mode by holding the button or selecting a specific setting.
8. **Move to Next Area:** After each flash, lift the device slightly and move it to the next adjacent area. Avoid overlapping flashes excessively.
9. **Treatment Frequency:** For optimal hair removal results, use the device 2-3 times per week for the first 4-6 weeks, then reduce frequency as hair growth diminishes. For anti-aging, follow the specific regimen outlined in the detailed manual.
10. **Power Off:** After your session, long-press the main power button to turn off the device.

PAIN-FREE & COMFORTABLE TREATMENTS

CryoSoothe™ cooling technology



Figure 4: The CryoSoothe cooling technology in action, providing comfort during IPL treatment.

visible hair reduction after 1st use

BEFORE



AFTER



- o ingrown hairs
- o bumps + irritation
- o same day stubble
- o expensive laser treatments



Figure 5: Example of visible hair reduction after initial use of the Nood Flasher Pro.

visibly younger skin

JUST 10 MINUTES A WEEK

BEFORE



AFTER



- o reduce the look of fine lines and wrinkles
- o improve skin texture
- o stimulate collagen production
- o reduce appearance of dark spots and hyperpigmentation

Figure 6: Example of visibly younger skin results after using the Ageless Lens for skin rejuvenation.

5. MAINTENANCE

Cleaning and Storage:

- **Cleaning the Device:** After each use, ensure the device is powered off and unplugged. Gently wipe the treatment window and the device's exterior with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Lens Care:** If necessary, carefully clean the lens surface with a lint-free cloth. Avoid scratching the lens.
- **Storage:** Store the Nood Flasher Pro in its original packaging or a protective case in a cool, dry place, away from direct sunlight and out of reach of children.

6. TROUBLESHOOTING

If you encounter any issues with your Nood Flasher Pro, please refer to the following common troubleshooting steps:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power button not long-pressed.	Ensure the power adapter is securely connected and plugged into a working outlet. Long-press the power button for a few seconds.
Device does not flash.	Insufficient skin contact; incorrect skin tone/hair color; device malfunction.	Ensure the treatment window is fully and firmly pressed against the skin. Verify your skin tone and hair color are compatible with the device (refer to safety section). If issues persist, contact customer support.
Skin feels too warm or uncomfortable.	Intensity level is too high for your skin sensitivity.	Reduce the intensity level. Ensure the CryoSoothe cooling is active. If discomfort continues, discontinue use and consult a dermatologist.

For further assistance, please refer to the full user manual or contact Nood customer support.

7. SPECIFICATIONS

Feature	Detail
Brand	Nood
Model Number	T16A
UPC	850071754124
Package Dimensions	9.88 x 7.6 x 3.9 inches
Item Weight	3.13 Pounds
Manufacturer	Nood
Country of Origin	China
Technology	IPL (Intense Pulsed Light) with CryoSoothe Cooling
Lenses Included	Hairless Lens, Ageless Lens

8. WARRANTY & SUPPORT

Warranty Information:

The Nood Flasher Pro comes with a **lifetime flash replacement guarantee**, ensuring the longevity of your device's effectiveness. Additionally, it is backed by a **90-day satisfaction promise**.

For specific return policy details, please refer to your purchase documentation or the retailer's website. Typically, products are eligible for a refund or replacement within 30 days of purchase.

Customer Support:

For any questions, technical assistance, or warranty claims, please contact Nood customer support. Contact information can usually be found on the official Nood website or within your product packaging.

Sold by: Nood Beauty

© 2025 Nood. All rights reserved.