

## Morelax M893

# Morelax Mini Massage Gun M893 User Manual

Model: M893 | Brand: Morelax

## INTRODUCTION

---

Thank you for choosing the Morelax Mini Massage Gun. This portable deep tissue muscle percussion massager is designed to provide effective muscle relief and recovery. Please read this manual carefully before use to ensure safe and optimal operation of your device. Keep this manual for future reference.

## SAFETY INFORMATION

---

**WARNING: To reduce the risk of electric shock, fire, and personal injury, this product must be used in accordance with the following instructions:**

- Use for no more than 15 minutes at a time to prevent overheating. Allow the device to cool down before resuming use.
- Keep the massage gun away from ears to avoid potential hearing damage.
- Do not immerse the device in water or use it in wet environments.
- Do not operate the device if it has been damaged, dropped, or exposed to water.
- Keep out of reach of children and pets.
- Consult a physician before use if you have any medical conditions, are pregnant, or have a pacemaker.
- Do not use on open wounds, swollen areas, or areas with skin irritation.
- This device is not intended to diagnose, treat, cure, or prevent any disease or health condition.

## PRODUCT OVERVIEW

---

The Morelax Mini Massage Gun is a compact and powerful percussion massager designed for deep tissue

relief. It features a brushless motor, multiple speed settings, and interchangeable massage heads to target various muscle groups.



Image: The Morelax Mini Massage Gun, showcasing its compact design.

## Key Features

- **Versatile Deep Tissue Relief:** Equipped with 6 interchangeable massage heads and 4 adjustable speed levels.
- **Long Battery Life & Quick Charge:** 2500mAh battery, full recharge in 2.5 hours via USB-C, up to 5 hours of use.
- **Compact & Portable Design:** Ultra-lightweight, easily fits into bags for on-the-go muscle recovery.
- **Smart 15-Minute Auto Shut-Off:** Built-in safety timer prevents overheating and conserves battery life.
- **Ultra-Quiet Operation:** Brushless motor operates at less than 45dB.
- **7mm Amplitude:** Effective percussion for deep muscle penetration.

## 6 Personalized Massage Heads

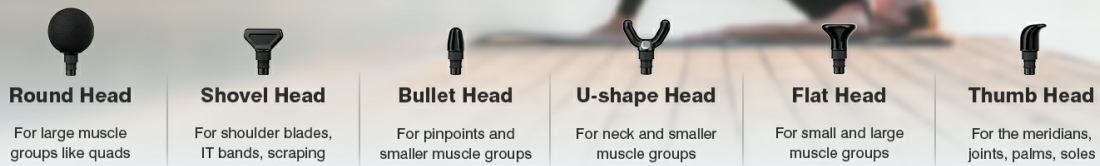


Image: Visual representation of the massage gun's quiet operation, 15-minute auto-off, and 7mm amplitude.

### Package Contents

- Morelax Mini Massage Gun
- 6 Interchangeable Massage Heads
- USB-C Charging Cable
- User Manual
- Storage Bag (may vary)

## Powerful and quiet Brushless Motor

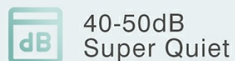


Image: A user sitting with the massage gun, illustrating the six different massage heads and their uses.

## SETUP

### Charging the Device

1. Before first use, fully charge the massage gun.
2. Connect the USB-C charging cable to the charging port at the bottom of the device.
3. Connect the other end of the cable to a compatible USB power adapter (not included) or a computer USB port.
4. The indicator lights will show charging status. A full charge takes approximately 2.5 hours.
5. Once fully charged, disconnect the charging cable.

# Conquer Muscle Cramps

The Portable Deep Tissue Tool!



Image: The massage gun being charged via a USB-C cable, highlighting its fast charging capability.

## Attaching Massage Heads

1. Ensure the device is turned off before attaching or removing massage heads.
2. Select the desired massage head.
3. Gently push the stem of the massage head into the opening at the front of the massage gun until it is securely in place.
4. To remove, gently pull the massage head straight out.

## OPERATING INSTRUCTIONS

---

### Power On/Off and Speed Adjustment

1. **Power On:** Press and hold the power button at the bottom of the device for 2 seconds to turn it on. The device will start at the lowest speed setting.
2. **Adjust Speed:** Briefly press the power button to cycle through the 4 available speed levels. The indicator lights will show the current speed setting.
3. **Power Off:** Press and hold the power button for 2 seconds to turn off the device.
4. **Auto Shut-Off:** The device will automatically turn off after 15 minutes of continuous use to prevent overheating and conserve battery.

### Using the Massage Heads

Select the appropriate massage head for the area you want to target:

- **Round Head:** For large muscle groups like quads, glutes, and hamstrings.
- **Shovel Head:** For shoulder blades, IT bands, and scraping.
- **Bullet Head:** For pinpoint and smaller muscle groups.
- **U-shape Head:** For neck and smaller muscle groups.
- **Flat Head:** For small and large muscle groups.
- **Thumb Head:** For impacting deep tissue, meridians, joints, palms, and soles.

Apply the massage gun to the desired muscle area. Start with the lowest speed setting and gradually increase the intensity as needed. Move the device slowly over the muscle, applying gentle pressure. Do not apply excessive pressure or use on bony areas.

Your browser does not support the video tag.

Video: An official product video demonstrating the portable design and various features of the Morelax Mini Massage Gun.

## MAINTENANCE

---

### Cleaning

- Always turn off and unplug the device before cleaning.
- Wipe the surface of the massage gun with a slightly damp cloth. Do not use abrasive cleaners or solvents.
- Clean massage heads with a mild soap and water solution, then rinse and air dry completely before reattaching.

### Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the device and all accessories in the provided storage bag to protect them from dust and damage.

## TROUBLESHOOTING

---

Problem	Possible Cause	Solution
Device does not turn on	Low battery; Device not powered on correctly	Charge the device fully; Press and hold the power button for 2 seconds
Device stops during use	15-minute auto shut-off activated; Low battery	Restart the device if desired; Recharge the device
Weak percussion force	Low battery; Incorrect speed setting	Recharge the device; Increase the speed level
Massage head falls off	Head not securely attached	Ensure the massage head is pushed firmly into place

## SPECIFICATIONS

---

Feature	Detail
Model Number	M893
Product Dimensions	3.15 x 1.77 x 5.91 inches
Item Weight	1.21 Pounds (0.55 Kilograms)
Battery Type	Nonstandard Battery (2500mAh)
Charging Time	2.5 Hours
Battery Average Life	Up to 5 Hours (depending on speed)

Feature	Detail
Power Source	Battery Powered (USB-C Charging)
Speed Levels	4
Amplitude	7mm
Noise Level	<45dB
Auto Shut-Off	15 Minutes
Material Type	Metal, Plastic, Silicone

## WARRANTY AND SUPPORT

---

### Warranty Information

Morelax products are manufactured to the highest quality standards. This product comes with a limited warranty against defects in materials and workmanship from the date of purchase. Please retain your proof of purchase for warranty claims.

The warranty does not cover damage caused by improper use, accidents, unauthorized repairs, or normal wear and tear. For specific warranty terms and conditions, please refer to the warranty card included with your product or contact customer support.

### Customer Support

If you have any questions, concerns, or require assistance with your Morelax Mini Massage Gun, please contact our customer support team. Our representatives are available to help you with product inquiries, troubleshooting, and warranty services.

For the most up-to-date contact information, please visit the official Morelax website or refer to the contact details provided on your product packaging.