

H2ofloss HF-9

H2ofloss HF-9 Water Flosser Instruction Manual

Model: HF-9 | Brand: H2ofloss

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using this appliance.

- Do not use while plugged in or charging.
- Keep away from children.
- Do not use if you have certain medical conditions without consulting a healthcare professional.
- Always unplug after use.
- Do not submerge the main unit in water.
- Allow the unit to cool before storing if it has been in continuous operation.

2. PRODUCT OVERVIEW

The H2ofloss HF-9 Water Flosser is designed to provide effective oral hygiene through high-frequency pulse water flow. It features an 800ML water tank for extended use, 5 adjustable pressure settings, and a variety of nozzles to cater to different oral care needs. The magnetic handle design ensures convenient storage, and the unit includes intelligent overheat protection for safe operation.



Figure 2.1: H2ofloss HF-9 Water Flosser with included nozzles.

Product Details



Figure 2.2: Detailed view of the H2ofloss HF-9 Water Flosser components and dimensions.

Key Features:

- **High-Frequency Pulse Technology:** Delivers up to 1200 pulses per minute for effective cleaning.
- **Large 800ML Water Tank:** Provides continuous flossing without frequent refills.
- **5 Adjustable Water Pressure Settings:** Ranging from 5-110 PSI for personalized comfort.
- **Hygienic Nozzle Storage:** Integrated compartment holds up to 20 tips, keeping them clean and organized.
- **Magnetic Handle Design:** Allows for seamless docking and easy storage of the handle.
- **Anti-leak Design:** Features an IPX7 waterproof structure with double-sealing rings.
- **Smart Overheat Protection:** Automatically stops operation if overheated, restarting after 2 hours.

3. PACKAGE CONTENTS

Your H2ofloss HF-9 Water Flosser package includes the following items:

- H2ofloss HF-9 Water Flosser Main Unit
- Standard Tips (x4) for general cleaning
- Periodontal Tips (x2) for gum pockets
- Orthodontic Tips (x2) for braces
- Tongue Scrapers (x2) for tongue cleaning
- Nasal Cleaning Tips (x2)
- Whitening Tip (x1)
- Power Adapter
- Instruction Manual



Figure 3.1: Assortment of multifunctional nozzles included.

4. SETUP

Follow these steps to set up your H2ofloss HF-9 Water Flosser for first use:

1. **Fill the Water Tank:** Remove the 800ML water tank from the base. Fill it with lukewarm water. Ensure the water temperature is comfortable for oral use.
2. **Attach Nozzle:** Select the appropriate nozzle for your needs and insert it firmly into the handle until it clicks into place.
3. **Place Tank:** Return the filled water tank to the base, ensuring it is seated correctly.
4. **Connect Power:** Plug the power adapter into a suitable electrical outlet.



Figure 4.1: Visual guide for initial setup and usage steps.

5. OPERATING INSTRUCTIONS

To use your H2ofloss HF-9 Water Flosser:

1. **Adjust Pressure:** Use the control dial on the handle to select your desired water pressure setting (5-110 PSI). It is recommended to start with a lower setting and gradually increase.
2. **Position Nozzle:** Lean over a sink and place the tip of the nozzle into your mouth, directing it towards your gumline at a 90-degree angle.
3. **Turn On:** Press the ON/OFF button on the main unit to start the water flow.
4. **Floss:** Guide the tip along your gumline and between your teeth, pausing briefly between each tooth. Allow the water to flow into the sink.
5. **Pause Function:** The handle features a pause button to temporarily stop the water flow without turning off the unit.
6. **360° Rotatable Nozzle:** The nozzle can be rotated 360 degrees to easily reach all areas of your mouth.
7. **Turn Off:** Once finished, press the ON/OFF button on the main unit to stop the water flow.

Smart Handle Design

with adjustable pressure settings



Figure 5.1: Smart handle design with adjustable pressure settings.

Multifunction nozzles tackle a variety of oral issues



Deeply clean the gumline

Clean Orthodontic Gums

Prevent Bad Breath

Figure 5.2: Recommended water pressure levels for various cleaning needs.



Figure 5.3: Pause function and 360° rotatable nozzle for enhanced control.

6. MAINTENANCE AND CLEANING

Regular cleaning and maintenance will ensure the longevity and hygienic operation of your water flosser.

- **Empty and Dry Tank:** After each use, empty any remaining water from the tank. It is recommended to remove the tank and allow it to air dry to prevent mold and mildew buildup.
- **Clean Nozzles:** Remove the nozzle from the handle. Rinse the nozzle under running water. For deeper cleaning, nozzles can be soaked in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes, then rinsed thoroughly.
- **Store Nozzles:** Utilize the hygienic nozzle storage compartment to keep your nozzles clean and protected from dust.
- **Wipe Unit:** Wipe the exterior of the main unit with a soft, damp cloth. Do not use abrasive cleaners or submerge the unit.
- **Drainage Hole:** The unit features a drainage hole at the bottom to help keep internal parts dry. Ensure it is not obstructed.

Two Reservoir Design

Built-in Nozzle
storage tank

800 ML Large
Water Tank



Figure 6.1: The 800ML water tank and integrated nozzle storage compartment.



Figure 6.2: Hygienic tips storage compartment for up to 20 nozzles.

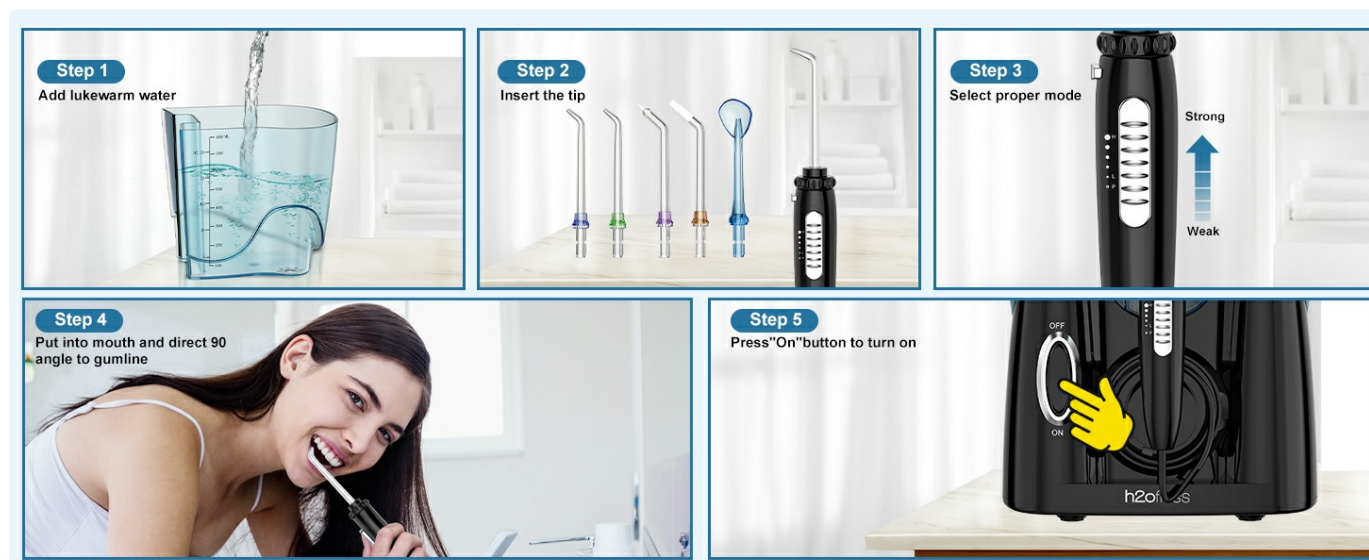


Figure 6.3: Anti-leakage design with double-sealing rings.

7. TROUBLESHOOTING

If you encounter any issues with your H2ofloss HF-9 Water Flosser, please refer to the following common problems and solutions:

- **Unit stops working / Overheating:** The unit is equipped with smart overheat protection. If it stops working, it may have overheated. Allow the unit to cool down for at least 2 hours before attempting to restart it.
- **Reduced Water Pressure or No Water Flow:**
 - Ensure the water tank is filled and properly seated.
 - Check if the nozzle is securely attached to the handle.
 - Inspect the nozzle for any blockages and clean if necessary.
 - Verify that the water hose is not kinked or obstructed.
- **Water Leakage:**
 - Ensure the water tank is correctly positioned and sealed.
 - Check the nozzle connection for a secure fit.
 - Inspect the O-rings and seals for any damage or displacement.
- **Unit is Noisy:** Some operational noise is normal. Ensure the unit is placed on a stable, flat surface to minimize vibration.

If the problem persists after attempting these solutions, please contact customer support.

8. SPECIFICATIONS

Feature	Detail
Product Dimensions	6.5 x 4.3 x 8.27 inches

Feature	Detail
Weight	4.03 Pounds
Model Number	HF-9
Power Source	Corded Electric
Water Tank Capacity	800ML
Water Pressure Range	5-110 PSI
Pulse Rate	Up to 1200 pulses/minute
Wattage	39 watts

9. WARRANTY AND CUSTOMER SUPPORT

The H2ofloss HF-9 Water Flosser typically comes with a **2-year warranty** from the date of purchase, covering manufacturing defects. Please retain your proof of purchase for warranty claims.

For technical assistance, troubleshooting beyond this manual, or warranty inquiries, please contact H2ofloss customer support through their official website or the retailer where the product was purchased.

Manufacturer: H2ofloss

© 2025 H2ofloss. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question