



Manuals.plus /

› Tao Clean /

› TAO Clean UV Sanitizing Sonic Toothbrush and Cleaning Station User Manual

## Tao Clean BA-0130-LIL-US

# TAO Clean UV Sanitizing Sonic Toothbrush and Cleaning Station User Manual

Model: BA-0130-LIL-US | Brand: Tao Clean

## 1. INTRODUCTION

This manual provides essential information for the safe and effective use of your TAO Clean UV Sanitizing Sonic Toothbrush and Cleaning Station. Please read it thoroughly before first use and retain it for future reference.



Image 1: TAO Clean UV Sanitizing Sonic Toothbrush and Cleaning Station.

## 2. PRODUCT FEATURES

- **Sonic Technology:** Utilizes 40,000 tiny brush strokes per minute for effective plaque removal and gentle polishing of surface stains.
- **UV Sanitizing Station:** The patented docking station dries, sanitizes, and protects your toothbrush. Powerful UV-C light eliminates 99.9% of bacteria, ensuring a clean, sanitized toothbrush for every use.
- **Dual Speed Settings:** Choose between 'Super Clean' for a thorough clean and 'Sensitive Clean' for gentle brushing, allowing for a customized experience.
- **2-Minute Timer with Quadrant Pacer:** Features a dentist-recommended 2-minute timer. A brief pause every 30 seconds reminds you to switch to another section of your mouth, ensuring comprehensive cleaning.
- **Ergonomic Design:** The handle is designed to fit comfortably in your hand for ease of use.

**UV-C** LIGHT KILLS  
99.9% OF GERMS  
AND VIRUSES.



*including but not limited to*

HUMAN CORONAVIRUS • FLU  
COMMON COLD • STREP

*\*independently tested*

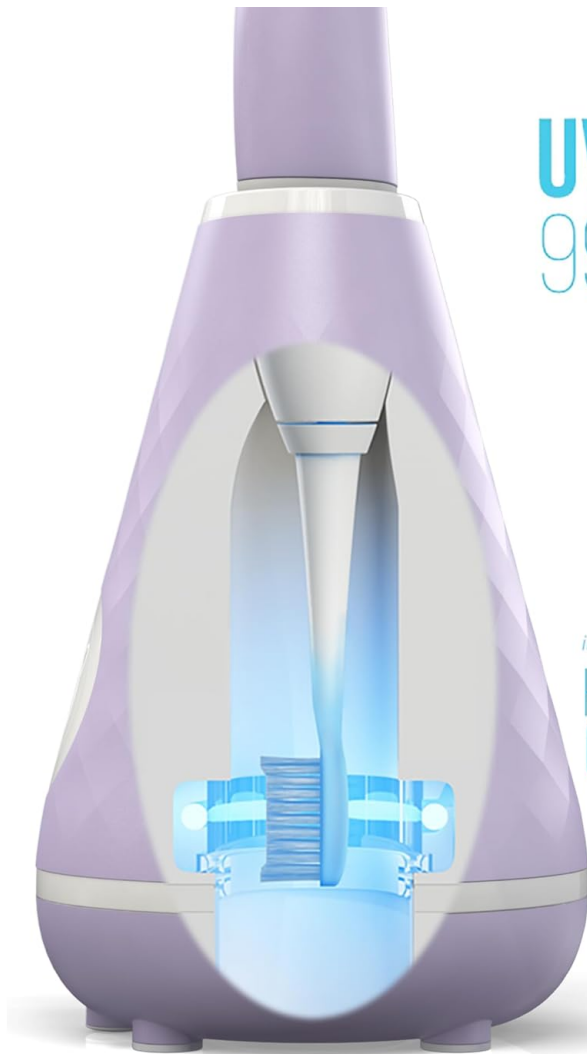


Image 2: Illustration of UV-C light sanitizing the toothbrush head within the station.



Image 3: Visual representation of the dual speed settings for customized brushing.

### 3. SETUP

1. **Unpack:** Carefully remove all components from the packaging.
2. **Charge:** Connect the charging station to a power outlet. Place the toothbrush handle into the station. The station will begin charging the toothbrush and initiate the UV sanitization cycle.
3. **Initial Sanitize:** Allow the toothbrush to fully charge and complete its first sanitization cycle before initial use.

### 4. OPERATING INSTRUCTIONS

#### 4.1 Brushing Your Teeth

1. **Apply Toothpaste:** Moisten the brush head and apply a small amount of toothpaste.

2. **Power On:** Press the power button on the toothbrush handle to turn it on.
3. **Select Mode:** Press the mode button to cycle between 'Super Clean' and 'Sensitive Clean' settings.
4. **Brush:** Place the brush head against your teeth at a slight angle towards the gum line. Move the brush slowly across your teeth, allowing the sonic vibrations to clean each surface.
5. **Quadrant Pacer:** The brush will briefly pause every 30 seconds. This is your cue to move to the next quadrant of your mouth (upper right, upper left, lower right, lower left).
6. **Auto Shut-off:** The toothbrush will automatically turn off after 2 minutes, indicating the completion of the recommended brushing time.



Image 4: Proper brushing technique with the sonic toothbrush.

#### 4.2 Using the UV Sanitizing Station

1. **Rinse Brush Head:** After brushing, rinse the brush head thoroughly under running water to remove any remaining toothpaste.
2. **Place in Station:** Insert the toothbrush handle, brush head down, into the UV Sanitizing Station.
3. **Sanitization Cycle:** The station will automatically initiate a UV-C sanitization cycle, indicated by a blue light. This cycle also includes a drying function to prevent mildew buildup.



Image 5: Placing the toothbrush into the station for sanitization and charging.

## 5. MAINTENANCE

- **Replace Brush Head:** Replace the brush head every 3 months or sooner if bristles become frayed.
- **Clean Station:** Periodically clean the inside of the UV Sanitizing Station to ensure optimal performance and hygiene.
- **Clean Handle:** Wipe the toothbrush handle with a damp cloth as needed.

## 6. SPECIFICATIONS

Feature	Detail
Model Number	BA-0130-LIL-US
UPC	810099660012
Package Dimensions	7.8 x 6.89 x 4.88 inches
Weight	1.3 Pounds
Batteries	1 Nonstandard Battery (included)
Power Source	Battery Powered
Item Firmness	Soft
Age Range	Adult

## 7. SAFETY INFORMATION

- Do not use this product for purposes other than oral hygiene.
- Keep out of reach of children.
- Replace brush head every 3 months or when bristles become frayed.
- Do not use the toothbrush if the charging cord or plug is damaged.
- Always supervise children under 12 years old while brushing.

## 8. WARRANTY AND SUPPORT

While specific warranty details are not provided in the product information, some customer feedback indicates a potential lifetime warranty. For official warranty information and customer support, please refer to the manufacturer's website or contact Tao Clean directly.