

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [MySmile](#) /

› [MySmile LP211 Cordless Water Flosser User Manual](#)

MySmile LP211

MySmile LP211 Cordless Water Flosser User Manual

Model: LP211

1. INTRODUCTION

Thank you for choosing the MySmile LP211 Cordless Water Flosser. This device is designed to enhance your oral hygiene routine by effectively removing plaque, food debris, and bacteria from between teeth and below the gumline. Please read this manual thoroughly before first use to ensure proper operation and maintenance.



Image 1: MySmile LP211 Cordless Water Flosser, showing the main unit, 8 different replacement tips, a storage pouch, and product packaging.

2. SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- **Charging:** Only charge the device using the provided Type-C USB charging cord (5V 1A standard). Using an incorrect charger may damage the battery or device.
- **Water Temperature:** Do not fill the reservoir with water hotter than 40°C (104°F).
- **Water Additives:** Do not add iodine, saline, or water-insoluble concentrated essential oils to the reservoir, as these can reduce product performance and shorten its lifespan.
- **Children:** This device is not a toy. Supervise children during use, especially when using the Child Mode.
- **Medical Conditions:** Consult your dentist or physician before using if you have any oral medical conditions or have undergone recent oral surgery.

- **Disassembly:** Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- MySmile LP211 Cordless Water Flosser Unit
- 4x Classic Jet Tips
- 1x Plaque Seek Tip
- 1x Orthodontic Tip
- 1x Periodontal Pocket Tip
- 1x Tongue Cleaner Tip
- USB-C Charging Cable
- PU Storage Bag
- User Guide (this document)



Image 2: Visual representation of the MySmile LP211 package contents.

4. PRODUCT FEATURES AND COMPONENTS

The MySmile LP211 Water Flosser offers advanced features for effective oral care:

- **Powerful Cleaning:** Utilizes 40-150 PSI water pressure and 2000-2800 pulses per minute to remove plaque and debris.
- **Long-lasting Battery:** Equipped with a 2000mAh rechargeable battery, providing up to 30 days of use on a 3-hour charge.
- **5 Cleaning Modes:** Strong, Normal, Soft, Pulse, and Child modes for personalized oral care.
- **IPX7 Waterproof:** Allows for safe use in the shower and easy cleaning.
- **8 Replacement Tips:** Includes various tips for different dental needs.

The advertisement features a pink MySmile LP211 Cordless Water Flosser as the central product. It is positioned on a light-colored wooden tray, next to a stack of three white folded towels. To the right of the device, eight different replacement tips are displayed in a row. The background is a dark blue gradient. Text on the right side lists features: 'Family-Friendly', '8 Replacement tips', '5 cleaning modes', and 'cater to everyone's needs-kids, adults, and seniors alike'. Below this, it says 'Trusted by Thousands of Families!'. At the bottom, a white banner contains the text 'AS SEEN IN:' followed by logos for People, Health, TikTok, ABC, and CNN.

MySmile® | LP211 Cordless Water Flosser

Family-Friendly
8 Replacement tips
5 cleaning modes

cater to everyone's
needs-kids, adults, and
seniors alike

Trusted by Thousands of Families!

AS SEEN IN:

People Health TikTok abc CNN

Image 3: Overview of the MySmile LP211's key features, including multiple tips and cleaning modes.

4.1. Nozzle Types and Their Uses

- **Standard Tip:** For general cleaning between teeth and along the gumline.
- **Plaque Jet Tip:** Designed for areas prone to plaque accumulation.
- **Orthodontic Tip:** Features soft bristles for cleaning around braces and dental work.
- **Periodontal Pocket Tip:** A soft, rubberized tip for gently cleaning deep pockets below the gumline.
- **Tongue Scraper Tip:** Shaped to effectively clean the tongue for fresh breath.

5 Types of Tips

Standard Tip

creates an ultra-fine jet stream for cleaning between teeth and along the gumline

Plaque Jet Tip

Use for dental areas that are susceptible to plaque accumulation

Orthodontic Tip

features soft bristles for dislodging plaque and debris from braces and dental work

Periodontal Tip

features a soft, rubberized tip for gently cleaning deep pockets below the gumline

Tongue Scraper Tip

is ideally shaped to clean the tongue

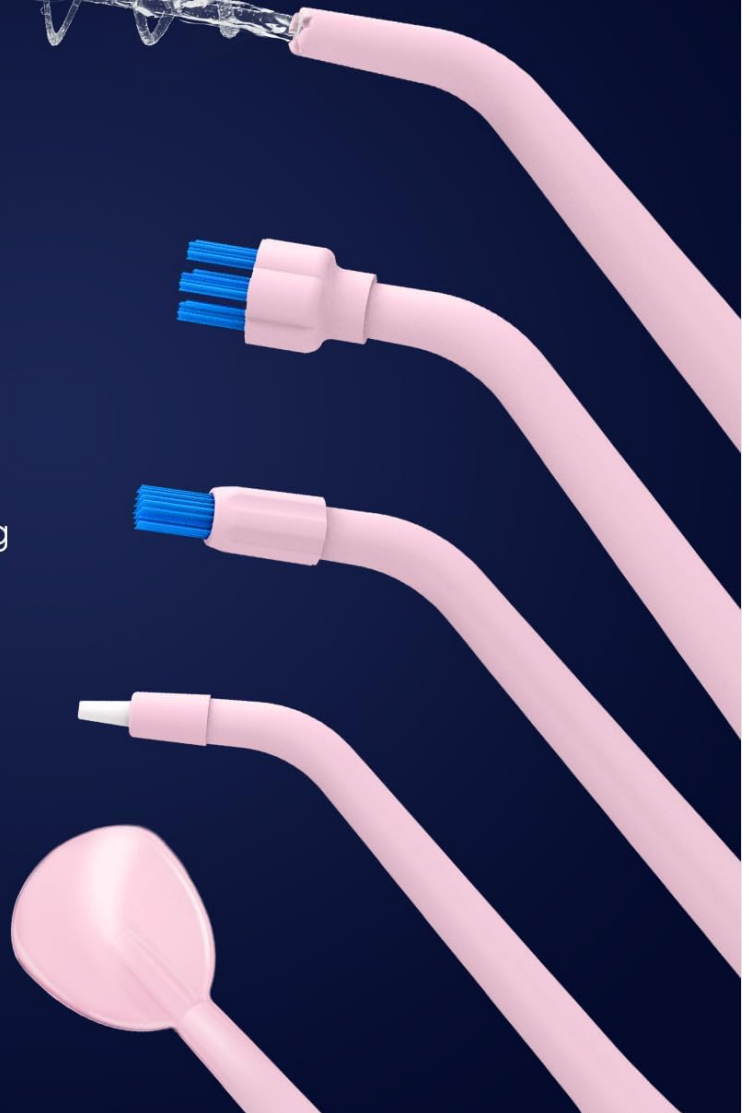


Image 4: Detailed view of the five distinct nozzle types and their specific applications.

5. CHARGING THE DEVICE

Before first use, fully charge the water flosser. A full charge takes approximately 3 hours and provides up to 30 days of cordless use.

1. Locate the charging port on the device, usually covered by a rubber lid.
2. Open the rubber lid and insert the USB-C charging cable.
3. Connect the other end of the USB-C cable to a 5V 1A USB power adapter (not included) or a compatible USB port.
4. The indicator lights on the device will show charging progress. All lights will illuminate when fully charged.
5. After charging, ensure the rubber lid is securely closed to maintain the IPX7 waterproof rating.

2,000mAh Rechargeable Battery

Built in long-lasting battery



3 Hours

Fully Charge



Up to 30 Days

Flossing/Charge

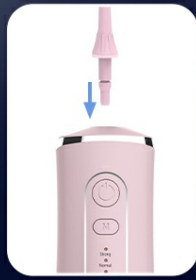


Image 5: Illustration of the internal battery and charging capabilities of the MySmile LP211.

6. SETUP AND FIRST USE

1. **Attach Nozzle:** Insert your desired nozzle into the hole on the top of the flosser until it clicks into place.
2. **Fill Reservoir:** Open the water inlet cover on the side of the device. Fill the reservoir with lukewarm water. Close the cover securely.
3. **Positioning:** Lean over a sink and place the nozzle tip into your mouth, aiming it at your gumline at a 90-degree angle.
4. **Power On:** Press the Power button to turn on the device.
5. **Select Mode:** Press the Mode button to cycle through the 5 cleaning modes (Strong, Normal, Soft, Pulse, Child) until you find your preferred setting. For first-time users, start with the "Soft" or "Child" mode.
6. **Begin Flossing:** Keep your mouth slightly open to allow water to flow into the sink. Guide the water stream along your gumline, pausing briefly between teeth.
7. **Power Off:** Press the Power button again to turn off the device when finished.

Easy to Use



Installation of
nozzles



Open the water
inlet



Clean your teeth



Switch the mode
that suits you



Timed cleaning

Image 6: Easy-to-follow steps for setting up and operating the water flosser.

7. OPERATING MODES

The MySmile LP211 offers five distinct cleaning modes to suit various oral care needs:

- **Strong Mode:** Provides a powerful water stream for a deep, thorough clean. Ideal for users accustomed to water flossing.
- **Normal Mode:** A balanced pressure setting for everyday use and effective plaque removal.
- **Soft Mode:** Gentle pressure for sensitive teeth and gums, or for new users.
- **Pulse Mode:** Delivers intermittent bursts of water to massage gums and improve circulation.
- **Child Mode:** An extra-gentle setting specifically designed for children or individuals with very sensitive gums, helping to build healthy oral care habits.

You can adjust the water pressure freely within a range of 40-150 PSI to match your comfort and oral care requirements.

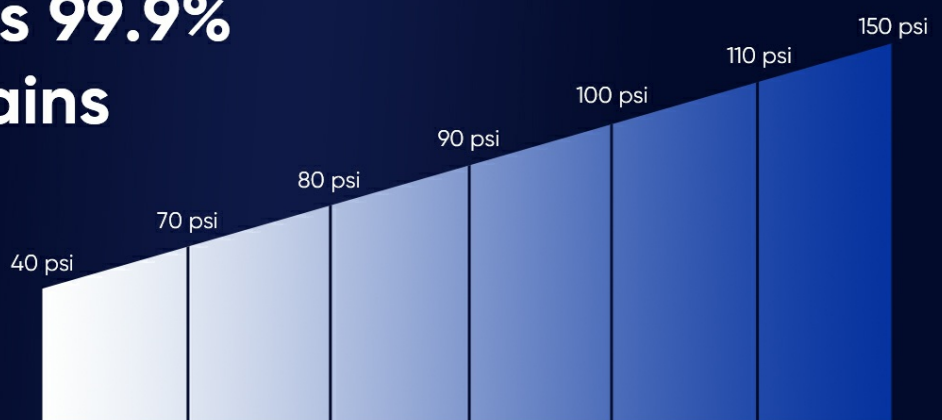
Deep cleans 99.9% of Tooth Stains

Best choice for

people with

sensitive

teeth



Freely adjust the water pressure (40 psi-150 psi) as your oral care need

Image 7: Illustration of the adjustable water pressure settings.

Gentle Startup

Why MySmile Water Flosser Stands Out!

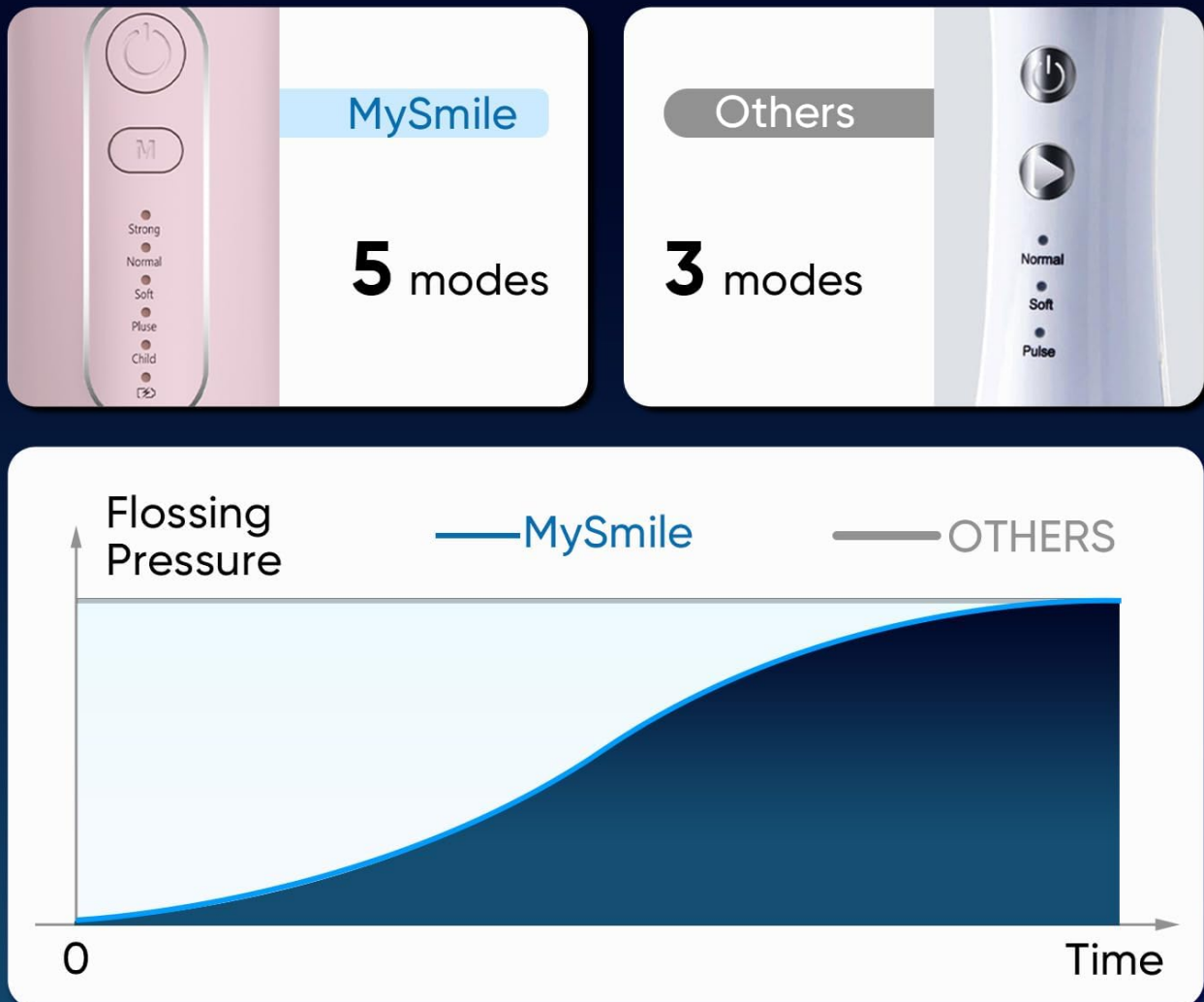


Image 8: Comparison of MySmile's 5 cleaning modes and gentle startup feature.

8. MAINTENANCE AND CLEANING

Regular cleaning and proper maintenance will ensure the longevity and hygiene of your water flosser.

1. **Empty Reservoir:** After each use, empty any remaining water from the reservoir.
2. **Clean Reservoir:** The bottom of the flosser can be completely removed to allow the entire water reservoir to air out and prevent mold or mildew. Press the release button and pull the bottom off. Wash the reservoir with warm, soapy water.
3. **Clean Nozzle:** Remove the nozzle by pressing the release button on the side of the unit. Rinse the nozzle under running water.
4. **Clean Unit Exterior:** Wipe the main unit with a soft, damp cloth. Do not immerse the main unit in water, except for rinsing the exterior.
5. **Drying:** Allow all parts to air dry completely before reassembling or storing.
6. **Storage:** Store the water flosser in a cool, dry place, ideally in the provided PU storage bag, when not in use.

The device is IPX7 waterproof, meaning it can withstand immersion in water up to 1 meter for 30 minutes. However, always ensure the charging port's rubber lid is securely closed to prevent water ingress during use or cleaning.



Image 9: The MySmile LP211's IPX7 waterproof feature.

9. TROUBLESHOOTING

If you encounter any issues with your MySmile LP211 Water Flosser, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery; device not charged.	Charge the device fully (approx. 3 hours).
Weak water pressure.	Low battery; clogged nozzle; incorrect mode selected.	Charge the device; clean the nozzle; select a stronger mode (e.g., Normal or Strong).
Water leaks from the device.	Reservoir not properly sealed; nozzle not securely attached; charging port cover open.	Ensure reservoir cover is closed tightly; re-attach nozzle until it clicks; close charging port rubber lid.
Device stops working after a short time.	Low battery.	Recharge the device.

If the problem persists, please contact MySmile customer support for further assistance.

10. SPECIFICATIONS

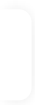
- **Model Number:** LP211
- **Power Source:** Battery Powered
- **Battery:** 2000mAh, 12V (included)
- **Charging Time:** Approximately 3 hours
- **Operating Time:** Up to 30 days per charge
- **Water Pressure:** 40-150 PSI

- **Water Pulses:** 2000-2800 times per minute
- **Waterproof Rating:** IPX7
- **Package Dimensions:** 9.96 x 3.82 x 3.66 inches
- **Item Weight:** 15.52 ounces
- **Manufacturer:** MySmile

11. WARRANTY AND SUPPORT


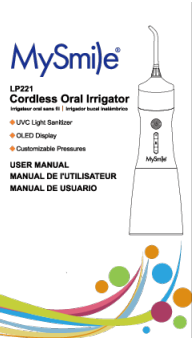


MySmile is committed to providing high-quality products. For any questions, concerns, or warranty claims regarding your LP211 Cordless Water Flosser, please contact MySmile customer support. Refer to the product packaging or the official MySmile website for contact information.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.



© 2025 MySmile. All rights reserved.

Related Documents - LP211

	<p>MySmile LP211 Portable Oral Irrigator User Manual</p> <p>User manual for the MySmile LP211 Portable Oral Irrigator, detailing specifications, operating guide, troubleshooting, and cleaning instructions. Learn how to use and maintain your oral irrigator for optimal dental hygiene.</p>
	<p>MySmile LP221 Cordless Oral Irrigator User Manual</p> <p>Comprehensive user manual for the MySmile LP221 Cordless Oral Irrigator, detailing features, operation, maintenance, and troubleshooting for optimal oral hygiene.</p>
	<p>MySmile Powerful Cordless Water Dental Flosser User Manual</p> <p>User manual for the MySmile Powerful Cordless Water Dental Flosser, featuring a 350ML capacity, OLED display, 5 pressure modes, and 8 replaceable jet tips for home and travel use.</p>
	<p>MySmile LP233 Pro Electric Power Flosser User Manual</p> <p>Comprehensive user manual for the MySmile LP233 Pro Electric Power Flosser, detailing setup, usage, safety precautions, and troubleshooting. Features UVC sanitizer, water heating, and customizable pressure.</p>

[\[pdf\]](#) Guide Warranty

LP211 BK MS IN User Guide MySmile Water Dental Flosser for Teeth Cordless Oral Irrigator 5 Cleaning Modes 4 Replaceable Jet Tips IPX 7 Waterproof USB Rechargeable Picks with PU Bag Black Health Household C1ZUyApFPnL m media amazon images I |||



1. If you have oral or dental surgery in the past two months, please consult your dentist before using this product, and follow your dentist s advice. 2. You may experience minor gum bleeding when using this product for the first time because your gums are being stimulated by the jet s pressure. You...

lang:en score:22 filesize: 5.95 M page_count: 1 document date: 2023-06-30

[MySmile LP211 Portable Oral Irrigator User Manual](#)

User manual for the MySmile LP211 Portable Oral Irrigator, detailing specifications, operating guide, troubleshooting, and cleaning instructions. Learn how to use and maintain your oral irrigator for optimal dental hygiene.

lang:en score:21 filesize: 2.83 M page_count: 1 document date: 2023-11-30

