

[Manuals.plus](#) /

> [APXPF](#) /

> APXPF Shiatsu Neck and Back Massager with Heat (Model A-Gray-1N) - Instruction Manual

APXPF A-Gray-1N

APXPF Shiatsu Neck and Back Massager with Heat

MODEL: A-GRAY-1N

Introduction

The APXPF Shiatsu Neck and Back Massager with Heat is designed to provide deep tissue massage for various body parts. This versatile massager targets areas such as the neck, back, shoulders, legs, and feet, offering relief from tension and discomfort. Its ergonomic design ensures comfortable use in multiple environments, including home, office, and vehicle.

Key features include:

- **Multi-Purpose Design:** Effective for neck, back, shoulders, legs, and feet.
- **Heat Function:** Soothing heat enhances relaxation and pain relief.
- **Premium Quality Materials:** Durable and skin-friendly synthetic fabric and plastic/metal alloy construction.
- **Ergonomic Shape:** Contours comfortably to the body for optimal massage.
- **Portable:** Compact and lightweight for easy transport and use anywhere.

Safety Information

- Avoid use on broken skin.
- Do not use during pregnancy.
- Ensure the device is used on a stable surface.
- Consult a physician before use if you have any pre-existing medical conditions.
- Do not immerse the device in water or use it in excessively humid environments.
- Keep out of reach of children.

Package Contents

- APXPF Shiatsu Neck and Back Massager (Model A-Gray-1N)
- AC Power Adapter (for home/office use)

- Car Adapter (for in-car use)
- User Manual



Image: The massager unit, along with its AC power adapter, car adapter, and user manual, as included in the package.

Setup

1. **Connect Power:** Insert the appropriate power adapter (AC for home/office or car adapter for vehicle) into the massager's power port. Ensure a secure connection.
2. **Positioning:** Place the massager on the desired body area. For neck and back, use the integrated elastic strap to secure it to a chair or car seat for stability.



Image: The massager highlighting its adjustable Velcro straps for securing to a chair and the included car adapter for portable use.

Operating Instructions

1. **Power On/Off:** Press and hold the power button (usually marked with a power symbol) to turn the massager on or off. The heat function typically activates automatically upon power-on.
2. **Adjust Massage Direction:** Press the direction button (often a circular arrow symbol) to reverse the rotation direction of the massage nodes.
3. **Adjust Intensity/Speed:** Use the intensity button (usually marked with bars) to cycle through different massage speeds (low, medium, high).
4. **Toggle Heat Function:** Press the heat button (often a wavy lines symbol) to turn the heat function on or off.

Easy to use with multi-function button



Push firmly until you hear a 'click' or 'snap' sound. This means the connection is successful.

Image: A detailed view of the massager's control panel, illustrating the power, direction, intensity, and heat buttons for easy operation.

Usage Areas

- **Neck:** Position the massager behind your neck while seated or lying down.
- **Back:** Place it against your lower, middle, or upper back, securing with the strap if needed.
- **Shoulders:** Adjust to target shoulder muscles.
- **Legs & Feet:** Rest your calves or feet on the massager for targeted relief.

Relax any part of your body

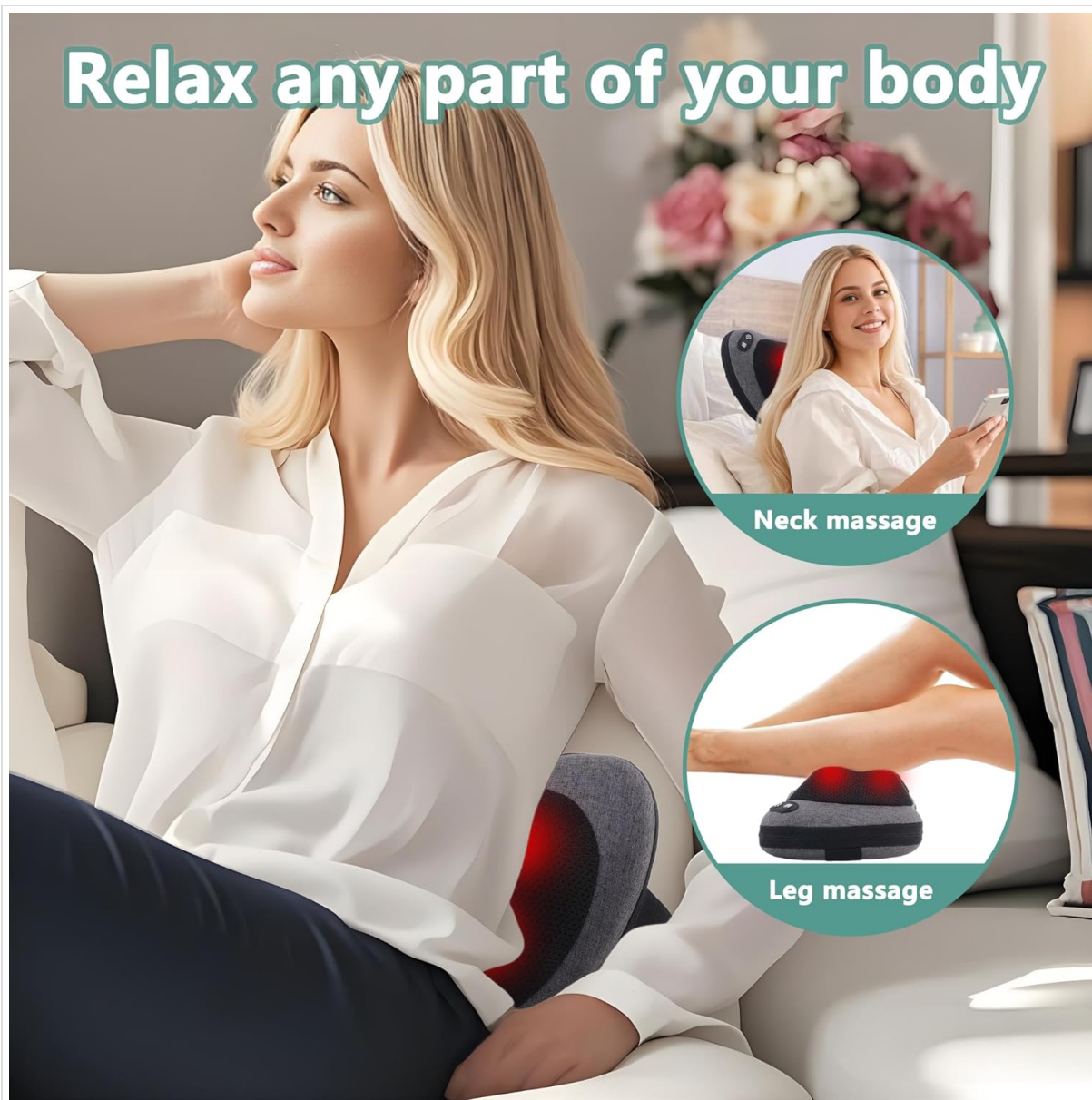


Image: A woman demonstrating the versatile use of the massager on her neck while relaxing and on her leg for calf massage.

Official Product Video: APXPF Neck and Back Massager

Video: A short demonstration by APXPF showcasing the neck and back massager in action, highlighting its features and ease of use.

Maintenance

- **Cleaning:** Disconnect the massager from the power source before cleaning. Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or immerse in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and moisture.
- **Care:** Avoid placing heavy objects on the massager. Do not attempt to disassemble or repair the unit yourself.

Troubleshooting

- **Device not turning on:** Ensure the power adapter is securely connected to both the massager and a working power outlet/car socket.

- **Heat function not working:** Verify that the heat button has been pressed and the indicator light (if any) is on. Allow a few minutes for the heat to build up.
- **Massage nodes stop or change direction unexpectedly:** The massager may have an automatic shut-off feature after a certain period (e.g., 15 minutes) to prevent overheating. Allow it to cool down before restarting.

Specifications

Specification	Detail
Product Dimensions	12.3 x 7.1 x 4.2 inches
Item Model Number	A-Gray-1N
Item Weight	3.22 Pounds (1.46 Kilograms)
Brand	APXPF
Power Source	Corded Electric
Material	Synthetic Fabric, Plastic/Metal Alloy
Use For	Back, Neck

Warranty and Support

For any product inquiries, warranty claims, or technical support, please contact APXPF customer service through the retailer where the product was purchased. Please refer to the included user manual for specific warranty details and contact information.

Legal Disclaimer

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for personal, non-medical use only.
