

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

› [Nekteck](#) /

› [Nekteck Shiatsu Foot & Neck/Shoulder Massager User Manual](#)

Nekteck Shiatsu Foot & Neck/Shoulder Massagers

Nekteck Shiatsu Foot & Neck/Shoulder Massager User Manual

Comprehensive Guide for Safe and Effective Operation

INTRODUCTION

This manual provides essential information for the proper use, setup, and maintenance of your Nekteck Shiatsu Foot Massager and Nekteck Shiatsu Neck/Shoulder Massager. Please read all instructions carefully before operating the devices to ensure safety and optimal performance. Keep this manual for future reference.

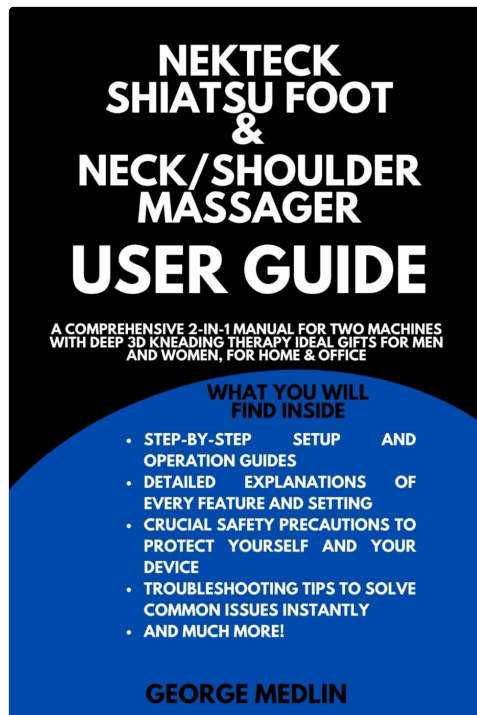


Image: Front cover of the Nekteck Shiatsu Foot & Neck/Shoulder Massager User Guide. This image displays the title and branding of the manual itself.

IMPORTANT SAFETY INFORMATION

To reduce the risk of electric shock, burns, fire, or injury, always follow these basic safety precautions:

- Always unplug the massager from the electrical outlet immediately after use and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store the massager where it can fall or be pulled into a tub or sink.
- Do not immerse in water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on sensitive skin areas or where circulation is poor.
- Do not use on infants, invalids, or sleeping or unconscious persons.
- Do not crush. Avoid sharp folds.

SETUP INSTRUCTIONS

Unpacking

1. Carefully remove the massager from its packaging.
2. Inspect the device for any signs of damage. If damaged, do not use and contact customer support.
3. Ensure all components are present: massager unit, power adapter.

Power Connection

- Connect the power adapter cable to the massager's power input port.
- Plug the power adapter into a suitable electrical outlet (110-240V AC, 50/60Hz).
- Ensure the connection is secure.

OPERATING INSTRUCTIONS

General Operation (Foot Massager)

1. Place the foot massager on a flat, stable surface.
2. Sit comfortably and place your feet into the designated footrests.
3. Press the Power button to turn on the device. The default massage mode will activate.

4. Use the Mode button to cycle through different massage programs (e.g., Shiatsu, rolling, air compression).
5. Press the Heat button to activate or deactivate the heating function.
6. Adjust intensity levels if available, using the Intensity button.
7. The massager typically has an auto-shutoff timer (e.g., 15-20 minutes). To stop manually, press the Power button again.

General Operation (Neck/Shoulder Massager)

1. Drape the massager over your neck and shoulders, positioning the massage nodes as desired.
2. Insert your arms into the arm loops for better control and to apply gentle pressure.
3. Press the Power button to turn on the device. The massage nodes will begin rotating.
4. Press the Direction button to change the rotation direction of the massage nodes.
5. Press the Heat button to activate or deactivate the heating function.
6. Adjust speed levels if available, using the Speed button.
7. The massager typically has an auto-shutoff timer. To stop manually, press the Power button again.

CARE AND MAINTENANCE

Cleaning

- Always unplug the massager before cleaning.
- Wipe the surface of the massager with a soft, slightly damp cloth.
- Do not use abrasive cleaners, solvents, or strong chemicals.
- Ensure no liquid enters the device.

Storage

- Store the massager in a cool, dry place when not in use.
- Avoid direct sunlight and extreme temperatures.
- Do not wrap the power cord tightly around the device, as this can damage the cord.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Massager does not turn on.	Not plugged in correctly; Power outlet issue; Damaged power cord.	Check power connection; Try a different outlet; Inspect cord for damage (if damaged, discontinue use and contact support).
Heat function is not working.	Heat function not activated; Malfunction.	Ensure Heat button is pressed; Allow a few minutes for heat to build; If still not working, contact support.

Problem	Possible Cause	Solution
Massage feels weak or inconsistent.	Incorrect mode selected; Overuse protection activated; Internal malfunction.	Cycle through massage modes; Turn off and let cool for 30 minutes before restarting; Contact support if problem persists.
Unusual noise during operation.	Internal component issue.	Discontinue use and contact customer support.

PRODUCT SPECIFICATIONS

Please note that specific technical specifications for the Nekteck Shiatsu Foot Massager and Neck/Shoulder Massager are typically provided with the individual product packaging or on the manufacturer's official website. This manual serves as a general guide for operation.

- **Product Type:** Shiatsu Foot Massager, Shiatsu Neck/Shoulder Massager
- **Massage Modes:** Deep Kneading Shiatsu, Heat Function (specific modes may vary by model)
- **Power Input:** AC 110-240V, 50/60Hz (typical for household use)
- **Auto Shut-off:** Approximately 15-20 minutes (for safety and energy saving)
- **Material:** Durable fabric and plastic components

For precise specifications, refer to the product label or the official Nekteck website.

WARRANTY AND CUSTOMER SUPPORT

For information regarding product warranty, returns, or technical support, please refer to the warranty card included with your original product purchase or visit the official Nekteck website. You may also contact Nekteck customer service directly for assistance with any product-related inquiries.

Nekteck Official Website: www.nekteck.com (Note: This is a placeholder URL, please verify the actual manufacturer's website.)