

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Medcursor](#) /
- › [Medcursor Cordless Heated Knee Massager User Manual](#)

Medcursor B0DZX5SR1C

Medcursor Cordless Heated Knee Massager User Manual

Model: B0DZX5SR1C

INTRODUCTION

Thank you for choosing the Medcursor Cordless Heated Knee Massager. This device is designed to provide targeted massage and heat therapy to your knees, shoulders, and elbows, promoting comfort and relaxation. Please read this manual thoroughly before use to ensure safe and effective operation.



Image: The Medcursor Cordless Heated Knee Massager, illustrating its design and areas of heat and vibration application.

SAFETY INFORMATION

To reduce the risk of electric shock, burns, fire, or injury, please observe the following precautions:

- Do not use this device if you have a pacemaker, heart condition, or any other medical implant.
- Consult a physician before use if you are pregnant, have diabetes, phlebitis, thrombosis, or are at increased risk of developing blood clots.
- Do not use on swollen, inflamed, or injured skin areas.
- Avoid using the massager for more than 20 minutes at a time.
- Keep out of reach of children.
- Do not disassemble or attempt to repair the device yourself.
- Use only the provided charging cable.

PACKAGE CONTENTS

- Medcursor Cordless Heated Knee Massager
- USB-C Charging Cable
- Extension Strap (19 inches)
- User Manual

PRODUCT DIAGRAM

Familiarize yourself with the components of your Medcursor Knee Massager:



Image: Close-up of the Medcursor Knee Massager's LED screen, highlighting the control panel for various functions.

1. **Power Button:** Turns the device on/off.
2. **Heating Button:** Cycles through heating levels (Low, Medium, High).
3. **Vibration Button:** Cycles through vibration modes (Auto, Low, High).
4. **LED Display:** Shows current settings and battery level.
5. **Adjustable Straps:** Secure the massager to your body.
6. **USB-C Charging Port:** For recharging the internal battery.

SETUP

Charging the Device

Before first use, fully charge the massager. A full charge takes approximately 2-3 hours.

- Connect the USB-C charging cable to the charging port on the device.
- Connect the other end of the cable to a compatible USB power adapter (not included) or a computer USB port.
- The LED display will indicate charging status. The battery icon will show full when charging is complete.



Image: The Medcursor Knee Massager displaying battery capacity, working time, charging port type, and timer information.

Wearing the Massager

The massager is designed to fit most knee sizes. Use the included extension strap if needed.

- Place the massager on your knee, shoulder, or elbow, ensuring the massage nodes are positioned over the desired area.
- Wrap the adjustable hook-and-loop straps securely around your limb. Ensure a snug but comfortable

fit.

- If the standard straps are too short, attach the 19-inch extension strap to achieve a proper fit.



Image: A diagram illustrating the maximum dimensions of the Medcursor Knee Massager and the length of the adjustable strap, indicating its suitability for various body sizes.

OPERATING INSTRUCTIONS

Once the massager is securely in place and charged, you can begin your massage session.



3 Heating Modes & 3 Vibration Modes



Heating Modes

- Low: 104°F-113°F
- Medium: 115°F-122°F
- High: 124°F-136°F



Vibration Modes

- Auto: 
- Low: 
- High: 

Image: A user comfortably using the Medcursor Knee Massager while relaxing on a sofa, with an overlay detailing the available heating and vibration modes.

Power On/Off

- Press and hold the **Power Button** for 2 seconds to turn the device on. The LED screen will illuminate.
- Press and hold the **Power Button** again for 2 seconds to turn the device off.

Adjusting Heating Levels

The massager offers three heating levels:

- Press the **Heating Button** to cycle through the heat settings:
- **Low:** Approximately 104-113°F (40-45°C)
- **Medium:** Approximately 115-122°F (46-50°C)
- **High:** Approximately 124-136°F (51-58°C)
- The selected heat level will be displayed on the LED screen.

Adjusting Vibration Modes

The massager features three vibration modes:

- Press the **Vibration Button** to cycle through the vibration modes:
- **Auto:** Alternating vibration patterns.
- **Low:** Gentle, continuous vibration.
- **High:** Stronger, continuous vibration.
- The selected vibration mode will be displayed on the LED screen.

Automatic Timer

The device is equipped with an automatic 10-minute timer. After 10 minutes of continuous operation, the massager will automatically shut off to prevent overuse and conserve battery life. You can restart the device for another session if desired.

MAINTENANCE AND CARE

- **Cleaning:** Disconnect the device from power before cleaning. Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** To prolong battery life, charge the device regularly, even if not in frequent use. Avoid fully discharging the battery for extended periods.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Battery is low or depleted.	Charge the device fully using the provided USB-C cable.
No heat or vibration.	Function not activated or set to lowest setting.	Press the Heating/Vibration button to cycle through modes. Ensure the device is powered on.
Device turns off unexpectedly.	Automatic 10-minute timer activated.	This is normal operation. Restart the device if you wish to continue the session.
Straps do not fit.	Standard straps are too short.	Use the included 19-inch extension strap for a more comfortable fit.

SPECIFICATIONS

Model	B0DZX5SR1C
Brand	Medcursor
Power Source	Battery Powered (3000 mAh)
Material	Acrylonitrile Butadiene Styrene

Heating Levels	3 (104-136°F / 40-58°C)
Vibration Modes	3 (Auto, Low, High)
Automatic Timer	10 minutes
Item Weight	1.83 Pounds
Applicable Areas	Knee, Shoulder, Elbow



Image: The Medcursor Cordless Knee Massager demonstrating its versatility for use on the knee, shoulder, and elbow, addressing muscle stiffness, joint swelling, and aiding in injury recovery.

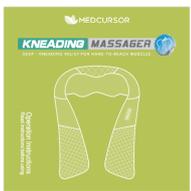
WARRANTY AND SUPPORT

Medcursor products come with a standard manufacturer's warranty. For detailed warranty information, product support, or to register your product, please visit the official Medcursor website or contact customer service.

Medcursor Official Website: [Visit Medcursor Store on Amazon](#)

Please retain your proof of purchase for warranty claims.

Related Documents - B0DZX5SR1C

	<p>Medcursor MD-KP02 Knee Massager with Heat User Manual - Operation, Care, Troubleshooting</p> <p>Comprehensive user manual for the Medcursor MD-KP02 Knee Massager with Heat. Learn about operation, safety precautions, care, maintenance, troubleshooting, and warranty information for this cordless massager.</p>
	<p>Medcursor Kneading Massager User Manual and Safety Guide</p> <p>Comprehensive guide for the Medcursor Kneading Massager, covering safe operation, features, specifications, maintenance, and warranty. Learn how to get deep kneading relief for hard-to-reach muscles.</p>
	<p>Medcursor MD-LM01 Air Compression Leg Massager with Heat - User Manual</p> <p>Comprehensive user manual for the Medcursor MD-LM01 Air Compression Leg Massager with Heat. Learn about product features, safety instructions, operation, troubleshooting, and care.</p>
	<p>Medcursor Kneading Massager: Operation Manual and Safety Guide</p> <p>User manual for the Medcursor Kneading Massager, providing detailed operation instructions, important safety warnings, product specifications, care guidelines, warranty information, and customer service contact details. Features deep tissue relief, shiatsu nodes, and heat functionality.</p>
	<p>Medcursor Shiatsu Foot Massager MDFM400 User Manual Deep Kneading & Heat Therapy</p> <p>Comprehensive user manual for the Medcursor Shiatsu Foot Massager (Model MDFM400). Learn about safety, specifications, operation, troubleshooting, and maintenance for this deep kneading and heat therapy device.</p>
	<p>Medcursor Smart Shiatsu Foot Massager User Manual</p> <p>User manual for the Medcursor Smart Shiatsu Foot Massager (Model: MDFM500), covering important safeguards, setup, operation, care, troubleshooting, FCC statement, and warranty information.</p>

