

[Manuals.plus](#) /

> [INIA](#) /

> INIA Red Light Therapy Mask for Face & Neck Instruction Manual

INIA B0F181GKGZ

INIA Red Light Therapy Mask for Face & Neck Instruction Manual

Model: B0F181GKGZ

1. PRODUCT OVERVIEW

The INIA Red Light Therapy Set provides comprehensive skincare for the face, neck, and chest. This dual-zone system includes a face mask and a neck & chest panel, designed to deliver synchronized light therapy across these areas. It targets common skin concerns such as uneven tone, lack of firmness, and dullness. The device is non-invasive, gentle on all skin types, and suitable for regular home use without discomfort or recovery time.

From Doubts to Visible Results

"I was skeptical, but WOW! My sunspots are fading, and my skin looks years younger!"

— Melissa P., 52



Image: The INIA GLOW 2 + GLOW Neck Dual-Zone Red Light Therapy Set, illustrating the face mask and the neck/chest panel with their respective controllers.

Video: A demonstration of the INIA Red Light Therapy Mask for Face & Neck, highlighting the dual components and their combined use.

2. KEY FEATURES

- **Dual-Zone Coverage:** Combines a face mask and a neck & chest panel for comprehensive treatment of the entire upper body.
- **Multiple Light Modes:** Offers Red-L (620nm), Red-H (630nm), Blue (460–470nm), Yellow (590nm), and Red+Blue (460–630nm) light modes.
- **Near-Infrared (NIR) Customization:** Each mode can be enhanced by activating 850nm near-infrared support for deeper penetration.
- **High-Coverage LED Panel:** Features 228 high-energy LED chips, delivering up to 105mW/cm² output for effective treatment.
- **Ergonomic Design:** Ensures full contact and even light diffusion across treated areas.
- **Tailored Timing:** High-capacity 2600mAh rechargeable controller with intuitive timer functions for 10, 20, or 30-minute sessions.

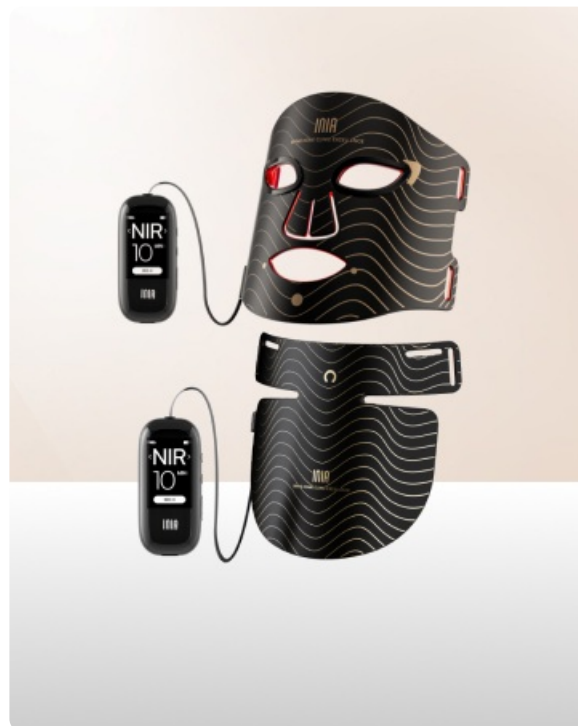


Image: A visual representation of the four light modes (Red + NIR, Blue + NIR, Red + Blue + NIR Combo, Yellow + NIR) and their targeted effects on skin layers.



Image: A cross-section diagram demonstrating the deep penetration of Near-Infrared (NIR) light from INIA GLOW 2

into the epidermis, dermis, and muscles, contrasting it with typical light penetration.

Video: A user demonstrating the INIA Red Light Therapy Mask, highlighting the red, amber, and NIR light modes and their benefits for collagen, firmness, and tone.

3. SETUP INSTRUCTIONS

1. **Unpack Components:** Carefully remove the face mask, neck & chest panel, and controller from the packaging.
2. **Charge Controller:** Ensure the 2600mAh rechargeable controller is fully charged before first use. Connect the charging cable to the controller and a suitable power source.
3. **Connect Masks:** Attach the face mask and neck & chest panel to the controller via their respective cables. Ensure a secure connection.
4. **Position Masks:** Place the face mask comfortably on your face, securing it with the adjustable straps. Position the neck & chest panel over your neck and décolleté area. The ergonomic design ensures a snug fit.

Video: A user demonstrating the wireless nature of the INIA Red Light Therapy Mask, showing how it is worn and adjusted.

4. OPERATING INSTRUCTIONS

1. **Power On:** Press and hold the power button on the controller to turn on the device.
2. **Select Light Mode:** Use the mode selection button to cycle through the available light modes: Red-L (620nm), Red-H (630nm), Blue (460–470nm), Yellow (590nm), and Red+Blue (460–630nm).
3. **Activate NIR (Optional):** Press the dedicated NIR button to activate near-infrared (850nm) support for deeper skin penetration, if desired.
4. **Set Timer:** Press the "TIME" button to select your desired session duration: 10, 20, or 30 minutes. The device will automatically shut off after the selected time.
5. **Begin Treatment:** Relax and allow the light therapy to work. You can perform other activities as the device is hands-free.
6. **Frequency:** For optimal results, use the system 3 to 7 times a week.
7. **Power Off:** After your session, the device will automatically turn off. You can also manually power it off by pressing and holding the power button.

5. MAINTENANCE

- **Cleaning:** Gently wipe the surface of the face mask and neck & chest panel with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.

- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the controller frequently. Charge it regularly, even if not in constant use.

6. TROUBLESHOOTING

- **Device Not Turning On:** Ensure the controller is fully charged. Check that all cables are securely connected to the masks and the controller.
- **Lights Not Activating:** Verify that the device is powered on and the desired mode has been selected. Ensure the masks are properly connected to the controller.
- **Discomfort During Use:** The device is designed to be gentle. If you experience any discomfort, discontinue use and consult the user manual or a healthcare professional.
- **Inconsistent Results:** Consistent use (3-7 times a week) is crucial for visible results. Ensure you are following the recommended session durations.

7. SPECIFICATIONS

Feature	Detail
Package Dimensions	13.43 x 11.46 x 4.17 inches; 4.41 Pounds
Batteries	2 Lithium Ion batteries required (included)
UPC	199284678712
Manufacturer	INIA
ASIN	B0F181GKGZ
Brand	INIA
Item Form	Sheet
Product Benefits	Anti Aging
Special Feature	Lightweight
Unit Count	1.0 Count
GTIN	199284678712

8. WARRANTY AND SUPPORT

INIA offers a 2-year protection plan for this product upon registration. For any support or warranty claims, please refer to the official INIA website or contact customer service through your purchase platform. The product is Verified by Transparency, ensuring authenticity and quality.

© 2026 INIA. All rights reserved.