

arboleaf MW10D

arboleaf MW10D Back Massager User Manual

MODEL: MW10D

Safety Information

Always unplug the massager when not in use. Do not use the massager while charging. Keep the massager away from water. Do not use the massager if you have any medical conditions or are pregnant without consulting a doctor first. Always follow the instructions provided with the massager for proper use and maintenance.

Product Overview

The arboleaf MW10D Back Massager provides deep tissue shiatsu massage with heat therapy. It features a cordless design for portability and Bluetooth app control for personalized settings. The ergonomic design is suitable for various body parts including the neck, shoulders, back, lumbar, and legs.



Image: The arboleaf MW10D Back Massager in ivory color, shown with its packaging and a smartphone displaying the control app.

Video: An official arboleaf video demonstrating the Shiatsu Back Massager in use, highlighting its features and portability.

Package Contents

- arboleaf MW10D Back Massager Unit
- USB-C Charging Cable
- User Manual
- Removable Cover



Image: The arboleaf MW10D Back Massager, USB cable, user manual, and removable cover neatly arranged in the product box.

Setup

1. **Charge the Device:** Connect the provided USB-C charging cable to the massager's charging port and a suitable power source. The charging indicator will show its status. Do not use the massager while charging.
2. **Attach to a Chair (Optional):** Use the adjustable straps on the back of the massager to securely fasten it to an office chair, car seat, or other seating.



Image: The arboleaf MW10D Back Massager connected to a laptop via its USB-C charging cable, demonstrating its cordless charging capability.



Image: Close-up view of the adjustable straps on the arboleaf MW10D Back Massager, showing how they can be used to secure the device to a chair.

Operating Instructions

1. **Power On/Off:** Press and hold the power button on the massager's control panel to turn the device on or off.
2. **Select Massage Mode:** Press the 'M' button to cycle through the 4 preset shiatsu massage modes. These include clockwise, counter-clockwise, and alternating directions.
3. **Adjust Heat Level:** Press the heat button (indicated by wavy lines) to select between two adjustable heat settings (low or high) or to turn the heat off.

4. **App Control:** For advanced customization, download the arboleaf app. Connect via Bluetooth to adjust massage duration (10-30 minutes), speed, and specific massage modes.

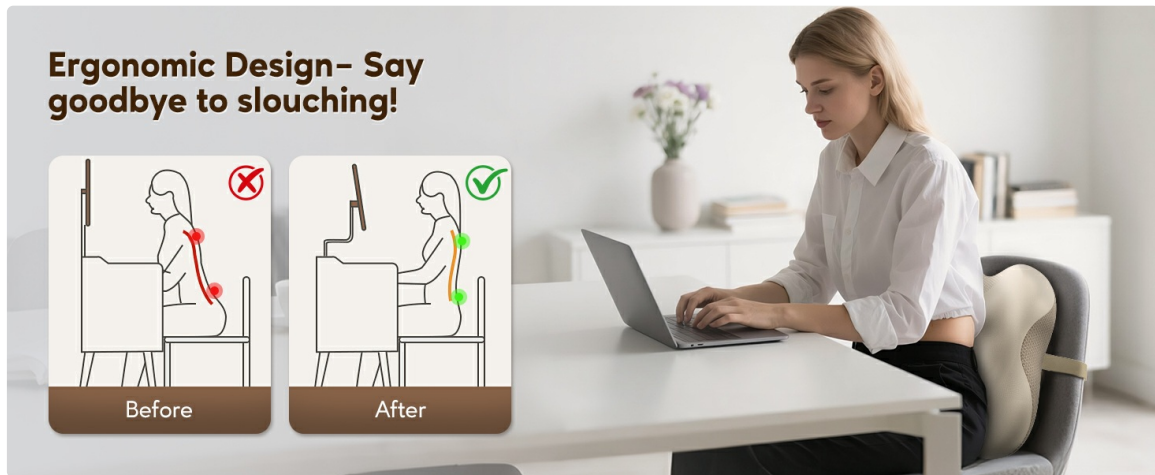


Image: The control panel of the arboleaf MW10D Back Massager with buttons for power, mode, heat, and intensity, alongside a smartphone showing the Bluetooth app interface for advanced control.

Video: An official arboleaf video demonstrating the Bluetooth app control for the Back Massager, showing how to customize settings and modes.

The massager features 4-node kneading shiatsu massage technology, providing bidirectional massage. It automatically pauses after 15 minutes to prevent overuse, a duration that can be adjusted via the app.



Image: Diagram illustrating the 4-node human-like, bi-directional shiatsu massage provided by the arboleaf MW10D Back Massager.

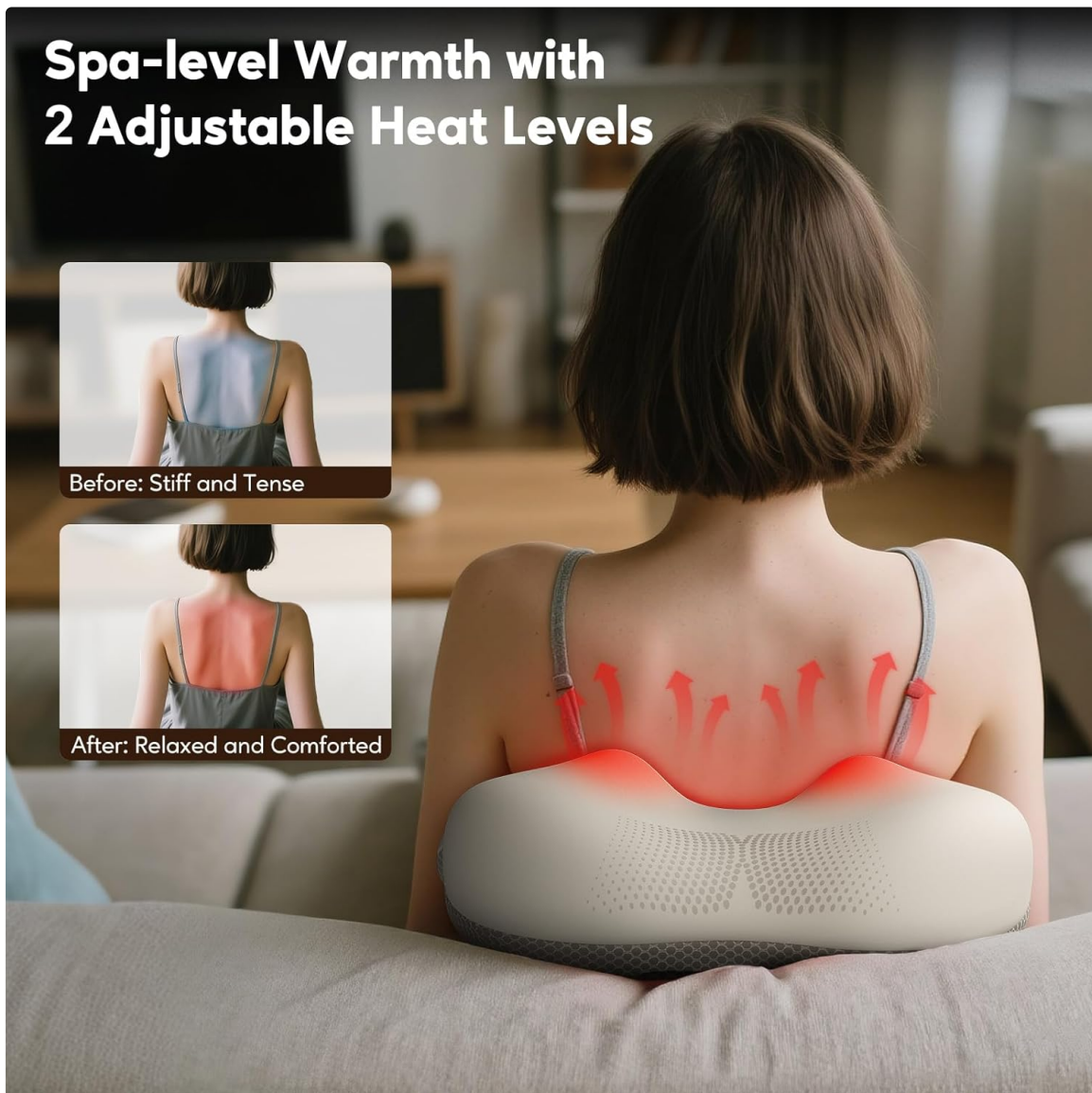


Image: A person using the arboleaf MW10D Back Massager, highlighting the dual-temperature heat therapy for muscle relaxation.

Maintenance

The massager comes with a thoughtfully designed removable cover. To clean, simply unzip and remove the cover. Hand wash or machine wash on a gentle cycle with mild detergent. Air dry completely before reattaching.



Image: The arboleaf MW10D Back Massager with its removable, dust-resistant, and washable cover being detached for cleaning.

Troubleshooting

- **Device not turning on:** Ensure the massager is fully charged. If the issue persists, try using a different USB-C cable or power adapter.
- **Massage intensity too low/high:** Adjust the intensity via the arboleaf app. For optimal comfort, start at the lowest intensity and gradually increase.
- **Heat not working:** Verify that the heat function is activated via the control panel or app. Allow a few minutes for the heat to build up.
- **App connectivity issues:** Ensure Bluetooth is enabled on your smartphone and the massager is within range. Restart both the app and the massager if necessary.

Specifications

Product Dimensions	12 x 12 x 3.5 inches
Item Weight	3.53 Pounds (1.6 Kilograms)
Item Model Number	850059066812
Batteries	2 AAA batteries required (included)
Material	Nylon
Power Source	Battery Powered
Manufacturer	arboleaf

Warranty and Support

The arboleaf MW10D Back Massager comes with a standard manufacturer's warranty. For detailed warranty information, product support, or assistance with any issues, please refer to the full user manual included in the package or visit the official arboleaf website.