

ZZP ZZP-NFD-2025

ZZP Nail Fungus Light Device for Toenails and Fingernails, Model 2025

USER INSTRUCTION MANUAL

[Introduction](#)[Safety Information](#)[Product Overview](#)[Setup](#)[Operating Instructions](#)[Maintenance](#)[Troubleshooting](#)[Specifications](#)[Warranty & Support](#)

1. Introduction

This manual provides essential information for the safe and effective use of your ZZP Nail Fungus Light Device. This device is designed for home use to assist in the repair of damaged and discolored toenails and fingernails through light therapy. Please read this manual thoroughly before operation and retain it for future reference.

2. Safety Information

- Read all instructions carefully before using the device.
- Use only the standard DC5V/2000mA adapter. Using an incorrect adapter may cause device damage.
- Do not expose the device to water or other liquids to prevent product damage.
- Disconnect the adapter if the device is not in use for an extended period.
- Avoid using the wrong adapter to prevent potential hazards.
- Do not operate the device continuously for more than 600 seconds (10 minutes) to prolong its lifespan.
- Do not look directly at the LED lights to avoid eye strain.
- This device is not intended to diagnose, treat, cure, or prevent any disease. Consult a healthcare professional for any medical conditions.

3. Product Overview

The ZZP Nail Fungus Light Device utilizes a combination of 470nm blue light and 905nm invisible infrared light, emitted by 18 LEDs, to interact with nail fungus without affecting surrounding skin. Its design allows for simultaneous treatment of multiple fingers or toes.



Figure 1: The ZPP Nail Fungus Light Device, a white, U-shaped unit with an LED display and control buttons, designed for treating multiple nails simultaneously.

Key Features:

- **Dual Light Therapy:** Combines 470nm blue light and 905nm invisible infrared light for optimal results.
- **Large Treatment Area:** Accommodates multiple fingers or toes at once.
- **Adjustable Design:** The device can be adjusted to fit different hand and foot sizes.
- **Automatic Sensing:** Features an intelligent sensor that activates the light when a hand or foot is inserted.
- **Timer Settings:** Includes 30-second, 60-second, and 99-second timer options.
- **Portable:** Lightweight and compact for easy transport and home use.

A BOON TO FINGERNAIL AND TOENAIL PROBLEMS



No Side Effects



Light



Safe & Painless



Effective



2 Times
a day



10 Mins
once

Figure 2: An illustration demonstrating the interaction of 470nm blue light and 905nm infrared light with nails, targeting fungus at its core.

4. Setup

4.1 Unboxing and Initial Inspection

Carefully remove the device and all accessories from the packaging. Inspect for any signs of damage. The package should contain the Nail Fungus Light Device and a USB charging cable.

Your browser does not support the video tag.

Video 1: An unboxing video demonstrating the contents of the ZZP Nail Fungus Light Device package, including the device and charging cable.

4.2 Charging the Device

Connect the provided USB charging cable to the device's charging port and plug the other end into a compatible DC5V/2A power adapter (not included) or a USB power source. The device must be plugged in to operate.

4.3 Nail Preparation

Before each use, ensure your hands or feet are clean and dry. Trim your nails to a suitable length to allow the light to effectively reach the affected areas.



Figure 3: A visual guide illustrating the steps for nail preparation, including soaking, trimming, and placing nails into the device.

5. Operating Instructions

5.1 Device Activation and Placement

Once plugged in, the device features an automatic sensing function. Simply place your fingers or toes into the designated area of the device. The light will automatically turn on.

Your browser does not support the video tag.

Video 2: A demonstration of the Nail Fungus Light Device in use, showing how to place a foot into the device and the automatic activation of the LED lights.

5.2 Setting the Timer

The device offers preset timer options: 30 seconds, 60 seconds, and 99 seconds. Press the corresponding button to select your

desired treatment duration. If no timer is set, the device will automatically turn off after 99 seconds when a hand or foot is inserted.

5.3 Recommended Usage

For optimal results, use the device twice a day for 10 minutes each time. The device will automatically extinguish once the set time is complete. Consistent use is crucial for effective nail repair.



Figure 4: A visual timeline showing expected improvements in nail appearance over 1 week, 2-3 weeks, 4-6 weeks, and 8 weeks of consistent device use.

5.4 Expected Results

Noticeable improvements in nail appearance typically occur within 4-6 weeks. For effective addressing of weak and damaged nails, continuous use for 3-4 months is recommended to achieve long-lasting results and restore nails to their original beauty.

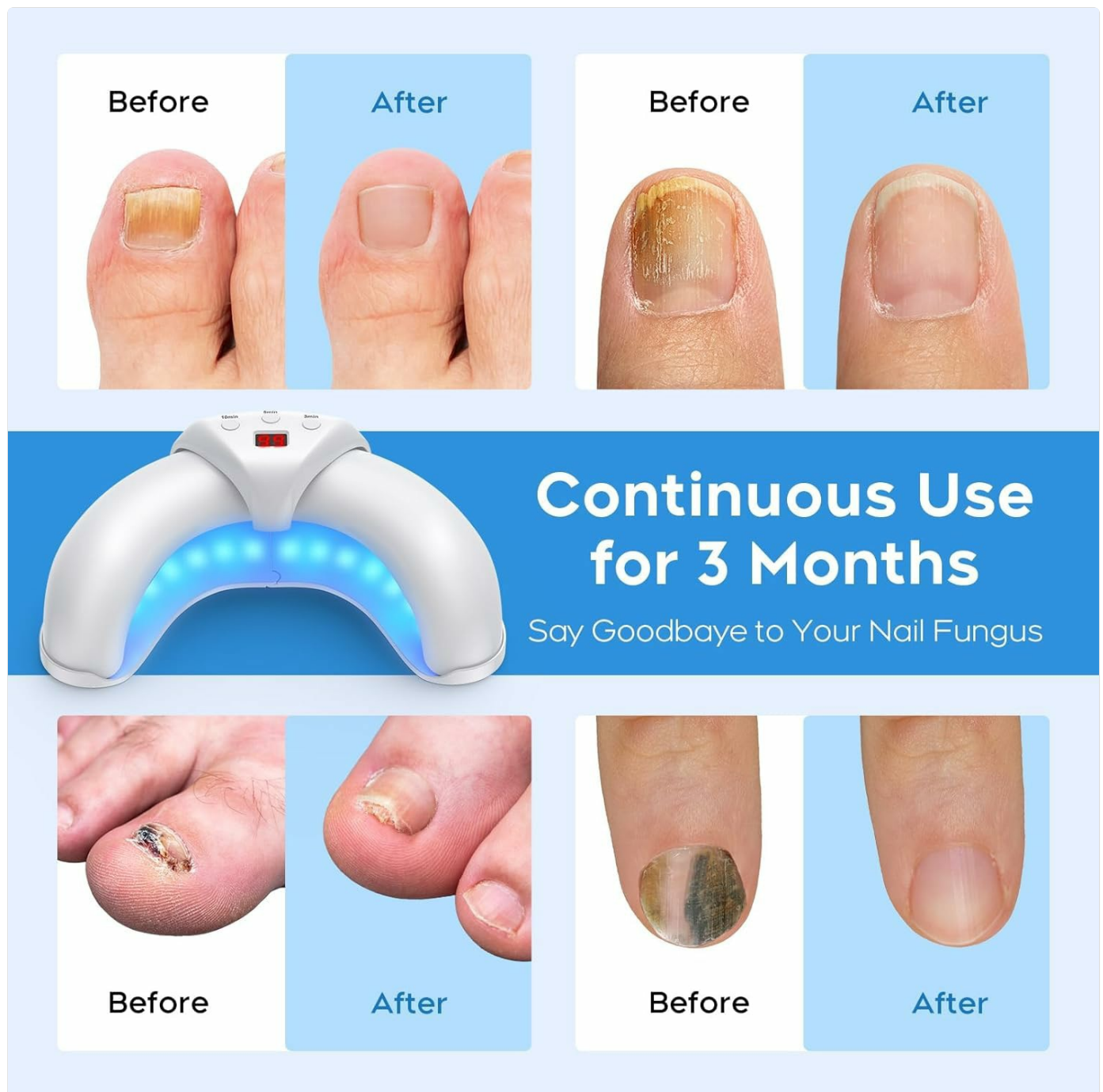


Figure 5: Before and after comparison images illustrating the improvement in nails after continuous use of the device for 3 months.

6. Maintenance

6.1 Cleaning the Device

To maintain hygiene and device longevity, clean the device regularly with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.

6.2 Storage

Store the device in a cool, dry place away from direct sunlight and moisture. Its foldable design allows for compact storage.

7. Troubleshooting

7.1 Device Not Turning On

- Ensure the device is properly plugged into a working power source using the provided USB cable.

- Verify that the power adapter meets the DC5V/2A specification.
- Check if your hand or foot is fully inserted into the device to activate the automatic sensor.

7.2 Lights Not Emitting Correctly

- Confirm the device is receiving adequate power.
- If the issue persists, contact customer support.

8. Specifications

Brand	ZZP
Model	ZZP-NFD-2025
Light Source	18 LEDs (470nm Blue Light, 905nm Invisible Infrared Light)
Power	45W (Max Power)
Operation Voltage	DC5V/2000mA
Timer Settings	30s, 60s, 99s (Automatic Sensing)
Dimensions	6.65 x 6.22 x 3.11 inches
Weight	8.64 ounces (0.54 Pounds)
UPC	759071763672
ASIN	B0F1T2VSQ4

9. Warranty & Support

For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official ZZP website. Keep your purchase receipt as proof of purchase.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question