

COMEBY J03

User Manual: Body Sculpting Machine

Model: J03 | Brand: COMEBY

1. INTRODUCTION

Thank you for choosing the COMEBY Body Sculpting Machine. This professional cordless cellulite massager is designed to assist in reducing the appearance of stretch marks, tightening skin, and promoting a more sculpted body. Featuring an upgraded thermostat function and multiple modes, this device offers a user-friendly experience for home or travel use. Please read this manual thoroughly before operation to ensure proper use and maintenance.

2. PRODUCT FEATURES

- **Multi-functional Design:** Helps reduce stretch marks, tighten skin, and promote a sculpted body.
- **Upgraded Thermostat Function:** Features a temperature range of 37°C - 42°C (98.6°F - 107.6°F) to quickly warm the body, soothe tight muscles, and enhance results.
- **Skin Firming Technology:** Incorporates unique repairing technology to assist in burning cellulite, promoting metabolism, and making the skin firmer.
- **User-Friendly & Portable:** Cordless design with Type-C charging, compact, and lightweight for easy carrying and use at home or while traveling.
- **Adjustable Strength Levels:** Offers 1 to 10 strength levels for customized intensity.
- **Multiple Modes:** Includes D-Line (Reduce Stretch Marks), B-Fat (Get Rid of Fat), and L-Skin (Tighten Skin) modes.



Figure 2.1: Front and back view of the COMEBY Body Sculpting Machine, showcasing its ergonomic design and control panel.



3 LEVELS OF SUCTION

3 Different Functions in 1 Body Sculpting Machine You Deserve

D-Line
Reduce Stretch Marks

B-Fat
Get Rid of Fat

L-Skin
Tighten Skin

Figure 2.2: Illustration detailing the three primary functions: D-Line for stretch mark reduction, B-Fat for fat reduction, and L-Skin for skin tightening.

UPGRADE CONSTANT TEMPERATURE FUNCTION (37°C-42°C)

Which quickly warms up, Soothes tight muscles, and Promotes better results



New upgraded Technology



Which quickly Warms up



Soothe tight Muscles

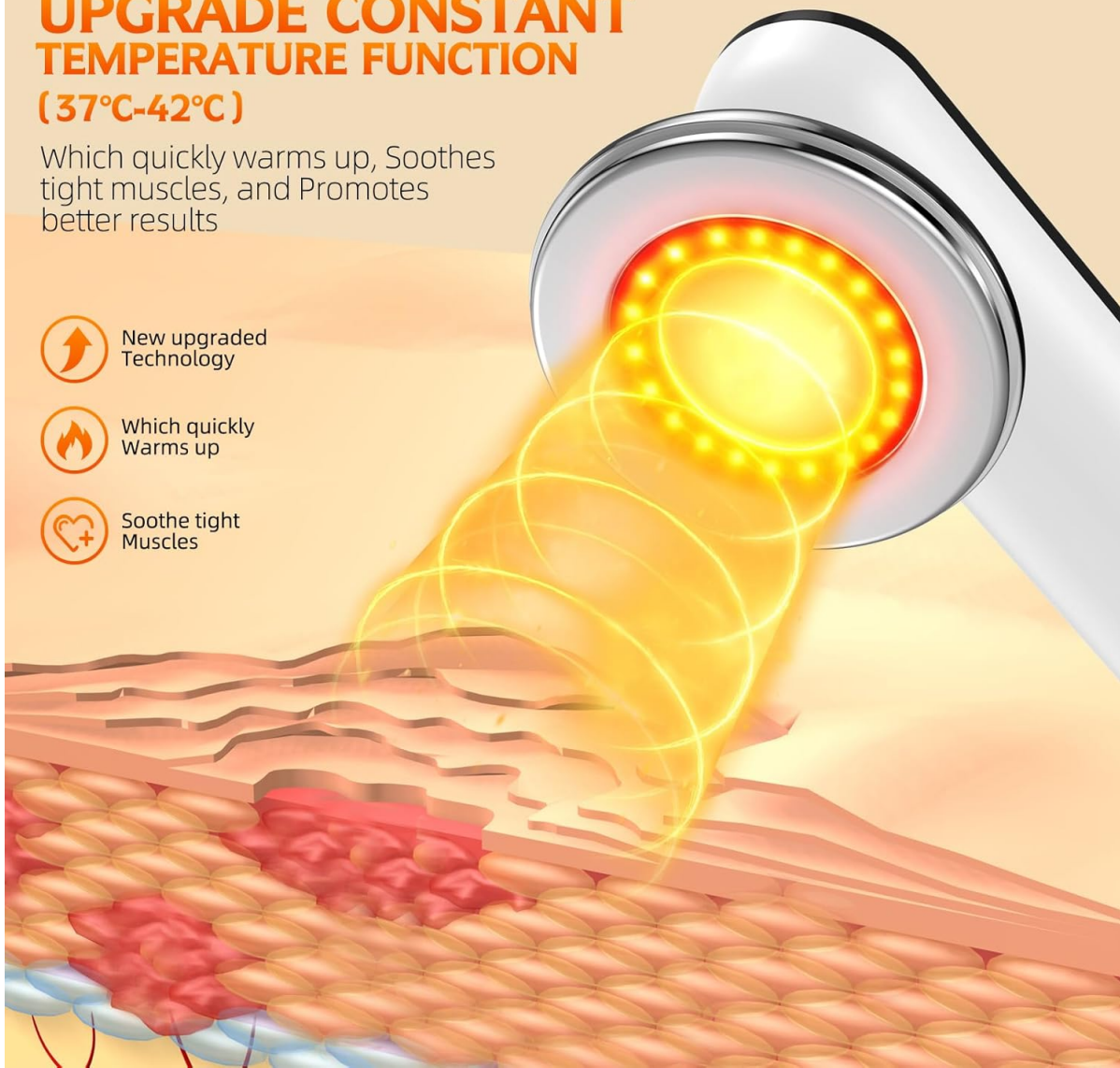


Figure 2.3: Visual representation of the upgraded constant temperature function, highlighting its warming and muscle-soothing benefits.

3. SAFETY INFORMATION AND IMPORTANT SAFEGUARDS

- **Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.
- **Do Not Use on Face:** This device is not intended for facial use.
- **Consult a Physician:** If you have any pre-existing medical conditions, are pregnant, or have sensitive skin, consult your physician before using this device.
- **Keep Out of Reach of Children:** Store the device in a safe place away from children.
- **Proper Charging:** Only use the provided Type-C charging cable. Do not use the device while charging.
- **Avoid Water:** Do not immerse the device in water or use it in wet environments.
- **Discontinue Use If Irritation Occurs:** If you experience any discomfort, redness, or irritation, stop using the device immediately.
- **Handle with Care:** Avoid dropping or subjecting the device to strong impacts.

4. SETUP

1. **Initial Charge:** Before first use, fully charge the device using the provided Type-C charging cable. The display will light up to indicate charging status.
2. **Prepare Skin:** Ensure the skin area where you intend to use the massager is clean and dry.
3. **Apply Gel/Cream:** For optimal results and smooth gliding, apply a suitable body gel or face cream evenly to the target area. (Gels and creams are not included).



Figure 4.1: Overview of the device's controls, including the strength display, heating switch, strength adjustment buttons, and power/mode switch. Also shows compact dimensions.

5. OPERATING INSTRUCTIONS

1. **Power On:** Long press the power button () to turn on the device. The display will light up.
2. **Select Mode:** Short press the power button () to cycle through the available modes: D-Line, B-Fat, and L-Skin. Select the mode appropriate for your needs.
3. **Adjust Strength:** Use the '+' and '-' buttons to adjust the strength level from 1 to 10. Start with a lower setting and gradually increase as comfortable.
4. **Activate Heating (Optional):** Press the flame icon button to activate the thermostat function. The

device will warm up to the optimal temperature range (37°C - 42°C).

5. **Application:** Gently place the massage head on the prepared skin area. Move the device in slow, circular, or upward motions. Ensure continuous contact with the skin.
6. **Recommended Usage:** Use for approximately 15-20 minutes per target area daily. The device has an automatic shut-off feature after 15 minutes of continuous operation. If you need to continue, simply restart the device.
7. **Target Areas:** Suitable for use on the belly, neck, legs, hips, and arms. (Refer to Figure 5.1 for visual guidance).
8. **Power Off:** Long press the power button () to turn off the device after use.



Figure 5.1: Demonstrates suitable application areas for the massager, including the neck, arm, belly, waist, legs, and buttocks. Note: Do not use on face.

SHAPING A GOOD BODY

Reshape Your Body for 15 Minutes Every Day

- Stretch Marks Removal
- Fat Burning
- Body Sculpting
- Body Relaxation
- Tightening Skin
- Make Skin Delicate

Figure 5.2: Highlights the benefits of consistent use, such as stretch marks removal, fat burning, body sculpting, body relaxation, tightening skin, and making skin delicate.

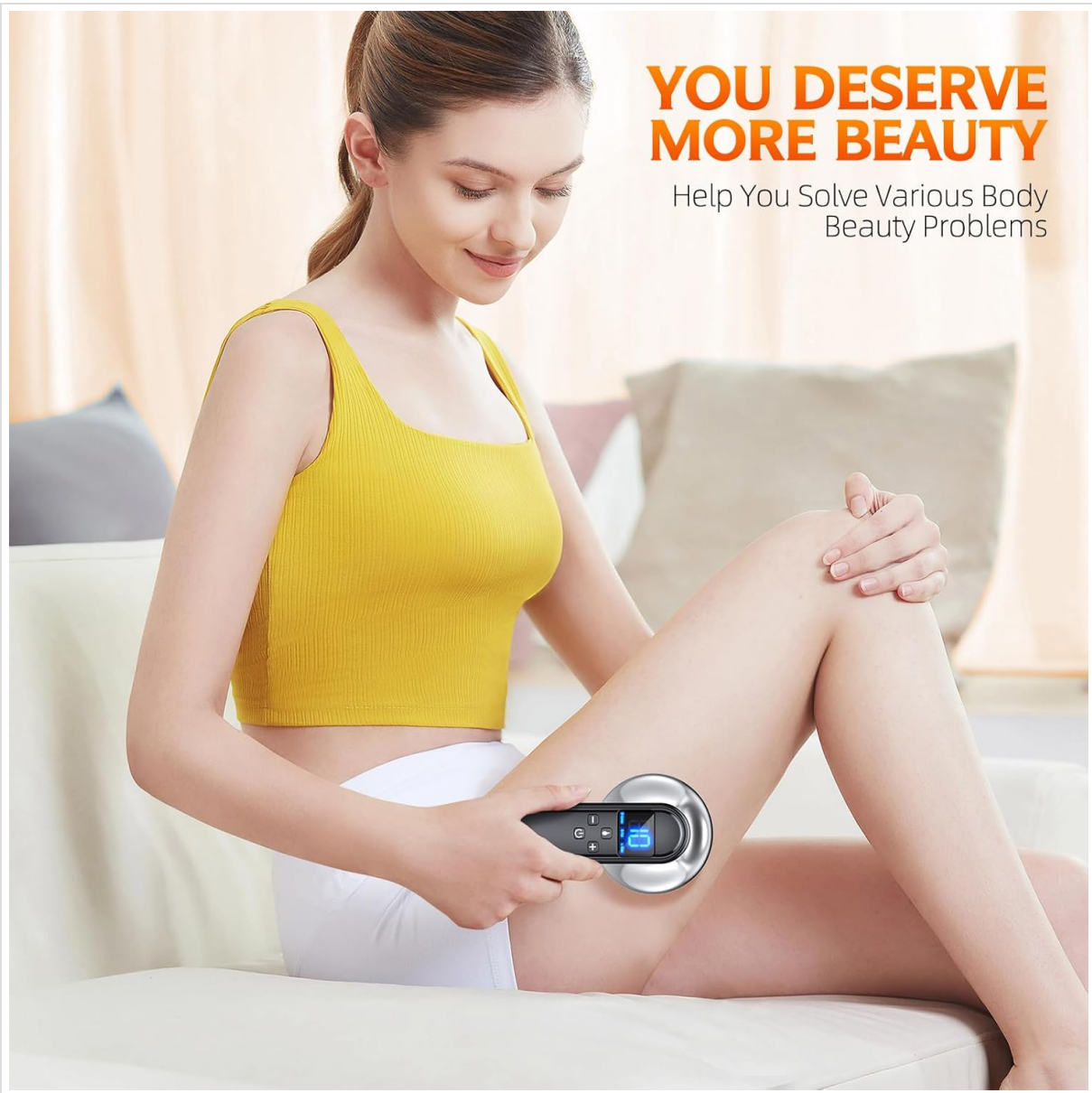


Figure 5.3: A user demonstrating the application of the massager on the leg, emphasizing ease of use.

6. MAINTENANCE AND CARE

- **Cleaning:** After each use, wipe the massage head with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.
- **Charging:** Recharge the device when the battery indicator is low to ensure optimal performance.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Battery is low or depleted.	Charge the device fully using the Type-C cable. Ensure the power button is long-pressed.

Problem	Possible Cause	Solution
No heat or weak heat.	Heating function not activated or device not fully warmed up.	Press the flame icon button to activate heating. Allow a few moments for the device to reach temperature.
Device stops during use.	Automatic 15-minute shut-off activated.	Restart the device by long-pressing the power button.
Skin irritation after use.	Strength level too high or sensitive skin.	Reduce the strength level. If irritation persists, discontinue use and consult a physician.

If you encounter any issues not listed here, please contact COMEBY customer support for assistance.

8. SPECIFICATIONS

- **Model Number:** J03
- **Brand:** COMEBY
- **Dimensions:** 7.8 x 3.94 x 2.52 inches (approximate package dimensions)
- **Weight:** 11.68 ounces (approximate package weight)
- **Charging Interface:** Type-C
- **Temperature Range (Thermostat):** 37°C - 42°C (98.6°F - 107.6°F)
- **Strength Levels:** 1-10
- **Auto Shut-off:** 15 minutes

9. WARRANTY AND SUPPORT

For warranty information, product support, or any inquiries, please refer to the contact details provided with your purchase or visit the official COMEBY website. Keep your purchase receipt as proof of purchase for warranty claims.

This manual is subject to change without prior notice.