

XAUCY XAUCY

XAUCY EMS Microcurrent Face Neck Beauty Device

Model: XAUCY

INTRODUCTION

The XAUCY EMS Microcurrent Face Neck Beauty Device is designed to provide comprehensive skin care for the face and neck. Utilizing EMS microcurrent, LED photon therapy, constant temperature heating, and vibration massage, this device aims to promote skin firming, rejuvenation, and improved elasticity. It features a 160° curve massage head for ergonomic application and offers multiple modes to address various skin concerns.

PRODUCT FEATURES

- **EMS Microcurrent Technology:** Stimulates facial and neck muscles for a lifting and firming effect.
- **LED Photon Therapy:** Features 3 color lights with different functions to enhance skin tone and texture. Blue LED specifically aids in collagen production and improving skin elasticity.
- **45°C Constant Heating:** Provides a gentle warming sensation to stimulate blood circulation and promote nutrient absorption.
- **Vibration Massage:** Enhances product penetration and relaxes facial muscles.
- **160° Curve Massage Head:** Ergonomically designed to fit the contours of the face and neck.
- **4-Speed Adjustable Modes:** Offers customizable intensity levels for personalized treatment.
- **Rechargeable Battery:** Equipped with a 400 mAh lithium polymer battery, fully charged in approximately 3 hours.

YOU DESERVE TO HAVE IT!



Image: Overview of the XAUCY EMS Microcurrent Face Neck Beauty Device, illustrating its key features such as 3 color light functions, 160° curve massage head, 45°C heating, 4-speed adjustment, and benefits like deep cleaning and skin rejuvenation.

PACKAGE CONTENTS

- XAUCY EMS Microcurrent Face Neck Beauty Device
- USB Charging Cable
- User Manual



Image: The XAUCY EMS Microcurrent Face Neck Beauty Device shown with its retail packaging and included instruction manual.

SETUP

Charging the Device

1. Connect the USB charging cable to the device's charging port.
2. Plug the other end of the USB cable into a 5V/1A USB power adapter (not included).
3. The device will indicate charging status. A full charge takes approximately 3 hours.
4. Ensure the device is fully charged before first use for optimal performance.

OPERATING INSTRUCTIONS

Preparation for Use

1. Cleanse your face and neck thoroughly to remove any makeup or impurities.
2. Apply your preferred skincare products, such as serum, cream, or gel, evenly to the treatment area. This helps the device glide smoothly and enhances product absorption.

Power On/Off and Mode Selection

1. **To Power On:** Long press the switch button on the device.
2. **To Select Mode:** Each subsequent short press of the switch button will cycle through the available modes (three gears, seven modes). Choose the mode that best suits your individual needs.
3. **To Power Off:** Press and hold the switch button again to turn off the device.

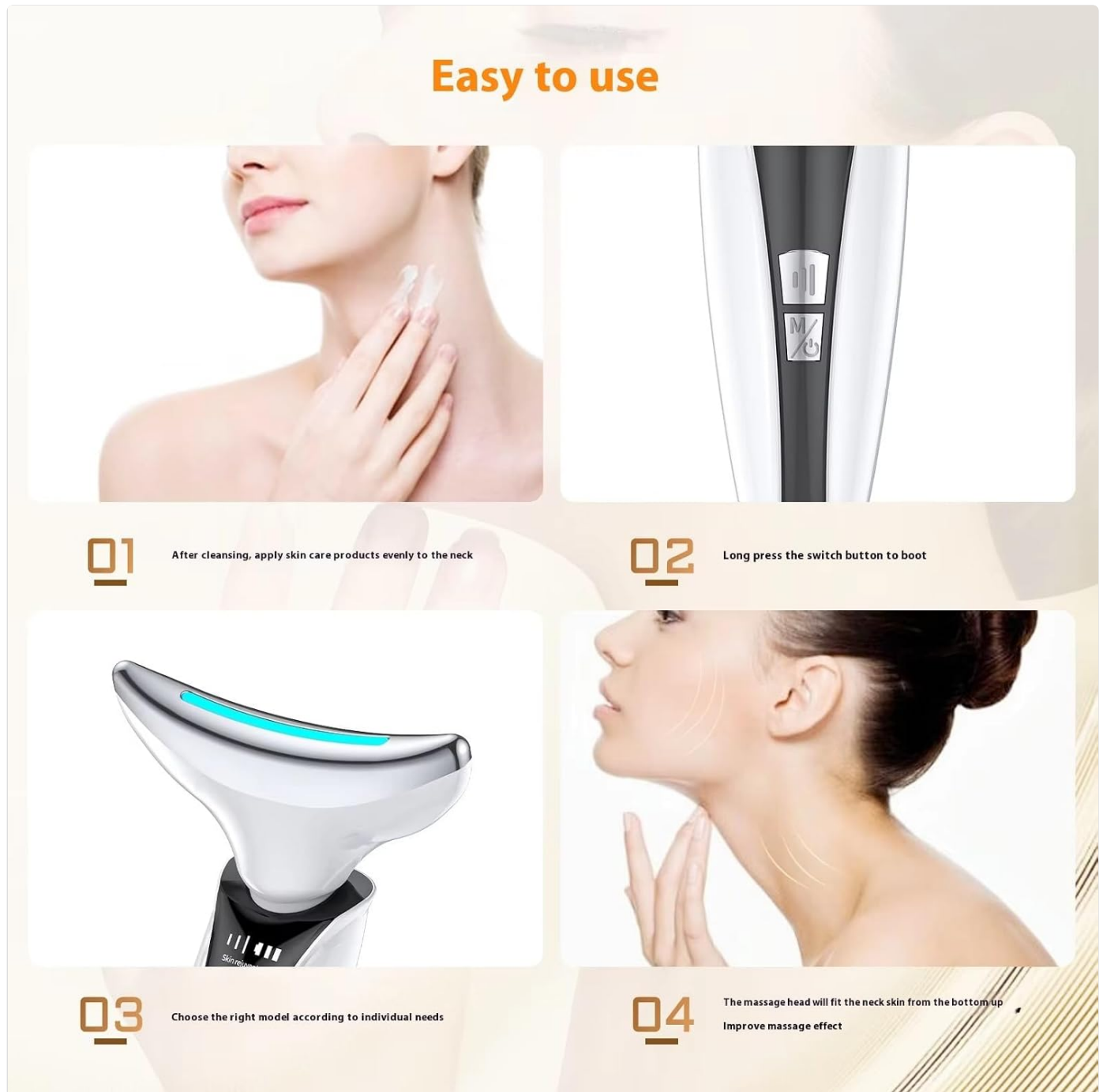


Image: A visual guide showing the four easy steps to use the device: 1. Apply skincare products to the neck after cleansing. 2. Long press the switch button to power on. 3. Choose the appropriate mode. 4. Apply the massage head to the neck for improved massage effect.

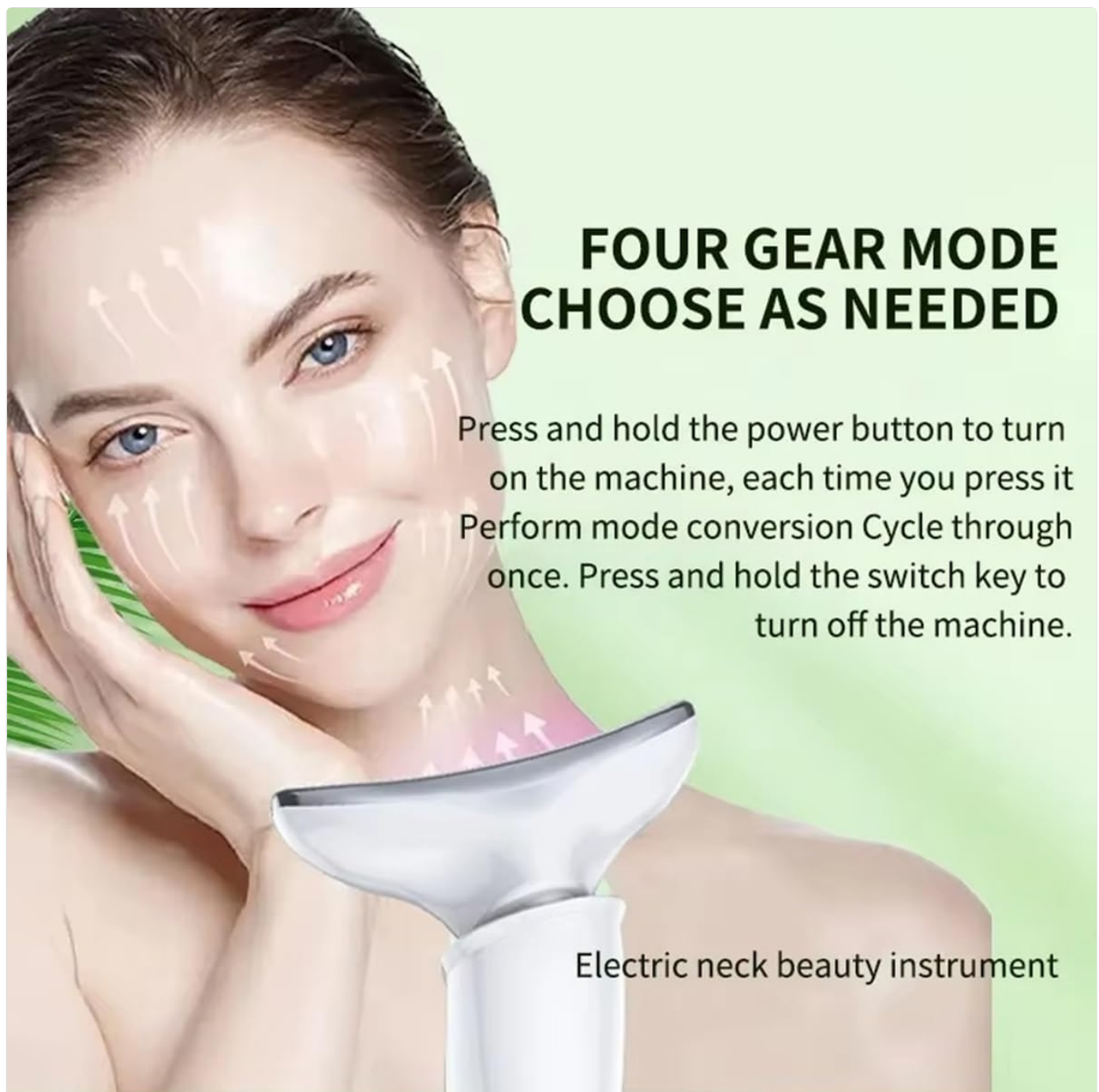


Image: A woman demonstrating the use of the XAUCY device on her neck, with text indicating to press and hold the power button to turn on and cycle through four gear modes as needed.

Application Areas and Techniques

The device can be used on various areas of the face and neck. Follow the recommended movements for optimal results:

- **Neck:** Glide the device upwards from the base of the neck towards the chin.
- **Face (Chin to Cheek):** Lift from the chin towards the cheekbones.
- **Face (Cheek to Outer Canthus):** Lift from the cheek towards the outer corner of the eye.
- **Eye Area:** Gently massage around the eye area, avoiding direct contact with the eyeball.
- **Forehead:** Move the device upwards from the eyebrows towards the hairline.

RECOMMENDED USING METHODS & TARGETS AREAS OF CONCERN

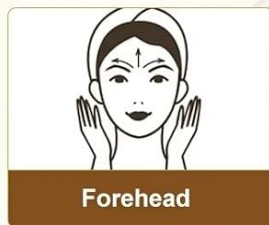


Image: Diagrams illustrating recommended massage techniques and target areas of concern, including lifting from chin to cheek, neck, face, eye, and forehead.

Recommended Usage Frequency

For best results, it is recommended to use the device 2-3 times a week. Each session should last approximately 5-10 minutes, ideally once in the morning and once in the evening.

PERSISTENT USE FOR VISIBLE RESULTS

It is recommended to use it once in the morning and once in the evening for 5-10 minutes each time

BEFORE



AFTER



Image: A collage of 'Before' and 'After' images demonstrating visible results from persistent use of the device, with text recommending use once in the morning and once in the evening for 5-10 minutes each time.

MODES AND FUNCTIONS EXPLAINED

- **LED Photon Therapy:** Harnesses LED technology to improve skin tone and texture. The blue LED light stimulates collagen production and enhances skin elasticity.
- **Constant Temperature Heat (45°C):** Provides a gentle warming sensation that stimulates blood circulation, promoting microcirculation, and delivering oxygen and nutrients to the skin for enhanced health and radiance. This also helps in fading neck lines by tightening and smoothing the skin.
- **EMS Microcurrent:** Targets fine lines and wrinkles, particularly around the neck area, to diminish their appearance for a smoother, more contoured profile. Regular use helps reduce the visibility of neck lines.
- **Vibration Massage:** Aids in deep cleaning, nutrient import, and toxin flush out, contributing to overall skin rejuvenation.

45°C CONSTANT HEATING

Constant heat combined with coloured light technology accelerates blood circulation, rejuvenates the neck and face tightens and smoothes it.

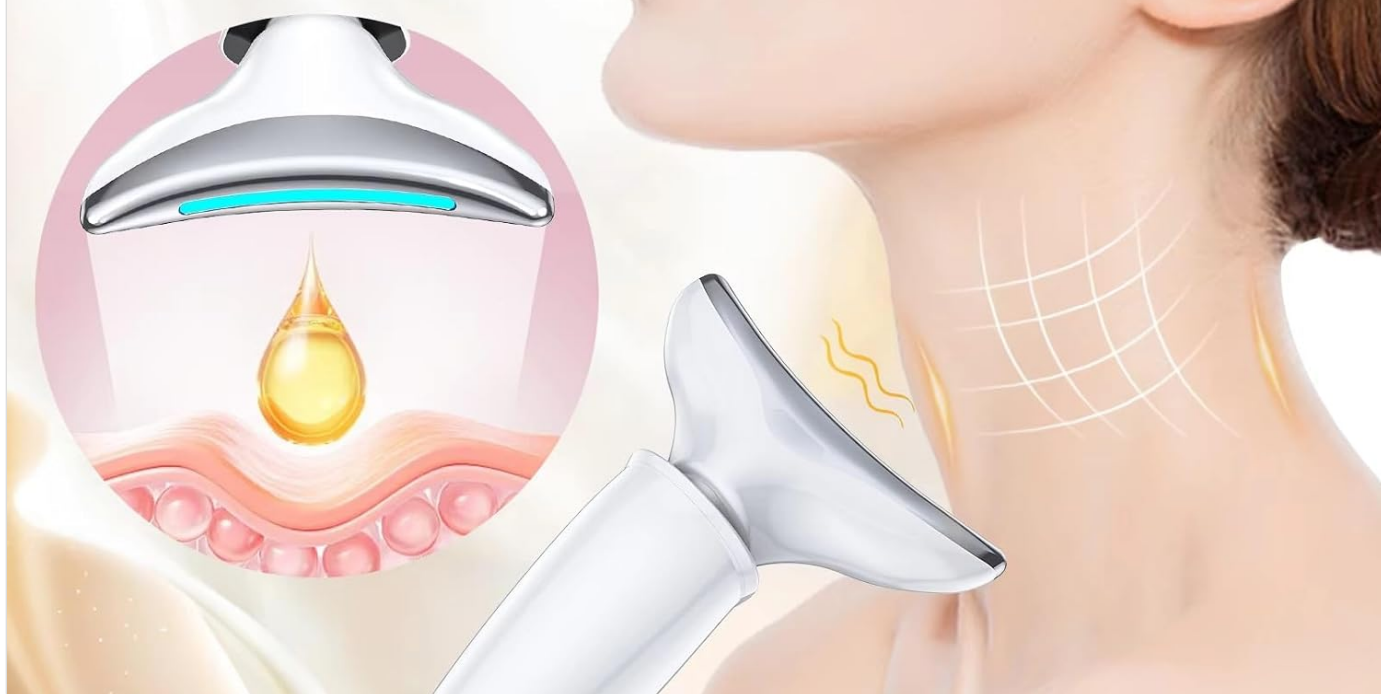


Image: A diagram explaining the 45°C constant heating function, showing how constant heat combined with colored light technology accelerates blood circulation, rejuvenating and tightening the neck and face.

MAINTENANCE

Cleaning the Device

1. Ensure the device is powered off and disconnected from the charging cable before cleaning.
2. Wipe the massage head and body of the device with a soft, damp cloth.
3. Do not use abrasive cleaners, alcohol, or harsh chemicals, as these may damage the device.
4. Allow the device to air dry completely before storing.

Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep out of reach of children.

SPECIFICATIONS

Product Name	7 Colour Light Neck Beauty Instrument
Model Number	XAUCY
Material	ABS + Electroplating Technology
Battery Capacity	400 mAh
Battery Type	Lithium Polymer Battery
Charging Voltage	5V/1A
Rated Voltage	3.7V
Rated Power	3W
Functions	Three Gears, Seven Modes, Constant Temperature Heat, Vibration Massage
Usage Duration (per charge)	1 to 2 hours
Product Size	139 x 90 x 38 mm
Power Source	USB Charger

© 2025 XAUCY. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question