

## COMFIER CF-2025 Massage System

# COMFIER CF-2025 Massage System User Manual

Massage Seat Cushion & Shiatsu Foot Massager

## 1. INTRODUCTION

Thank you for choosing the COMFIER CF-2025 Massage System. This product combines a vibrating massage seat cushion with a shiatsu foot massager to provide a comprehensive relaxation experience. This manual provides important information regarding the safe operation, setup, maintenance, and troubleshooting of your new massage system. Please read this manual thoroughly before use and retain it for future reference.



**Figure 1.1:** The COMFIER CF-2025 Massage System, featuring both the vibrating massage seat cushion and the shiatsu foot massager.

## 2. IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, basic precautions should always be followed to reduce the risk of electric shock, fire, and injury to persons, including the following:

- Read all instructions before using this appliance.
- Do not use while bathing or in a shower.
- Do not use near water.
- Never insert pins or other metallic fasteners into the appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- This appliance is designed for household use only.
- Do not attempt to repair the appliance. All servicing must be performed by authorized service personnel.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. Consult your physician before using if you have any medical conditions.

### 3. PACKAGE CONTENTS

---

Please check the package contents upon unboxing to ensure all items are present and undamaged:

- 1 x COMFIER Vibrating Massage Seat Cushion
- 1 x COMFIER Shiatsu Foot Massager
- 1 x AC Power Adapter (for Seat Cushion)
- 1 x AC Power Adapter (for Foot Massager)
- 1 x User Manual

### 4. PRODUCT FEATURES

---

#### 4.1. Vibrating Massage Seat Cushion

The COMFIER Massage Seat Cushion is designed to provide soothing vibration massage and gentle heat therapy for your back, waist, and buttocks. It features:

- **10 Vibration Motors:** Strategically placed to target neck/shoulder, mid-back, lower back, and hips/thighs.
- **Heat Function:** Gentle warmth for the upper back, waist, and buttocks, enhancing comfort and relaxation.
- **Adjustable Intensity:** Multiple intensity levels for customized massage.
- **Comfortable Design:** Made with short plush fabric for a softer and more comfortable experience.
- **Versatile Fit:** Longer length suitable for various chair types and body shapes.

**Please note: This seat cushion is a vibrating massager only and does not feature Shiatsu or roller balls.**



## Others



Heating for upper back, waist, buttocks



Larger size vibration motor with wide vibration coverage



Made of short plush fabric, softer and more comfortable



Longer length, suitable for more types of chairs & different body shapes



Only heat one of the back or buttocks



Small size vibration motor



Ordinary cotton fabric



Short size does not fit most chairs

**Figure 4.1.1:** Key advantages of the COMFIER seat cushion, including heating for multiple zones, larger vibration motors, plush fabric, and extended length.

**Please note:** this is a **vibrating massager only** and does not feature Shiatsu or roller balls.



***Vibration Massage with 10 Motors***  
*Help relieve stress, aches, tension and eliminate fatigue*

**Figure 4.1.2:** Illustration of the 10 vibration motors providing targeted massage for neck/shoulder, mid-back, lower back, and hips/thighs.





Gentle Warmth,  
Enhances Comfort, Ensures Safety

Feel Warmth in 3 Mins

Reach 140°F for about 18 Mins

**Figure 4.1.3:** The seat cushion's gentle warmth feature, which reaches 140°F in approximately 18 minutes, enhancing comfort and safety.

## 4.2. Shiatsu Foot Massager

The COMFIER Shiatsu Foot Massager delivers deep kneading, rolling, and heat therapy to your feet, designed to accommodate a wide range of foot sizes. Key features include:

- **Shiatsu & Deep Kneading:** Provides a professional 360° full wrap massage experience.
- **3 Intensity Levels:** Choose from Gentle, Relaxing, or Deep massage.
- **3 Heating Levels:** Adjustable heat from L-104°F to H-130°F for soothing warmth.
- **3 Air Intensity Levels:** For customizable compression massage.
- **Timer Function:** Set massage duration for 15, 20, or 30 minutes.
- **Extra-Large Design:** Comfortably fits feet up to men's size 14.

# Suitable for Whole Family

Fits up to **size 14**



**Figure 4.2.1:** The foot massager's extra-large design, accommodating feet up to size 14, making it suitable for the entire family.



# Larger Shiatsu Foot Massager



 **COMFIER**

**VS**

**Others**

Up to men size 14	↑	×
3 levels heat	↑	×
15/20/30 Mins Timer	↑	×
131F Heating Temperature	↑	×

**Figure 4.2.2:** A comparison showing the COMFIER foot massager's advantages, including larger size, 3 heat levels, multiple timer options, and higher heating temperature.



## 3 Adjustable Levels Soothing Heat



Feel Warmth in 3 Mins  
Reach **130°F** for about 6 Mins



Figure 4.2.3: The foot massager offers three adjustable heat levels (L-104°F, M-122°F, H-130°F), providing soothing warmth that can be felt within 3 minutes.

# Professional 360°full warp massage design



**3 Knead Intensity**  
L / M / H



**3 Heating level**  
L / M / H



**3 Air Intensity**  
L / M / H



**3 Timer**  
15 / 20 / 30Mins



**Figure 4.2.4:** The foot massager's comprehensive 360° full wrap massage design, offering 3 knead intensities, 3 heating levels, 3 air intensities, and timer options.



# Shiatsu & Deep Kneading & Rolling Massager

## 3 Intensity levels

L


Gentle massage

M

Relaxing massage

H

Deep massage

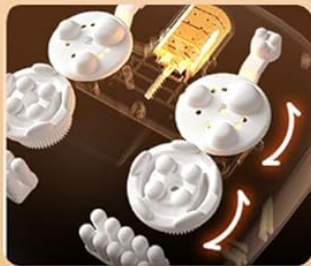

**Figure 4.2.5:** The foot massager provides three intensity levels: Gentle, Relaxing, and Deep massage, for a customized Shiatsu and deep kneading experience.

## 5. SETUP

### 5.1. Massage Seat Cushion Setup



1. Unpack the seat cushion and its power adapter from the packaging.
2. Place the seat cushion on a chair, sofa, or bed. Ensure it is securely positioned and will not slip during use.
3. Connect the AC power adapter to the input jack on the side of the seat cushion.
4. Plug the power adapter into a standard electrical outlet (110-240V).
5. Ensure the remote control is easily accessible.

## 5.2. Shiatsu Foot Massager Setup

1. Unpack the foot massager and its power adapter.
2. Place the foot massager on a flat, stable surface, such as the floor, in front of your chair or sofa.
3. Connect the AC power adapter to the input jack on the back or side of the foot massager.
4. Plug the power adapter into a standard electrical outlet (110-240V).

## 6. OPERATING INSTRUCTIONS

### 6.1. Massage Seat Cushion Operation

The seat cushion is operated via an attached wired remote control.

1. **Power On/Off:** Press the **POWER** button to turn the unit on or off.
2. **Vibration Modes:** Use the **MODE** button to cycle through different vibration patterns (e.g., continuous, pulsing, wave).
3. **Intensity Adjustment:** Press the **INTENSITY** button (or +/- buttons) to increase or decrease the vibration intensity.
4. **Heat Function:** Press the **HEAT** button to activate or deactivate the warming function. Allow a few minutes for the heat to reach its optimal temperature.
5. **Targeted Zones:** Some models may have buttons to select specific massage zones (e.g., upper back, lower back, thighs). Press the corresponding button to activate or deactivate a zone.
6. **Timer:** The unit may have an auto-shutoff timer (e.g., 15, 20, 30 minutes). Refer to the remote for specific timer settings if available.

### 6.2. Shiatsu Foot Massager Operation

The foot massager features an intuitive control panel on its top surface.

1. **Power On/Off:** Press the **POWER** button to turn the unit on or off.
2. **Knead Intensity:** Press the **KNEAD** button to cycle through Low (L), Medium (M), and High (H) kneading intensity levels.
3. **Heat Function:** Press the **HEAT** button to select desired heat level (L, M, H) or turn off heat.
4. **Air Intensity:** Press the **AIR** button to adjust the air compression intensity (L, M, H).
5. **Timer Setting:** Press the **TIMER** button to set the massage duration to 15, 20, or 30 minutes. The unit will automatically shut off after the selected time.

## 7. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and optimal performance of your massage system.

- **Cleaning:** Always unplug the unit from the power outlet before cleaning. Wipe the surface of both the seat cushion

and foot massager with a soft, slightly damp cloth. Do not use abrasive cleaners, brushes, gasoline, kerosene, glass cleaner, furniture polish, or paint thinner to clean the massager.

- **Storage:** When not in use, store the massager in a cool, dry place. Avoid direct sunlight and extreme temperatures. Do not place heavy objects on the massager.
- **Avoid Sharp Objects:** Keep sharp objects away from the fabric and internal components to prevent damage.
- **Power Cord:** Regularly inspect the power cord for any signs of damage. If damaged, do not use the unit and contact customer support.

## 8. TROUBLESHOOTING

If you encounter issues with your COMFIER Massage System, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in; Power outlet malfunction; Damaged power cord.	Ensure power adapter is securely plugged into both the unit and a working electrical outlet. Try a different outlet. Check power cord for damage. If damaged, discontinue use and contact support.
No heat or insufficient heat.	Heat function not activated; Unit not warmed up yet.	Ensure the HEAT button is pressed. Allow 5-10 minutes for the heat to fully develop.
Vibration/Kneading is weak or inconsistent.	Intensity setting is too low; Unit overloaded.	Increase the intensity level. Ensure the unit is not being pressed down too hard or used on an uneven surface.
Unit suddenly stops during use.	Auto-shutoff timer activated; Overheating protection.	This is normal if the timer has expired. If not, the unit may have overheated. Turn off and unplug the unit, allow it to cool down for 30-60 minutes before restarting.

If the problem persists after attempting these solutions, please contact COMFIER customer support.

## 9. SPECIFICATIONS

Feature	Detail
Product Name	COMFIER CF-2025 Massage System
Model Number	CF-2025
ASIN	B0F2MY18C5
Power Input	AC 100-240V ~ 50/60Hz
Rated Power	Varies by unit (refer to power adapter label)
Seat Cushion Dimensions	Approx. 43.3 x 18.9 inches (L x W)

Feature	Detail
Foot Massager Dimensions	Approx. 14.65 x 12.40 x 6.69 inches (L x W x H)
Foot Massager Foot Size Accommodated	Up to Men's Size 14
Auto Shut-off Timer	15/20/30 minutes (Foot Massager), Varies (Seat Cushion)
Heat Temperature (Foot Massager)	Up to 130°F (54°C)
Heat Temperature (Seat Cushion)	Up to 140°F (60°C)

## 10. WARRANTY AND SUPPORT

COMFIER products are designed for durability and performance. This product comes with a limited warranty against defects in materials and workmanship. Please refer to the warranty card included in your package for specific terms and conditions, or visit the official COMFIER website for detailed warranty information.

For technical support, troubleshooting assistance, or warranty claims, please contact COMFIER Customer Service:

- **Website:** [Visit the COMFIER Store on Amazon](#)
- **Email:** Refer to your product packaging or official website for contact email.
- **Phone:** Refer to your product packaging or official website for contact number.

Please have your product model number (CF-2025) and ASIN (B0F2MY18C5) ready when contacting support.