

## COMFIER CF-6505

# COMFIER Cordless Neck and Shoulder Massager CF-6505 User Manual

## 1. INTRODUCTION

The COMFIER Cordless Neck and Shoulder Massager (Model CF-6505) is designed to provide deep tissue shiatsu massage with an optional heating function for various body areas. This portable and rechargeable device offers relief from muscle tension and stiffness, suitable for use at home, in the office, or while traveling. It features 3D kneading nodes, adjustable intensity, and a gentle heating function to enhance relaxation.



*Image 1.1: The COMFIER Cordless Neck and Shoulder Massager in use, highlighting its ergonomic design and heating function.*

## 2. SAFETY INFORMATION

**Please read all instructions carefully before use and retain for future reference.**

- Always unplug the massager when not in use.
- Do not use while sleeping or taking a bath. Keep away from water.

- Do not use if you have a pacemaker or other medical device.
- Consult your doctor before use if you are pregnant, have a heart condition, or any other medical condition.
- Do not use on broken skin or areas with inflammation.
- Always follow the instructions provided with the massager.

### 3. PRODUCT OVERVIEW

The massager features an ergonomic U-shape design, constructed with premium PU leather and breathable mesh for comfort and durability. It includes an integrated control panel for easy operation.

#### 3.1 Components

- **Main Massager Unit:** U-shaped body with internal kneading nodes.
- **Control Panel:** Located on one of the arm straps, featuring buttons for operation.
- **Arm Straps:** Used to hold and adjust the massager's position and pressure.
- **Charging Port:** USB-C port for recharging the internal battery.



Image 3.1: Front view of the COMFIER Cordless Neck and Shoulder Massager.

### 3.2 Control Panel

The control panel provides intuitive access to all functions:




Image 3.2: Detailed view of the control panel buttons.

- **Power Button:** Turns the massager on or off.
- **Direction Button:** Changes the rotation direction of the massage nodes.
- **Intensity Button:** Adjusts the massage intensity (3 levels).
- **Heating Button:** Activates or deactivates the heat function.

## 4. CHARGING THE DEVICE

The massager is equipped with a 2200mAh rechargeable battery, allowing for cordless operation.

1. Locate the USB-C charging port on the massager.
2. Connect the provided USB-C cable to the massager and a compatible USB power adapter (not included).
3. The power indicator light will pulse white during charging.
4. A full charge takes approximately 4 hours and provides up to 1 hour of cordless operation.
5. The power light flashes red when the battery level is low, indicating it's time to recharge.

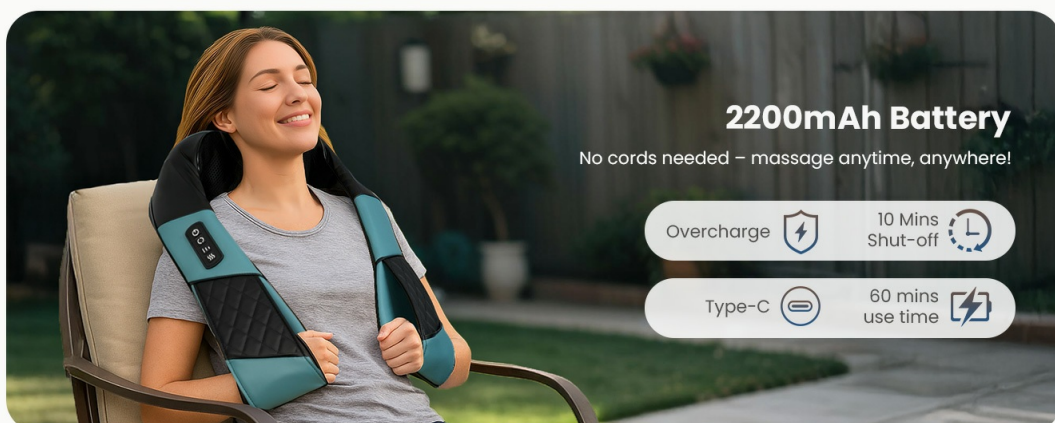

## 5. OPERATING INSTRUCTIONS

Follow these steps to operate your COMFIER Cordless Neck and Shoulder Massager:

1. **Power On/Off:** Press the **Power button** once to turn the massager on. Press it again to turn the massager off. The device has a 10-minute auto shut-off feature for safety and to prevent overuse.
2. **Adjust Massage Direction:** Press the **Direction button** to change the rotation direction of the kneading nodes. This allows for varied massage sensations.
3. **Select Intensity Level:** Press the **Intensity button** to cycle through the 3 adjustable intensity levels (low, medium, high). Choose the level that best suits your comfort and needs.
4. **Activate Heat Function:** Press the **Heating button** to activate or deactivate the gentle warming function. The heat can help soothe muscles and enhance relaxation.

The quiet motor ensures an undisturbed massage experience, allowing you to relax while reading, watching TV, or working.

# Why Choose Comfier?

## Others



- ✗ Tangled Cords, Limited Use
- ✗ Uncomfortable, Harsh Massage

- ✓ Wireless Freedom, Use Anywhere
- ✓ Deep Kneading, Soothing Comfort

Image 5.1: The control panel for adjusting settings during use.

## 6. APPLICATION AREAS

The ergonomic U-shape design allows for versatile use on various body parts, providing targeted relief.

- **Neck and Shoulders:** Position the massager around your neck and shoulders. Use the arm straps to adjust the pressure and find the optimal massage spot.
- **Lower Back:** Place the massager against your lower back while seated in a chair or reclining on a sofa.
- **Abdomen:** Apply gently to the abdominal area for a soothing massage.
- **Calves:** Position the massager around your calves for leg muscle relief after exercise or a long day.



## Flexible Massage Solution




**Neck/Shoulder**



**Lower Back**



**Abdomen**

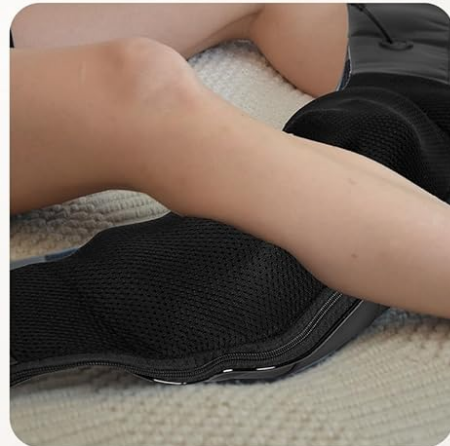

**Calf**

*Image 6.1: Various application areas for the COMFIER massager.*

## Convenient to Use in Any Scenario



*Image 6.2: The massager can be used conveniently in various settings, including a car.*

## 7. CARE AND MAINTENANCE

### 7.1 Cleaning

- Always unplug the device before cleaning.
- Wipe the surface of the massager with a soft, slightly damp cloth.
- Do not use abrasive cleaners, solvents, or harsh chemicals.
- Never immerse the device in water or any other liquid.
- The premium PU leather and breathable mesh materials are designed for easy spot cleaning.

## 7.2 Storage

- Store the massager in a cool, dry place when not in use.
- Keep it away from direct sunlight and extreme temperatures.
- Ensure the device is out of reach of children.
- Do not place heavy objects on the massager.

## 8. TROUBLESHOOTING

If you encounter any issues with your massager, please refer to the following common solutions:

- **Device not turning on:** Ensure the battery is fully charged. Connect the massager to a power source and allow it to charge for at least 30 minutes before attempting to power on again.
- **No heat:** Verify the heat function is activated via the Heating button. Allow a few minutes for the heat to become noticeable, as it is designed to be gentle.
- **Massage feels too intense/weak:** Adjust the intensity level using the Intensity button. You can also modify the pressure by holding the arm straps more loosely or firmly against your body.
- **Device stops unexpectedly:** The massager has a built-in 10-minute auto shut-off feature. If you wish to continue your massage session, simply restart the device. Overheating protection may also cause the device to shut off; allow it to cool down before restarting.

## 9. SPECIFICATIONS

Feature	Detail
Model Number	CF-6505
Product Dimensions	7.28 x 6.46 x 15.59 inches
Item Weight	3.5 Pounds
Batteries	1 Nonstandard Battery (included)
Power Source	Battery Powered
Material	Plastic, Polyester, Premium PU Leather, Breathable Mesh
Manufacturer	COMFIER
Battery Capacity	2200mAh
Charging Time	Approximately 4 hours (USB-C)
Cordless Use Time	Up to 1 hour
Auto Shut-off	10 minutes










## 10. WARRANTY AND SUPPORT

For warranty information, product registration, or customer support, please refer to the documentation included with your purchase or visit the official COMFIER website. Please retain your purchase receipt as proof of purchase for any warranty claims.

© 2025 COMFIER. All rights reserved.

### Related Documents - CF-6505

 User Manual  CF-6701 DEEP-KNEADING SHIATSU NECK & SHOULDER MASSAGER	<a href="#">Comfier CF-6701 Deep-Kneading Shiatsu Neck &amp; Shoulder Massager User Manual</a> User manual for the Comfier CF-6701 Deep-Kneading Shiatsu Neck & Shoulder Massager, covering product features, technical data, setup, operation, safety instructions, care, troubleshooting, and warranty information.
 COMFIER shiatsu neck and shoulder massager with heat ES neck & shoulder shiatsu massager with heat GS massager shiatsu pour le cou et les épaules avec chaleur FI massaggiatore shiatsu per collo e spalle con calore IT masajeador de cuello y hombros shiatsu con calor ES User Manual CF-6302	<a href="#">COMFIER CF-6302 Shiatsu Nacken- und Schultermassagegerät mit Wärmefunktion – Bedienungsanleitung</a> Entdecken Sie die COMFIER CF-6302 Bedienungsanleitung für das Shiatsu Nacken- und Schultermassagegerät mit Wärmefunktion. Erfahren Sie mehr über die Funktionen, Anwendung und Pflege dieses entspannenden Massagegeräts für Nacken und Schultern.
 COMFIER Shiatsu Neck and Back Massager User Manual CF-2309A	<a href="#">Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</a> User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.
 COMFIER Neck & Shoulder Shiatsu Massager with Heat User Manual CF-6302	<a href="#">Comfier Neck &amp; Shoulder Shiatsu Massager with Heat User Manual (CF-6302)</a> Comprehensive user manual for the Comfier CF-6302 Neck & Shoulder Shiatsu Massager with Heat, covering safety instructions, features, setup, operation, care, troubleshooting, and warranty information.

<div><div><div>CONFIER</div><div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div><div>NECK &amp; SHOULDER SHIATSU MASSAGER WITH HEAT</div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div></div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-6302</div></div></div>	<div><div><a href="#">Comfier CF-6302 Neck &amp; Shoulder Shiatsu Massager with Heat User Manual</a></div><div>Comprehensive user manual for the Comfier CF-6302 Neck &amp; Shoulder Shiatsu Massager with Heat. Learn about safety instructions, features, operation, maintenance, and warranty.</div></div>
<div><div><div>CONFIER</div><div><div>SHIATSU MASSAGE PILLOW WITH HEAT</div></div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-1228</div></div></div>	<div><div><div><a href="#">Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</a></div><div>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</div></div></div>