

GDCAEEEGP-NATURAL-150X40CM

Over The Bed and Recliner Table User Manual

Model: GDCAEEEGP-NATURAL-150X40CM | Brand: Generic

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your new Over The Bed and Recliner Table. Designed for convenience and versatility, this mobile laptop cart and workstation is ideal for use in homes and hospitals, offering an adjustable height and sturdy construction. Please read this manual thoroughly before initial use and retain it for future reference.



Image 1: The Over The Bed and Recliner Table positioned over a bed, showcasing its natural wood top and white frame.

2. PRODUCT FEATURES

Your Over The Bed and Recliner Table is equipped with several features designed for comfort, safety, and durability:

- **Wheel Design:** Equipped with smooth-rolling wheels for easy mobility. A simple foot-step mechanism locks the wheels, ensuring stability during use.
- **Adjustable Height:** The table's height can be easily adjusted to suit various activities, from dining to working on a laptop.
- **Rounded Corners:** Features carefully polished rounded corners to minimize the impact of accidental bumps and enhance safety.
- **Sturdy Build:** Constructed with thickened tubing and a stable framework, ensuring high durability and the ability to hold heavy items without wobbling.
- **Sleek Design:** Its minimalist design with clean lines and a natural wood finish blends seamlessly with various room decors.

- **Quality Material:** Made from a combination of a steel frame and solid wood ecological board for reliability and longevity.

Product parameters

Material: Steel frame+solid wood ecological board

Note: Due to manual measurement, there may be errors. Please understand



Image 2: Detailed product dimensions, showing adjustable height ranges and table top sizes (150x40 CM, 150x50 CM, 120x40 CM, 120x50 CM).

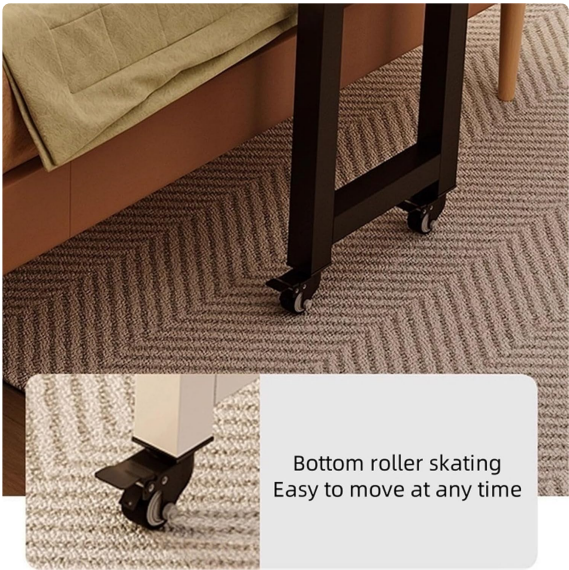


Image 3: Close-up view of the bottom roller wheels, highlighting their ease of movement and the foot-step locking mechanism.



Stable triangle design
Stable and durable
Stronger resistance to pressure



Height adjustable
Multiple adjustable levels
customizable height

Image 4: Illustrates the stable triangle design for enhanced pressure resistance and the multiple adjustable levels for customizable height.



Safe rounded corner design
Carefully polished edges and corners
Reduce collisions



Square tube 5 * 5mm
Bold pipe material
sturdy and durable

Image 5: Shows the safe rounded corner design to reduce collisions and the sturdy 5x5mm square tube material used in construction.

Solid wood ecological board
Rest assured to use



Image 6: A close-up of the solid wood ecological board, emphasizing its quality and suitability for use.



Image 7: The stable framework of the table, demonstrating its load-bearing capability and robust construction.

3. SPECIFICATIONS

Attribute	Detail
Brand	Generic
Model Name	GDCAEEEGP-NATURAL-150X40CM
Shape	Rectangular
Color	Natural
Size (Current Model)	150X40 CM (59.1X15.7 IN)

Adjustable Height Range	67-92 CM (26.4-36.2 IN)
Material	Steel frame + Solid wood ecological board
Item Weight	16000 Grams
Assembly Required	No
Usage	Home, Hospital

4. SETUP AND INITIAL USE

Your Over The Bed and Recliner Table comes pre-assembled, requiring no complex setup. Follow these simple steps for initial use:

- 1. Unpacking:** Carefully remove the table from its packaging. Inspect for any signs of damage during transit.
- 2. Positioning:** Roll the table to your desired location, such as over a bed or recliner. Ensure there is adequate clearance for movement.
- 3. Locking Wheels:** Once positioned, press down on the foot-step locks located on each wheel to secure the table in place. This prevents unintended movement.
- 4. Adjusting Height:** To adjust the table's height, locate the adjustment mechanisms on the legs. Typically, these involve a pin-and-hole system or a lever. Release the lock, raise or lower the tabletop to your preferred height, and then re-engage the lock securely. Ensure both sides are adjusted to the same height for stability.

5. OPERATING INSTRUCTIONS

Using your Over The Bed and Recliner Table is straightforward:

- Moving the Table:** To move the table, first unlock all wheels by lifting the foot-step locks. Gently push or pull the table to its new location.
- Securing the Table:** Always engage the wheel locks when the table is in its desired position to prevent accidental rolling, especially when placing items on it or leaning on it.
- Adjusting for Use:** Adjust the height as needed for different activities. For example, a lower height for reading in bed, and a higher height for using a laptop while sitting upright.
- Weight Capacity:** While sturdy, avoid placing excessively heavy items on the table.

6. MAINTENANCE AND CARE

Proper maintenance will ensure the longevity and optimal performance of your table:

- Cleaning:** Wipe the tabletop and frame with a soft, damp cloth. For stubborn stains, a mild soap solution can be used, followed by drying with a clean cloth. Avoid abrasive cleaners or harsh chemicals that could damage the finish.
- Wheel Care:** Periodically check the wheels for any accumulated dust or debris that might hinder smooth rolling. Clean as necessary.
- Frame Inspection:** Regularly inspect the steel frame and connections for any signs of loosening or

damage.

- **Avoid Extreme Conditions:** Do not expose the table to extreme temperatures, direct prolonged sunlight, or excessive moisture, as this may affect the wood and metal components.

7. TROUBLESHOOTING

If you encounter any issues with your table, refer to the following common solutions:

- **Table Wobbles:**
 - Ensure all wheel locks are fully engaged.
 - Verify that the height adjustment mechanisms on both legs are securely locked at the same level.
 - Check if the table is on an even surface.
- **Difficulty Moving Table:**
 - Ensure all wheel locks are disengaged.
 - Check wheels for debris (hair, lint) and clean if necessary.
- **Height Adjustment Issues:**
 - Ensure no weight is on the table when attempting to adjust height.
 - Check if the locking mechanism is fully disengaged before attempting to move the table up or down.

If problems persist, please contact customer support.

8. SAFETY INFORMATION

- Always ensure the table's wheels are locked when in use to prevent accidental movement.
- Do not overload the table beyond its intended weight capacity.
- Keep children and pets away from the table during height adjustment to prevent pinching.
- Avoid standing or sitting on the table.
- The rounded corners are designed for safety, but exercise caution to avoid bumps.

9. WARRANTY AND SUPPORT

For any questions, concerns, or support needs regarding your Over The Bed and Recliner Table, please contact the seller or manufacturer directly. While specific warranty details are not provided in this manual, standard consumer rights apply. Please retain your purchase receipt for any warranty claims.

For further assistance, refer to the contact information provided at the point of purchase or on the product packaging.

