



[Manuals.plus](#) /

> [CallToU](#) /

> CallToU Wireless Caregiver Pager Smart Call System User Manual

CallToU CC01-01A-06W

CallToU Wireless Caregiver Pager Smart Call System

Model: CC01-01A-06W

1. INTRODUCTION

The CallToU Wireless Caregiver Pager Smart Call System is designed to provide reliable and convenient assistance for individuals requiring care, such as the elderly, sick, or disabled. This system allows caregivers to be alerted instantly when help is needed, ensuring timely response and enhanced safety. Its user-friendly design and robust features make it an essential tool for home and personal attention.



Image: Overview of the CallToU Wireless Caregiver Pager Smart Call System, showing two white plug-in receivers, two red and white call buttons with lanyards, and one orange and white wrist call button.

2. PACKAGE CONTENTS

Please check the package to ensure all items are present:

- 2 x CallToU Plug-in Receivers
- 2 x Call Buttons (with lanyards)
- 1 x Wrist Call Button
- User Manual

3. PRODUCT FEATURES

- **Wireless Coverage:** Offers a wide range of up to 500 feet (150 meters) in open areas, ensuring reliable signal transmission throughout your home or facility.

FULLY EXPANDABLE SYSTEM

500ft distance coverage, Add extra receivers and transmitters for a wider range of security

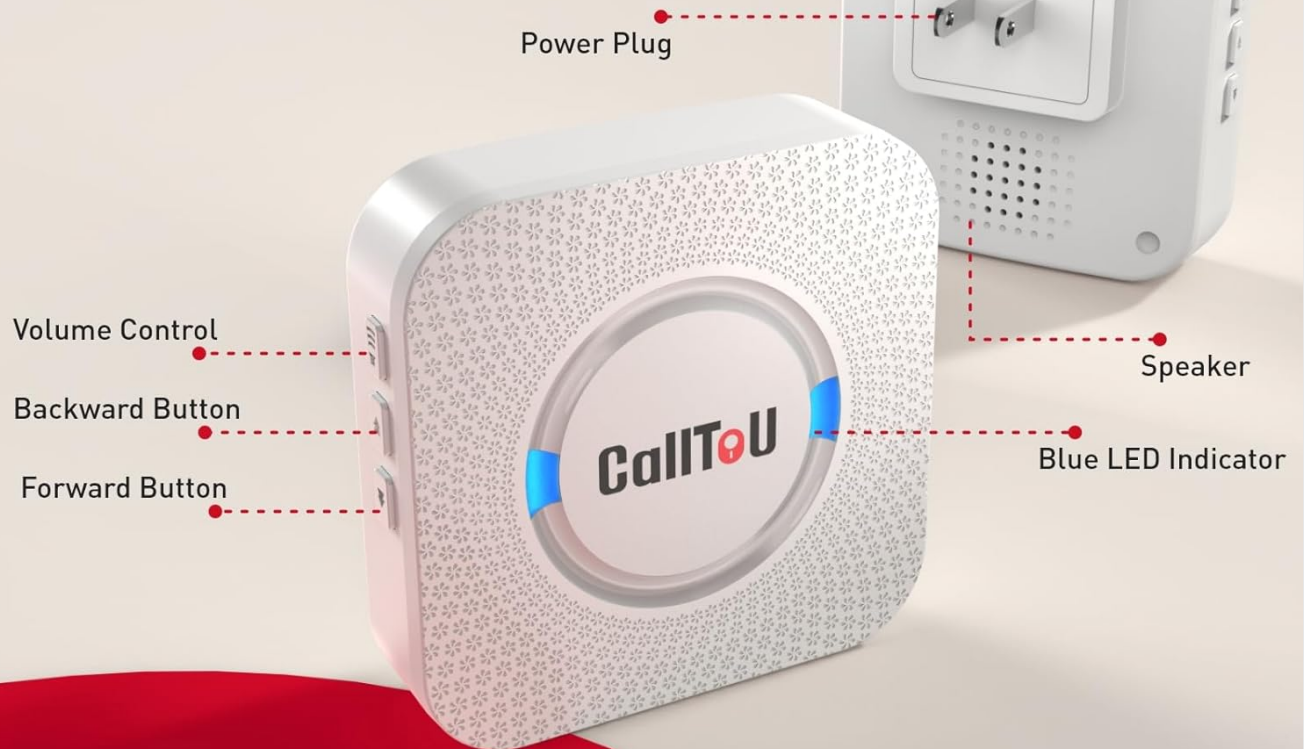


Image: An illustration demonstrating the 500ft wireless coverage of the CallToU system, showing receivers and call buttons placed around a house.

- **Customizable Alerts:** The receiver features 55 selectable ringtones and 5 adjustable volume levels, ranging from 0 dB (mute) to 110 dB, allowing you to choose the perfect alert for your environment.

Easy of Use

Opt for simple, intuitive designs that are easy for seniors to operate.



 **5 Volume Levels**

 **55 Ringtone Options**

 **Support 30 Call Buttons**

Image: A close-up view of the CallToU receiver, highlighting its power plug, speaker, blue LED indicator, volume control, backward button, and forward button, along with text indicating 5 volume levels, 55 ringtone options, and support for 30 call buttons.

- **Durable and Water-Resistant:** The call buttons are rated IP55 waterproof, and the wrist call button is IP66 waterproof, making them suitable for use in bathrooms or during light rain.

IP55 Waterproof



Portable



Replaceable Battery



Long Battery Life

Image: A red and white call button being splashed with water, demonstrating its IP55 waterproof rating.


IP66 Waterproof





Image: A person washing their hands while wearing the orange and white wrist call button, illustrating its IP66 waterproof rating.


- **Portable and Comfortable:** The call buttons can be worn around the neck or secured with a bracket/sticker. The wrist call button features a lightweight, adjustable rubber strap for comfortable wear.

Wrist Call Button

 Easy to Carry

 Rubber Strap

 Light in Quality

 Removable Strap


 Easy to Press



Image: A collage showing the wrist call button being worn, its rubber strap, light quality, removable strap, and ease of pressing, with an elderly person and caregiver in the background.

- **Expandable System:** The receiver supports pairing with up to 30 call buttons, allowing for a fully customizable and expandable alert system.

4. SETUP GUIDE

4.1 Initial Setup

1. Unpack all components from the box.
2. Plug the CallToU Receivers into standard 120V AC electrical outlets in desired locations within your home or care facility. Ensure they are within range of where the call buttons will be used.
3. The receivers will power on and emit a default ringtone.

4.2 Pairing Call Buttons with Receiver

To ensure your call buttons communicate with the receiver, follow these pairing steps:



Image: A four-step diagram illustrating the pairing process for the CallToU receiver and call buttons, showing button presses and LED indicators.

1. **Select Ringtone:** On the receiver, press the 'Forward' (▲) or 'Backward' (▼) button to cycle through the available melodies until you find your preferred ringtone.
2. **Enter Pairing Mode:** Hold the 'Volume Control' button on the receiver until the blue LED indicator flashes and you hear a 'ding' sound. This indicates the receiver is in pairing mode.
3. **Activate Call Button:** Quickly press the desired Call Button (either the pendant or wrist button). The receiver will emit a 'ding' sound, and its LED will flash, confirming successful pairing.
4. **Verify Pairing:** Press the call button/watch again. If the receiver plays the selected ringtone, the pairing is complete. If not, repeat steps 2 and 3.

Repeat these steps for each additional call button you wish to pair with the receiver.

5. OPERATING INSTRUCTIONS

5.1 Making a Call

- To request assistance, simply press the large 'CALL' button on the pendant call button or the 'SOS' button on the wrist call button.
- Upon activation, the paired receiver(s) will immediately sound the selected ringtone and flash its blue LED indicator, alerting the caregiver.



Image: An elderly man in bed pressing his wrist call button, with a dotted red line indicating the signal transmission to a caregiver in another room who

is alerted by the receiver's flashing light and sound.



Image: An elderly woman wearing a pendant call button, pressing it, with a dotted red line showing the signal reaching a caregiver in another room who is alerted by the receiver.

5.2 Adjusting Volume and Ringtones

- **Volume:** Press the 'Volume Control' button on the receiver to cycle through the 5 available volume levels (0 dB to 110 dB).
- **Ringtone:** Use the 'Forward' (▲) or 'Backward' (▼) buttons on the receiver to select a different ringtone from the 55 options.

6. MAINTENANCE

6.1 Battery Replacement (Wrist Call Button)

The wrist call button is powered by a replaceable CR2023 Lithium Coin Battery. A low battery reminder will activate when the battery needs replacement.

Low Battery Reminder



How to Change the Battery?



Image: A three-step diagram showing how to replace the battery in the CallToU wrist call button: 1. Use a tool to lift the battery cover. 2. Replace with a 3V CR2023 Lithium Coin Battery. 3. Put the cover back to complete the replacement.

1. Use a small tool (e.g., a flathead screwdriver) to gently lift the battery cover on the back of the wrist call button.
2. Remove the old battery and insert a new 3V CR2023 Lithium Coin Battery, ensuring correct polarity.
3. Carefully place the battery cover back and press firmly until it clicks into place, ensuring a secure seal for water resistance.

6.2 Cleaning

- Wipe the receivers and call buttons with a soft, dry cloth.
- Do not use abrasive cleaners or solvents.
- Ensure the devices are dry before plugging in or using, especially after exposure to water.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Receiver does not sound when button is pressed.	Button not paired; Receiver unplugged; Out of range; Low battery in call button.	Ensure receiver is plugged in. Re-pair the call button (see Section 4.2). Move receiver closer to call button. Replace battery in call button (see Section 6.1).
Receiver volume is too low/high.	Volume setting needs adjustment.	Press the 'Volume Control' button on the receiver to adjust the volume level.
Receiver plays an unexpected ringtone.	Ringtone setting needs adjustment.	Use the 'Forward' (▲) or 'Backward' (▼) buttons on the receiver to select the desired ringtone.
Call button seems unresponsive.	Battery depleted; Button damaged.	Replace the battery in the call button. If issue persists, contact customer support.

8. SPECIFICATIONS

Feature	Detail
Model	CC01-01A-06W
Brand	CallToU
Power Source (Receiver)	Corded Electric
Voltage (Receiver)	120 Volts (AC)
Control Method	Touch
Ringtones	55 options
Volume Levels	5 levels (0 dB to 110 dB)
Wireless Range	Up to 500 feet (open area)
Call Button Waterproof Rating	IP55
Wrist Call Button Waterproof Rating	IP66
Wrist Call Button Battery	CR2023 Lithium Coin Battery
Max Paired Buttons	30

9. WARRANTY AND SUPPORT

CallToU products are designed for reliability and performance. For specific warranty information, please refer to the documentation included with your purchase or visit the official CallToU website.

If you encounter any issues or have questions regarding your CallToU Wireless Caregiver Pager Smart Call System, please contact CallToU customer support for assistance. Contact details can typically be found on the product packaging or the official brand website.

