

Ya-man YJEA4NAM

Ya-man Reibote Hyper ZERO Light Beauty Device (YJEA4NAM) Instruction Manual

Model: YJEA4NAM

1. INTRODUCTION

Thank you for purchasing the Ya-man Reibote Hyper ZERO Light Beauty Device. This manual provides essential information for the safe and effective use of your device. Please read it thoroughly before first use and keep it for future reference.

The Ya-man Reibote Hyper ZERO is a high-power light beauty device designed for home use, offering efficient hair removal and skin care benefits. It features a long-life lamp, multiple attachments for various body areas, and intelligent functions for optimal results.

2. SAFETY INFORMATION AND PRECAUTIONS

WARNING: This product is not a medical device. Do not use if any of the following apply to you or on the specified body parts.

2.1 Contraindications (Do Not Use If):

- Using medical electrical devices (e.g., implanted pacemakers, life-support systems, wearable electrocardiographs).
- Individuals with heart disease, disabilities, or suspected conditions.
- Infants or individuals unable to express their intentions.
- Individuals unable to perceive temperature or stimulation.
- Individuals with autoimmune diseases.
- Individuals in early pregnancy to immediately postpartum.
- Individuals with hyperpigmentation.
- Individuals with epilepsy.
- Individuals with photosensitivity.

- Individuals using photosensitizing agents for photodynamic therapy.

2.2 Consult a Doctor Before Use If:

- Acute illness.
- Infectious disease.
- Malignant tumor.
- Febrile illness.
- High blood pressure or blood disorders.
- Fever.
- Individuals with sensory impairment due to severe peripheral circulatory dysfunction (e.g., diabetes).
- Nerve disorders.
- Skin diseases or atopic dermatitis.
- Under medical treatment or taking medication.
- Allergic constitution, contact dermatitis, etc., or skin discomfort during use.

2.3 Do Not Use On The Following Areas:

- Areas with bone fractures.
- Skin with wounds, rashes, inflammation, or other abnormalities.
- While intoxicated.
- Head, inside the mouth, intimate areas*, mucous membranes.
- Eyeballs, eyelids (eye sockets).
- Throat.
- Swollen moles.
- Viral warts.
- Scars, blemishes, or dark areas.
- Eyebrows, forehead, scalp, nape of the neck.
- Sideburns.
- Ears.
- Nipples, areola.
- Inside the navel.
- Areas where blood vessels are visible, such as large veins.
- Tattoos, body paint.
- Areas with internal bleeding, swelling, or scabs after cosmetic medical procedures.

*For intimate areas, please consult a doctor after 2 weeks post-procedure if there are no issues. (Examples: Hyaluronic acid injection, silicone injection, Botox, thread lift, double eyelid surgery, HIFU, etc.)

If you feel unwell or notice any skin abnormalities during use, stop immediately. If symptoms persist, consult a doctor.

3. PRODUCT OVERVIEW AND FEATURES

The Ya-man Reibote Hyper ZERO is a compact yet powerful light beauty device designed for comprehensive hair removal and skin care. It incorporates Ya-man's advanced technology to deliver

efficient results.

- **High Power:** Equipped with a W-lamp for powerful light output, effectively targeting even thick hair and facial hair.
- **Speed Care:** Achieve full-body care in approximately 3 minutes with an irradiation interval of about 0.2 seconds in roller mode.
- **Skin Beauty Benefits:** Continuous use helps make unwanted hair less noticeable, and the light treatment promotes moisturized and firm skin.
- **Long-Life Lamp:** Features a durable lamp with approximately 1.2 million flashes, suitable for family use.
- **Intelligent Roller:** A uniquely developed intelligent roller helps prevent missed spots and over-treatment, ensuring speedy and uniform care.
- **Automatic Level Adjustment:** A built-in skin tone sensor reads your skin color density and automatically adjusts the light level to the optimal setting from 5 stages.



Image: The Ya-man Reibote Hyper ZERO device in use, highlighting its sleek design and application on the leg.

4. COMPONENTS

The Ya-man Reibote Hyper ZERO package includes the following items:

- Main Unit
- AC Adapter
- Roller Attachment
- LED Attachment
- Spot Attachment
- V-Design Tool
- Pearl Skin Lotion (100ml)



Image: All components included with the Ya-man Reibote Hyper ZERO device, including the main unit, various attachments, the V-Design Tool, and the accompanying Pearl Skin Lotion.

5. SETUP

1. **Unpack the Device:** Carefully remove all components from the packaging.
2. **Connect the AC Adapter:** Insert the AC adapter into the device's power port and then plug it into a wall outlet. The device operates with a wired connection.
3. **Attach the Desired Head:** Select the appropriate attachment (Roller, LED, or Spot) for the area you intend to treat and securely attach it to the main unit.
4. **Prepare the Skin:** Ensure the treatment area is clean, dry, and free of any lotions or oils. Shave the area thoroughly before use to maximize effectiveness and prevent discomfort.

6. OPERATING INSTRUCTIONS

6.1 Power On/Off

- To turn on, press and hold the power button.
- To turn off, press and hold the power button again.

6.2 Selecting Modes

The device offers different modes for various treatment needs:

- **Roller Mode:** For continuous, high-speed treatment over large areas. The intelligent roller ensures even coverage.
- **Quick Mode:** Provides rapid continuous flashes, ideal for efficient treatment.
- **Single Mode:** For precise, single-flash treatment, recommended for sensitive areas like the face and underarms.



クイックモードは10回、
ローラーモードは
7回連続照射で
全身を高速ムダ毛ケア



連続照射モードで時短ケア
単発照射は顔・脇にオススメ

Image: Demonstrates the Quick Mode and Roller Mode for continuous irradiation, and Single Mode for precise application on the face and underarms.

6.3 Adjusting Irradiation Levels

The device has 5 adjustable irradiation levels (MIN to MAX) to suit your skin sensitivity and hair type.

- Press the "LEVEL" button to cycle through the intensity levels.
- Start with a lower level and gradually increase if comfortable.

ハイパワーでしっかりケアしたい方はもちろん、
痛みをなるべく抑えてお手入れしたい方も
心地よくお使いいただけます



**照射レベルを
5段階から選べます**

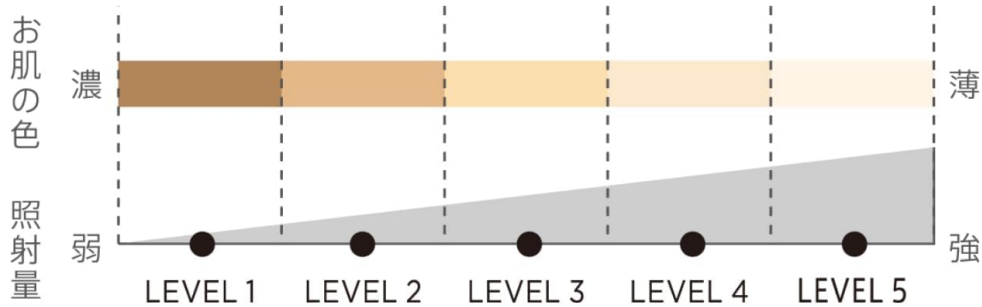
Image: Close-up of the device's control panel, illustrating the 5-stage irradiation level selection from MIN to MAX.

6.4 Skin Tone Sensor and Auto Mode

The device features an "AUTO" mode that uses a skin tone sensor to detect your skin color and automatically select the optimal irradiation level from the 5 available stages. This ensures safe and effective treatment tailored to your skin.

- To activate, select "AUTO" mode.
- Place the device on your skin; the sensor will detect your skin tone and adjust the flash intensity.

オートレベルでは、
肌の色濃度をセンサーで読み取り、
最適な光レベルを5段階で
自動調整します



肌色センサーで 最適レベルを自動調整

Image: A diagram illustrating how the skin tone sensor automatically adjusts the optimal light level based on skin color density, from light to dark.

6.5 Treatment Frequency

For optimal results, it is generally recommended to use the device once a week. Hair growth cycles vary, and consistent use over several months is necessary to see significant reduction in hair growth. Do not treat the same area more than once per day.

7. ATTACHMENTS AND USAGE AREAS

The Ya-man Reibote Hyper ZERO comes with three specialized attachments for different body parts:

- **Roller Attachment:** Designed for large body areas like legs and arms. Its wide surface and high output ensure efficient treatment.
- **LED Attachment:** Specifically for facial use, this attachment includes LEDs to enhance skin care benefits during treatment.
- **Spot Attachment:** Ideal for precise treatment on smaller or sensitive areas, such as the VIO line (bikini area) and underarms.

ケアにあったアタッチメントを選べます



**ローラー
アタッチメント**

大面積高出力の
ボディ向け
アタッチメント



**LED
アタッチメント**

LED付きの
フェイス特化
アタッチメント



**スポット
アタッチメント**

スポット照射や
VIOラインに特化した
アタッチメント

適切な出力でのケアが可能な 3つのアタッチメント

Image: The three distinct attachments: Roller Attachment for large areas, LED Attachment for facial care, and Spot Attachment for precise VIO line treatment.

8. POST-TREATMENT CARE

After using the light beauty device, it is highly recommended to moisturize your skin. The included Pearl Skin Lotion is formulated to provide hydration and promote skin elasticity and radiance.

8.1 Pearl Skin Lotion Ingredients:

Water, DPG, Glycerin, Pentylene Glycol, Diglycerin, Triethylhexanoin, Caprylic/Capric Triglyceride, Betaine, BG, Sodium Hyaluronate, Hydrolyzed Hyaluronic Acid, Dipotassium Glycyrrhizate, Papain, Laminaria Ochroleuca Extract, Centifolia Rose Flower Extract, Pomegranate Flower Extract, Camellia Japonica Seed Extract, Daisy Flower Extract, Hydrogenated Lecithin, Sucrose Cocoate, Polysorbate 80, Acrylates/Hydroxyethyl Acrylate/Sodium Acryloyldimethyl Taurate Copolymer, Sorbitan Stearate, Dextrin, Ethylhexylglycerin, Rosemary Leaf Oil, Meadowfoam Oil.

光美容器によるムダ毛ケアの後に 使うのがおすすめです



レイボーテシリーズを使ってムダ毛ケアをした後は...

パーリースキンローションで保湿ケア♪

しっかり保湿して、
うるおいとハリつやに溢れたお肌へ

Image: The Ya-man Reibote Hyper ZERO device alongside the Pearl Skin Lotion, emphasizing the importance of moisturizing after light hair removal.

9. MAINTENANCE

Proper maintenance ensures the longevity and performance of your device.

- **Cleaning the Device:** After each use, wipe the light window and the device body with a soft, dry cloth. Ensure no hair or debris remains on the light window.
- **Cleaning Attachments:** Remove attachments and clean them separately. Do not use abrasive cleaners or immerse the main unit in water.
- **Storage:** Store the device and its accessories in a cool, dry place, away from direct sunlight and out of reach of children.

10. TROUBLESHOOTING

If you encounter any issues with your Ya-man Reibote Hyper ZERO, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	AC adapter not properly connected or power outlet issue.	Ensure the AC adapter is securely plugged into both the device and a working power outlet.
No flash or irregular flashing.	Device not making full contact with skin; skin tone sensor issue; device is wet.	Ensure the light window is fully pressed against the skin. Make sure the skin is dry. Clean the light window if dirty.
Discomfort or excessive heat during use.	Irradiation level is too high for your skin type or area.	Reduce the irradiation level. Ensure the skin is clean and dry. Do not treat the same area multiple times in quick succession.
Hair does not seem to be reducing.	Insufficient treatment duration or frequency; hair is not in the active growth phase.	Consistent use over several months is required. Ensure you are using the correct level and technique. Shave the area thoroughly before each treatment.

If the problem persists after trying these solutions, please contact customer support.

11. SPECIFICATIONS

Feature	Detail
Brand	Ya-man
Model Number	YJEA4NAM
Product Dimensions	16.9 cm (Length) x 5.4 cm (Width)
Package Dimensions	27.7 x 23.6 x 12.9 cm
Package Weight	1.27 kg
Power Source	AC Adapter (Wired)
Irradiation Count (Max)	Approx. 1,200,000 flashes
Irradiation Interval (Min)	Approx. 0.2 seconds (in roller mode, level 1)
Color	Gold
Included Components	AC Adapter, LED Attachment, V-Design Tool, Spot Attachment, Roller Attachment

12. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the warranty card included with your product or visit the official Ya-man website. Keep your proof of purchase for warranty claims.



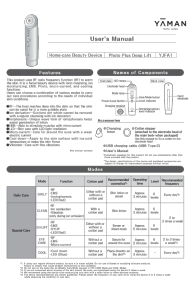
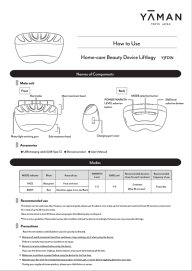
If you have any questions or require assistance with your device, please contact Ya-man customer service through their official channels.



Ya-man Official Website: www.ya-man.com

© 2023 Ya-man Ltd. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.

Related Documents - YJEA4NAM

	<p>YA-MAN Photo PLUS Hyper HRF-11 Instruction Manual</p> <p>Comprehensive instruction manual for the YA-MAN Photo PLUS Hyper HRF-11 home beauty device, covering safety precautions, features, treatment modes, usage instructions, troubleshooting, and after-service information.</p>
	<p>YA-MAN Bloom 6 User Manual: Operation, Safety, and Maintenance</p> <p>Comprehensive user manual for the YA-MAN Bloom 6 home-care beauty device. Learn about safe operation, features, troubleshooting, and maintenance for optimal results.</p>
	<p>YA-MAN Photo Plus Deep Lift YJFA1 User Manual - Home Care Beauty Device</p> <p>Comprehensive user manual for the YA-MAN Photo Plus Deep Lift YJFA1, a home-care beauty device featuring RF, EMS, Iontophoresis, LED, and Micro-current technologies for facial skin treatment.</p>
	<p>YA-MAN Liftlogy Home-care Beauty Device: How to Use Guide</p> <p>A comprehensive guide on how to use the YA-MAN Liftlogy home-care beauty device, covering its features, modes, and recommended usage for face and body treatments.</p>

	<p>YA-MAN ACETINO SMART IB-7 Facial and Body Massager User Manual</p> <p>User manual for the YA-MAN ACETINO SMART IB-7, a home-use beauty device for facial and body treatments. Includes safety precautions, operating instructions, maintenance, and after-sales service information.</p>
	<p>YA-MAN YJMF4 Blue Green Mask Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN YJMF4 Blue Green Mask Lift home-care beauty device, covering setup, usage instructions, safety precautions, troubleshooting, and product specifications.</p>