

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [SALLOOT](#) /

› [SALLOOT Nova P5 IPL Laser Hair Removal Device User Manual](#)

SALLOOT D-T003B

SALLOOT Nova P5 IPL Laser Hair Removal Device User Manual

Model: D-T003B

1. PRODUCT OVERVIEW

The SALLOOT Nova P5 is a compact IPL (Intense Pulsed Light) device designed for at-home hair removal. It incorporates advanced sapphire ice-cooling technology and dual pulse functionality to provide a comfortable and effective treatment experience. This device is suitable for various body areas, including the bikini line, face, legs, and armpits, aiming for long-term hair reduction.



Image 1: The SALLOOT Nova P5 IPL Laser Hair Removal Device.

2. SAFETY INFORMATION

Before using the SALLOOT Nova P5, please read all safety instructions carefully. Improper use can lead to skin irritation or injury.

- **Skin Tone and Hair Color Compatibility:** IPL technology is most effective on certain skin tones and hair colors. Refer to the 'Skin and Hair Suitability Chart' in Section 6 to ensure the device is appropriate for your complexion and hair. Using the device on incompatible skin or hair types may result in burns or ineffective treatment.
- **Eye Protection:** Always wear the provided protective eyewear during treatment to shield your eyes from the intense pulsed light. Do not look directly into the light window.
- **Skin Preparation:** Ensure the treatment area is clean, dry, and shaved before use. Do not use on unshaved skin, as this can cause discomfort and reduce effectiveness.

- **Avoid Sensitive Areas:** Do not use the device on tattoos, permanent makeup, moles, freckles, or areas with skin conditions such as eczema, psoriasis, or open wounds.
- **Pregnancy and Medical Conditions:** Consult a healthcare professional before use if you are pregnant, breastfeeding, or have any underlying medical conditions, especially those related to skin sensitivity or light.
- **Children:** Keep the device out of reach of children. This device is not intended for use by individuals under 18 years of age.
- **Power Supply:** Use only the original power adapter and cable provided with the device.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- SALLOOT Nova P5 IPL Hair Removal Device
- Power Adapter
- Power Cable
- Protective Eyewear
- Shaver
- Storage Pouch
- User Manual
- App Download/Connection Guide

4. PRODUCT FEATURES

The SALLOOT Nova P5 is equipped with several features designed for effective and comfortable hair removal:

- **Full Spectrum Long Lightwave (560-1200nm):** Utilizes a broad spectrum of light for safe and permanent hair removal.
- **Super Clean Dual Pulses (22J Energy):** Targets deep hair follicles, contributing to up to 93% hair reduction within 4 weeks. The 78W powerful adapter ensures consistent energy delivery.
- **Advanced Silkcryo Cooling Technology:** Features a 360° sapphire ice-cooling window that instantly cools the skin contact surface to 39°F (4°C) and maintains 59°F (15°C) for up to 35 minutes of continuous use, enhancing comfort, especially for sensitive areas.
- **Compact and Lightweight Design:** Weighing only 188g, the Nova P5 is highly portable.
- **Rapid Auto-Glide Flashing Mode:** Delivers 120 flashes per minute, allowing for faster full-body sessions in approximately 8-10 minutes.
- **AI Coaching App:** A smart application that provides personalized treatment plans, recommends energy levels, and suggests suitable flash settings based on your skin tone, hair color, thickness, and treatment area. It also offers real-time AI analysis reports.



Image 2: Full Spectrum Lightwave for Hair Removal.

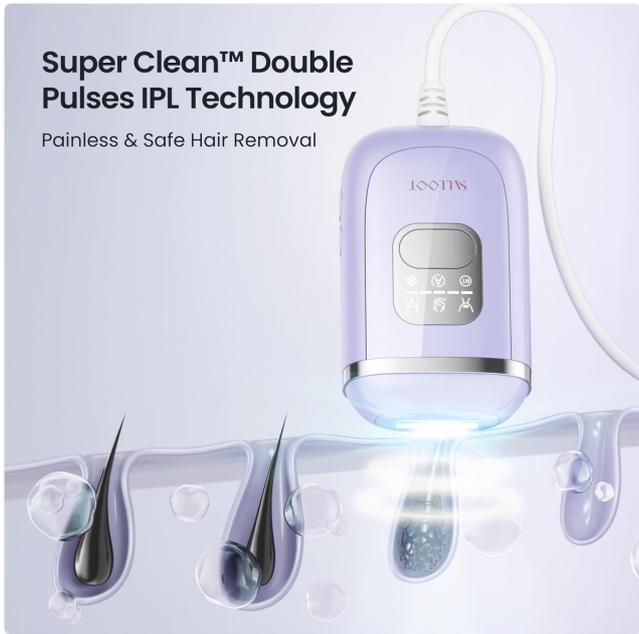


Image 3: Super Clean Dual Pulses IPL Technology.



Image 4: Advanced Silkryo Cooling Technology.



Image 5: Mini Tech Yet Powerful Design.

AI Smart Coaching Personalized Treatments

The simplest hair removal, Like a beautician at home

- ✓ skin tone and hair color match
- ✓ personalized treatment
- ✓ step-by-step tutorials
- ✓ Precision reminder
- ✓ Results tracking & AI reports
- ✓ Smart email reminder

The bottom section of the image features three panels. The first panel shows a woman using the device on her arm, with a smartphone displaying a 'Step-by-Step Tutorials' app screen. The second panel shows the woman holding the device, with a smartphone displaying a 'Personalized Treatments' app screen. The third panel shows a close-up of the device being used on an arm, with a smartphone displaying a 'Results Tracking' app screen.

Image 6: AI Smart Coaching Personalized Treatments.

Compact, Powerful And Travel Friendly

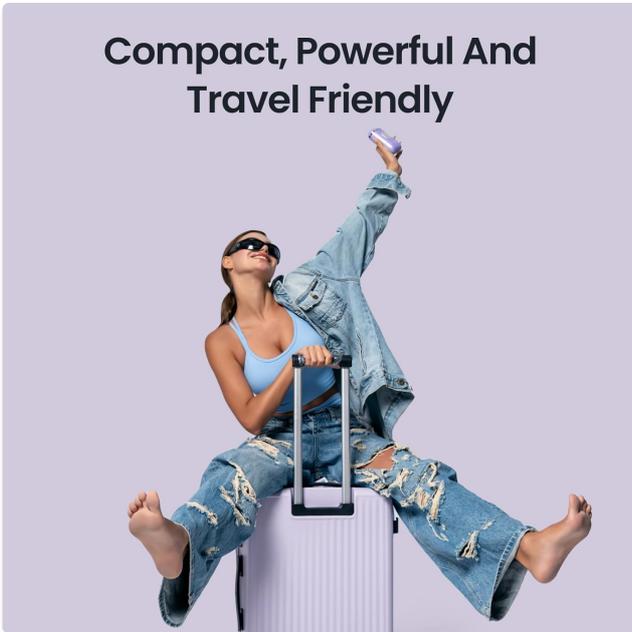


Image 7: Compact and Travel-Friendly Design.

5. SETUP

Follow these steps to set up your SALLOOT Nova P5 device:

1. **Unpack the Device:** Carefully remove all components from the packaging.
2. **Connect Power:** Connect the power cable to the adapter, and then plug the adapter into the device. Finally, plug the power cable into a suitable electrical outlet.
3. **Download the App (Optional but Recommended):** Scan the QR code provided in the 'App Download/Connection Guide' or search for the 'CTONE' app in your device's app store (Google Play or App Store). The app offers personalized treatment plans and guidance.
4. **Connect to App:** Follow the in-app instructions to connect your device via Bluetooth. This typically involves opening the app, ensuring Bluetooth is enabled, and powering on the device.

6. OPERATION

To use your SALLOOT Nova P5 for hair removal:

1. **Prepare the Skin:** Shave the area to be treated. Ensure the skin is clean, dry, and free from any lotions, creams, or deodorants.
2. **Wear Eye Protection:** Always put on the protective eyewear before starting the treatment.
3. **Power On:** Press and hold the power button on the device to turn it on.
4. **Select Energy Level:** The device typically has multiple energy levels. Start with the lowest setting and gradually increase it as your skin tolerates. The app can help recommend an appropriate level.
5. **Choose Mode:** Select between manual flash mode (for small, precise areas) or auto-glide mode (for larger areas like legs).
6. **Apply to Skin:** Place the device's light window firmly against your skin, ensuring full contact. The device will not flash unless proper contact is made.
7. **Initiate Flash:** In manual mode, press the flash button. In auto-glide mode, glide the device slowly over the skin while it automatically flashes.
8. **Treat the Area:** Move the device systematically across the treatment area, ensuring each section is covered once. Avoid overlapping flashes excessively.

9. **Power Off:** After completing the treatment, press and hold the power button to turn off the device.
10. **Post-Treatment Care:** You may apply a soothing lotion or gel to the treated area if desired. Avoid sun exposure and use sunscreen on treated areas.

7. TREATMENT SCHEDULE

For optimal results, follow the recommended treatment schedule:

- **Weeks 1-2:** 3-4 times per week
- **Weeks 3-4:** 2-3 times per week
- **After 4 Weeks:** 2 times per month for long-lasting results

Note: For best results, complete the full 8-week treatment cycle. Up to 93% hair reduction can be achieved within 4 weeks.

How Often Should I Use It?

weeks 1-2: 3-4 times/week
weeks 3-4: 2-3 times/week
After 4 weeks: 2 times/monthly for long-lasting results

Results: Up To 93% Hair Reduction Within 4 Weeks.*
 Note : For Best Results, Complete The Full 8-Week Treatment Cycle.



Image 8: Recommended Treatment Schedule.

8. SKIN AND HAIR SUITABILITY

The effectiveness and safety of IPL hair removal depend on your skin tone and hair color. IPL primarily targets melanin in the hair follicle. Darker skin tones have higher melanin levels, which increases the risk of burns. Light hair colors lack sufficient melanin for IPL to be effective.

Identify Your Skin Tone and Hair Color:



PCOS Friendly

Is Nova P5 Right for My Skin?

IPL primarily targets melanin, which means darker skin tones with higher melanin levels increase the risk of burns. Light hair color lacks sufficient melanin for IPL to be effective.

Step 1 Identify your skin tone

<input checked="" type="checkbox"/> White	<input checked="" type="checkbox"/> Beige	<input checked="" type="checkbox"/> Light Brown	<input checked="" type="checkbox"/> Brown	<input checked="" type="checkbox"/> Dark Brown
---	---	---	---	--

Step 2 Identify the color of your unwanted hair

<input checked="" type="checkbox"/> Black	<input checked="" type="checkbox"/> Dark Brown	<input checked="" type="checkbox"/> Brown	<input checked="" type="checkbox"/> Light Brown	<input checked="" type="checkbox"/> Red	<input checked="" type="checkbox"/> Light Blonde	<input checked="" type="checkbox"/> White/Gray
---	--	---	---	---	--	--

Skin Tone + Hair Color = Applicable
 Either Skin Tone or Hair Color = Not Applicable

Image 9: Skin Tone and Hair Color Suitability Chart.

Skin Tone + Hair Color = Applicable: The device is suitable for your combination.

Either Skin Tone or Hair Color = Not Applicable: The device is not recommended for your combination.

If you are unsure about your suitability, please consult the app or contact customer support.

9. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and effectiveness of your device:

- **Cleaning:** After each use, wipe the light window and the device surface with a soft, dry cloth. Do not use abrasive cleaners or solvents. Ensure the device is unplugged before cleaning.
- **Storage:** Store the device in its original pouch or a clean, dry place away from direct sunlight and extreme temperatures.
- **Inspection:** Regularly inspect the device, power adapter, and cable for any signs of damage. Do not use if any part is damaged.

10. TROUBLESHOOTING

If you encounter any issues with your SALLOOT Nova P5, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power button not held long enough.	Ensure the power adapter and cable are securely connected and plugged into a working outlet. Press and hold the power button for 2-3 seconds.
Device does not flash.	Light window not in full contact with skin; skin tone/hair color incompatibility; device malfunction.	Ensure the light window is pressed firmly and completely against the skin. Check the skin tone/hair color chart. If the issue persists, contact customer support.
Skin feels too hot or uncomfortable.	Energy level is too high for your skin sensitivity.	Reduce the energy level. Ensure the sapphire ice-cooling function is active. If discomfort continues, discontinue use.
App connection issues.	Bluetooth not enabled; device not in pairing mode; app outdated.	Ensure Bluetooth is on and the device is powered on. Restart the app and device. Check for app updates.

For further assistance, please refer to the 'Warranty and Support' section.

11. SPECIFICATIONS

- **Model Number:** D-T003B
- **UPC:** 762949721384
- **ASIN:** B0F5MRM54
- **Package Dimensions:** 8.15 x 7.01 x 4.21 inches
- **Item Weight:** 2.98 Pounds
- **Manufacturer:** Shenzhen Beauty Every Moment intelligent electric Co., Ltd.

- **Energy Output:** Up to 22J
- **Cooling Temperature:** Down to 39°F (4°C) on skin contact, maintains 59°F (15°C)
- **Lightwave Spectrum:** 560-1200nm
- **Flash Speed:** 120 flashes/minute

12. WARRANTY AND SUPPORT

SALLOOT is committed to providing high-quality products and customer satisfaction. Your Nova P5 device comes with a standard warranty. For specific warranty details, please refer to the documentation included in your package or visit the official SALLOOT website.

If you encounter any issues, have questions about product usage, or require technical assistance, please contact our customer support team. Contact information can be found in the 'App Download/Connection Guide' or on the SALLOOT official website.

Please kindly check your skin tone and hair color for compatibility before use. Any issues or concerns, please reach out to our customer support team for guidance.