

KIMTAR B0F7FB1GDG

KIMTAR LED Face Mask User Manual

Model: B0F7FB1GDG

1. INTRODUCTION

Thank you for choosing the KIMTAR LED Face Mask. This device is designed to provide professional-grade light therapy for facial skin care in the comfort of your home. It utilizes advanced LED technology with 7 distinct light modes to address various skin concerns, promoting skin rejuvenation and overall skin health. Please read this manual thoroughly before first use to ensure safe and effective operation.



Image 1.1: The KIMTAR LED Face Mask with its remote control unit.

2. SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- Do not use the device if you have any open wounds, skin infections, or severe skin conditions. Consult a dermatologist if you have concerns.
- Avoid direct eye contact with the LED lights. The mask includes eye protection; ensure it is properly positioned during use.
- Keep the device away from water or excessive moisture.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.
- Store the device out of reach of children.
- If you experience any discomfort or adverse reactions, discontinue use immediately and consult a medical professional.
- The infrared light is invisible to the naked eye. Do not assume the device is not working if you do not see a red glow for infrared mode.

3. PRODUCT COMPONENTS

The KIMTAR LED Face Mask package typically includes:

- KIMTAR LED Face Mask
- Remote Control Unit
- USB Charging Cable
- User Manual



Image 3.1: Multi-function remote control diagram.

4. SETUP AND FIRST USE

Follow these steps to set up your KIMTAR LED Face Mask for the first time:

1. **Charge the Device:** Connect the remote control unit to a power source using the provided USB charging cable. The red light on the remote control will indicate charging. Ensure the device is fully charged before initial use.
2. **Prepare Your Skin:** Cleanse your face thoroughly to remove any makeup, dirt, or oils. Pat your skin dry.
3. **Position the Mask:** Place the LED face mask comfortably on your face. Adjust the straps to ensure a snug but not overly tight fit. The mask features soft silicone eye pads for comfort and protection.
4. **Connect Remote:** Ensure the mask is securely connected to the remote control unit via its USB connection port.



Image 4.1: Visual guide for usage steps.

5. OPERATING INSTRUCTIONS

Operating your KIMTAR LED Face Mask is straightforward:

1. **Power On:** Press and hold the Power On/Off/Time Selection key on the remote control to turn on the device.
2. **Select Light Mode:** Press the Light Switch button to cycle through the 7 available light colors. Each press will change the color.
3. **Set Timer:** Use the dedicated timer buttons (5 minutes, 10 minutes, 15 minutes) to select your desired treatment duration. The device will automatically shut off after the selected time.
4. **Infrared Light:** The infrared light can be toggled on or off using the Infrared Light On/Off key. Note that infrared light is invisible.
5. **Relax During Treatment:** Once the light mode and timer are set, relax for the duration of the treatment. The hands-free design allows for various activities during use.
6. **Power Off:** After the session, the device will automatically turn off. You can also manually power it off by pressing and holding the Power On/Off/Time Selection key.
7. **Post-Treatment:** Follow with your regular skincare routine, such as applying serums or moisturizers.

Recommended Usage:

- Use for 15 minutes per day.
- Consistency is key for optimal results.

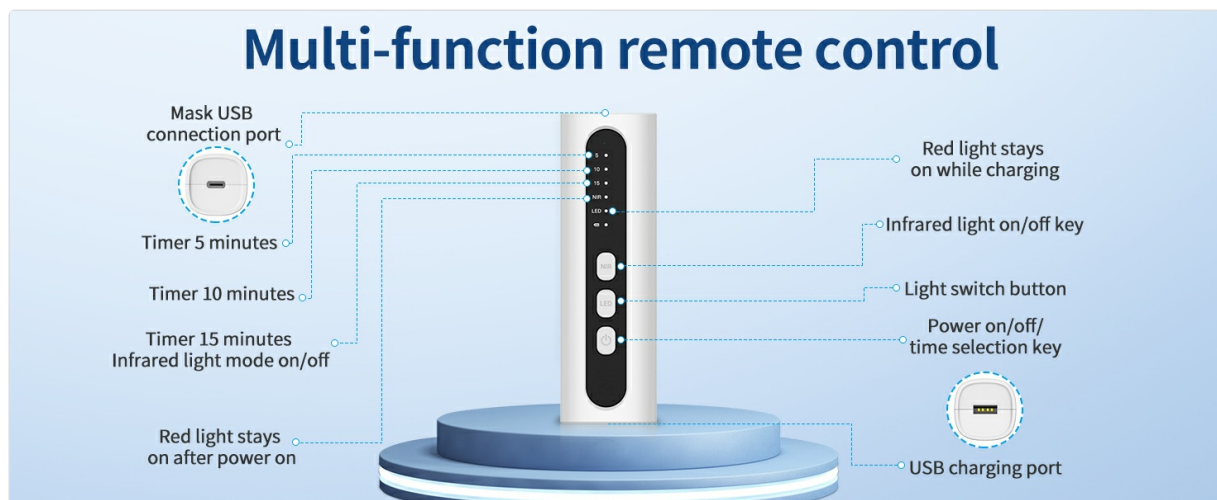


Image 5.1: The mask's hands-free design allows for use during various activities.

6. LIGHT THERAPY MODES

The KIMTAR LED Face Mask offers 7 distinct light modes, each targeting specific skin concerns:

- **Infrared Light (850nm):** Invisible light that penetrates deeply to stimulate collagen production and promote skin healing.
- **Red Light (620-625nm):** Helps reduce fine lines and wrinkles, promotes collagen synthesis, and improves blood circulation.
- **Green Light (520-525nm):** Balances skin tone, reduces hyperpigmentation, and calms sensitive skin.
- **Blue Light (390nm):** Targets acne-causing bacteria, reduces inflammation, and helps minimize pores.
- **Yellow Light (Dual frequency of red and green light):** Improves oxygen exchange in cells and brightens the complexion.

- **Purple Light (Dual frequency of red and blue light):** Combines the benefits of red and blue light, aiding in skin repair and acne reduction.
- **Orange Light:** Helps alleviate skin and color problems, promoting a vibrant glow.

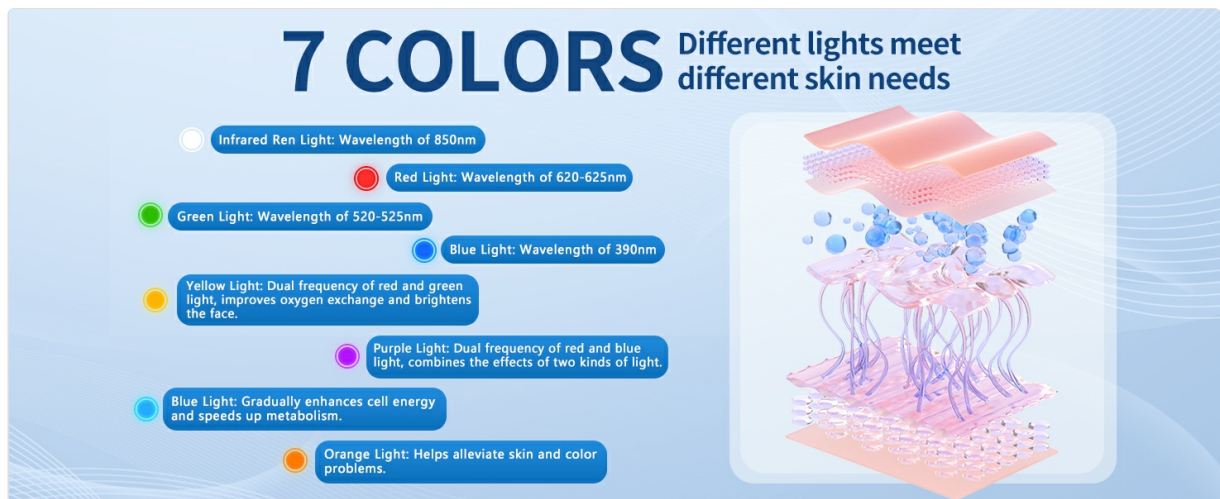


Image 6.1: Overview of 7 light colors and their skin benefits.

KIMTAR

LED 7 COLOR LIGHT SKIN BEAUTY CONTOUR FACE

- Silicone Material
- Soft Eye Pads
- 288 LEDs
- 576 Beams of Light
- 7 Colors Light
- Adjustable Strap

Image 6.2: Comparison of light penetration depth.

7. MAINTENANCE

Proper care ensures the longevity and effectiveness of your device:

- **Cleaning:** After each use, gently wipe the mask's interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** To maintain battery health, avoid fully discharging the device frequently. Charge it regularly, even if not in use for extended periods.

8. TROUBLESHOOTING

If you encounter any issues with your KIMTAR LED Face Mask, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully using the provided USB cable.
Lights are not visible.	Infrared light is invisible. Incorrect light mode selected.	Ensure you have selected a visible light mode (Red, Green, Blue, Yellow, Purple, Orange). Infrared light is not visible.
Remote control not responding.	Remote control battery low or connection issue.	Charge the remote control. Ensure the mask is properly connected to the remote control unit.
Mask feels uncomfortable.	Straps are too tight or mask is improperly positioned.	Adjust the straps for a comfortable fit. Reposition the mask to ensure eye protection is correctly aligned.

If the problem persists after attempting these solutions, please contact customer support for further assistance.

9. SPECIFICATIONS

- **Brand:** KIMTAR
- **Model:** B0F7FB1GDG
- **Color:** White
- **Item Form:** Stick (referring to the mask's flexible material)
- **Skin Type:** All skin types
- **Special Features:** Rechargeable, Reusable
- **Number of Light Modes:** 7 (Red, Green, Blue, Yellow, Purple, Orange, Infrared)
- **Power Source:** 1 Lithium-ion battery (included)
- **Package Dimensions:** 35.1 x 24.6 x 4.7 cm
- **Item Weight:** 770 grams
- **Country of Origin:** China

10. WARRANTY AND SUPPORT

Information regarding product warranty and customer support was not provided in the available product data. Please refer to the product packaging or the seller's website for specific warranty details and contact information for support.