

KINGDO B0F8VHV4HT

KINGDO Red Light Therapy Face Mask Instruction Manual

Model: B0F8VHV4HT

			Introduction	Safety Information	Product	
Overview	Setup	Operation	Maintenance	Troubleshooting	Specifications	Warranty & Support

1. INTRODUCTION

Thank you for choosing the KINGDO Red Light Therapy Face Mask. This manual provides essential information for the safe and effective use of your device. Designed for home use, this portable and rechargeable LED face mask offers four distinct light modes to support your daily skincare routine.

The mask features a comfortable and flexible design, made from soft, skin-friendly silicone with a dual-sided jelly texture. Its ergonomic shape ensures full-face coverage, and an adjustable strap provides a secure fit for various face types. Built-in eye shields help reduce light exposure for a relaxing experience.

Enjoy Beauty Care Anytime, Anywhere



Silicone Materials

Soft, Safe and Durable
Complete Facial Coverage



Portable Power

No Worries and
Convenience



Adjustable Strap

Fit Various Facial
Shapes



Protective Eye Pads

Block out Light



Exercising



Travelling



Working



Reading

Image: The KINGDO Red Light Therapy Face Mask worn by a user, showcasing its ergonomic design and red light emission.

2. SAFETY INFORMATION

Important Safety Precautions:

- Always ensure the battery is charged before use.
- Do not connect the controller to both the mask and the charging cable simultaneously, as this will prevent the mask from operating.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is designed for cosmetic use and general well-being.
- Consult a healthcare professional if you have any pre-existing skin conditions, light sensitivity, or are taking medications that may increase light sensitivity.

- Keep out of reach of children.
- Do not use if the device is damaged.

3. PRODUCT OVERVIEW

3.1 Key Features

- **Comfortable & Flexible Design:** Made from soft, skin-friendly silicone with a dual-sided jelly texture. Ergonomic shape for full-face coverage.
- **Adjustable Strap:** Ensures a secure and comfortable fit for various face types.
- **Built-in Eye Shields:** Reduces light exposure to the eyes for a more relaxing experience.
- **Advanced LED Technology:** Features 100 high-intensity LEDs for effective light therapy.
- **4 Light Modes:** Red, Blue, Infrared, and Mixed light options.
- **Built-in 10-Minute Timer:** For convenient and efficient skincare sessions.
- **Portable and Rechargeable:** Allows for use at home or on the go without being tethered to a power outlet.

3.2 Light Modes Explained

The KINGDO LED Face Mask offers four versatile light modes, each designed to address specific skin concerns:

1. **Red Light (620nm):** Helps improve skin smoothness and reduce the appearance of fine lines.
2. **Blue Light (460nm):** Supports clearer-looking skin by targeting blemishes.
3. **Infrared Light (850nm):** Gently warms the skin surface, providing a soothing care experience. Note: This light is invisible to the human eye, but it is functioning normally.
4. **Mixed Light:** Combines all three wavelengths (Red, Blue, and Infrared) for a balanced and comprehensive skincare experience.

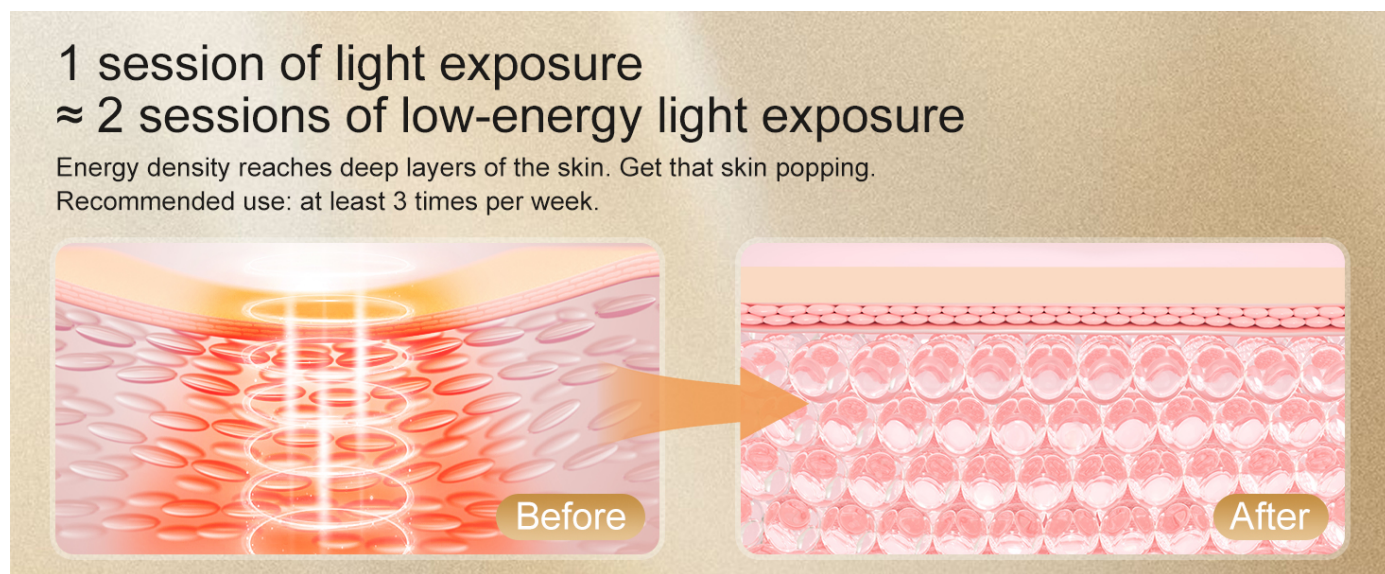


Image: Visual representation of the four distinct light modes and their targeted skin benefits.

4. SETUP

4.1 Charging the Mask

Before first use, or when the battery is low, charge the mask. Disconnect the controller from the mask body and connect

the charging cable directly to the mask. A full charge typically takes about two hours.

4.2 Preparing for Use

1. Ensure the mask's battery is fully charged.
2. Cleanse and thoroughly dry your face before applying the mask.
3. Place the mask on your face, ensuring the eye shields align with your eyes.
4. Adjust the elastic strap at the back of your head to achieve a comfortable and secure fit.



Image: A visual guide demonstrating the preparation and application steps for the LED face mask.

5. OPERATING INSTRUCTIONS

5.1 Turning On and Mode Selection

1. Connect the remote controller to the mask.
2. Press and hold the power button on the remote controller for 1 second to turn on the device.
3. Select your preferred light mode (M1-M4) by pressing the corresponding button on the remote.
4. The mask has a built-in 10-minute timer. It will automatically turn off after 10 minutes in manual mode, or cycle through modes in automatic mode.



Effortless Operation

Equipped with an automatic timer that shuts off after 10 minutes



Boosts Blood Circulation

Works on the skin's surface to improve blood circulation



Relaxation Benefits

Relieve muscle tension, stress, and fatigue



Image: The remote controller for the KINGDO LED Face Mask, showing mode selection and power buttons.

5.2 Recommended Usage

- Wear the mask for 15 minutes per session.
- Use the mask at least 3 times a week for optimal results.
- For blue light therapy, it is recommended not to exceed 2 times a week. For other light modes, 3-4 times a week is generally recommended.

6. MAINTENANCE

6.1 Cleaning

After each use, gently clean the mask with alcohol wipes or a soft, damp cloth. Ensure the mask is completely dry before storing.

6.2 Storage

Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging to protect it from dust and damage.

7. TROUBLESHOOTING

7.1 Controller Not Starting

If the controller does not start, ensure it is only connected to the mask itself, not to the charging cable at the same time. The mask will not operate if both are connected simultaneously.

7.2 M3 Mode (Infrared Light) Appears Dim or Off

The M3 setting uses 850nm near-infrared light, which is invisible to the human eye. Although the LED beads may not visibly light up, they are still functioning normally and emitting therapeutic light.

8. SPECIFICATIONS

Feature	Detail
Model	B0F8VHV4HT
Brand	KINGDO
Light Modes	Red (620nm), Blue (460nm), Infrared (850nm), Mixed Light
LED Count	100 high-intensity LEDs
Timer	Built-in 10-minute automatic shut-off
Power Source	Rechargeable Lithium Ion Battery (1 included)
Charging Time	Approximately 2 hours for full charge
Material	Soft, skin-friendly silicone
Package Dimensions	13.7 x 9.88 x 2.52 inches
Weight	2.34 Pounds
UPC	762849536385

9. WARRANTY & SUPPORT

The KINGDO LED face mask is developed based on standards recognized by trusted U.S. institutions for consumer electronics, designed to deliver safe, consistent performance. For any questions, technical support, or warranty inquiries, please contact KINGDO customer service through the retailer where the product was purchased or visit the official KINGDO website for contact information.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question