

## ANLAN B0F9D6SDY5

# ANLAN V-Face Lifting Facial Massager User Manual

Model: B0F9D6SDY5

## INTRODUCTION

Welcome to the ANLAN V-Face Lifting Facial Massager user manual. This innovative device is designed to help you achieve a more defined facial contour, improve skin elasticity, and address various skin concerns through advanced technologies such as EMS microcurrent, LED phototherapy, adjustable vibration, and constant temperature hot compress. Please read this manual thoroughly before first use to ensure safe and effective operation.

## IMPORTANT SAFETY INFORMATION

**For external use only. Avoid contact with eyes. Keep away from fire or direct heat sources.**

- Do not use if you have sensitive skin or a thin stratum corneum without consulting a dermatologist.
- Keep out of reach of children.
- Do not disassemble or modify the device.
- Ensure the device is fully charged before use to avoid interruptions.
- Clean the device regularly as instructed in the maintenance section.

## PRODUCT COMPONENTS

The ANLAN V-Face Lifting Facial Massager system includes the main facial device and a convenient remote control for easy operation.



Figure 1: The ANLAN V-Face Lifting Facial Massager device and its accompanying remote control.



Figure 2: The remote control allows for easy adjustment of modes and intensity during use.

## KEY FEATURES

### 1. LED Phototherapy (Red & Blue Light)

The device is equipped with 36 bi-color LED lights, offering both red and blue light therapy. Red light promotes collagen production, improves skin tone, and reduces wrinkles, restoring skin elasticity. Blue light helps cleanse the skin, combat acne, inhibit inflammation, and soothe the skin.



Figure 3: Illustration of Red and Blue Light Therapy targeting different skin concerns.

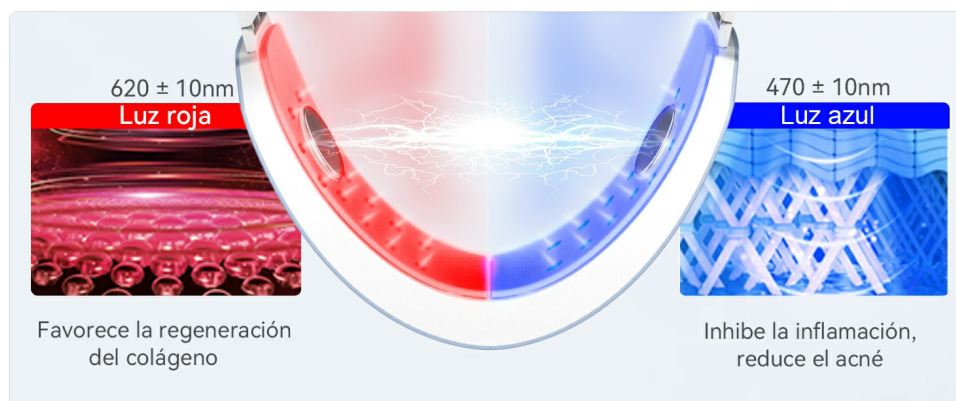


Figure 4: Red light ( $620 \pm 10\text{nm}$ ) promotes collagen regeneration, while blue light ( $470 \pm 10\text{nm}$ ) inhibits inflammation and reduces acne.

## 2. EMS Microcurrent Lifting

This V-shape facial lifting instrument utilizes microcurrent stimulation to firm the skin. Energy is evenly distributed across the face, restoring muscle elasticity and helping to reduce the appearance of a double chin. The device offers various massage techniques including kneading, rolling, acupuncture, and tapping.



Figure 5: EMS microcurrent technology works to tighten and firm the skin, restoring muscle elasticity.

### 3. Adjustable Vibration Massage

The facial lifting instrument provides 9 adjustable vibration levels. Levels 1-3 are suitable for sensitive skin, 4-5 for prominent masticatory muscles and double chin (recommended 3-5 weeks), and 6-8 for wide jaws and tense muscles (recommended 3 months). Vibration helps relax muscles and alleviate fatigue.



Figure 6: The device offers multiple vibration intensity levels to customize your massage experience.

### 4. 42°C Constant Temperature Hot Compress

The device features an intelligent 42°C constant temperature heating design. This hot compress provides a comfortable experience, similar to a warm towel, promoting pore opening and enhancing the absorption rate of skincare products. Heat conduction at 42°C activates collagen molecules in the dermis, helping to alleviate facial swelling, firm contours, smooth expression lines, and protect the skin.



Figure 7: The constant temperature heating function enhances product absorption and promotes circulation.

## 5. Ergonomic Design & Remote Control

The double chin reducer features a 180° ring design that conforms to the chin line, along with adjustable and extendable shoulder straps, making it suitable for various face shapes and providing a comfortable, hands-free experience. It comes with a practical remote control for easy operation, offering 9 vibration levels and 4 adjustable massage intensity levels for a personalized and comprehensive care experience.

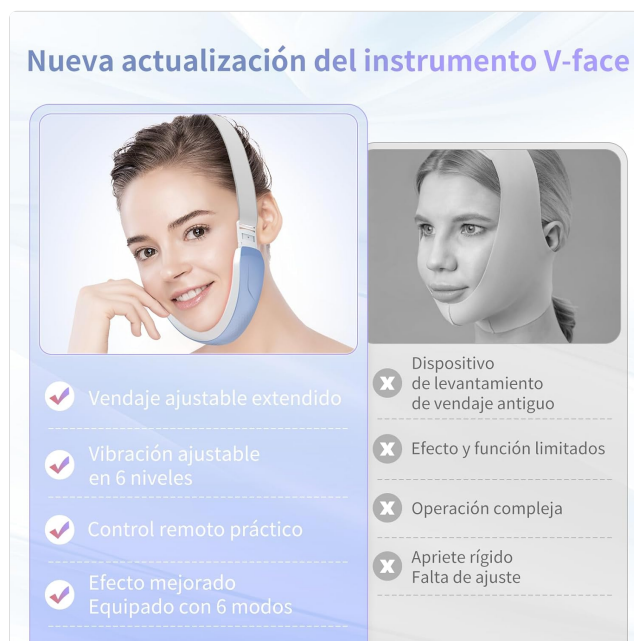


Figure 8: The updated V-face instrument offers an adjustable bandage, 6-level vibration, practical remote control, and improved effects with 6 modes, unlike older, rigid devices.

## SETUP GUIDE

### Charging the Device

Before first use, fully charge the device. Connect the provided USB cable to the charging port on the device



and a suitable USB power adapter. A full charge typically takes approximately 3 hours and provides up to 90 minutes of continuous use, or about six days of 15-minute daily sessions.

## Adjusting the Fit

Place the device around your chin and jawline. Adjust the extendable straps to ensure a snug yet comfortable fit. The 180° ring design is made to adapt to various face shapes.

## OPERATING YOUR DEVICE

### Basic Operation

To power on the device, press and hold the power button on the remote control. Use the remote to select your desired mode and adjust intensity levels. The device will automatically turn off after approximately 10-15 minutes of use.

### Six Facial Care Modes

The device offers six distinct facial care modes, each combining different technologies for specific benefits:



Figure 9: Overview of the six distinct facial care modes available on the device.

1. **Lifting:** Combines EMS and Blue Light.
2. **V-face:** Utilizes EMS only.
3. **Physiotherapy:** Combines EMS and Red Light.
4. **Sterilization:** Combines EMS, Heat, and Red Light.
5. **Acne Treatment:** Uses Red Light only.
6. **Skin Rejuvenation:** Combines EMS, Blue Light, and Vibration.

### Adjusting Intensity

You can adjust the vibration intensity across 9 levels and massage intensity across 4 levels using the remote control. Start with lower levels and gradually increase as your skin adapts and tolerates the sensation.

### Recommended Usage Steps

Follow these four steps for optimal results:



Figure 10: The four essential steps for achieving a V-shaped face with the device.

1. **Step 1: Thoroughly Cleanse Your Skin.** Ensure your face is clean and free of makeup or impurities.
2. **Step 2: Apply the Double Chin Reducer.** Position the device comfortably around your chin and jawline.
3. **Step 3: Adjust Intensity and Treatment Mode.** Use the remote control to select your preferred mode and intensity level.
4. **Step 4: Use Hydrating Products.** After treatment, apply moisturizing products to further nourish your skin.

## CARE AND MAINTENANCE

### Cleaning

- After each use, gently wipe the device with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the device in water.
- Ensure the device is completely dry before storing.

### Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep it in its original packaging or a protective case to prevent damage.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or device not charged.	Charge the device fully. Ensure the charging cable is properly connected.
No vibration or weak sensation.	Intensity level is too low or improper contact with skin.	Increase the vibration intensity using the remote control. Ensure the device makes full contact with your skin.

Problem	Possible Cause	Solution
Device feels too hot.	Hot compress mode is active.	The device is designed for a constant 42°C. If it feels uncomfortably hot, discontinue use and check for skin irritation.
Remote control not responding.	Remote battery low or interference.	Ensure the remote control has working batteries. Point the remote directly at the device.

## PRODUCT SPECIFICATIONS

Feature	Specification
Model Number	B0F9D6SDY5
Dimensions (Package)	15.9 x 15.5 x 5.4 cm
Weight	300 g
Color	White
LED Lights	Red (620 ± 10nm), Blue (470 ± 10nm)
Vibration Levels	9 Levels
Massage Intensity	4 Levels
Hot Compress Temperature	42°C (Constant)
Charging Time	Approx. 3 hours
Battery Life	Approx. 90 minutes (continuous use)
Country of Origin	China

## WARRANTY AND SUPPORT

ANLAN is committed to providing high-quality products. For any questions, concerns, or technical support regarding your ANLAN V-Face Lifting Facial Massager, please contact ANLAN customer service through the retailer's platform or the official ANLAN website. Please retain your proof of purchase for warranty claims.



