

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [seago](#) /
- › [Seago SG2607 Electric Toothbrush User Manual](#)

seago SG-2607

Seago SG2607 Electric Toothbrush User Manual

Your comprehensive guide to operating and maintaining your Seago SG2607 Electric Toothbrush.

PRODUCT OVERVIEW

The Seago SG2607 Electric Toothbrush is designed to provide an effective and comfortable brushing experience. It features multiple cleaning modes, a long-lasting battery, and a waterproof design for convenient use.



Figure 1: Seago SG2607 Electric Toothbrush, showing the main unit, travel case, and eight replacement brush heads, along with its packaging.

This package includes the main toothbrush unit, a travel case, a USB Type-C charging cable, and eight round replacement brush heads, ensuring extended use and convenience.

SETUP AND INITIAL USE

1. Charging the Toothbrush

Before first use, fully charge the toothbrush. The device is equipped with a 700mAh rechargeable battery.

- Connect the provided USB Type-C cable to the charging port at the bottom of the toothbrush.
- Plug the other end of the cable into a compatible USB power adapter (not included).
- The charging indicator light on the toothbrush will illuminate. A full charge typically takes approximately 3 hours.
- Once fully charged, the toothbrush can provide up to 30 days of regular use.

ULTRA-LANGE AKKULAUFZEIT



≤ 3 Stunden
Vollständige Ladung in



700mAh



30 Tage
Nutzung

Figure 2: The toothbrush connected to a USB-C charging cable, illustrating its rechargeable battery and long battery life.

2. Attaching a Brush Head

Select one of the eight provided round brush heads. Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.

- The brush heads feature color-coded rings for easy identification if multiple users share the same handle.
- The soft bristles are designed for effective plaque removal and gum care.

5 Personalisierte Modi

Ein/Aus-Schalter

Modus-Taste

Weiß
Intensität
☺☺☺

Reinigen
Intensität
☺☺☺

Morgens
Intensität
☺☺☺

Abends
Intensität
☺☺☺☺☺

Empfindlich
Intensität
☺☺

Ladeanzeige

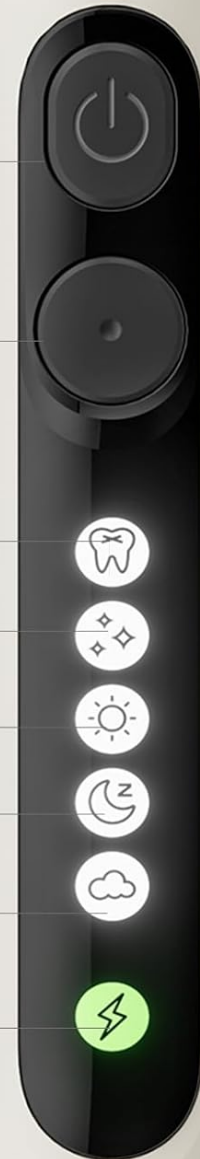


Figure 3: A set of eight replacement brush heads, highlighting their color-coded rings, soft bristles, and advanced design for thorough cleaning.

OPERATING INSTRUCTIONS

Power On/Off and Mode Selection

The Seago SG2607 features a simple two-button interface for power and mode selection.

Wert-Pack

8 Köpfe = 2 Jahre Zähneputzen

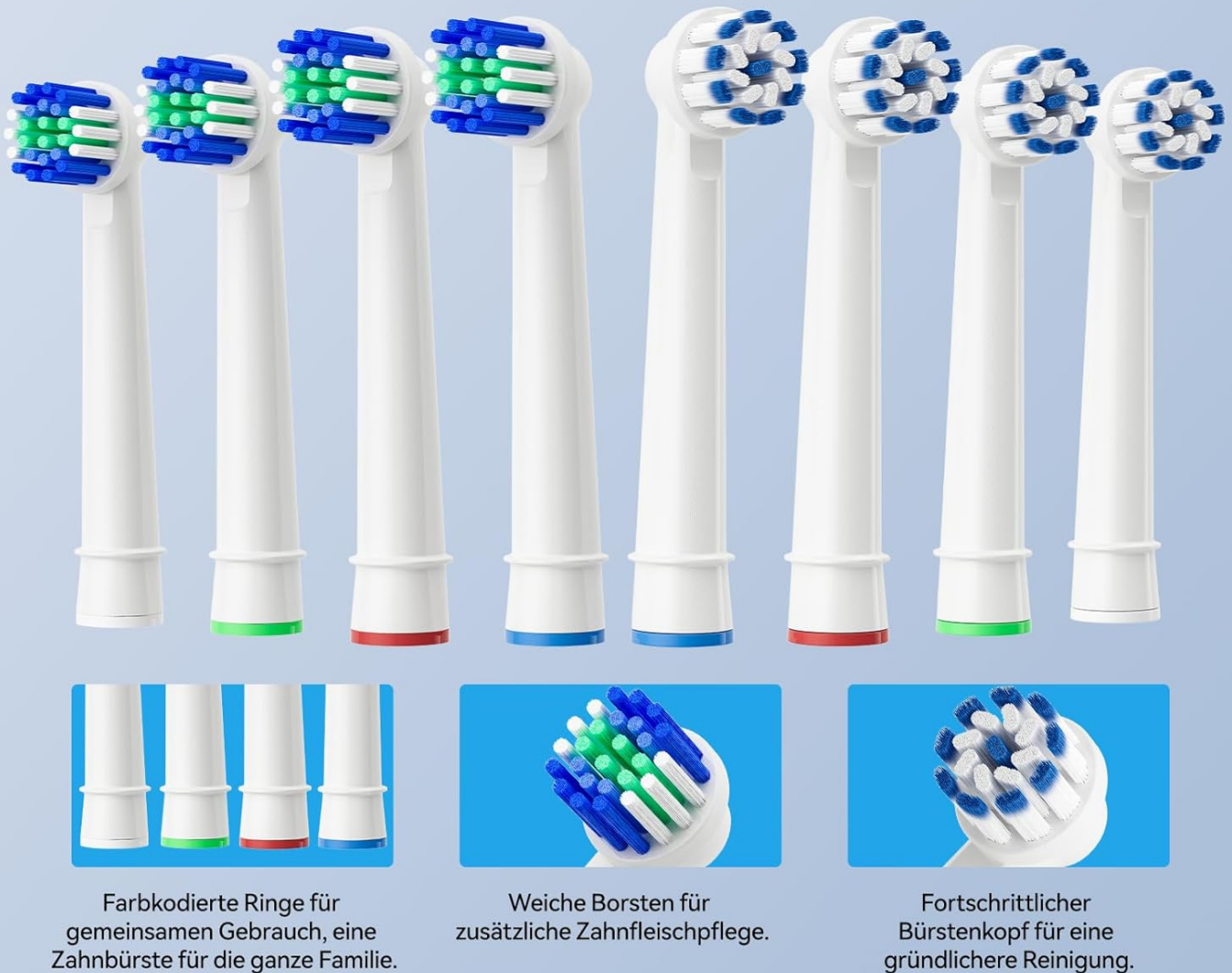


Figure 4: Detailed view of the toothbrush control panel, indicating the power button, mode button, and the five distinct cleaning modes with their respective intensity levels.

- To turn on the toothbrush, press the **Power On/Off button** (top button).
- To cycle through the cleaning modes, press the **Mode button** (bottom button) repeatedly while the toothbrush is on.
- To turn off the toothbrush, press the Power On/Off button again.

Cleaning Modes

The toothbrush offers 5 personalized cleaning modes to suit various oral care needs:

1. **White Mode:** For effective stain removal and whitening.
2. **Clean Mode:** Standard mode for daily comprehensive cleaning.
3. **Morning Mode:** Gentle yet effective cleaning for morning use.
4. **Evening Mode:** Thorough cleaning for evening use.
5. **Sensitive Mode:** Gentle cleaning for sensitive teeth and gums.

Brushing Technique and Smart Features

For optimal results, apply a small amount of toothpaste to the brush head. Place the brush head against your teeth at a slight angle towards the gumline. Move the brush slowly from tooth to tooth, allowing the rotating bristles to clean each surface.



Figure 5: Illustration of the 360° rotating cleaning power, demonstrating effective plaque removal and comprehensive tooth coverage.

- **2-Minute Timer:** The toothbrush features a built-in 2-minute timer, which is the dentist-recommended brushing time. The brush will automatically pause briefly every 30 seconds to signal you to move to a different quadrant of your mouth.
- **Quiet Design:** Engineered for a low-noise operation, enhancing the brushing experience.



Figure 6: An internal view highlighting the advanced technology, including the 2-minute timer and quiet operational design.

MAINTENANCE AND CARE

Cleaning the Toothbrush

- After each use, rinse the brush head thoroughly under running water while the toothbrush is briefly turned on to remove any remaining toothpaste or debris.
- Turn off the toothbrush and remove the brush head. Rinse both the brush head and the handle separately under running water.
- Wipe the handle with a soft cloth to dry it.
- The Seago SG2607 is **IPX7 waterproof**, meaning it can be safely rinsed under running water and used in the shower or bath. Do not submerge the entire unit for extended periods.



Figure 7: The toothbrush submerged in water, illustrating its IPX7 waterproof rating, making it suitable for use in the shower or bath.

Brush Head Replacement

It is recommended to replace your brush head every 3 months, or sooner if the bristles become frayed or worn. Regular replacement ensures optimal cleaning performance.

Storage and Travel

Store the toothbrush and brush heads in a clean, dry place. For travel, use the included portable travel case to protect the toothbrush and brush heads.

Reisekonform

Mit tragbarem Etui für eine einfache und komfortable Reise.



Figure 8: The toothbrush and brush heads neatly stored within the portable travel case, designed for easy and comfortable travel.

TROUBLESHOOTING

If you encounter any issues with your Seago SG2607 Electric Toothbrush, please refer to the following common solutions:

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Battery is depleted.	Charge the toothbrush fully (approx. 3 hours).
Weak vibration or reduced performance.	Low battery or worn brush head.	Recharge the toothbrush. Replace the brush head if it is old or frayed.

Problem	Possible Cause	Solution
Brush head does not stay on.	Brush head not properly attached or damaged.	Ensure the brush head is pushed down firmly until it clicks. Inspect for damage and replace if necessary.
Charging indicator not lighting up.	Charging cable or adapter issue.	Ensure the USB cable is securely connected to both the toothbrush and the power source. Try a different USB adapter or port.

SPECIFICATIONS

Feature	Detail
Model Number	SG-2607
Brand	Seago
Power Source	Battery Powered (Rechargeable)
Battery Capacity	700mAh
Full Charge Time	Approx. 3 hours
Usage Time per Charge	Up to 30 days
Cleaning Modes	5 (White, Clean, Morning, Evening, Sensitive)
Brush Head Type	Round, Soft Bristles
Included Brush Heads	8
Waterproof Rating	IPX7
Color	Beige
Dimensions (Package)	23.9 x 15.2 x 4.2 cm
Weight (Package)	390 g
Special Features	2-Minute Smart Timer, 30-Second Interval Reminder, Quiet Design, Travel Case Included

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official Seago website. Keep your proof of purchase for any warranty claims.

For further assistance, please contact customer service.