

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Zima Dental](#) /

› [Zima Dental Pod PRO Ultrasonic Cleaner Instruction Manual](#)

Zima Dental PRO1

Zima Dental Pod PRO Ultrasonic Cleaner Instruction Manual

1. INTRODUCTION

The Zima Dental Pod PRO is an advanced ultrasonic cleaning device designed for comprehensive cleaning of various oral appliances. This manual provides essential information for the safe and effective operation, maintenance, and troubleshooting of your device. Please read these instructions thoroughly before first use and retain them for future reference.



Image 1: The Zima Dental Pod PRO, its retail packaging, and a clear dental aligner, showcasing the product's design and purpose.

2. PRODUCT OVERVIEW AND FEATURES

The Zima Dental Pod PRO incorporates patented British engineering to deliver enhanced cleaning performance for oral appliances. Key features include:

- **Increased Cleaning Power:** Utilizes next-generation technology for effective cleaning.
- **Award-Winning Innovation:** Features intelligent, patented, dynamic 42k-47k Hz ultrasonic technology, recognized as the 2024 British Dental Association Product of the Year. This system automatically adjusts to create optimal ultrasonic conditions for any oral appliance.
- **Removable Water Tank:** A modular design facilitates easier water refilling and cleaning, offering flexibility in placement.
- **Advanced Cleaning Mode:** An additional 'Max Clean' setting combines ultrasonic programs to address stubborn stains and plaque buildup.

Patented Industry-Leading Cleaning Technology



Max Clean Mode
for stubborn build-up



Modular Design
for easy refilling & cleaning



Dynamic Frequency
for optimal cleaning power

Image 2: An illustration highlighting key features of the Zima Dental Pod PRO: Max Clean Mode for stubborn buildup, Modular Design for easy refilling and cleaning, and Dynamic Frequency for optimal cleaning power.

3. SETUP

Before using your Zima Dental Pod PRO, follow these steps:

1. **Unpack:** Carefully remove all components from the packaging.
2. **Placement:** Place the unit on a stable, flat, and dry surface near a power outlet.
3. **Connect Power:** Plug the power adapter into the unit and then into a suitable electrical outlet.



Image 3: A close-up image showing hands filling the removable water tank of the Zima Dental Pod PRO directly from a faucet, illustrating the ease of refilling.

4. OPERATING INSTRUCTIONS

The Zima Dental Pod PRO offers two cleaning modes:

1. **Fill with Water:** Remove the water tank from the base. Fill the tank with clean tap water up to the indicated fill line. Do not overfill.
2. **Place Appliance:** Carefully place your dental appliance (aligner, retainer, nightguard, mouthguard, denture) into the water tank. Ensure the appliance is fully submerged.
3. **Select Cleaning Mode:**
 - **Standard Clean:** Tap the power button once for a regular cleaning cycle.
 - **Max Clean:** Tap the power button twice for the advanced Max Clean setting, designed for tougher stains and buildup.
4. **Start Cycle:** The device will begin the selected cleaning cycle automatically.
5. **Cycle Completion:** The device will automatically shut off once the cleaning cycle is complete.
6. **Retrieve Appliance:** Carefully remove your clean dental appliance from the tank. Rinse it under

running water before use.

- Optional:** For enhanced cleaning power, a Zima Dental tablet can be added to the water before starting a cycle.

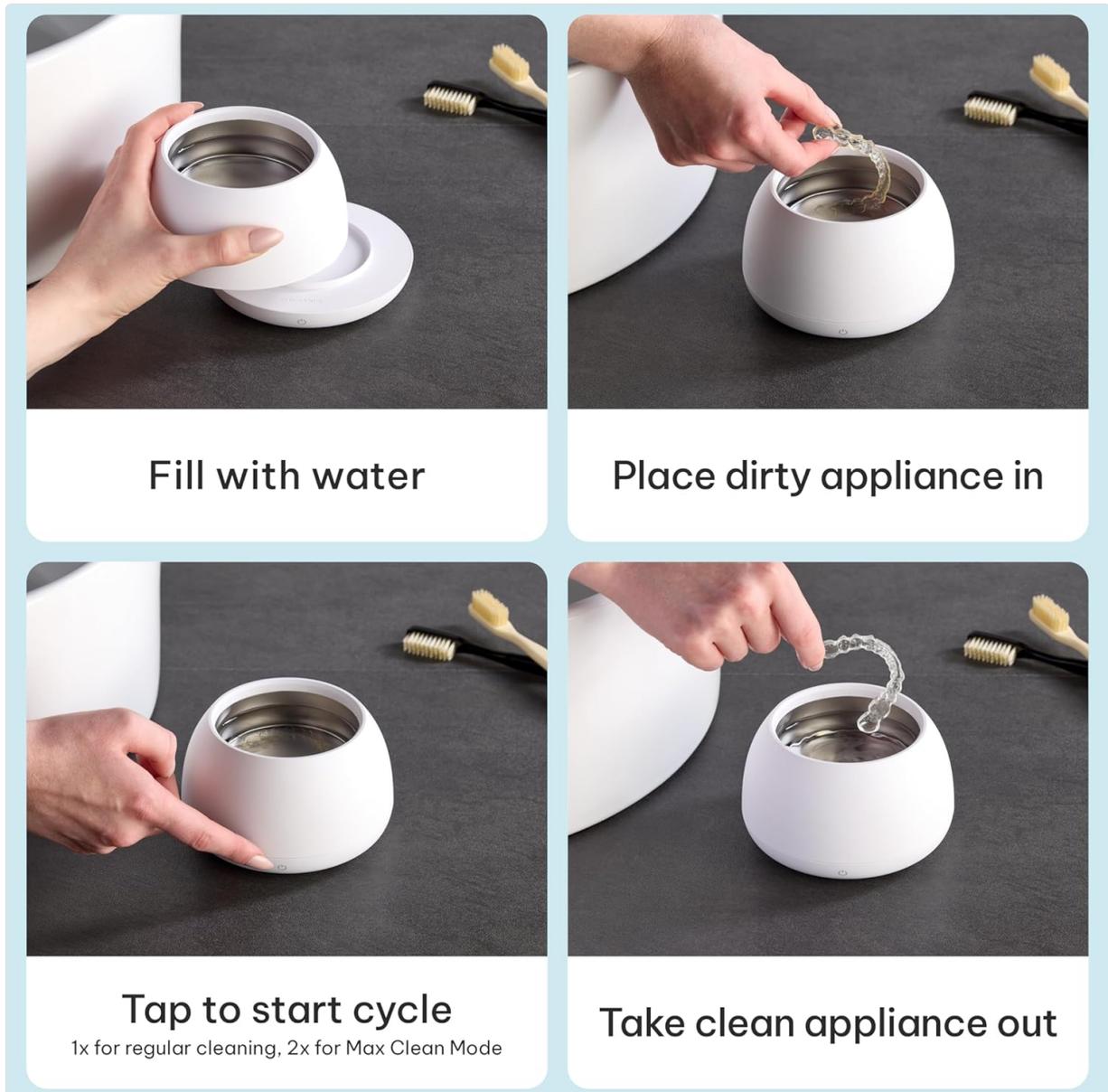


Image 4: A visual guide demonstrating the four simple steps to use the Zima Dental Pod PRO: 1. Fill the pod with water. 2. Place your dental appliance inside. 3. Tap the power button to start a cleaning cycle. 4. Remove your clean appliance after the cycle completes.

5. MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your Zima Dental Pod PRO:

- **Empty and Rinse:** After each use, empty the water from the tank and rinse it thoroughly with clean water.
- **Wipe Down:** Wipe the exterior of the unit with a soft, damp cloth.
- **Do Not Submerge:** Never submerge the main unit (base with electronics) in water. Only the removable tank is designed to hold water.
- **Storage:** Store the unit in a cool, dry place when not in use.

6. TROUBLESHOOTING

If you encounter issues with your Zima Dental Pod PRO, refer to the following common solutions:

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in correctly; Power outlet issue.	Ensure power adapter is securely connected to the unit and a working electrical outlet.
Cleaning is not effective.	Insufficient water; Appliance not fully submerged; Heavy buildup.	Ensure water level is adequate and appliance is submerged. Use 'Max Clean' mode or add a Zima Dental tablet for stubborn buildup.
Unit stops mid-cycle.	Overheating; Power interruption.	Unplug the unit, allow it to cool for 10-15 minutes, then try again. Check power connection.

If the problem persists after attempting these solutions, please contact Zima Dental customer support.

7. SPECIFICATIONS

Detailed technical specifications for the Zima Dental Pod PRO:

- **Model Number:** PRO1
- **Ultrasonic Frequency:** Dynamic 42,000 - 47,000 Hz
- **Body Weight:** 525g (approximately 1.16 lbs)
- **Dimensions:** 123mm x 123mm x 90mm (approximately 4.84 x 4.84 x 3.54 inches)
- **Tank Material:** Medical-grade stainless steel and silver alloy
- **Tank Volume:** 190ml
- **Power Source:** Corded Electric
- **Power Supply:** DC 12V/2A, 239V/50Hz
- **Special Feature:** Multiple Operation Modes (Standard Clean, Max Clean)
- **Product Benefits:** Removes plaque, grime, and biofilm; Extends appliance lifetime.



Image 5: An infographic detailing the technical specifications of the Zima Dental Pod PRO, including its dynamic

8. WARRANTY AND SUPPORT

The Zima Dental Pod PRO comes with a full **12-month warranty** from the date of purchase. For any warranty claims, technical assistance, or product inquiries, please contact Zima Dental customer support via the Amazon Buyer-Seller Messaging platform. Our team aims to respond within 24 hours.

9. SAFETY INFORMATION

Please observe the following safety precautions:

- Keep the device out of reach of children.
- Do not operate the device if the power cord or plug is damaged.
- Do not disassemble or attempt to repair the unit yourself. Refer all servicing to qualified personnel.
- Ensure the unit is unplugged before cleaning or when not in use for extended periods.
- Use only the provided power adapter.
- Avoid placing the unit near heat sources or in direct sunlight.
- This product is intended for cleaning oral appliances only. Do not use it for other purposes.