

Manuals.plus /

› amilibeaity /

› amilibeaity 9-in-1 Multifunction Facial & Body Care Device User Manual

amilibeauty 9-in-1 Multifunction Facial & Body Care Device

amilibeauty 9-in-1 Multifunction Facial & Body Care Device User Manual

Model: 9-in-1 Multifunction Facial & Body Care Device

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective operation of your amilibeaity 9-in-1 Multifunction Facial & Body Care Device. This device is designed for home use, offering multiple modes and adjustable intensity settings for various skin care needs.

Please read this manual thoroughly before first use and retain it for future reference.



Figure 1: The amilibeaauty 9-in-1 Multifunction Facial & Body Care Device with various handles and accessories.

2. SAFETY INFORMATION

To ensure safe operation and prevent injury or damage to the device, please observe the following safety precautions:

- Read all instructions before using the device.
- Do not use the device on broken, irritated, or infected skin.
- Consult a physician before use if you have any medical conditions, are pregnant, or have electronic implants.
- Keep the device out of reach of children.
- Do not immerse the main unit in water or any other liquid.
- Use only the provided power adapter and accessories.
- Unplug the device when not in use and before cleaning.
- If the device appears damaged or malfunctions, discontinue use immediately and contact customer support.

- Always use the device with appropriate gels, creams, or milk as instructed for specific functions.

3. PRODUCT COMPONENTS

The amilibeauty 9-in-1 Multifunction Facial & Body Care Device includes a main unit and various interchangeable handles and accessories. Familiarize yourself with each component:



Figure 2: Overview of the main unit and various treatment handles.

3.1 Main Unit

The main unit houses the control panel, power supply, and connection ports for all handles. It features an intelligent operation system with multiple modes and intensity adjustments.

PRODUCT PARAMETER

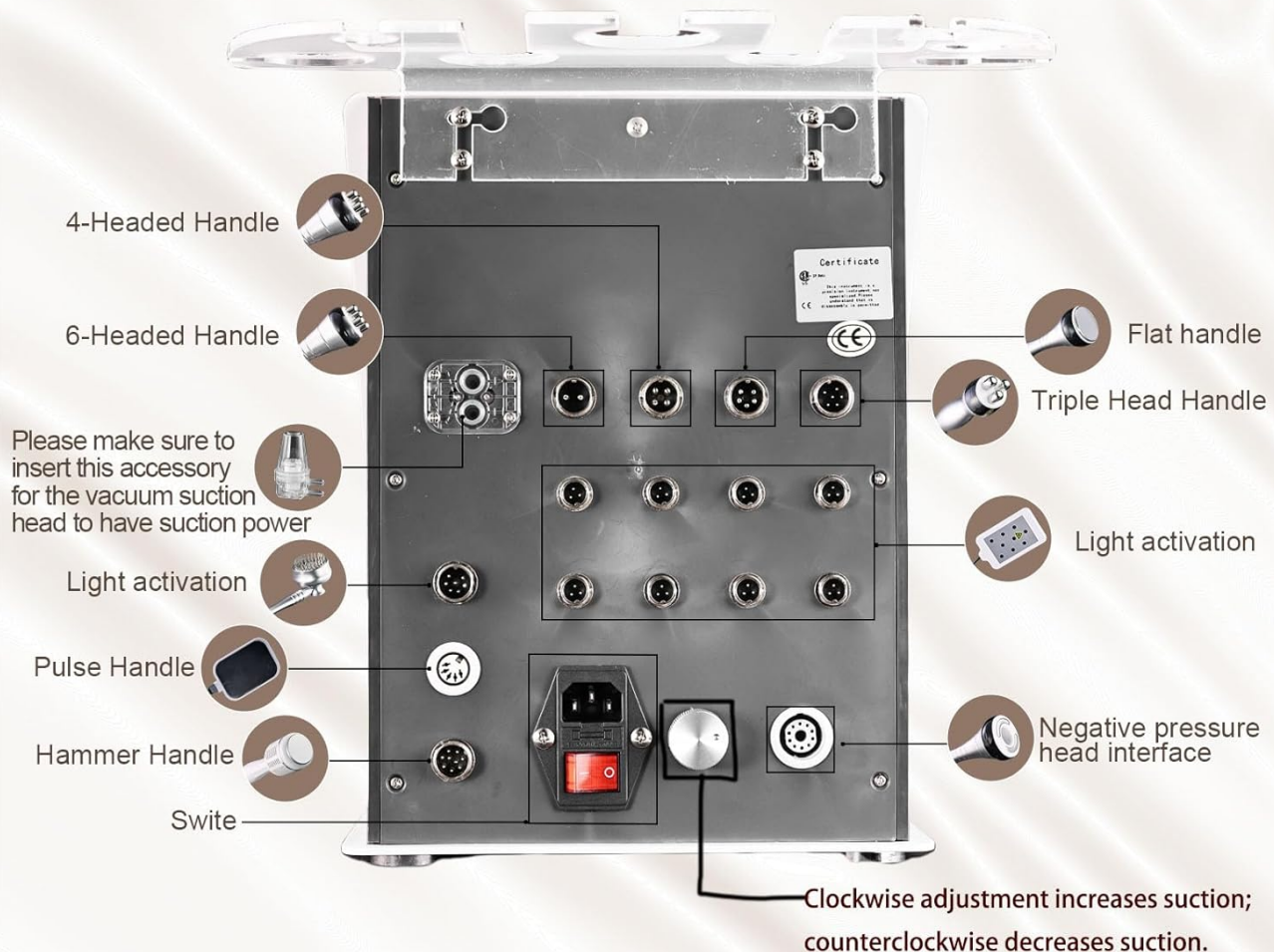


Figure 3: Rear view of the main unit showing connection ports and controls. Note the suction adjustment knob: clockwise increases suction, counter-clockwise decreases suction.

3.2 Treatment Handles and Accessories

- **4-Headed Handle:** Typically used for specific facial areas.
- **6-Headed Handle:** Designed for broader facial or smaller body areas.
- **Flat Handle:** General purpose handle.
- **Triple Head Handle:** Often used for targeted treatments.
- **Pulse Handle:** For specific pulse-based treatments.
- **Hammer Handle:** Used for specific massage or treatment applications.
- **Light Activation Pads:** For light-based therapies.
- **Negative Pressure Head Interface:** Connects to vacuum suction accessories. Ensure the correct accessory is inserted for proper suction.
- **Other Accessories:** Includes various pads and connectors for different applications (e.g., for waist, whole body, face/eyes, face/neck, neck/arm/back).

4. SETUP

Follow these steps to set up your device for first use:

1. **Unpack:** Carefully remove all components from the packaging.
2. **Placement:** Place the main unit on a stable, flat surface.
3. **Connect Power:** Plug the power cord into the main unit's power input, then plug the other end into a suitable electrical outlet.
4. **Attach Handles:** Select the desired treatment handle and securely connect it to the corresponding port on the main unit. Ensure all connections are firm. For vacuum suction, ensure the correct accessory is inserted into the negative pressure head interface.
5. **Power On:** Flip the main power switch (labeled "Swite" on the diagram) to the "ON" position. The display screen should illuminate.

5. OPERATING INSTRUCTIONS

The device offers intelligent operation with multiple modes and adjustable intensity. Always ensure your skin is clean and prepared before treatment.



Figure 4: General usage steps: 1. Cleanse skin, 2. Apply cream/gel, 3. Use facial handle, 4. Use body handle.

5.1 General Usage Steps

1. **Prepare Skin:** Cleanse the treatment area thoroughly.
2. **Apply Medium:** Apply a suitable amount of milk, cream, or gel to the skin area to be treated. This is crucial for effective and comfortable operation.
3. **Select Handle:** Choose the appropriate handle for the desired treatment area (e.g., 4-headed for face, 6-headed for body).
4. **Select Mode & Intensity:** On the main unit's display, select the desired mode and adjust the intensity level according to your comfort and skin type. Start with lower intensity and gradually increase if needed.
5. **Begin Treatment:** Gently move the handle over the prepared skin area. Follow specific instructions for each handle type as detailed in the device's on-screen guide or a more detailed product-specific manual if available.
6. **Monitor Progress:** The device may offer time progress tracking. Adhere to recommended treatment durations.
7. **Clean After Use:** Refer to the "Maintenance" section for cleaning instructions.

5.2 Application Areas

The device is suitable for various body parts:



Figure 5: Illustrates the various body areas where the device can be applied, including face, arms, waist, back, buttocks, and legs.

- **Face:** Use appropriate handles for facial contours and delicate areas.
- **Arms:** For upper and lower arm areas.
- **Waist & Abdomen:** For targeted treatments around the midsection.
- **Back:** For broader back areas.
- **Buttocks:** For specific treatments in this region.
- **Legs:** For thigh and calf areas.

6. MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your device.

6.1 Cleaning

- **After Each Use:** Disconnect the device from power. Wipe all handles and the main unit's surface with a soft, damp cloth. For stubborn residue, use a mild, non-abrasive cleaner.
- **Do Not Immerse:** Never immerse the main unit or any electrical components in water.
- **Dry Thoroughly:** Ensure all parts are completely dry before storing or next use.

6.2 Storage

- Store the device and its accessories in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the device in its original packaging or a protective case to prevent damage.

7. TROUBLESHOOTING

If you encounter any issues with your device, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not power on.	Power cord not properly connected; Power switch is off; No power from outlet.	Check power cord connections; Ensure power switch is ON; Test outlet with another device.
Handle not functioning.	Handle not securely connected; Incorrect mode selected; Intensity set too low.	Ensure handle is firmly connected; Select correct mode on display; Increase intensity gradually.
No suction from vacuum head.	Accessory not properly inserted into negative pressure head interface; Suction adjustment knob set too low.	Verify accessory connection; Turn suction adjustment knob clockwise to increase suction.
Skin irritation after use.	Intensity too high; Insufficient application of cream/gel; Allergic reaction.	Reduce intensity; Apply more cream/gel; Discontinue use and consult a physician if irritation persists.

If the problem persists after attempting these solutions, please contact customer support.

8. SPECIFICATIONS

- **Model:** 9-in-1 Multifunction Facial & Body Care Device
- **Brand:** amilibeaauty
- **ASIN:** B0FD9XCJR7
- **GTIN/UPC:** 607603218773
- **Package Dimensions:** 16.38 x 15.94 x 12.24 inches
- **Item Weight:** 18.67 Pounds
- **Manufacturer:** amilibeaauty

9. WARRANTY AND SUPPORT

amilibeauty provides the following warranty for this product:

- **Host Unit Warranty:** One (1) year from the date of purchase.
- **Accessories Warranty:** Thirty (30) days from the date of purchase.

If you have any questions about your product, require technical assistance, or need to claim warranty service, please contact amilibeaauty customer support. Your satisfaction is our priority.

Contact information for support can typically be found on the product packaging or the official amilibeaauty website.