

[Manuals.plus](#) /

> [Kaoudt](#) /

> Kaoudt 3D Knee Massager with Heat and Red Light Therapy User Manual (Model 701551400836)

Kaoudt 701551400836

Kaoudt 3D Knee Massager with Heat and Red Light Therapy

MODEL: 701551400836

INTRODUCTION

The Kaoudt 3D Knee Massager is designed to provide targeted relief and support for joint discomfort. This wireless device integrates three unique vibration modes, three heat intensity settings, and five red light intensity settings to stimulate muscles and deliver deep tissue warmth. Its ergonomic design ensures a comfortable and secure fit, making it suitable for various body parts beyond the knee. The built-in 5000mAh rechargeable battery offers extended wireless use, allowing for convenient application at home, in the office, or outdoors. A 30-minute auto shut-off feature ensures safe and worry-free operation.



The Kaoudt 3D Knee Massager offers multiple modes for comprehensive joint care.

SAFETY INFORMATION

Please read and understand all safety instructions before using this device. Failure to do so may result in injury or damage to the product.

- This device is intended for home wellness use only and is not a substitute for professional medical advice, diagnosis, or treatment.
- Do not use on open wounds, swollen or inflamed areas, or directly over the eyes.
- Individuals who are pregnant, have a pacemaker, or suffer from serious medical conditions should consult a healthcare professional before use.
- Discontinue use immediately if you experience discomfort or irritation.

PRODUCT COMPONENTS

- Main Massager Unit with integrated control panel
- Adjustable Straps with Velcro fasteners
- USB-C Charging Cable
- User Manual

Three Modes - Let your joints get triple the care!

3D Massage Plus & Heating & Red Light Therapy



Reduce Inflammation



Enhance Blood Circulation



Relieve muscle strain



Boost Cellular Energy



Speed Up Recovery



Deep Soothe



Pulse Circulation




Intensive Relief

The main massager unit securely attaches to the knee with adjustable straps.

SETUP

- 1. Charge the Device:** Before first use, fully charge the massager using the provided USB-C cable. Connect the cable to the Type-C charging port on the device and a suitable power source. The battery indicator on the control panel will show charging status.
- 2. Position the Device:** Securely fasten the massager around your knee or other desired treatment area using the adjustable straps. Ensure the massager and red light panel are properly aligned with the joint for optimal contact. The ergonomic design allows for comfortable wear.
- 3. Adjust for Comfort:** Pull the Velcro straps as tight or loose as needed to achieve a snug yet


comfortable fit. The breathable medical-grade neoprene material is designed to minimize skin irritation.

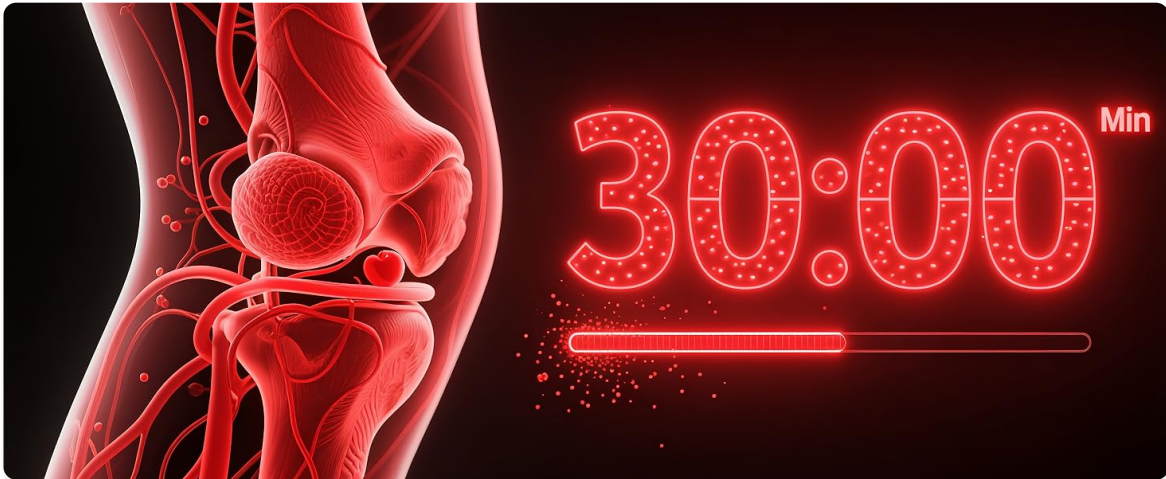

The wireless design and 5000mAh battery allow for use in various settings.

OPERATING INSTRUCTIONS

The device features an intelligent touchscreen for easy control of all functions.

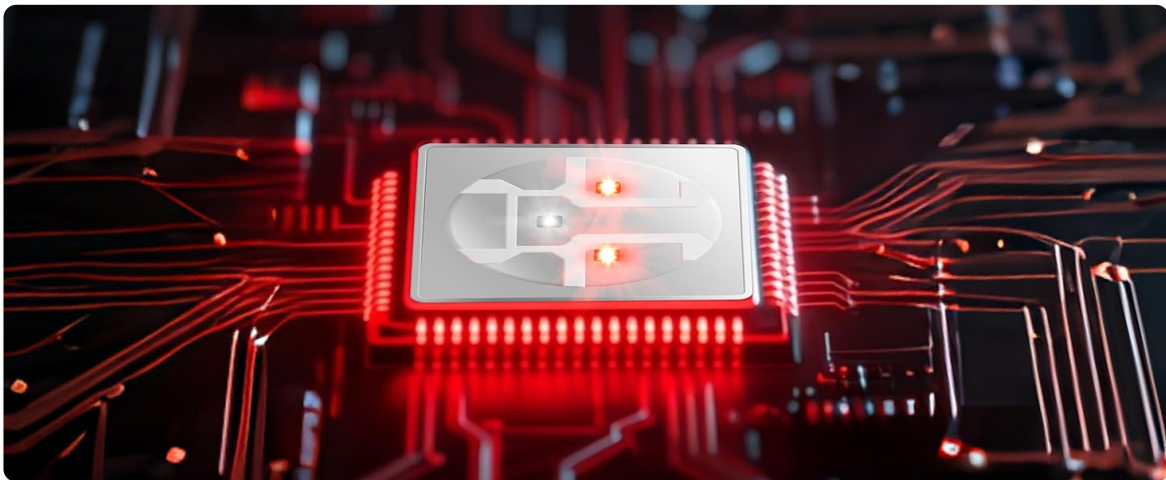

The intelligent touchscreen allows for easy customization of your therapy session.

1. **Power On/Off:** Press and hold the power button (usually located near the charging port) for 2 seconds to turn the device ON. To turn OFF, long press the power button again. The device will automatically power off after 30 minutes of continuous use.
2. **Select Vibration Mode:** Tap the Vibration button (often represented by wavy lines) to cycle through the three unique vibration modes (M1, M2, M3). Each mode offers a different massage pattern. Tap again to turn off vibration.
3. **Adjust Red Light Intensity:** Tap the Red Light button (often represented by a sun or light bulb icon) to cycle through five red light intensity levels (P1-P5). Tap again to turn off red light therapy.
4. **Select Heat Setting:** Tap the Heat button (often represented by a flame icon) to cycle through three heat intensity levels (L1, L2, L3). Tap again to turn off heat therapy.
5. **Activate Lock Mode (Anti-Mistouch):** After setting your desired modes and intensities, briefly press the power button to activate the anti-mistouch lock mode. This prevents accidental changes

during use. Briefly press the power button again to deactivate.



Visual guide for operating the Kaoudt 3D Knee Massager.

Versatile Application

While primarily designed for the knee, the adjustable ergonomic embrace allows the massager to be used on other body parts such as the elbow, shoulder, thigh, or calf, providing targeted relief wherever needed.



The massager's design allows for use on various body areas.

MAINTENANCE

- **Cleaning:** Disconnect the device from power before cleaning. Wipe the surface of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures. Ensure the device is fully charged before long-term storage to preserve battery life.
- **Battery Care:** To prolong battery life, avoid completely draining the battery frequently. Recharge the device when the battery indicator is low.

TROUBLESHOOTING

- **Device Not Turning On:** Ensure the device is fully charged. Press and hold the power button for at least 2 seconds.
- **No Vibration/Heat/Red Light:** Check if the specific function is activated on the control panel. Ensure the device is not in lock mode.
- **Weak Vibration/Heat:** The battery might be low. Recharge the device.
- **Device Shuts Off Unexpectedly:** The 30-minute auto shut-off feature may have activated. If it shuts off before 30 minutes, ensure the battery is sufficiently charged.

SPECIFICATIONS

Model Number	701551400836
Product Dimensions	10.2 x 8.6 x 1 inches
Item Weight	11.2 ounces
Battery Type	1 Lithium Ion battery (included)
Battery Capacity	5000mAh
Manufacturer	Kaoudt

WARRANTY AND SUPPORT

For any questions, concerns, or support needs regarding your Kaoudt 3D Knee Massager, please contact the manufacturer directly. Refer to the product packaging or the seller's information for specific warranty details and contact methods.



Contact customer support for assistance.

