



Manuals.plus /

› SDKWDH /

› Red Light Therapy Face Mask User Manual

SDKWDH White A1

SDKWDH 7-Color LED Face Mask User Manual

Model: White A1 | Brand: SDKWDH

1. INTRODUCTION AND OVERVIEW

Thank you for choosing the SDKWDH 7-Color LED Face Mask. This device is designed to provide at-home light therapy for various skin concerns, offering a convenient and non-invasive skincare method. It features 7 different LED colors, each targeting specific skin needs, and is made with comfortable silicone materials for extended use.

This manual provides essential information on how to safely and effectively use your LED face mask, including setup, operation, maintenance, and troubleshooting.

SDKWDH

7-COLOR LED LIGHTS

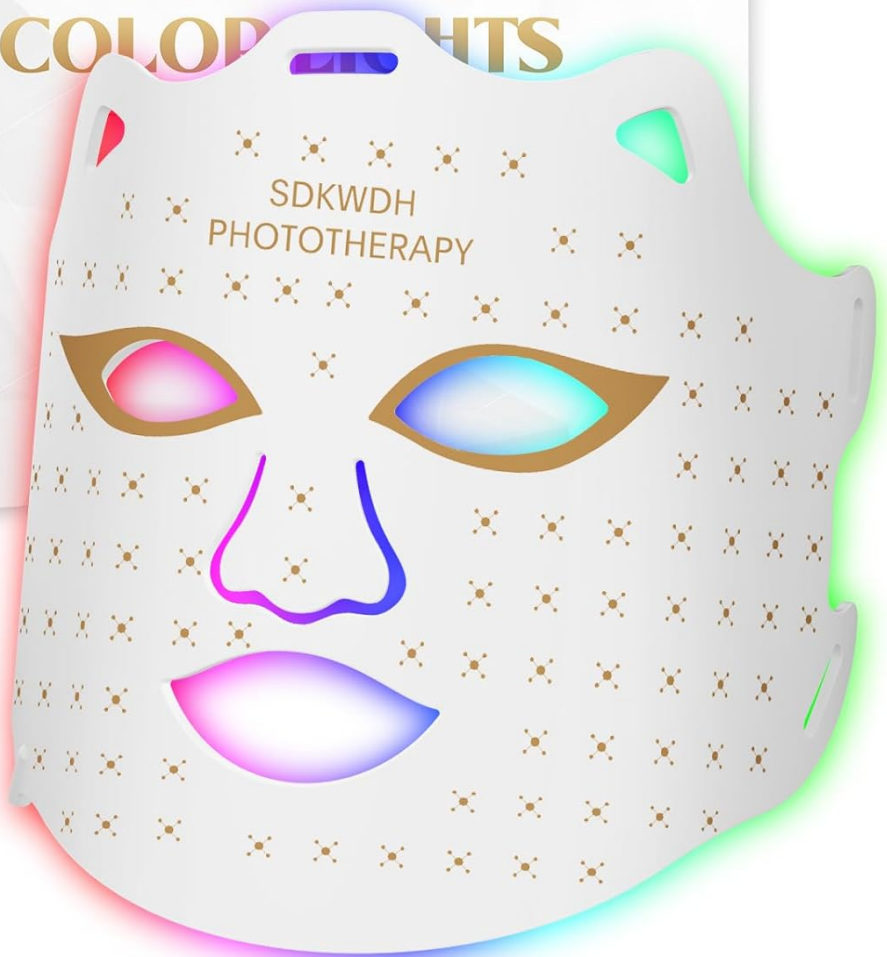


Image: The SDKWDH 7-Color LED Face Mask, showcasing its design and accompanying product box.

2. SAFETY INFORMATION

While no specific safety warnings were provided by the manufacturer, it is important to follow general safety guidelines when using any electronic skincare device:

- Do not use the device if you have any open wounds, skin infections, or severe skin conditions. Consult a dermatologist if you have concerns.
- Avoid direct eye contact with the LED lights. The mask is designed with soft eye pads for protection.
- Keep the device dry and do not immerse it in water.
- Store the device out of reach of children.
- If you experience any discomfort or adverse reactions, discontinue use immediately and consult a healthcare professional.

- Use only the provided charging cable and power adapter.

3. PRODUCT COMPONENTS

The SDKWDH 7-Color LED Face Mask package typically includes the following components:

- SDKWDH LED Face Mask
- Remote Control
- USB Charging Cable
- User Manual (this document)

7 COLORS LED LIGHTS



Red Light 630nm



Blue Light 415nm



Green Light 534nm



Yellow Light 590nm



Cyan Light 470nm



Purple Light 1050nm



White Light 1410nm

Image: Detailed diagram illustrating the various components of the LED face mask, including LED lamp beads, bandage connection lock, medical grade pads, VPS material, and charging port.

4. SETUP

Before your first use, please ensure the remote control is adequately charged.

1. **Charge the Remote Control:** Connect the USB charging cable to the remote control's electrical port (DC5V-1A~2A plug) and plug it into a suitable power source. The remote control has an internal battery that needs

charging.

2. **Connect Remote to Mask:** Once the remote control is charged, attach it to the LED face mask. The remote control must be connected to the mask for the mask to function.
3. **Adjust Strap:** Adjust the adjustable strap on the mask to ensure a snug yet comfortable fit on your face.

EASY TO USE, CONVENIENT AND COMFORTABLE

Use after cleansing, once a day, 15~30min each time



Wash your face, and apply
a little serum



Put on the mask and
choose the color



Lie flat with your eyes
closed and enjoy 15 minutes
of light therapy



Apply lotion to celebrate
the resolution

Image: A detailed view of the remote control, highlighting its various functions such as Mode Display, Time Display, Power Display, Power Switch Key, Color Cycle Button, Brightness Adjustment, Time Setting, and Electrical Port.

5. OPERATING INSTRUCTIONS

Follow these steps for an effective light therapy session:

1. **Prepare Your Skin:** Wash your face thoroughly and apply a small amount of your preferred serum.
2. **Put on the Mask:** Carefully place the LED face mask on your face, ensuring it fits comfortably.
3. **Select Color and Settings:** Use the remote control to switch on the mask and select your desired LED color, treatment time (15-30 minutes per session), and light intensity.
4. **Relax During Treatment:** Lie flat with your eyes closed and enjoy the light therapy session.
5. **Post-Treatment Care:** After the session, remove the mask and apply lotion or essence to accelerate absorption and enhance results.

Understanding the 7 LED Colors:

- **Red Light (630nm):** Promotes collagen regeneration, improves skin firmness, and reduces fine lines.
- **Blue Light (415nm):** Targets acne-causing bacteria, reduces inflammation, and calms the skin.
- **Green Light (534nm):** Balances skin tone, reduces hyperpigmentation, and calms sensitive skin.
- **Yellow Light (590nm):** Improves lymphatic drainage, reduces redness, and enhances skin radiance.
- **Cyan Light (470nm):** Soothes and calms irritated skin, reduces inflammation, and promotes healing.
- **Purple Light (1050nm):** A combination of red and blue light, beneficial for anti-aging and anti-acne effects.
- **White Light (1410nm):** Penetrates deeply to accelerate metabolism, improve fine lines, and tighten skin.

LED LIGHT THERAPY TECHNOLOGY

Take care of your skin and keep it young



Image: A visual guide demonstrating the four simple steps for using the LED face mask: washing face, putting on mask, enjoying light therapy, and applying lotion.

7 LED COLORS BEAUTY MASK

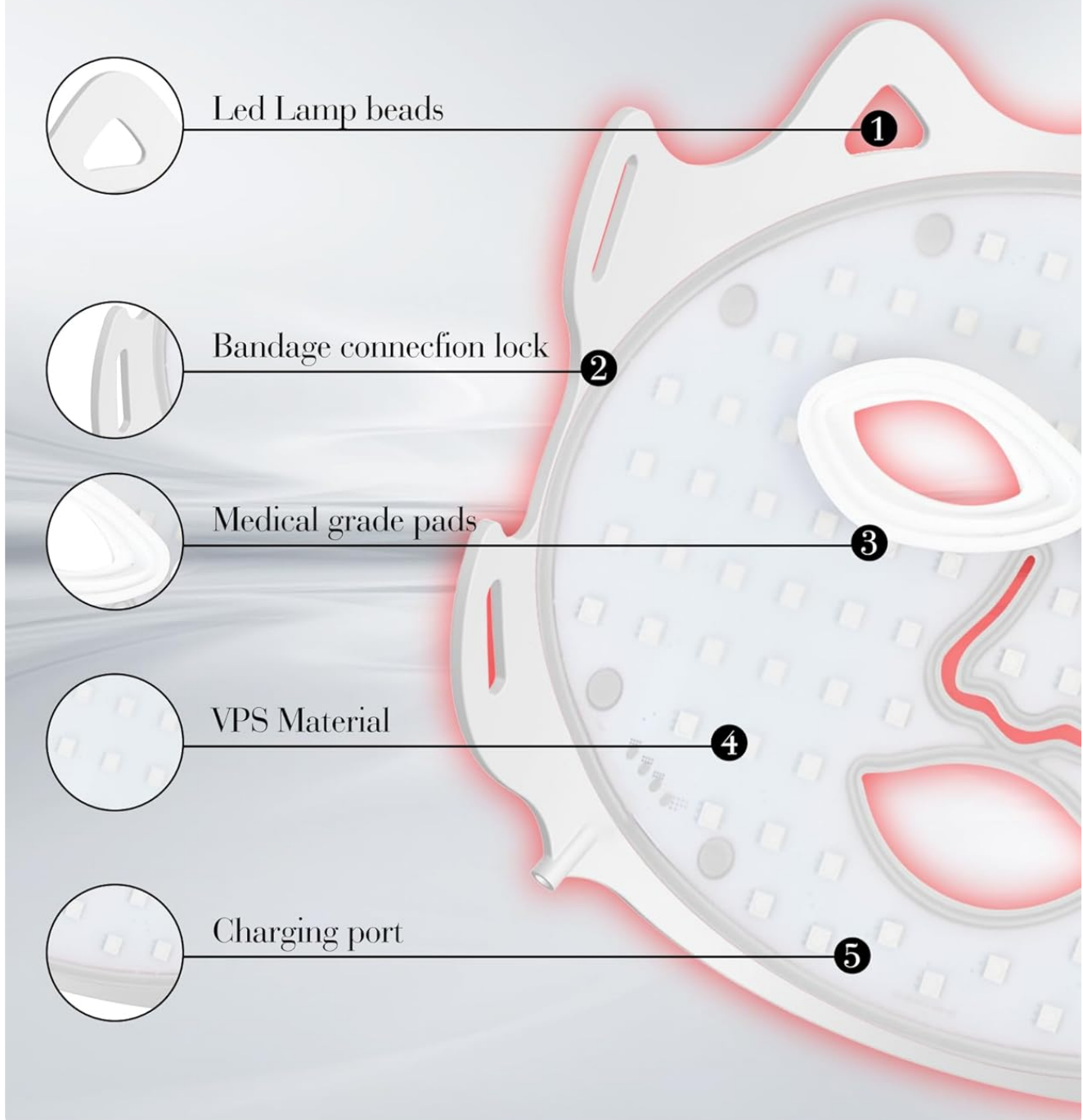


Image: An illustration detailing the 7 distinct LED light colors (Red, Blue, Green, Yellow, Cyan, Purple, White) and their corresponding wavelengths, indicating their specific applications.

6. MAINTENANCE

Proper maintenance ensures the longevity and effectiveness of your SDKWDH LED Face Mask:

- **Cleaning:** After each use, gently wipe the mask with a soft, damp cloth. Do not use harsh chemicals or abrasive cleaners. Ensure the mask is completely dry before storing.

- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** For optimal battery life, charge the remote control regularly and avoid completely draining the battery frequently.

7. TROUBLESHOOTING

If you encounter any issues with your SDKWDH LED Face Mask, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Mask does not turn on.	Remote control battery is low or not charged. Remote control not properly connected to mask.	Charge the remote control fully. Ensure the remote control is securely attached to the mask.
LED lights are dim or flickering.	Low battery on remote control. Incorrect brightness setting.	Recharge the remote control. Adjust the brightness setting using the 'E+' button on the remote.
Cannot change LED colors.	Incorrect button pressed.	Use the 'M+' or 'M-' (Color Cycle Button) on the remote control to cycle through colors.
Mask feels uncomfortable.	Strap is too tight or mask is not positioned correctly.	Adjust the adjustable strap for a looser fit. Reposition the mask on your face to ensure the eye pads align properly.

If the problem persists after trying these solutions, please contact customer support.

8. SPECIFICATIONS

Key specifications for the SDKWDH 7-Color LED Face Mask:

Feature	Detail
Model	White A1
Brand	SDKWDH
Number of LED Colors	7 (Red, Blue, Green, Yellow, Cyan, Purple, White)
Power Source	1 Lithium Ion battery (for remote control)
Input Voltage	DC5V-1A~2A (for remote control charging)
Dimensions (Package)	13.07 x 8.94 x 2.13 inches
Weight (Package)	1.48 Pounds
Material	Special silicone materials
Special Feature	Clinically Tested
UPC	761483199987

9. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact the seller directly. The manufacturer is Sunpint international beauty technology (guangzhou) co., LTD.

Please retain your proof of purchase for any warranty claims.