

AUVON TX6224

AUVON Wireless 24 Modes TENS Unit Muscle Stimulator (Model TX6224) User Manual

Your guide to effective pain relief and muscle stimulation.

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your AUVON Wireless 24 Modes TENS Unit Muscle Stimulator, Model TX6224. Please read this manual thoroughly before first use and retain it for future reference.

The AUVON Wireless TENS Unit is designed to provide targeted pain relief and muscle stimulation using Transcutaneous Electrical Nerve Stimulation (TENS) and Electrical Muscle Stimulation (EMS) technologies. Its compact, wireless design allows for convenient use at home, in the office, or while traveling.

2. IMPORTANT SAFETY INFORMATION AND LEGAL DISCLAIMER

WARNING: This product is not suitable for individuals with certain medical conditions or circumstances. Please read the following points carefully before using this product.

- Do not use this product if you have an implanted defibrillator, cardiac pacemaker, spinal cord stimulators or any other implanted metallic or electronic device.
- Do not use this product if you have heart rhythm problems, epilepsy, brain diseases, or neurological disorders.
- Do not use this product if you have undiagnosed chronic pain.
- Do not use this product if you are pregnant or in labor.
- Do not use this product if you have hypertension (high blood pressure) or thrombosis (blood clots).
- Do not use this product if you have cancer.
- Do not use this product if you are sensitive to electrotherapy.
- Do not use this product excessively.
- Do not use this product above the neck, especially on the head or neck area.
- Do not use this product if you are under 22 years of age.

Individual symptoms, situations, circumstances, and results may vary. The information provided is not intended to be used for medical diagnosis or treatment or as a substitute for professional medical advice.

Please consult your physician or qualified health care provider regarding your condition and appropriate medical treatment.

3. WHAT'S IN THE BOX

Upon unpacking your AUVON Wireless TENS Unit, please ensure all components are present:

- 1 Wireless Charging Case
- 2 Wireless Pods
- 4 Small Electrode Pads (2.4" x 1.7")
- 2 Large Electrode Pads (2.2" x 5.1")
- 1 User Manual (this document)
- 1 Quick Guide
- 1 Type-C Charging Cable
- 1 Pads Use Tips Card
- 1 Pouch
- 1 Gift Box



Image: All components included in the AUVON Wireless TENS Unit package, including the charging case, pods, various

4. SETUP

4.1 Charging the Device

1. Ensure the charging case has sufficient power. The case contains a 790 mAh rechargeable battery.
2. Place the wireless pods into the charging case. They will begin charging automatically.
3. To charge the case, connect the included Type-C cable to the charging port on the case and plug it into a wall adapter, power bank, or laptop.

A 10-minute charge in the case provides approximately 1 hour of therapy. A full charge of the case can power a single pod for up to 40 hours or two pods for 20 hours. Each pod runs up to 8 hours on a full charge.



Image: Diagram illustrating the charging process and battery life of the AUVON Wireless TENS Unit, showing 8 hours per pod and 40 hours total with the charging case.

4.2 Attaching Electrode Pads

1. Select the appropriate size of electrode pad for the target area (small for ankles, waist; large for shoulders, back).

2. Ensure your skin is clean and dry before applying the pads.
3. Remove the protective film from the gel side of the electrode pad.
4. Snap the wireless pod onto the electrode pad using the magnetic snap button.
5. Firmly press the electrode pad onto the desired area of your skin.

Medical-Grade & Latex-Free Pads



Multiple sizes of pads with industry-leading adhesive, individually packaged for safe, reusable use.

Soft PU Leather


Conductive Layer

Insulating Layer


Gel Layer

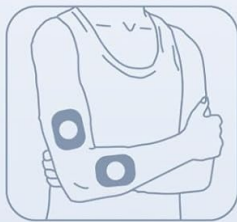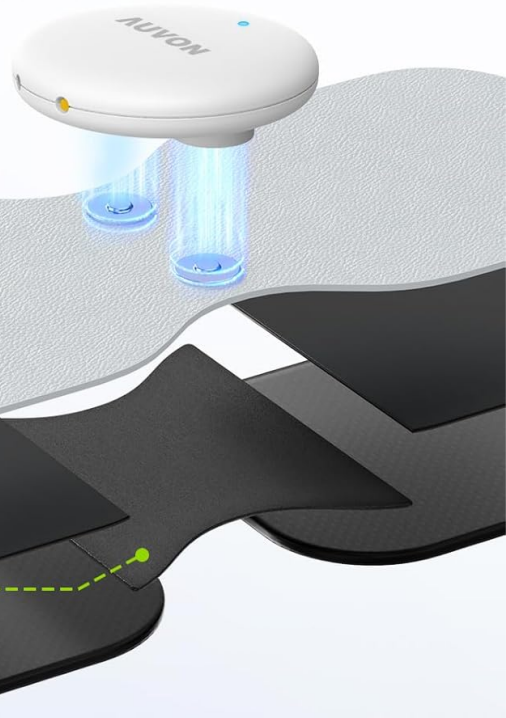

2.4" x 1.7" Electrodes
for arms, legs, etc.

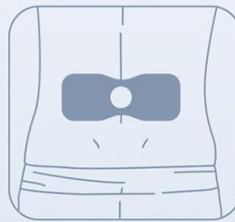

2.2" x 5.1" Electrodes
for back, shoulder, etc.

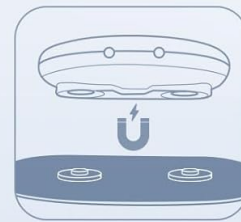

Magnetic Snap Button

Image: An exploded view of the medical-grade, latex-free electrode pad layers, showing the PU leather, conductive, insulating, and gel layers. Also depicted are the two sizes of electrode pads (2.4" x 1.7" and 2.2" x 5.1") and the magnetic snap button connection.

4.3 Initial Connection

The device follows Bluetooth headset settings and comes pre-paired from the factory. No manual pairing is needed under normal circumstances. If a Bluetooth disconnection occurs, place the pods back in the charging case for about 30 seconds, then remove them to automatically reconnect. If issues persist, contact customer service for assistance.

5. OPERATING INSTRUCTIONS

5.1 Basic Operation Flow

1. **Step 1:** Take out the Wireless Pods from the charging case. The device will power on automatically.
2. **Step 2:** Press the 'P' button on the charging case to select a program (mode). There are 24 available

modes.

3. **Step 3:** Press the '+' or '-' buttons to adjust the current intensity level. There are 30 adjustable intensity levels.
4. **Step 4:** Press the clock button to set the timer. Session times can be set from 10 to 90 minutes (default is 30 minutes).

Tips:

- There is no need to manually switch the charging case on or off.
- Press and hold the '-' button to quickly return the current intensity level to zero.



Image: A diagram showing the charging case display and buttons, including the time indicator, program indicator, current level indicator, timing button, program adjustment button, and intensity adjustment buttons.

5.2 Pad Placement and Usage

The AUVON TENS Unit can be used on various body parts for targeted relief. Refer to the quick guide for specific pad placement recommendations for different areas like shoulders, back, neck, arms, legs, abdomen, and hips.



Relieve Your Pain

A safe, non-invasive and drug-free method for pain relief and pain management.



Image: Examples of the AUVON Wireless TENS Unit being applied to different body areas, including the neck, arm, lower back, and leg, demonstrating its versatility for pain relief.

5.3 Before, During, and After Use

- **Before Use:** Ensure the charging case and pods have sufficient power. Check that the pods are securely attached to the pads and the pads are properly adhered to your skin.
- **During Use:** You can use one pod or both pods simultaneously. The pods will share the same mode and intensity. Normal connection status is indicated by the pod's Bluetooth light either flashing or staying steady.
- **After Use:** Return the pods to the charging case to recharge when the battery is low.

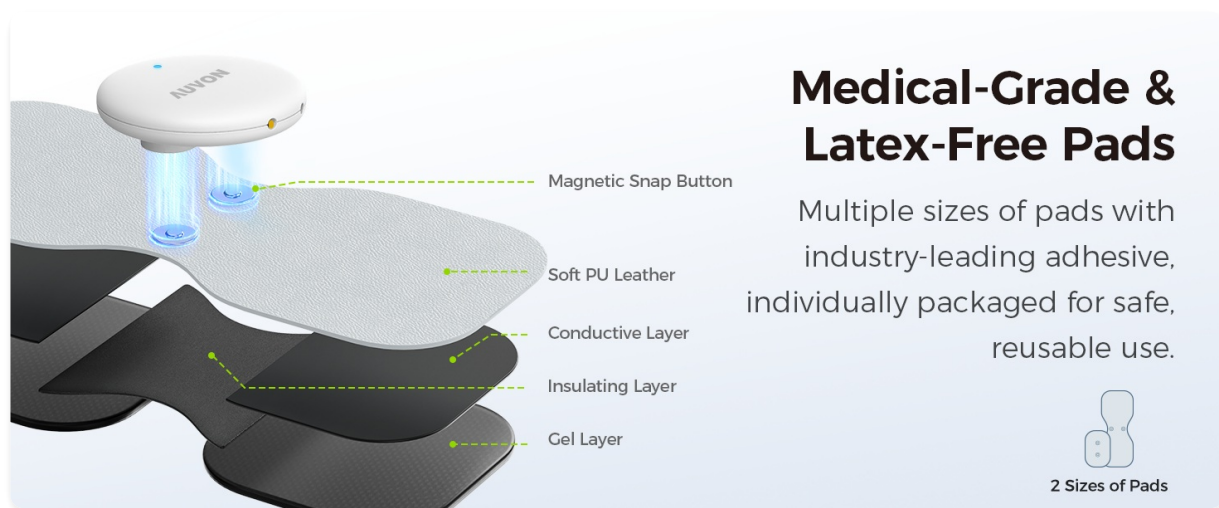


Image: A visual guide detailing the steps for using the AUVON TENS Unit: ensuring power and proper pad adhesion before use, simultaneous use of pods during therapy, and returning pods to the case for recharging after use.

6. MAINTENANCE

6.1 Electrode Pad Care

- After each use, place the protective film back onto the gel side of the electrode pads to prevent drying and maintain adhesion.
- If the pads lose their stickiness, gently wipe the gel surface with a damp cloth (water only) and let them air dry. Do not use soap or alcohol.
- Store pads in a cool, dry place, away from direct sunlight.
- Replace pads when they no longer adhere well to the skin or if they show signs of wear and tear. High-quality pads can last up to 45 uses with proper care.

6.2 Device Cleaning

- Wipe the charging case and pods with a soft, dry cloth.
- Do not use abrasive cleaners or immerse the device in water.

7. TROUBLESHOOTING

- **Device not turning on:** Ensure the pods are properly seated in the charging case and the case itself is charged. If the pods are out of the case, they should turn on automatically.
- **No sensation or weak sensation:**
 - Increase the intensity level using the '+' button.
 - Check that the electrode pads are firmly attached to your skin and are not dried out.
 - Ensure the pods are securely snapped onto the pads.
 - Verify the battery level of the pods and charging case.
- **Pads not sticking:**
 - Clean your skin thoroughly before application.
 - Clean the gel surface of the pads as described in the Maintenance section.
 - Replace old or worn-out pads with new ones.
- **Bluetooth disconnection:** Place the pods back into the charging case for approximately 30






seconds, then remove them. They should automatically reconnect. If the issue persists, contact AUVON customer support.





8. SPECIFICATIONS

Feature	Specification
Model Number	TX6224
Modes	24
Intensity Levels	30
Session Timer	10-90 minutes (default 30 min)
Pod Battery Life	Up to 8 hours per pod
Charging Case Battery	790 mAh (up to 40 hours total use)
Charging Time (10 min)	1 hour of therapy
Charging Port	Type-C
Electrode Pad Sizes	2.4" x 1.7" (small), 2.2" x 5.1" (large)
Pad Material	Medical-grade, latex-free
Package Dimensions	7.48 x 5.12 x 2.21 inches
Item Weight	12 ounces
Batteries Required	1 Lithium Ion (included)

9. WARRANTY AND SUPPORT

AUVON provides a **24-month quality warranty** for this product. Additionally, lifetime customer support is available to assist you with any questions or issues you may encounter. For support, please refer to the contact information provided in your Quick Guide or visit the official AUVON website.

<div></div> <div><p>USER MANUAL</p><p>AUVON® TENS Massager</p><p>Relaxation, Anywhere, Anytime</p><p>Model: SM9126</p></div> <div><p>Device Name: TENS & TMS</p></div>	<p>AUVON SM9126 TENS Massager User Manual: Features, Operation, and Troubleshooting</p> <p>Comprehensive user manual for the AUVON SM9126 TENS Massager. Learn about its features, operating instructions, safety precautions, troubleshooting tips, and maintenance for effective pain relief and muscle stimulation.</p>
<div></div> <div><p>USER MANUAL</p><p>Dual Channel TENS Massager</p><p>Relaxation, Anywhere, Anytime</p><p>Model: AS8012</p></div> <div><p>Device Name: TENS and Powerful Muscle Stimulator</p></div>	<p>AUVON AS8012 Dual Channel TENS Massager User Manual</p> <p>Comprehensive user manual for the AUVON AS8012 Dual Channel TENS Massager, detailing features, operating instructions, safety precautions, and troubleshooting for effective pain relief and muscle stimulation.</p>
<div><p>Quick Guide</p><p>Attention</p><p>The purpose of this Quick Guide is to help you understand the operation of the TENS unit much more quickly, and the contents of the quick guide are not a substitute for the user manual. Please read the user manual carefully before use to understand the contraindications and warnings. Please consult your doctor before use if you're not sure whether you can use this device.</p><p>Operation</p><div></div><p>1. Clean the device with water.</p><p>2. Connect the lead wire to a power source and the device.</p></div>	<p>AUVON AS8019 TENS and EMS Muscle Stimulator User Manual & Quick Guide</p> <p>Comprehensive user manual and quick guide for the AUVON AS8019 TENS and EMS muscle stimulator, covering operation, pad placement, features, troubleshooting, and safety.</p>
<div></div> <div><p>MANUAL DE USUARIO</p><p>Masajeador TENS</p><p>Relajación en cualquier lugar en cualquier momento</p><p>Modelo: AS8016</p></div> <div><p>Nombre del dispositivo: TENS y Estimulador Muscular Automático</p></div>	<p>AUVON AS8016 TENS Massager User Manual - Pain Relief and Muscle Stimulation</p> <p>Comprehensive user manual for the AUVON AS8016 TENS massager. Learn about features, operation, safety precautions, troubleshooting, and maintenance for effective pain relief and muscle stimulation.</p>

<div><div>Quick Guide</div><div>Attention</div><p>The purpose of this Quick Guide is to help you understand the operation of the TENS unit much more quickly, and the contents of the quick guide are not a substitute for the user manual. Please read the user manual carefully before use to understand the contraindications and warnings. Please consult your doctor before use if you're not sure whether you can use this device.</p><div>Operation</div><p>Note: Do not use it while charging.</p><div><div><p>1. Clean and dry the application area.</p></div><div><p>2. Connect the lead wire to a pair of pads and the device.</p></div></div></div>	<p>AUVON TENS & EMS Unit Quick Guide and User Manual - Model AS8016</p> <p>Comprehensive quick guide and user manual for the AUVON TENS and Powered Muscle Stimulator Model AS8016. Learn how to operate the device, understand TENS and EMS therapies, proper pad placement for various body areas, battery charging, troubleshooting, and safety precautions for effective pain relief and muscle stimulation.</p>
<div><div>Quick Guide</div><div>Attention</div><p>The purpose of this Quick Guide is to help you understand the operation of the TENS unit much more quickly, and the contents of the quick guide are not a substitute for the user manual. Please read the user manual carefully before use to understand the contraindications and warnings. Please consult your doctor before use if you're not sure whether you can use this device.</p><div>Operation</div><p>Note: Do not use it while charging.</p><div><div><p>1. Clean and dry the application area.</p></div><div><p>2. Connect the lead wire to a pair of pads and the device.</p></div></div></div>	<p>AUVON TU2224-A TENS and Muscle Stimulator User Manual - Pain Relief Guide</p> <p>Comprehensive user manual for the AUVON TU2224-A TENS and Powered Muscle Stimulator. Learn how to use, troubleshoot, and care for your device for effective pain relief and muscle stimulation.</p>