



Manuals.plus /

- › Philips Sonicare /
- › Philips Sonicare S2 Sensitive Replacement Toothbrush Heads HX6056/77 User Manual

Philips Sonicare HX6056/77

Philips Sonicare S2 Sensitive Replacement Toothbrush Heads User Manual

Model: HX6056/77

INTRODUCTION TO PHILIPS SONICARE S2 SENSITIVE BRUSH HEADS

The Philips Sonicare S2 Sensitive replacement toothbrush heads are designed for individuals with sensitive teeth and gums. These genuine Philips Sonicare brush heads provide a gentle yet effective cleaning experience, removing plaque efficiently while being kind to oral tissues.

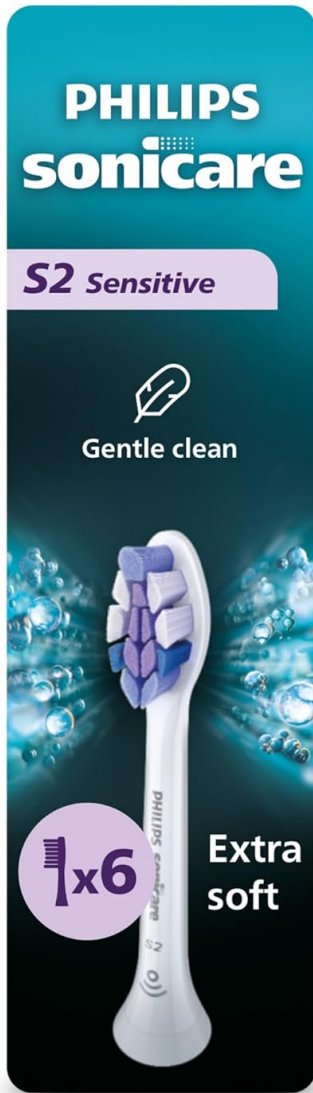


Image: Philips Sonicare S2 Sensitive replacement brush heads, showing a pack of six and individual heads.

Each S2 Sensitive brush head features long, thin, ultra-soft bristles, specifically engineered for a gentle touch. With over 3,000 densely packed bristles, these heads are capable of removing significantly more plaque than a manual toothbrush, even in areas that are typically difficult to reach.



Gentle yet effective cleaning of sensitive teeth with ultra-soft bristles

Image: A magnified view of the ultra-soft bristles on a Philips Sonicare S2 Sensitive brush head, highlighting their gentle design.

INSTALLATION AND REPLACEMENT

Philips Sonicare S2 Sensitive brush heads are designed for easy attachment and removal. They are compatible with all Philips Sonicare click-on electric toothbrush handles, excluding Philips One and e-series handles.

Attaching a New Brush Head:

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks into place. Ensure it is securely attached to prevent detachment during use.

Detaching a Brush Head:

1. Grip the brush head firmly and pull it straight up and off the handle.
2. Always handle the brush head and handle with care during attachment and detachment to avoid injury.



Image: Illustration demonstrating the click-on mechanism for attaching a Philips Sonicare brush head to the handle.

OPERATING AND USAGE

The S2 Sensitive brush heads utilize proven Sonicare technology to deliver a thorough and gentle clean. The sonic vibrations generate up to 62,000 bristle movements per minute, creating a dynamic fluid action that extends deep between teeth and along the gumline.

Brushing Technique:

- Wet the bristles and apply a small amount of toothpaste.
- Place the brush head against your teeth at a slight angle towards the gumline.
- Turn on the toothbrush. Gently move the brush head slowly across each tooth surface, allowing the sonic action to do the work. Avoid scrubbing.
- Brush all four quadrants of your mouth evenly.

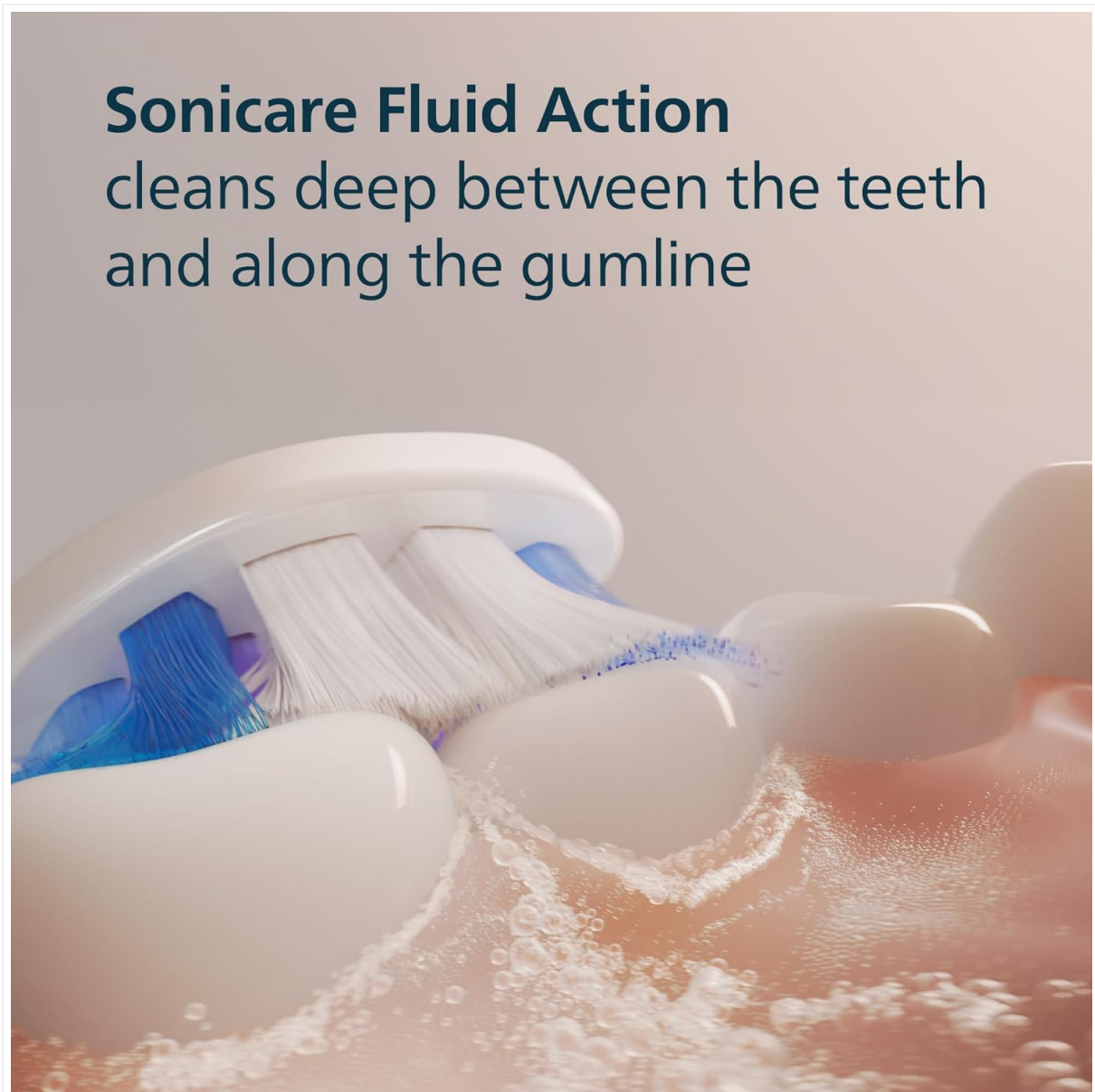


Image: A visual representation of Sonicare Fluid Action, showing water and toothpaste bubbles reaching deep between teeth and along the gumline for effective cleaning.

MAINTENANCE AND REPLACEMENT SCHEDULE

For optimal oral hygiene and performance, dental professionals recommend replacing your toothbrush head every 3 months. Over time, bristles can fray and become less effective at removing plaque.

BrushSync Replacement Reminder:

Many Philips Sonicare toothbrushes feature BrushSync technology. This system tracks how often and how hard you brush. When it's time to replace your brush head, a light on the toothbrush handle will illuminate. Additionally, the blue indicator bristles on the S2 Sensitive brush head will fade and turn white, providing a visual cue for replacement.

All genuine Philips Sonicare
brush heads fit all Philips
Sonicare click-on handles⁵



Image: A Philips Sonicare S2 Sensitive brush head showing the blue indicator bristles fading to white, signaling the need for replacement.

COMPATIBILITY

The Philips Sonicare S2 Sensitive replacement brush heads (Model HX6056/77) are designed to be compatible with all Philips Sonicare click-on electric toothbrush handles.

Note: These brush heads are not compatible with Philips One and e-series handles.



Image: An assortment of Philips Sonicare brush heads, including the S2 Sensitive, shown alongside a Philips Sonicare handle to illustrate universal click-on compatibility.

SPECIFICATIONS

Model Number	HX6056/77
Brand	Philips Sonicare
Color	White
Material	Nylon (bristles)
Age Range	Adult
UPC	075020114976
Package Dimensions	9.06 x 2.32 x 1.5 inches
Item Weight	2.47 ounces

SAFETY INFORMATION

- Always attach and detach the brush head from the electric toothbrush handle with care to avoid injury.
- Keep this product out of reach of young children to prevent choking hazards.

ENVIRONMENTAL INFORMATION

Philips is committed to reducing environmental impact. The S2 Sensitive sonic toothbrush replacement heads are made with 70% bio-based plastic (allocated to plastic of brush head on mass balance basis). The packaging is paper-based and recyclable where facilities exist.



Image: Graphic illustrating the use of 70% bio-based plastic in the brush head and the recyclable paper-based packaging, with an ISCC certified logo.