

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [VAAGHANM](#) /
- › [VAAGHANM VA-X-2 Electric Heating Pad Instruction Manual](#)

VAAGHANM VA-X-2

VAAGHANM VA-X-2 Electric Heating Pad Instruction Manual

Model: VA-X-2 | Brand: VAAGHANM

1. PRODUCT OVERVIEW

The VAAGHANM VA-X-2 Electric Heating Pad is designed to provide targeted heat therapy for the neck, shoulders, and back. Its weighted design ensures a snug fit, delivering consistent and soothing warmth. Crafted from soft flannel, it offers comfort and effective heat distribution.

Find Your Comfort Zone, with Weighted Comfort

Customized Contour for Perfect Support-Better fit for your neck & shoulder

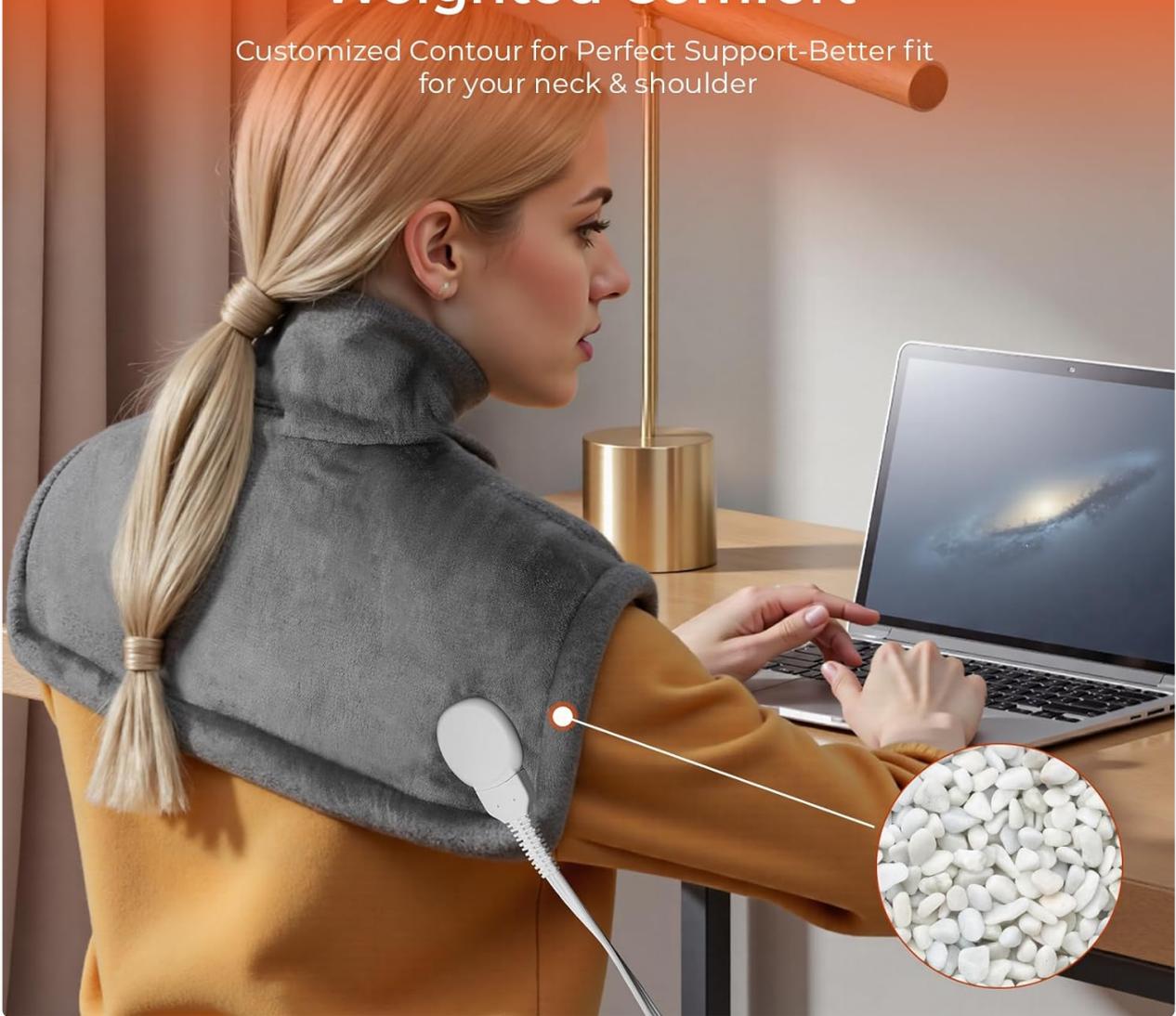


Image 1.1: VAAGHANM Electric Heating Pad in use, conforming to neck and shoulders.

2. SAFETY INFORMATION


WARNING: To avoid danger of suffocation, keep this bag away from babies and children. Do not use in cribs, beds, carriages, or playpens. This bag is not a toy. Thin film may cling to nose and mouth and prevent breathing.

- Always use on low heat settings for extended periods.
- Avoid direct contact with skin if sensitive to heat.
- Keep away from children and pets.
- Do not use while sleeping or bathing.
- Unplug and allow to cool before storing.
- Do not immerse the controller in water.

3. PRODUCT FEATURES

- **Weighted Design:** Evenly filled with micro-glass beads for gentle, consistent heat and a soft, conforming embrace.

- **Versatile Application:** Designed to naturally contour to your neck, shoulders, back, or waist.
- **Adjustable Heat Settings:** Features 6 adjustable heat levels (100°F to 150°F) to meet diverse needs safely.
- **Auto-Off Timer:** Includes 6 auto-off timers (0.5H to 3H) for personalized sessions, conserving energy and preventing overheating.
- **User-Friendly Controller:** Clear LED display with backlight for easy monitoring, ideal for night use.
- **Soft Flannel Fabric:** Crafted from luxurious soft flannel for even heat distribution and lasting coziness.
- **Machine-Washable:** Easy to care for (controller must be removed; not tumble-dried).

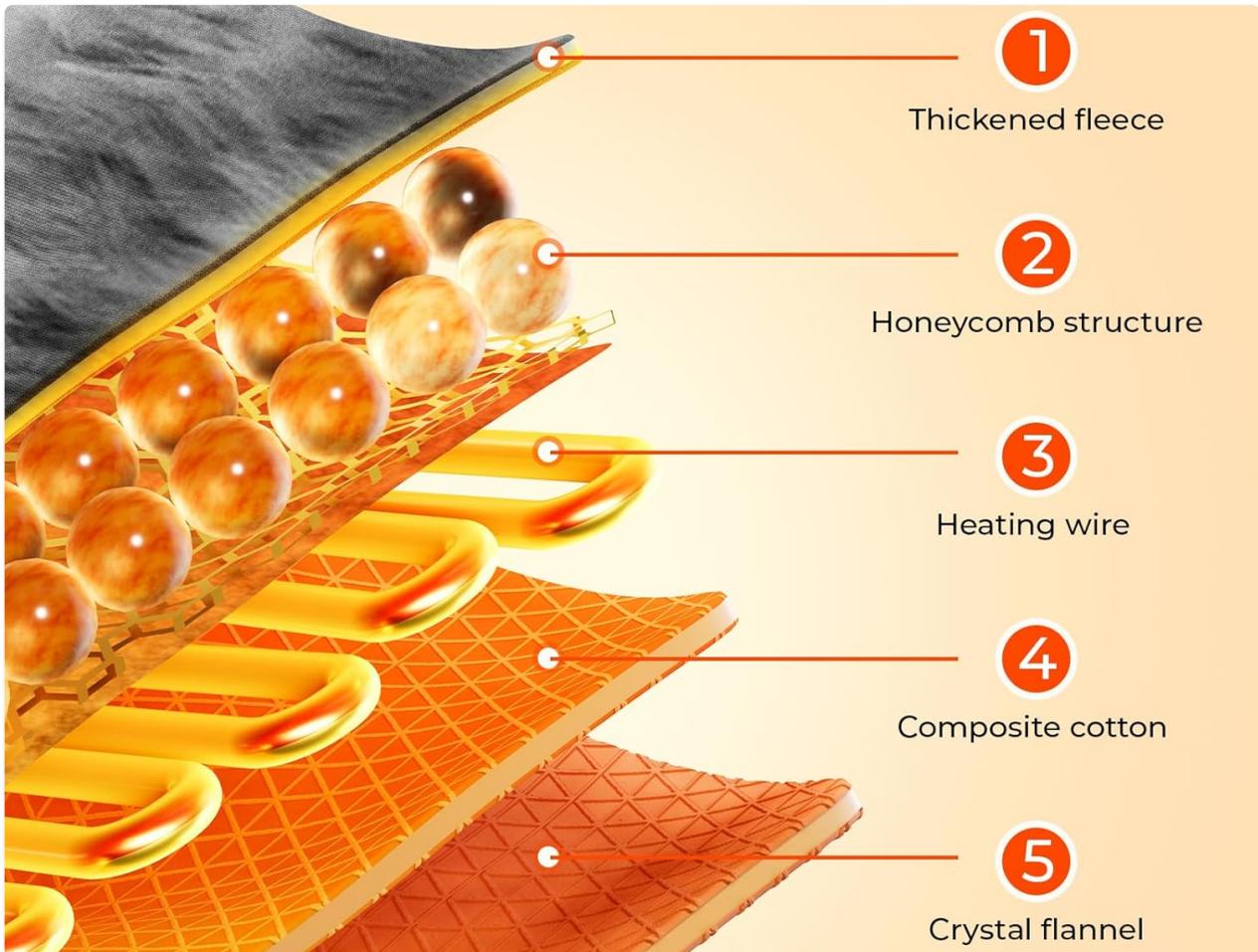

Safety-certified heating parts & controllers

Innovative technology for more safety and warmth

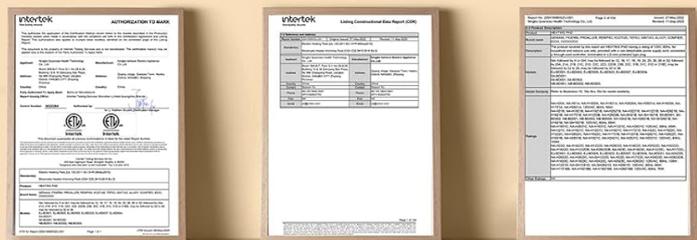


Image 3.1: Internal structure of the heating pad, highlighting material layers.

4. SETUP

1. Unpack the heating pad and controller from its packaging.
2. Ensure the controller is securely connected to the heating pad.
3. Place the heating pad on the desired area of your body (neck, shoulders, back, or waist).
4. Plug the power cord into a standard 120V~60Hz electrical outlet.

NECK & SHOULDER HEATING PAD

Customized Contour for Perfect Support-Better fit
for your neck & shoulder



Image 4.1: Proper placement of the heating pad for relaxation.

5. OPERATING INSTRUCTIONS

1. **Power On:** Press the power button on the controller to turn on the heating pad. The LED display will illuminate.
2. **Adjust Temperature:** Press the temperature button (indicated by a thermometer icon) to cycle through the 6 heat settings (100°F, 110°F, 120°F, 130°F, 140°F, 150°F). We recommend starting at a lower setting, as heat tolerance is personal.
3. **Set Timer:** Press the timer button (indicated by a clock icon) to select an auto-off duration from 0.5 hours to 3 hours (0.5H, 1.0H, 1.5H, 2.0H, 2.5H, 3.0H).
4. **Power Off:** Press the power button again to turn off the heating pad. Unplug from the outlet when not in use.

MORE TEMPERATURE & TIMER CHOICES



Image 5.1: Controller with temperature and timer settings.

6. MAINTENANCE

- **Cleaning:** The soft flannel heating pad is machine-washable. Ensure the controller is detached before washing. Do not tumble-dry.
- **Storage:** After use, unplug the heating pad and allow it to cool completely. Store it in a dry place.



Image 6.1: Detaching the controller for machine washing.

7. TROUBLESHOOTING

- **No Heat:** Ensure the power cord is securely plugged into a working outlet and the controller is firmly connected to the pad. Check if the power button has been pressed.
- **Inconsistent Heat:** Verify that the heat setting is appropriate for your preference. Ensure the pad is laid flat and not bunched up, which can affect heat distribution.
- **Controller Display Not Working:** Unplug the unit from the wall, wait a few minutes, and plug it back in. If the issue persists, contact customer support.

8. SPECIFICATIONS

Feature	Detail
Model Name	Heating Pad
Item Model Number	VA-X-2

Brand	VAAGHANM
Material	Polyester
Input	120V~60Hz
Power	95W
Temperature Settings	6 levels (100°F, 110°F, 120°F, 130°F, 140°F, 150°F)
Timer Modes	6 modes (0.5H, 1.0H, 1.5H, 2.0H, 2.5H, 3.0H)
Special Features	Automatic Shut Off, Overheat Protection, Portable, Timer
Package Dimensions	10.75 x 7.48 x 3.94 inches
Weight	2.16 Pounds
Neckline Circumference	22 Inches



Image 8.1: Heating pad dimensions and neckline circumference.

9. WARRANTY AND SUPPORT

VAAGHANM offers a 30-day return policy and 12-month limited warranty. Lifetime professional tech support is also provided.

For customer care, product inquiry, warranty claims, or replacement accessories, please contact us:

- **Email:** support@vaaghanm.net

You can also scan the QR code on the customer care card included with your product for more information.