

COSLUS C30

COSLUS C30 Sonic Electric Toothbrush User Manual

Model: C30 | Brand: COSLUS

1. INTRODUCTION

The COSLUS C30 Sonic Electric Toothbrush is designed to provide a superior and gentle cleaning experience for adults, teens, youth, kids, and seniors. Featuring advanced sonic technology, multiple cleaning modes, and intensity levels, this toothbrush offers a customizable solution for various oral care needs, including deep cleaning, whitening, and gum care. Its durable design and long-lasting battery make it ideal for both home and travel use.

2. PACKAGE CONTENTS

Upon opening your COSLUS C30 Sonic Electric Toothbrush package, you should find the following items:

- 1 × COSLUS C30 Electric Toothbrush Handle
- 1 × USB Charging Cable (adapter not included)
- 2 × U-shaped Soft Brush Heads with protective covers
- 2 × W-shaped Medium Brush Heads with protective covers

Real Value is The Lasting Motor Not The Head Count

3 YEARS
Metal motor, last longer

COSLUS Long-Term Partner		Others Wasted!	
Motor Lifespan	✔ >3-year	✘ <6 months	
Brush Heads	✔ 4 brush heads, enough for initial use.	✘ A trick of many heads.	
Result	✔ Lasting motor. Saves & easy. Just replace heads.	✘ Failed motor = wasted heads + new order.	
Safety	✔ Stable & safe, teeth protection	✘ Harsh vibration, harm teeth	
Noise	✔ <60dB quiet	✘ >80dB noisy	

*Based on the 2025 market research and analysis by the COSLUS electric toothbrush R&D team.

Image: The COSLUS C30 Electric Toothbrush handle, USB charging cable, and four individually packaged brush heads.

3. SETUP

- Charge the Toothbrush:** Before first use, fully charge the toothbrush for approximately 2 hours. Connect the USB charging cable to the charging port at the bottom of the handle and plug the USB end into a compatible 5V USB adapter (not included). The intensity indicator lights will illuminate during charging.
- Attach a Brush Head:** Select your desired brush head (U-shaped soft or W-shaped medium). Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.



Image: A quick guide illustrating the charging process and how to securely attach a brush head to the toothbrush handle.

4. OPERATING INSTRUCTIONS

Power On/Off and Mode Selection

1. **Power On:** Press the Power/Mode button (top button) to turn on the toothbrush. It will start in the last used mode and intensity.
2. **Change Mode:** Within 5 seconds of powering on, repeatedly press the Power/Mode button to cycle through the 5 cleaning modes: Clean, Soft, White, Massage, Polish. The corresponding indicator light will illuminate.
3. **Change Intensity:** Press the Intensity button (bottom button) to switch between 3 intensity levels (low, medium, high) for the selected mode. The intensity indicator lights will show the current level.
4. **Power Off:** To turn off the toothbrush, press and hold the Power/Mode button for more than 3 seconds, or let it complete its 2-minute cycle.

Cleaning Modes Explained

- **Clean:** Standard mode for superior daily teeth cleaning.
- **Soft:** Gentle care, suitable for sensitive teeth and gums.
- **White:** Teeth cleaning with polishing and whitening action.
- **Massage:** Pulsates for gentle gum stimulation.
- **Polish:** Polishes teeth to bring natural brightness.

Brushing Technique

Apply light pressure and move the brush head gently and slowly across the teeth in a small back-and-forth motion. To maximize cleaning effectiveness, divide your mouth into 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause briefly at 30-second intervals, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.

Advanced Sonic Technology Targets Hidden Stains Without the Damage

A more thorough clean contributes to
your daily oral care routine



Image: The COSLUS C30 toothbrush handle highlighting the 5 cleaning modes and 3 adjustable intensity levels.

Your browser does not support the video tag.

Video: An official COSLUS video demonstrating the C30 Electric Toothbrush's gum-friendly deep cleaning capabilities and various features.

5. MAINTENANCE

1. **Cleaning the Brush Head:** After each use, rinse the brush head thoroughly under running water.
2. **Cleaning the Handle:** Wipe the toothbrush handle with a damp cloth. Do not submerge the handle in water, as this may damage internal components. The COSLUS C30 is IPX7 waterproof, meaning it can withstand immersion in water up to 1 meter for 30 minutes, making it safe for shower use, but prolonged submersion is not recommended.
3. **Replacing Brush Heads:** For optimal hygiene and effectiveness, replace brush heads every 3 months, or sooner if bristles become frayed.



Image: A man brushing his teeth in the shower, demonstrating the IPX7 waterproof rating of the COSLUS C30 toothbrush.


6. CHARGING

The COSLUS C30 features fast charging, providing up to 30 days of battery life from a 2-hour charge (based on soft mode usage). The intensity indicator lights will show the charging status. Ensure the charging port is dry before connecting the cable.

1 Brush Head Outlasts 3 Ordinary Ones

Quality, not quantity

Other

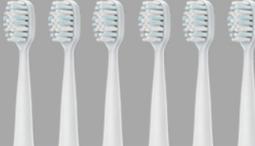

VS

COSLUS



 Ordinary Nylon, <1 Month Lifespan	 Low Density, Lower Cleaning Efficiency	 Enhanced Flexible Bristles, 3+ Months Durable	 20% Higher Bristle Density
 Conventional Plan-ting, Easy Shedding	 Rough Tips, Risk of Damage	 Advanced Implant Tech, <1% Fall-Out Rate	 End-Rounded, Gentle Safe Design

Image: Visual representation of the COSLUS C30's fast charging capability, achieving 30 days of battery life from a 2-hour charge.

7. TROUBLESHOOTING

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged. If the issue persists, contact customer support.
- **Weak vibration:** Check the battery level and recharge if necessary. Ensure the brush head is securely attached.
- **Pressure Sensor Lamp:** If you brush your teeth too hard, the pressure sensor lamp will light up, indicating you should reduce pressure to protect your gums and enamel.

8. SPECIFICATIONS

Feature	Detail
Model Number	C30
Dimensions	7.16 x 1.06 x 1.06 inches; 6.38 ounces
Power Source	Battery Powered (1 Lithium Ion battery included)
Waterproof Rating	IPX7
Cleaning Modes	5 (Clean, Soft, White, Massage, Polish)

Intensity Levels	3
Brush Head Types	U-shaped Soft, W-shaped Medium
Bristle Material	Nylon
Handle Material	Plastic
Smart Timer	2-minute auto-off, 30-second quadrant reminder

9. WARRANTY AND SUPPORT

The COSLUS C30 Sonic Electric Toothbrush comes with a **1-year manufacturer's warranty**. Additionally, COSLUS offers **24-hour after-sales service** and a **free 30-day refund/replacement** policy. For warranty activation or support, please refer to the warranty card included in your package or visit the official COSLUS website.

For customer service inquiries, you may contact COSLUS support. Please have your order number ready when contacting support.



Image: An illustration highlighting the durable magnetic levitation motor of the COSLUS C30, emphasizing its 3+ year lifespan.