

SKG YS100

SKG Foot Massager YS100 Instruction Manual

Brand: SKG | Model: YS100

INTRODUCTION

Thank you for choosing the SKG Foot Massager YS100. This electric foot massager is designed to provide comprehensive relief for your feet, targeting common discomforts such as plantar fasciitis, neuropathy, and general fatigue. Featuring advanced 360° air compression, deep shiatsu kneading, and optional heat therapy, the YS100 offers a personalized and relaxing massage experience. Please read this manual thoroughly before use to ensure safe and optimal operation.

SAFETY INFORMATION

- Do not stand on the device or place excessive weight on it.
- If you feel discomfort or pain, stop using immediately and consult your doctor.
- Use only with the provided power adapter (100–240V).
- Not suitable for children, pregnant women, or individuals with sensory impairment unless under supervision.
- Limit use to 15–30 minutes per session to avoid overuse.

PRODUCT FEATURES

- **Enhanced 360° Foot Massage:** Experience 20% stronger air compression and 30% deeper kneading with 12 rotating nodes on 4 plates, targeting arches, soles, and heels. Choose from 3 compression modes.

- **Personalized Massage Experience:** Customize your session with 3 compression levels and 3 shiatsu intensities, offering over 12 relaxing combinations.
- **Targeted Heat Therapy:** Optional 131°F heating gently warms arches, heels, and toes to soothe sore feet and improve circulation.
- **Dual Control Options:** Operate effortlessly with the included wireless remote control or the integrated touch panel.
- **Universal Fit:** Designed to comfortably accommodate feet up to US Men's size 13, suitable for both men and women.
- **Easy Maintenance:** Features removable, machine-washable liners for hygienic use. The sleek, supercar-inspired design is also easy to wipe clean.
- **Compact & Portable:** Weighing 5.51 lbs and measuring 6.97" tall, with a 72" cord, it is easy to carry and store.

UPGRADED +30% MASSAGE EXPERIENCE ANYWHERE



Image: The SKG Foot Massager YS100 shown in different environments such as an armchair, sofa, bedside, and office, highlighting its versatility for use anywhere.

SETUP

1. **Unpack:** Carefully remove the foot massager and all accessories from the packaging.
2. **Remote Activation:** Locate the remote control. Before first use, ensure the plastic tab on the remote control has been removed to activate the battery.
3. **Power Connection:** Connect the provided power adapter to the massager's power port and then plug it into a standard 100–240V electrical outlet. Use only the original adapter.
4. **Placement:** Place the foot massager on a flat, stable surface where you intend to use it.

OPERATING INSTRUCTIONS

The SKG Foot Massager YS100 can be operated using either the touch panel on the device or the wireless remote control.

Getting Started:

1. Sit comfortably and place your feet into the footwells of the massager.
2. Press the Power button on the touch panel or remote to turn on the device. The default settings will activate.

Adjusting Settings:

- **Message Modes (M):** Press the 'Mode' button to cycle through the 3 available shiatsu kneading massage modes.
- **Air Compression Intensity (+/-):** Press the 'Air Intensity' button to adjust between 3 levels of air compression (Low, Medium, High).
- **Heat Function (Heat):** Press the 'Heat' button to toggle the optional 131°F heat function on or off.
- **Timer (Time):** Press the 'Time' button to select a massage duration of 15 or 30 minutes. The device will automatically shut off after the selected time.

FITS UP TO SIZE 13

Perfect tool for nurses, teachers, athletes, and office workers to relax.

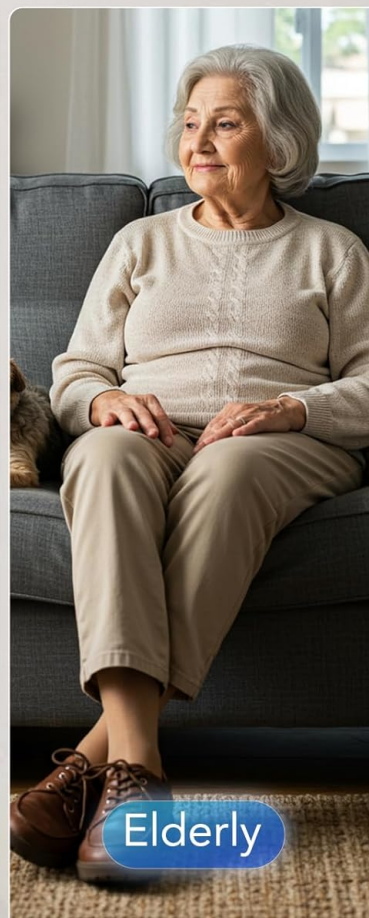


Image: Close-up of the SKG Foot Massager YS100's control panel and remote control, illustrating the various buttons for power, mode, air intensity, heat, and timer settings.

Remote Control Tips:

- Point the remote directly at the control panel of the massager for optimal signal reception.

Your browser does not support the video tag.

Video: A demonstration of how to change the remote control settings for the SKG Foot Massager YS100, showing the remote's functionality.

MAINTENANCE

- **Cleaning Liners:** The foot massager comes with removable, machine-washable liners. Unzip and remove the liners for cleaning to maintain hygiene.
- **Surface Cleaning:** The exterior of the device features a wipe-clean surface. Use a soft, damp cloth to clean the outer casing. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight when not in use. Its

compact design allows for convenient storage.

DUAL CONTROL WITH WIRELESS REMOTE

Personalize Your Foot Massage

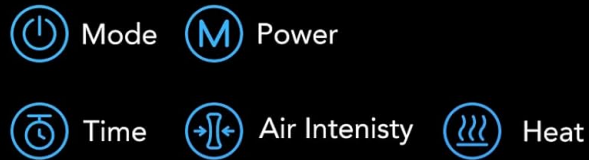


Image: The SKG Foot Massager YS100 highlighting its sleek design, removable liners for washing, and easy-to-clean surface, emphasizing its practical maintenance features.

TROUBLESHOOTING

- **Device Not Turning On:** Ensure the power adapter is securely plugged into both the massager and a working electrical outlet. Check if the remote's plastic tab has been removed.
- **Remote Control Unresponsive:** Verify the remote's battery is active (plastic tab removed). Ensure the remote is pointed directly at the massager's control panel.
- **Experiencing Discomfort:** Reduce the massage intensity or compression level. If discomfort persists, discontinue use and consult a healthcare professional.
- **Heat Function Not Perceptible:** The heat function provides a gentle warmth up to 131°F. Ensure the heat setting is activated. Individual perception of warmth may vary.

SPECIFICATIONS

Feature	Detail
Model Number	YS100
Product Dimensions	14.35 x 12.52 x 6.97 inches
Item Weight	5.51 Pounds
Power Source	Corded Electric
Material	Ethylene Vinyl Acetate, Polyurethane
Max Foot Size	US Men's Size 13

SUPPORT

For any questions, concerns, or assistance with your SKG Foot Massager YS100, please do not hesitate to contact our dedicated support team. We are always ready to help ensure your satisfaction and optimal product experience.